

AFRPL-TR-79-91

FOURIER TRANSFORM INFRARED SPECTROSCOPY

VOLUME II - USER'S MANUAL

PREPARED BY:

W. W. SCHWARZ  
O. C. SMITH  
THIOKOL CORPORATION/HUNTSVILLE DIVISION  
HUNTSVILLE, ALABAMA

DECEMBER 1979

Approved for Public Release, Distribution Unlimited

FINAL REPORT


PREPARED FOR:

AIR FORCE ROCKET PROPULSION LABORATORY  
DIRECTOR OF SCIENCE AND TECHNOLOGY  
AIR FORCE SYSTEMS COMMAND  
EDWARDS AIR FORCE BASE, CA 93523

DDC FILE COPY

ADA 079997

80 1 28 1980

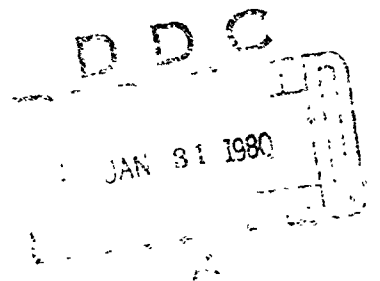

**Best  
Available  
Copy**

19 REPORT DOCUMENTATION PAGE		READ INSTRUCTIONS BEFORE COMPLETING FORM
1. REPORT NUMBER	2. GOVT ACCESSION NO.	3. RECIPIENT'S CATALOG NUMBER
13 AFRPL-TR-79-91-VOL-2		9
4. TITLE (and Subtitle)		5. TYPE OF REPORT & PERIOD COVERED
6 Fourier Transform Infrared Spectroscopy (FTIS). Final Report Volume II • User's Manual		Final Report • Apr 78 - Dec 79
7. AUTHOR(s)		8. PERFORMING ORG. REPORT NUMBER
W. W. Schwarz and O. C. Smith		14 U-79-18-VOL-2
9. PERFORMING ORGANIZATION NAME AND ADDRESS		10. PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS
Thiokol Corporation Huntsville Division Huntsville, AL 35807		16 JON 573013 LA 17 18
11. CONTROLLING OFFICE NAME AND ADDRESS		12. REPORT DATE
Air Force Rocket Propulsion Laboratory Edwards Air Force Base, CA 93523		11 December 1979
14. MONITORING AGENCY NAME & ADDRESS (if different from Controlling Office)		13. NUMBER OF PAGES
		222 12 220
		15. SECURITY CLASS. (of this report)
		Unclassified
		15a. DECLASSIFICATION/DOWNGRADING SCHEDULE
16. DISTRIBUTION STATEMENT (of this Report)		
Approved for Public Release. Distribution Unlimited.		
17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report)		
18. SUPPLEMENTARY NOTES		
19. KEY WORDS (Continue on reverse side if necessary and identify by block number)		
Fourier transform, infrared spectroscopy, solid rocket propellants, propellant aging, sample preparation, attenuated total reflectance, infrared measurement, computerized data reduction, physical property prediction, FTS-10 spectrophotometer, storage stability FORTRAN computer code.		
20. ABSTRACT (Continue on reverse side if necessary and identify by block number)		
This is the final report <del>of a contract</del> to establish the relationship between infrared spectra and propellant aging and to determine the applicability of this relationship to the prediction of service life of propellants. → next page Equipment for measurement of infrared spectra by both attenuated total reflectance and through transmission was identified, and procedures for sample preparation and testing were developed and demonstrated.		

Unclassified

SECURITY CLASSIFICATION OF THIS PAGE(When Data Entered)

The major accomplishments of this project are the establishment of methods for predicting propellant tensile properties from infrared spectral data during storage and the development of procedures and computer programs for rapidly and economically progressing large amounts of infrared data. The results of this project are a major step toward predicting the remaining shelf life of solid propellant rocket motors from infrared measurements using small amounts of propellant associated with them. Recommendations for application of these results and for improving and extending the developed methods are included in the report.

Cont.

This document (Volume II - User's Manual) describes the two computer programs E410 and E490 that were developed in order to reduce and statistically analyze the infrared spectral data. The version of these two codes described herein was made operational on the CDC 6600 computer at the Air Force Rocket Propulsion Laboratory. With minor modifications to these two FORTRAN codes themselves, along with the appropriate job control language, both E410 and E490 can also be made operational on an IBM 360/370 computer system.

Unclassified

SECURITY CLASSIFICATION OF THIS PAGE(When Data Entered)

## TABLE OF CONTENTS

	<u>Page</u>
1.0 ABSTRACT	2
2.0 OVERVIEW OF THE FTIS COMPUTER PROGRAMS	3
2.1 ORGANIZATION OF THE E410 AND E490 COMPUTER CODES	4
3.0 E490 - COMPUTER PROGRAM DESCRIPTION	5
3.0.1 MAIN PROGRAM - E490MAIN	5
3.0.2 SUBROUTINE APCALC	7
3.0.3 SUBROUTINE BASELN	8
3.0.4 SUBROUTINE COBASE	11
3.0.5 SUBROUTINE PKTABL	12
3.0.6 SUBROUTINE PLWAVE	16
3.0.7 SUBROUTINE REGRES	17
4.0 E490 - INPUT INSTRUCTIONS	18
4.1 JOB CONTROL LANGUAGE FOR CDC SYSTEMS	22
5.0 E490 - SAMPLE CASE	26
6.0 E410 - COMPUTER PROGRAM DESCRIPTION	89
7.0 E410 - INPUT INSTRUCTIONS	90
7.0.1 CREATE MODE	90
7.0.2 UPDATE MODE	91
8.0 E410 - SAMPLE CASE	93
APPENDIX A - FORTRAN SOURCE LISTINGS	96
APPENDIX B - PROGRAM FLOWCHARTS	146

## 1.0 ABSTRACT

In recent years, data have been obtained which show that changes in the chemical structure of the binder of a propellant can be correlated with changes in the mechanical properties of the propellant. It has also been demonstrated that these changes to the chemical structure of the binder can be detected using infrared spectroscopy. Therefore, the changes in the infrared absorption characteristics of the binder are directly related to changes in the chemical structure and, hence, to changes in mechanical properties of the propellant. The computer programs described in this manual were developed to demonstrate the feasibility of using Fourier transform infrared spectroscopy as a non-destructive tool to relate changes in binder structure to changes in mechanical properties of the propellant and eventually to use that correlation to predict the remaining service life of the propellant.

## 2.0 OVERVIEW OF THE FTIS COMPUTER PROGRAMS

The Fourier Transform Infrared Spectroscopy (FTIS) computer codes described in this manual, were written in FORTRAN IV - H EXTENDED for use on either CDC 6000 series or IBM 360/370 series digital computers. The primary purposes of these codes are to reduce infrared spectral data stored on a Digilab FTS-10 tape and correlate this spectral information with changes in propellant physical properties.

In order to accomplish this goal, the necessary two codes are:


- I. E490 - Performs the Infrared data reduction and statistical analysis.
- II. E410 - Generates a master file for Mechanical Property data.

## 2.1 ORGANIZATION OF THE E410 AND E490 COMPUTER CODES


The illustration below represents a functional flowchart for both the E410 and E490 programs.

### E490:

Logical Sequence in which  
Subroutines are called by  
E490MAIN.

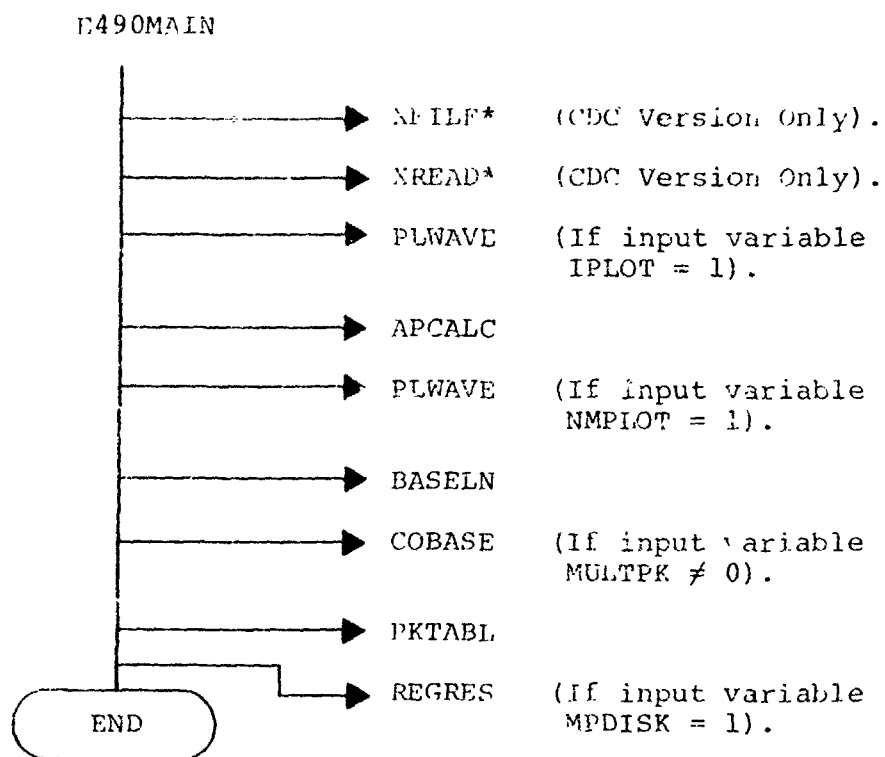

### E410:

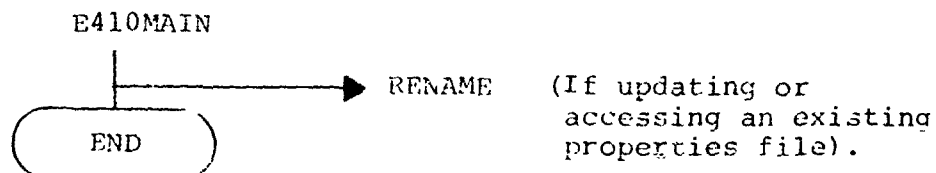

\*Please see description of  
E490MAIN (3.0.1).



### 3.0 E490 - COMPUTER PROGRAM DESCRIPTION

The following are general descriptions of both the main program plus the six (6) subroutines that comprise the E490 computer code.

#### 3.0.1 MAIN PROGRAM - E490MAIN

The MAIN program controls the flow of calculations for this FTIS spectral analysis utilizing both the NAMELIST input supplied by the user and a direct linkage to the six subroutines. Other key operations performed by E490MAIN are described below:

- a). Printing of the input variables supplied by the NAMELIST called /FILES/.
- b). Reads the Digilab FTS-10 data tape.  
An FTS tape file consists of two physical records...a header remarks record and a spectral data record. From the header record, the following variables are read:
  - IFILE() - The file number.
  - NWORDS() - File size in data words.
  - VSPR - Single / double precision designator.
  - NEXP - Exponent of two; (used to compute FTIS amplitude values).
  - II1 - Wave number (frequency) of the first point. (Integer portion).
  - II2 - Numerator of the decimal fraction for II1.
  - II3 - Wave number of the last point. (Integer portion).
  - II4 - Numerator of the decimal fraction for II3.

The spectral data record contains the raw spectral data in binary integer form.

NOTE: Because of the internal differences between IBM and CDC computer systems, the current method for reading this FTS-10 tape on the Air Force /RPL CDC 6600 computer is different. On the CDC system, E490MAIN calls two additional subroutines: XFILE and XREAD.

XFILE : Called at the beginning of E490MAIN, XFILE reads the entire FTS-10 input tape and copies it to a temporary file.

XREAD : By accessing the temporary file generated by XFILE, this routine will read a record from an FTS file.

- c). Calculation of the wave number (WAVENO(i)) and amplitude (AMPLTD(i)) for each spectral data point i within a file.
- d). Tabulation of these wave numbers and amplitudes.
- e). Normalization of these amplitudes, utilizing information computed in subroutine APCALC.
- f). Tabulation of these normalized amplitudes and their corresponding wave numbers.

### 3.0.2 SUBROUTINE APCALC

This subroutine calculates variables necessary for the normalization of the amplitude values (AMPLTD(i)) in all spectral files. The normalization of the array AMPLTD(i) is computed relative to the peak with the maximum amplitude nearest WAVNOR, a wave number value input by the user. For each individual FTIS file, the following normalization parameters are computed in APCALC then returned to the E490MAIN routine:

<u>Variable Name</u>	<u>Definition</u>
a) APMAX	Maximum amplitude nearest WAVNOR.
b) WE	Wave number at APMAX.
c) D2	Calculated amplitude at WE along the baseline of the peak.
d) AP2	Height of the peak located at WE.

### 3.0.3 SUBROUTINE BASELN

Subroutine BASELN searches for and records up to a maximum of fifty (50) peak heights within each individual spectrum. Due to the physical model of this project, the search for peaks is done only over certain wave number intervals, which are defined below:

- I.  $3200. \geq \text{Wave Numbers} \geq 2700.$
- II.  $1800. \geq \text{Wave Numbers} \geq 700.$


NOTE: For files with a last Wave Number greater than 700., the search will end at that final wave number.

#### Key Variables Used in BASELN:

DIFF(IF) : For the IF-th file,  
(Maximum normalized amplitude value) -  
(Minimum normalized amplitude value).

VALUE : Equal to the input variable VALID.  
(The "validity" factor).

FIGURE 1 - Example of Baseline Computations.

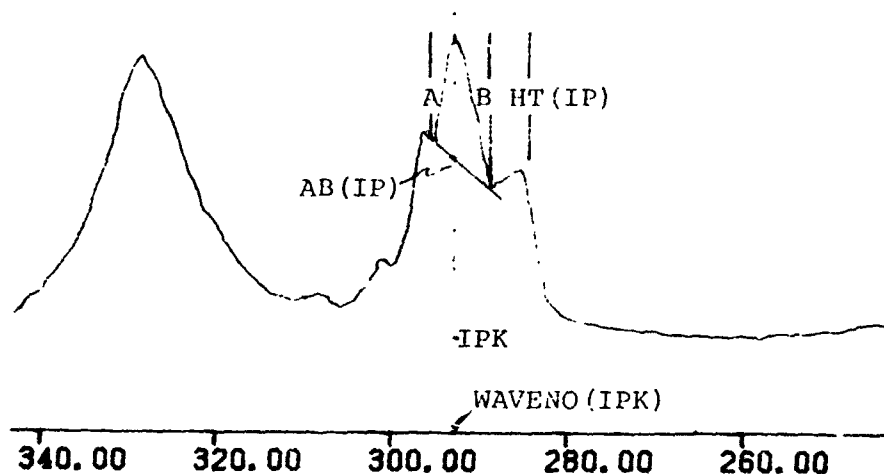

### Key Parameters Calculated by BASELN:

By referencing Figure 1, the following parameters are calculated in order to describe each individual spectral peak.

---

Variable Name(s)	Definition
<hr/>	
1) A and B	(Amplitude at the apex of the peak)  MINUS  (Amplitude at the minimum point(s) on the peak).
2) IP	The IP-th peak in the spectrum, reading from left to right.
3) AB(IP)	The baseline amplitude at wave number WAVENO(IPK).
4) HT(IP)	The actual peak height for the IP-th peak in the spectrum.  (Amplitude at the apex) - AB(IP).
5) IDWORD(IPK)	= 0 : No peak at data word IPK. = 1 : Invalid peak at data word IPK. = 2 : Valid peak at data word IPK. = 3 : Invalid peak at data word IPK, measured by the COBASE subroutine. = 4 : Valid peak at data word IPK, measured by the COBASE subroutine.

The measurement of peak "validity" is a technique utilized to try to eliminate random 'noise' within an FTIS spectrum.

The method for determining the "validity" of a peak is as follows:

- a) Calculate the values for A and B.
- b) From these two values, determine which variable has the maximum value, then set it equal to C.
- c) Multiply DIFF(IF) by VALUE: set product equal to D22.
- d) If C is greater than D22, then the peak at this wave number (WAVENO(IPK)) is a valid peak.

#### 3.0.4 SUBROUTINE COBASE

As an optional subroutine, COBASE will only be executed when the input variable MULTPK is greater than zero and less than seven. By utilizing the same general techniques described in subroutine BASELN, the logic in this routine will establish a "common" baseline for adjacent spectral peaks.

For example, to establish the "common" baseline for the n-th set of adjacent peaks, the COBASE algorithm first searches for NPHOTO(n), the file that contains the optimum picture (or example) of these adjoining peaks. After calculating this new "common" baseline, the data words locating the newly measured peak heights in this file also become the reference locations from which to measure this "common" baseline's peak heights in all the other FTIS files.

It should be noted that the maximum allowable number of peaks per "common" baseline area is twenty-five (25).

### 3.0.5 SUBROUTINE PKTABL

PKTABL is an output subroutine that generates two different types of tabulations listing FTIS peak height information. The first type of tabulation, as illustrated by Figure 2 , is a "Normalized Peak Height Information" listing that is printed for every FTIS file analyzed. As illustrated in Figure 3 , the second tabulation entitled "Peak Height Tabulation for all the Files" is strictly a listing, referenced by the Data Word location, of all the peak heights found in all the FTIS files examined.

A check is also made to determine whether the mechanical property correlation subroutine REGRES will be called. If so, then all the peak height values that are simultaneously located at the same Data Word location and valid in every FTIS file examined, will be stored on a temporary disk file (File 27) for use as input to the REGRES subroutine.



FIGURE 2

FOURIER TRANSFORM INFRARED SPECTROSCOPY -- NORMALIZED PEAK HEIGHT INFORMATION : FILE NUMBER 20

\*\*\* - DENOTES A VALID PEAK.  
CB. - DENOTES A PEAK MEASURED FROM A "COMMON BASELINE".

PEAK NUMBER	DATA WORD	WAVE NUMBER	PEAK HEIGHT	BASELINE AMPLITUDE
1	120	3084.1318	0.9578886	11.7398720
*** 2 CB.	130	3006.8142	3.4200430	11.4080601
*** 3 CB.	136	2960.4236	12.4571257	11.2908144
*** 4 CB.	141	2921.7649	18.6699829	11.1931105
*** 5 CB.	146	2883.1060	9.1964798	11.0954056
*** 5 CB.	150	2852.1790	10.7229776	11.0172415
7	166	2728.4707	0.1328259	10.6185064
9	170	2697.5437	0.2030725	10.3796721
*** 9	294	1738.8047	23.0211334	14.8958855
10	306	1646.0234	1.6773682	16.5301514
11	312	1599.6328	0.0498047	16.0350342
*** 12	323	1514.5835	0.5473433	19.5836029
*** 13	336	1414.0706	49.9136200	22.3494263
14	349	1313.5576	0.5114430	24.3073425
15	352	1290.3623	0.1609344	23.6483154
*** 16	359	1236.2400	2.5067139	24.3601227
*** 17	343	1050.6775	49.9581604	66.9271393
*** 18	394	965.6282	21.0735474	79.3897858
*** 19	398	934.7009	9.2646484	68.1747894
*** 20	401	911.5054	13.2904053	62.8501740
21	419	772.3340	1.0077057	37.4304047
*** 22	425	725.9434	5.3882509	35.2491509

FIGURE 3

***** PEAK HEIGHT TABULATION FOR ALL THE FILES *****									
DATA	FILE	FILE	FILE	FILE	FILE	FILE	FILE	FILE	FILE
WORD	16	20	24	28	0	0	0	0	0
113	0.0	0.0	0.0	0.191	0.0	0.0	0.0	0.0	0.0
114	0.0	0.0	0.241	0.0	0.0	0.0	0.0	0.0	0.0
120	1.058	0.958	1.031	0.809	0.0	0.0	0.0	0.0	0.0
130	3.603	3.420	3.505	3.351	0.0	0.0	0.0	0.0	0.0
136	11.895	12.457	12.289	13.109	0.0	0.0	0.0	0.0	0.0
141	18.336	18.670	18.798	19.183	0.0	0.0	0.0	0.0	0.0
146	9.024	9.196	9.182	9.418	0.0	0.0	0.0	0.0	0.0
150	10.723	10.723	10.723	10.723	0.0	0.0	0.0	0.0	0.0
166	0.176	0.133	0.211	0.222	0.0	0.0	0.0	0.0	0.0
169	0.146	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
170	0.0	0.203	0.0	0.0	0.0	0.0	0.0	0.0	0.0
288	0.0	0.0	0.081	0.0	0.0	0.0	0.0	0.0	0.0
294	21.621	23.021	23.274	25.219	0.0	0.0	0.0	0.0	0.0
303	0.0	0.0	0.200	0.0	0.0	0.0	0.0	0.0	0.0
306	1.972	1.677	1.402	1.792	0.0	0.0	0.0	0.0	0.0
311	0.0	0.0	0.084	0.0	0.0	0.0	0.0	0.0	0.0
312	0.0	0.050	0.0	0.0	0.0	0.0	0.0	0.0	0.0
314	0.143	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
316	0.368	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
323	0.796	0.547	0.602	0.0	0.0	0.0	0.0	0.0	0.0
336	45.982	49.214	48.401	50.002	0.0	0.0	0.0	0.0	0.0
349	0.513	0.513	0.449	0.570	0.0	0.0	0.0	0.0	0.0
352	0.154	0.161	0.170	0.150	0.0	0.0	0.0	0.0	0.0
358	2.465	2.507	2.522	2.515	0.0	0.0	0.0	0.0	0.0
393	44.327	49.959	49.027	44.155	0.0	0.0	0.0	0.0	0.0

FTIS \*\*\*\*\* OFAK HEIGHT TABULATION FOR ALL THE FILES \*\*\*\*\*

DATA WORD	FILE 16	FILE 20	FILE 24	FILE 28	FILE 0	FILE 0	FILE 0	FILE 0	FILE 0	FILE 0	FILE 0
394	20.645	23.074	22.071	21.456	0.0	0.0	0.0	0.0	0.0	0.0	0.0
399	6.477	8.265	7.339	7.983	0.0	0.0	0.0	0.0	0.0	0.0	0.0
401	11.813	13.290	12.457	13.037	0.0	0.0	0.0	0.0	0.0	0.0	0.0
419	1.252	1.009	1.303	1.212	0.0	0.0	0.0	0.0	0.0	0.0	0.0
425	5.446	5.389	5.425	3.581	0.0	0.0	0.0	0.0	0.0	0.0	0.0

NORMAL END OF PEAK HEIGHT TABLE :  
TOTAL NUMBER OF FILES LISTED = 4.

### 3.0.6 SUBROUTINE PLWAVE

When requested by the user, this routine plots the Wave Number versus Amplitude values for the individual FTIS spectral files. Depending upon the value of the input variables IPLOT and NMPLLOT, the option is available for plotting either "Non-Normalized" or "Normalized" Amplitude data values or if so requested, both types of spectral plots.

### 3.0.7 SUBROUTINE REGRES

Subroutine REGRES is also an option that, when requested, (input variable MPDISK = 1) will determine whether there are any correlations between FTIS spectral data and corresponding propellant mechanical properties tabulated by the E410 computer program.

The primary logic in REGRES provides a multiple linear regression model that defines as its independent variables the set of "valid" peak heights written on File 27 by subroutine PKTABL. Regarding the dependent variable(s), they are determined by the input variable NPHYSP(k), the array that specifically defines what mechanical properties the user wishes to examine.

For each dependent variable, a multiple correlation coefficient (RMULT) is computed. For any dependent variable where  $RMULT \geq .750$ , a non-linear regression attempt is made using that single dependent variable and those independent variables with non-zero linear regression coefficients. At the end of this attempt, a test is made comparing its multiple correlation coefficient with the linear one, thereby, determining the best "curve fit" for this particular dependent variable (or mechanical property).

#### 4.0 E490 INPUT INSTRUCTIONS

Input variables for E490 are contained in a single NAMELIST called /FILES/. By using NAMELIST, a free form input is available requiring only that the first column on each input card be blank.

The following table, TABLE 1, presents a complete listing and description of the variables necessary to execute the E490 computer program.

<u>VARIABLE NAME</u>	<u>DIMENSION</u>	<u>TYPE</u>	<u>DEFAULT VALUE</u>	<u>DEFINITION AND REMARKS</u>
<u>ITAPE</u>	1	Integer (Literal)	None	The name of the <u>FTIS</u> input tape.  Note: Maximum Hollerith character length for ITAPE is: (1) 4 for IBM 370 systems (2) 10 for CDC 6600 systems
<u>ITOTAL</u>	1	Integer	None	Total number of spectral files on the <u>FTIS</u> tape that will be examined.
<u>INFILE</u>	30	Integer	Ø	The set of <u>FTIS</u> file numbers whose spectral data will be read from the tape.  Note: Maximum of 30 files can be read per computer run.
<u>MULTPK</u>	1	Integer	Ø	Specific number of "Common Baseline" areas per <u>FTIS</u> spectrum.  Maximum Value = 6
<u>NPHOTO</u>	6	Integer	Ø	<u>NPHOTO(i)</u> = <u>FTIS</u> file number in which the i <sup>th</sup> "Common Baseline" areas peaks are most definitive.
<u>LLIMIT</u>	6	Real	Ø.Ø	<u>LLIMIT(i)</u> = A Wave Number value within file <u>NPHOTO(i)</u> , that defines an initial estimate for the left limit of "Common Baseline" Area i.

Note: Use a Wave Number value that is approx. 2/3  
the way down the left side of the left most peak

TABLE 1: F490 INPUT VARIABLES

<u>VARIABLE NAME</u>	<u>DIMENSION</u>	<u>TYPE</u>	<u>DEFAULT VALUE</u>	<u>DEFINITION AND REMARKS</u>
<u>RLIMIT</u>	6	Real	0.0	<u>RLIMIT(i)</u> = A Wave Number value within file <u>NPHOTO(i)</u> , that defines an initial estimate for the right limit of "Common Baseline" Area i.  Note: Use Wave Number value that is approx. 2/3 the way down the right side of the right most peak.
<u>VALID</u>	1	Real	.02	The Peak Height validity factor.
<u>WAVNOR</u>	1	Real	None	Normalizing Wave Number for all FTIS files.
<u>MPDISK</u>	1	Integer	0	Specifies whether statistical correlation of peak heights versus mechanical properties is requested.  If = 0, Correlation option is <u>not</u> requested. = 1, Correlation option is requested.
<u>MPRECS</u>	30	Integer	0	Record number(s) from the mechanical properties input tape that are in a one-to-one correspondence with the <u>INFILE(i)</u> array.  Example: If <u>INFILE(2)</u> = 3 and <u>MPRECS(2)</u> = 9, then the <u>Mechanical Properties</u> for FTIS File No. 3 are located at the 9th record of data on the Mechanical Properties tape.

TABLE 1: E490 INPUT VARIABLES (Cont.)



<u>VARIABLE NAME</u>	<u>DIMENSION</u>	<u>TYPE</u>	<u>DEFAULT VALUE</u>	<u>DEFINITION AND REMARKS</u>
NPHYSP	10	Integer	Ø	<p>Column numbers corresponding to the E410 physical property listing, which identify what properties to use as dependent variables in subroutine REGRES.</p> <p>In the current E410 configuration, the possible dependent variables are:</p> <ul style="list-style-type: none"> <li>1 - Modulus</li> <li>2 - Strain at Break</li> <li>3 - Strain at Maximum Stress</li> <li>4 - Maximum Stress</li> <li>5 - Strain Energy Density</li> <li>6 - Strain Endurance</li> <li>7 - TBD</li> <li>8 - TBD</li> <li>9 - TBD</li> <li>10 - TBD</li> </ul>

21

IPLOT	1	Integer	Ø	<p>"Non-normalized" plot designator.</p> <p>= 0 , No plots have been requested.</p> <p>= 1 , Plots have been requested.</p>
NMPLOT	1	Integer	Ø	<p>"Normalized" plot designator.</p> <p>= 0 , No plots have been requested.</p> <p>= 1 , Plots have been requested.</p>

TABLE 1: E490 INPUT VARIABLES (Cont.)

#### 4.1 JOB CONTROL LANGUAGE FOR CDC SYSTEMS

The following are examples of the Job Control Language needed to execute computer program E490 on a CDC-6000 series computer.

---

OPTION ONE: RUN REQUESTING NO PLOTS AND NO  
STATISTICAL CORRELATION.

```
SEQ,THA01.  
THA,T100,NT1,P4. LA,YOURNAME,30  
ATTACH,LGO,E490GO,ID=YOURNAME,MR=1.  
REQUEST,TAPE1,NT,HD,S,NORING,VSN=FTISTAPENAME  
FILE(TAPE1,BT=K,RT=U,RB=1,MBL=1257,MRL=1257,CM=NO,MNR=24,MNB=24)  
LCSET(FILES=TAPE1)  
LIBRARY,BIT8LIB.  
MAP,PART.  
LGO.  
UNLOAD,TAPE1.
```



Card Column 1

OPTION TWO: RUN REQUESTING PLOTS BUT NO  
STATISTICAL CORRELATIONS.

SEQ,THA01.  
THA,T100,NT1,P4. LA,YOURNAME,30  
ATTACH,LGO,E490GO,ID=YOURNAME,MR=1.  
REQUEST,TAPE1,NT,HD,S,NORING,VSN=FTISTAPENAME.  
FILE(TAPF1,BT=K,RT=U,RR=1,MBL=1257,MRL=1257,CM=NO,MNR=24,MNB=24)  
LDSET(FILEFS=TAPE1)  
LIPRARY,BITBLIB.  
MAP,PART.  
LGO.  
UNLOAD,TAPE1.  
REWIND,TAPE8.  
REQUEST,PLOT,HD,RING,VSN=PLOT. YOURNAME  
COPYBF,TAPE8,PLOT.  
RETURN,PLOT.

OPTION THREE: RUN REQUESTING STATISTICAL  
CORRELATIONS BUT NO PLOTS.

SEQ,THA01.  
THA,T100,NT1,P4. LA,YOURLNAME,30  
ATTACH,LGO,E490GO,ID=YOURLNAME,MR=1.  
ATTACH,TAPE10,MECHPROFFILE,ID=YOURLNAME,MR=1.  
REQUEST,TAPE1,IT,HD,S,NORING,VSN=FTISTAPENAME.  
FILE(TAPE1,BT=K,RT=U,RB=1,MBL=1257,MRL=1257,CM=NO,MNR=24,MNB=24)  
LDSET(FILE=TAPE1)  
LIBRARY,BIT8LIB.  
MAP,PART.  
LGO.  
UNLOAD,TAPE1.

OPTION FOUR: RUN REQUESTING PLOTS AND  
STATISTICAL CORRELATIONS.

```
SEQ,THA01.  
THA,T100,NT1,P4. LA,YOURNAME,30  
ATTACH,LGO,E490G0,ID=YOURNAME,MR=1.  
ATTACH,TAPE10,MFCHPROFILE,ID=YOURNAME,MR=1.  
REQUEST,TAPE1,NT,HD,S,NORING,VSN=FTISTAPENAME.  
FILE(TAPE1,BT=K,RT=U,RB=1,MBL=1257,MRL=1257,CM=NO,MNR=24,MNB=24)  
LOSET(FILE=TAPE1)  
LIBRARY,BIT8LIB.  
MAP,PART.  
LGO.  
UNLOAD,TAPE1.  
REWIND,TAPE8.  
REQUEST,PLOT,HD,RING,VSN=PLOT. YOURNAME  
COPY8F,TAPE8,PLOT.  
RETURN,PLOT.
```

## 5.0 E490 - SAMPLE CASE

The card deck shown in Figure 4 is the Job Control Language and input data that generated the sample case illustrated in Table 2.

### E490 SAMPLE CASE

```
SEQ,THA01.
THA,T100,NT1,P4.  LA,SMITHQ,30
ATTACH,LGO,E490GO,ID=SMITHQ,MR=1.
ATTACH,TAPE10,PROFFILETP25,ID=SMITHQ,MR=1.
REQUEST,TAPE1,NT,HQ,S,NORING,VSN=TAPE25.
FILE(TAPE1,RT=K,RT=U,PR=1,MRL=1257,MRL=1257,CM=NO,MNR=24,MNR=24)
LOSET(FILE$=TAPE1)
LIBRARY,BIT8LIB.
MAP,P'RT.
LGO.
UNLOAD,TAPE1.
$FILES
ITAPE='TP25',
ITOTAL=4,
INFILE=16,20,24,28,
MULTPK=1,
NPHOTO=24,
LLIMIT=3030.,
RLIMIT=2750.,
WAVNOR=2850.,
VALID = .02,
MPDISK=1,
MPRECS=17,21,25,29,
NPHYSP=1,2,3,4,
IPLUT = 0,
NMPLGT = 0,
$END
```

FIGURE 4. E490 SAMPLE CASE



16	468	0	1	401	13716
4004.7119		401.2092	-7.7152		
			4004	13902	



FTIS FILE NUMBER : 10

NON-NORMALIZED (PURE) SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
1	4004.2119	0.9686279E+01	53	3603.0212	0.9716797E+01	105	3201.8303	0.1785807E+02
2	3996.4966	0.9375000E+01	54	3595.3059	0.9564209E+01	106	3194.1150	0.1701050E+02
3	3988.7815	0.9161377E+01	55	3587.5908	0.9771729E+01	107	3186.3999	0.1655841E+02
4	3981.0662	0.9167480E+01	56	3579.8755	0.9649658E+01	108	3178.6846	0.1572876E+02
5	3973.3511	0.9600830E+01	57	3572.1604	0.9875488E+01	109	3170.9695	0.1505127E+02
6	3965.6357	0.8850098E+01	58	3564.4451	0.9826660E+01	110	3163.2542	0.1480103E+02
7	3957.9207	0.9155273E+01	59	3556.7300	0.1009777E+02	111	3155.5391	0.1444897E+02
8	3950.2053	0.8984375E+01	60	3549.0146	0.9826660E+01	112	3147.8237	0.1409912E+02
9	3942.4902	0.8728027E+01	61	3541.2996	0.9967041E+01	113	3140.1086	0.1368408E+02
10	3934.7749	0.8923340E+01	62	3533.5842	0.1033933E+02	114	3132.3933	0.1359253E+02
11	3927.0598	0.8612061E+01	63	3525.8691	0.985178E+01	115	3124.6782	0.1347048E+02
12	3919.3445	0.8537715E+01	64	3518.1538	0.1015625E+02	116	3116.9629	0.1318970E+02
13	3911.6294	0.8526611E+01	65	3510.4387	0.1074829E+02	117	3109.2478	0.1309814E+02
14	3903.9141	0.9191895E+01	66	3502.7234	0.1109009E+02	118	3101.5325	0.1324463E+02
15	3896.1990	0.8441162E+01	67	3495.0083	0.105298E+02	119	3093.8174	0.1345825E+02
16	3888.4836	0.8563232E+01	68	3487.2930	0.1008470E+02	120	3086.1021	0.1397329E+02
17	3880.7686	0.8978271E+01	69	3479.5778	0.1115112E+02	121	3078.3870	0.1353760E+02
18	3873.0532	0.8972168E+01	70	3471.8625	0.1130371E+02	122	3070.6716	0.1322632E+02
19	3865.3381	0.9326172E+01	71	3464.1475	0.1179199E+02	123	3062.9565	0.1270142E+02
20	3857.6228	0.9094238E+01	72	3456.4321	0.1231689E+02	124	3055.2412	0.1286011E+02
21	3849.9077	0.8990479E+01	73	3448.7170	0.1208055E+02	125	3047.5261	0.1318970E+02
22	3842.1924	0.9283447E+01	74	3441.0017	0.1274414E+02	126	3039.8108	0.1376953E+02
23	3834.4773	0.9375000E+01	75	3433.2866	0.1328403E+02	127	3032.0957	0.1387329E+02
24	3826.7620	0.9210205E+01	76	3425.5713	0.1329956E+02	128	3024.3804	0.1440710E+02
25	3819.0469	0.9228516E+01	77	3417.8562	0.1353149E+02	129	3016.6653	0.1586306E+02
26	3811.3315	0.9643355E+01	78	3410.1409	0.1411133E+02	130	3008.9502	0.1584473E+02
27	3803.6165	0.1052856E+02	79	3402.4258	0.1413574E+02	131	3001.2349	0.1560059E+02
28	3795.9011	0.9716797E+01	80	3394.7104	0.1496582E+02	132	2993.5198	0.1622013E+02
29	3788.1860	0.9594727E+01	81	3386.9954	0.1503296E+02	133	2985.8044	0.1732703E+02
30	3780.4707	0.9985352E+01	82	3379.2800	0.1608887E+02	134	2978.0894	0.1934214E+02
31	3772.7556	0.9948730E+01	83	3371.5649	0.1629028E+02	135	2970.3740	0.2249116E+02
32	3765.0403	0.9313965E+01	84	3363.8496	0.1697998E+02	136	2962.6589	0.2327527E+02
33	3757.3252	0.9716797E+01	85	3356.1345	0.1746663E+02	137	2954.9436	0.2260127E+02
34	3749.6099	0.1015015E+02	86	3348.4192	0.1867676E+02	138	2947.2285	0.2479853E+02
35	3741.8950	0.9600830E+01	87	3340.7041	0.1980591E+02	139	2939.5132	0.2789917E+02
36	3734.1797	0.9625244E+01	88	3332.9888	0.211815E+02	140	2931.7981	0.2982600E+02
37	3726.4646	0.9838867E+01	89	3325.2737	0.2285767E+02	141	2924.0828	0.2913213E+02
38	3718.7493	0.9759212E+01	90	3317.5583	0.2443237E+02	142	2916.3677	0.2732710E+02
39	3711.0342	0.9863281E+01	91	3309.8433	0.2566528E+02	143	2908.6525	0.2465820E+02
40	3703.3188	0.9875488E+01	92	3302.1279	0.2756349E+02	144	2900.9373	0.2225952E+02
41	3695.6038	0.9491455E+01	93	3294.4128	0.2789917E+02	145	2893.2219	0.2033091E+02
42	3687.8884	0.1024780E+02	94	3286.6975	0.285585E+02	146	2885.5068	0.2044067E+02
43	3680.1733	0.5716797E+01	95	3278.9824	0.2804565E+02	147	2877.7915	0.2027593E+02
44	3672.4580	0.9947695E+01	96	3271.2671	0.2682261E+02	148	2870.0764	0.2076415E+02
45	3664.7429	0.9564209E+01	97	3263.5520	0.2593384E+02	149	2862.3611	0.2111706E+02
46	3657.0276	0.9582520E+01	98	3255.8367	0.2440196E+02	150	2854.6460	0.2109522E+02
47	3649.3125	0.9759521E+01	99	3248.1216	0.2734747E+02	151	2846.9307	0.2039765E+02
48	3641.5972	0.9473828E+01	100	3240.4062	0.2707612E+02	152	2839.2156	0.1695557E+02
49	3633.8821	0.9655762E+01	101	3232.6912	0.2133179E+02	153	2831.5002	0.1435547E+02
50	3626.1667	0.9619141E+01	102	3224.9758	0.2017212E+02	154	2823.7852	0.1340337E+02
51	3618.4517	0.9744100E+01	103	3217.2607	0.1945901E+02	155	2816.0698	0.1291111E+02
52	3610.7366	0.9744100E+01	104	3209.5454	0.1810044E+02	156	2808.3544	0.1241111E+02

FTIS FILE NUMBER : 16

NON-NORMALIZED (PURE) SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
157	2800.6394	0.1240845E+02	209	2399.4485	0.1160889E+02	261	1998.2578	0.1259766E+02			
158	2792.9243	0.122145E+02	210	2391.7334	0.1153564E+02	262	1990.5425	0.1254272E+02			
159	2785.2090	0.1212769E+02	211	2384.0181	0.1152954E+02	263	1982.4274	0.1234131E+02			
160	2777.4939	0.1195679E+02	212	2376.3030	0.1165771E+02	264	1975.1121	0.1215210E+02			
161	2769.7786	0.1192017E+02	213	2368.5876	0.1160278E+02	265	1967.3970	0.1213379E+02			
162	2762.0635	0.1189575E+02	214	2360.8726	0.1160278E+02	266	1959.6816	0.1207275E+02			
163	2754.3481	0.1176756E+02	215	2353.1572	0.1156006E+02	267	1951.9666	0.1203613E+02			
164	2746.6331	0.1169434E+02	216	2345.4421	0.1162720E+02	268	1944.2512	0.1200562E+02			
165	2738.9177	0.1170044E+02	217	2337.7268	0.1169434E+02	269	1936.5361	0.1196899E+02			
166	2731.2026	0.1171265E+02	218	2330.0117	0.1157227E+02	270	1928.8208	0.1199341E+02			
167	2723.4873	0.1166382E+02	219	2322.2964	0.1156006E+02	271	1921.1057	0.1196289E+02			
168	2715.7722	0.1140747E+02	220	2314.5813	0.1148682E+02	272	1913.3904	0.1207886E+02			
169	2708.0569	0.1152344E+02	221	2306.8660	0.1147461E+02	273	1905.6753	0.1198120E+02			
170	2700.3418	0.1145020E+02	222	2299.1509	0.1149292E+02	274	1897.9600	0.1196899E+02			
171	2692.6265	0.1139526E+02	223	2291.4355	0.1154175E+02	275	1890.2449	0.1205444E+02			
172	2684.9114	0.1133423E+02	224	2283.7205	0.1145630E+02	276	1882.5295	0.1210327E+02			
173	2677.1960	0.1145020E+02	225	2276.0051	0.1153564E+02	277	1874.8145	0.1207275E+02			
174	2669.4810	0.1131592E+02	226	2268.2900	0.1140137E+02	278	1867.0991	0.1212769E+02			
175	2661.7656	0.1119385E+02	227	2260.5747	0.1138306E+02	279	1859.3840	0.1218267E+02			
176	2654.0505	0.1121826E+02	228	2252.8596	0.1135254E+02	280	1851.6687	0.1213989E+02			
177	2646.3352	0.1118164E+02	229	2245.1443	0.1141357E+02	281	1843.9536	0.1235352E+02			
178	2638.6201	0.1114502E+02	230	2237.4292	0.1133423E+02	282	1836.2383	0.1237183E+02			
179	2630.9048	0.1104736E+02	231	2229.7139	0.1137085E+02	283	1828.5232	0.1261597E+02			
180	2623.1897	0.1107178E+02	232	2221.9988	0.1141968E+02	284	1820.8079	0.1245728E+02			
181	2615.4744	0.1107798E+02	233	2214.2834	0.1143799E+02	285	1813.0928	0.1245728E+02			
182	2607.7593	0.1106567E+02	234	2206.5684	0.1143799E+02	286	1805.3774	0.1242065E+02			
183	2600.0439	0.1108022E+02	235	2198.8530	0.1139526E+02	287	1797.6624	0.1241455E+02			
184	2592.3289	0.1101685E+02	236	2191.1379	0.1141357E+02	288	1789.9470	0.1251831E+02			
185	2584.6135	0.1091919E+02	237	2183.4226	0.1143799E+02	289	1782.2319	0.1258545E+02			
186	2576.8984	0.1091919E+02	238	2175.7073	0.1145020E+02	290	1774.5166	0.1288540E+02			
187	2569.1831	0.1093750E+02	239	2167.9922	0.1146751E+02	291	1766.8015	0.1333618E+02			
188	2561.4680	0.1093140E+02	240	2160.2771	0.1149292E+02	292	1759.0862	0.1305737E+02			
189	2553.7527	0.1100464E+02	241	2152.5618	0.1154175E+02	293	1751.3711	0.1242497E+02			
190	2546.0376	0.1098022E+02	242	2144.8467	0.1161499E+02	294	1743.6558	0.1352661E+02			
191	2538.3223	0.1089478E+02	243	2137.1313	0.1153564E+02	295	1735.9407	0.1276855E+02			
192	2530.6072	0.1092745E+02	244	2129.4163	0.1154785E+02	296	1728.2253	0.1286987E+02			
193	2522.8918	0.1102295E+02	245	2121.7009	0.1149902E+02	297	1720.5103	0.1198486E+02			
194	2515.1768	0.1102295E+02	246	2114.9859	0.1162109E+02	298	1712.7949	0.1790161E+02			
195	2507.4614	0.1108398E+02	247	2106.2705	0.1166382E+02	299	1705.0798	0.1664429E+02			
196	2499.7463	0.1110940E+02	248	2098.5554	0.1179810E+02	300	1697.3645	0.1160217E+02			
197	2492.0310	0.1098633E+02	249	2090.8401	0.1179810E+02	301	1689.6494	0.1509399E+02			
198	2484.3159	0.1110840E+02	250	2083.1250	0.1198120E+02	302	1681.9341	0.1522217E+02			
199	2476.6006	0.1121216E+02	251	2075.4097	0.1195679E+02	303	1674.2190	0.1555176E+02			
200	2468.8855	0.1121826E+02	252	2067.6946	0.1208476E+02	304	1666.5037	0.1567383E+02			
201	2461.1702	0.1128540E+02	253	2059.9792	0.1209482E+02	305	1658.7886	0.1596069E+02			
202	2453.4551	0.1125488E+02	254	2052.2642	0.1223145E+02	306	1651.0732	0.1596069E+02			
203	2445.7397	0.1136033E+02	255	2044.5488	0.1229077E+02	307	1643.3582	0.1601564E+02			
204	2438.0247	0.1132813E+02	256	2036.8337	0.1222534E+02	308	1635.6424	0.1491039E+02			
205	2430.3093	0.1141968E+02	257	2029.1184	0.1231689E+02	309	1627.9277	0.1450806E+02			
206	2422.5942	0.1142578E+02	258	2021.4033	0.1246338E+02	310	1620.2124	0.1445319E+02			
207	2414.8789	0.1142578E+02	259	2013.6882	0.1254883E+02	311	1612.4973	0.1438599E+02			
208	2407.1638	0.1141133E+02	260	2005.9729	0.1258545E+02	312	1604.7870	0.1436768E+02			

FTIS FILE NUMBER : 16

NON-NORMALIZED (PURE) SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
313	1597.0669	0.1427002E+02	365	1195.8760	0.2277832E+02	417	794.6851	0.3062.1			
314	1589.3516	0.1441040E+02	366	1188.1606	0.2447510E+02	418	785.9700	0.3093267E+02			
315	1581.6365	0.1428833E+02	367	1180.4456	0.2645874E+02	419	779.2546	0.3118423E+02			
316	1573.9211	0.1478827E+02	368	1172.7302	0.2933960E+02	420	771.5396	0.3049927E+02			
317	1566.2061	0.1461187E+02	369	1165.0151	0.3143311E+02	421	763.8242	0.3005371E+02			
318	1558.4907	0.1563721E+02	370	1157.2998	0.3468018E+02	422	755.1091	0.3016967E+02			
319	1550.7756	0.1587524E+02	371	1149.5847	0.4074707E+02	423	748.3938	0.3112793E+02			
320	1543.0603	0.1665649E+02	372	1141.8694	0.4663696E+02	424	740.6787	0.3273315E+02			
321	1535.3452	0.1724243E+02	373	1134.1543	0.5092773E+02	425	732.9634	0.3391113E+02			
322	1527.6299	0.1760254E+02	374	1126.4390	0.5470581E+02	426	725.2483	0.3323364E+02			
323	1519.9148	0.1834717E+02	375	1118.7239	0.5917358E+02	427	717.5330	0.3283691E+02			
324	1512.1995	0.1811523E+02	376	1111.0085	0.6463622E+02	428	709.8179	0.3243408E+02			
325	1504.4844	0.1836548E+02	377	1103.2935	0.7027588E+02	429	702.1025	0.3210449E+02			
326	1496.7690	0.1858521E+02	378	1095.5781	0.7467651E+02	430	694.3875	0.3155518E+02			
327	1489.0540	0.1976929E+02	379	1087.8630	0.7954712E+02	431	686.6721	0.3126831E+02			
328	1481.3386	0.2127686E+02	380	1080.1477	0.8591309E+02	432	678.9570	0.3034058E+02			
329	1473.6235	0.2564087E+02	381	1072.4326	0.9144897E+02	433	671.2417	0.2899780E+02			
330	1465.9082	0.2857666E+02	382	1064.7173	0.9447632E+02	434	663.5266	0.2767334E+02			
331	1458.1931	0.3274536E+02	383	1057.0022	0.9584961E+02	435	655.8113	0.2601318E+02			
332	1450.4778	0.3718262E+02	384	1049.2869	0.9490356E+02	436	648.0962	0.2640341E+02			
333	1442.7627	0.4520874E+02	385	1041.5718	0.9298096E+02	437	640.3809	0.4874258E+02			
334	1435.0474	0.5422974E+02	386	1033.8567	0.9061279E+02	438	632.6658	0.8100586E+02			
335	1427.3323	0.6139526E+02	387	1026.1414	0.8741455E+02	439	624.9504	0.1044434E+03			
336	1419.6169	0.6190186E+02	388	1018.4263	0.8071289E+02	440	617.2354	0.9537355E+02			
337	1411.9019	0.5666504E+02	389	1010.7109	0.7901001E+02	441	609.5200	0.3342830E+02			
338	1404.1865	0.4954224E+02	390	1002.9958	0.7901001E+02	442	601.8049	0.3665723E+02			
339	1396.4714	0.4344882E+02	391	995.2805	0.7515259E+02	443	594.0896	0.6726076E+02			
340	1388.7561	0.3815308E+02	392	987.5654	0.7125244E+02	444	586.3745	0.4687477E+02			
341	1381.0410	0.3274536E+02	393	979.8501	0.7406616E+02	445	578.6592	0.4738413E+02			
342	1373.3257	0.2973633E+02	394	972.1350	0.8401499E+02	446	570.9441	0.4708477E+02			
343	1365.6106	0.2787476E+02	395	964.4197	0.7623291E+02	447	563.2288	0.6137971E+02			
344	1357.8953	0.2603149E+02	396	956.7046	0.6483154E+02	448	555.5137	0.5061401E+02			
345	1350.1802	0.2437134E+02	397	949.9893	0.5566406E+02	449	547.7983	0.4056519E+02			
346	1342.4648	0.2274780E+02	398	941.2742	0.6162109E+02	450	540.0813	0.6015345E+02			
347	1334.7498	0.2144165E+02	399	933.5588	0.5565796E+02	451	532.3679	0.5762933E+02			
348	1327.0344	0.2106323E+02	400	925.8437	0.5589379E+02	452	524.6528	0.7426636E+02			
349	1319.3193	0.2127192E+02	401	918.1284	0.4194453E+02	453	516.9375	0.5762933E+02			
350	1311.6040	0.2105101E+02	402	910.4133	0.5543721E+02	454	509.2224	0.4956444E+02			
351	1303.8889	0.2012324E+02	403	902.6980	0.4823608E+02	455	501.5071	0.4345441E+02			
352	1296.1736	0.2015991E+02	404	894.9829	0.4436033E+02	456	493.7920	0.4056519E+02			
353	1288.4585	0.1900947E+02	405	887.2676	0.4190674E+02	457	486.0767	0.5027466E+02			
354	1280.7432	0.1999345E+02	406	879.5525	0.3933533E+02	458	478.3616	0.5167505E+02			
355	1273.0281	0.2053223E+02	407	871.8372	0.3810423E+02	459	470.6462	0.5120859E+02			
356	1265.3127	0.2141724E+02	408	864.1221	0.3686523E+02	460	462.9312	0.3942627E+02			
357	1257.5977	0.2194824E+02	409	856.4067	0.3526001E+02	461	455.2158	0.4687500E+02			
358	1249.8823	0.2251597E+02	410	848.6917	0.3439331E+02	462	447.5007	0.4415946E+02			
359	1242.1672	0.2283325E+02	411	840.9763	0.3332520E+02	463	439.7855	0.4175415E+02			
360	1234.4519	0.2275952E+02	412	833.2612	0.3263550E+02	464	432.0703	0.4175415E+02			
361	1226.7368	0.2192383E+02	413	825.5459	0.3190308E+02	465	424.3550	0.4175415E+02			
362	1219.0215	0.2127686E+02	414	817.8308	0.3140869E+02	466	416.6399	0.4175415E+02			
363	1211.3064	0.2109375E+02	415	810.1155	0.3081665E+02	467	408.9246	0.4175415E+02			
364	1203.5913	0.2127686E+02	416	802.4003	0.3040027E+02	468	401.2095	0.4175415E+02			

20	468	0	1	13902	393	31295
4004.2119	393.4773	-7.7318	4004			

FTIS FILE NUMBER : 20

NON-NORMALIZED (PURE) SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
1	4004.2119	0.5373535E+02	53	3602.02	0.6213379E+01	105	3200.1084	0.1395.00
2	3996.4800	0.5373535E+02	54	3594.185	0.6353760E+01	106	3192.3767	0.1371.00
3	3988.7483	0.5373535E+02	55	3586.6968	0.6414795E+01	107	3184.6449	0.1282.00
4	3981.0166	0.5373535E+02	56	3578.9648	0.6738281E+01	108	3176.9131	0.1196.00
5	3973.2847	0.5373535E+02	57	3571.2332	0.6579500E+01	109	3169.1814	0.1156.00
6	3965.5530	0.6408691E+01	58	3563.5015	0.6549072E+01	110	3161.4495	0.111084.00
7	3957.8213	0.5120850E+01	59	3555.7695	0.6518555E+01	111	3153.7173	0.1095581E+02
8	3950.0894	0.5432129E+01	60	3548.0378	0.6530762E+01	112	3145.9861	0.1041870E+02
9	3942.3577	0.4943848E+01	61	3540.3062	0.6774902E+01	113	3138.2542	0.1009571E+02
10	3934.6260	0.5004883E+01	62	3532.5742	0.6787109E+01	114	3130.5225	0.967941E+01
11	3926.8940	0.5432129E+01	63	3524.8425	0.6378174E+01	115	3122.7908	0.9418713E+01
12	3919.1624	0.4949951E+01	64	3517.1108	0.6817627E+01	116	3115.0591	0.9478760E+01
13	3911.4307	0.4577637E+01	65	3509.3789	0.7397461E+01	117	3107.3271	0.9674072E+01
14	3903.6987	0.5407715E+01	66	3501.6472	0.7678223E+01	118	3099.5955	0.9722900E+01
15	3895.9670	0.4882813E+01	67	3493.9155	0.6805420E+01	119	3091.8638	0.1011353E+02
16	3888.2354	0.4248047E+01	68	3486.1838	0.6811523E+01	120	3084.1318	0.9965928E+01
17	3880.5037	0.4016113E+01	69	3478.4519	0.7617188E+01	121	3076.4001	0.9448242E+01
18	3872.7717	0.5285645E+01	70	3470.7202	0.7647705E+01	122	3068.6685	0.9271240E+01
19	3865.0400	0.5303955E+01	71	3462.9885	0.7910156E+01	123	3060.9365	0.9179688E+01
20	3857.3083	0.5497061E+01	72	3455.2566	0.8099365E+01	124	3053.2048	0.9527588E+01
21	3849.5764	0.5035400E+01	73	3447.5249	0.8471680E+01	125	3045.4731	0.1006470E+02
22	3841.8447	0.5145264E+01	74	3439.7932	0.8782959E+01	126	3037.7412	0.1034546E+02
23	3834.1130	0.4583720E+01	75	3432.0613	0.9216309E+01	127	3030.0095	0.1190185E+02
24	3826.3811	0.5358887E+01	76	3424.3296	0.9290045E+01	128	3022.2778	0.1181030E+02
25	3818.6494	0.5426025E+01	77	3416.5979	0.9759521E+01	129	3014.5459	0.1228027E+02
26	3810.9177	0.5578613E+01	78	3408.8660	0.1024780E+02	130	3006.8142	0.1345215E+02
27	3803.1858	0.5816650E+01	79	3401.1343	0.1047363E+02	131	2999.0825	0.1546631E+02
28	3795.4541	0.5694580E+01	80	3393.4026	0.1082768E+02	132	2991.3508	0.1627277E+02
29	3787.7224	0.5725098E+01	81	3385.6707	0.1181641E+02	133	2983.6189	0.1768188E+02
30	3779.9905	0.5474854E+01	82	3377.9390	0.1211548E+02	134	2975.8872	0.1817627E+02
31	3772.2588	0.5383301E+01	83	3370.2073	0.1357427E+02	135	2968.1555	0.2015381E+02
32	3764.5271	0.5627441E+01	84	3362.4756	0.1435547E+02	136	2960.4236	0.2297589E+02
33	3756.7954	0.5401611E+01	85	3354.7437	0.1493530E+02	137	2952.6919	0.2398652E+02
34	3749.0635	0.5804443E+01	86	3347.0120	0.1656494E+02	138	2944.9602	0.2738540E+02
35	3741.3320	0.6060791E+01	87	3339.2803	0.1806020E+02	139	2937.2283	0.2919139E+02
36	3733.6003	0.6256104E+01	88	3331.5483	0.1947021E+02	140	2929.4966	0.3177937E+02
37	3725.8684	0.6121826E+01	89	3323.8167	0.2105713E+02	141	2921.7649	0.3437937E+02
38	3718.1367	0.6616211E+01	90	3316.0850	0.2237891E+02	142	2914.0330	0.3709373E+02
39	3710.4050	0.6610107E+01	91	3308.3530	0.2378210E+02	143	2906.3013	0.4014.00
40	3702.6731	0.6561279E+01	92	3300.6213	0.2528210E+02	144	2898.5696	0.4312136E+02
41	3694.9414	0.6231659E+01	93	3292.8896	0.2692951E+02	145	2890.8376	0.4676636E+02
42	3687.2097	0.6451416E+01	94	3285.1577	0.2823029E+02	146	2883.1060	0.5162111E+02
43	3679.4778	0.6365947E+01	95	3277.4260	0.2920374E+02	147	2875.3743	0.5681071E+02
44	3671.7461	0.6127930E+01	96	3269.6943	0.3092725E+02	148	2867.6426	0.6274756E+02
45	3664.0144	0.6134033E+01	97	3261.9624	0.3209285E+02	149	2859.9106	0.6986308E+02
46	3656.2825	0.6357760E+01	98	3254.2307	0.3392791E+02	150	2852.1790	0.7731567E+02
47	3648.5508	0.6530760E+01	99	3246.4990	0.3577.4260	151	2844.4472	0.8590576E+02
48	3640.8191	0.6498691E+01	100	3238.7673	0.3731.0354	152	2836.7153	0.9579373E+02
49	3633.0872	0.6567383E+01	101	3231.0337	0.3885.1577	153	2828.9835	0.1079712E+02
50	3625.3555	0.6643210E+01	102	3223.3020	0.4016113E+01	154	2821.2520	0.9997549E+01
51	3617.6238	0.64414.55E+01	103	3215.5720	0.4246.4990	155	2813.5200	0.9609346E+01
52	3609.8921		104	3207.8401	0.4444.299E+02	156	2805.7893	0.105.00

FTIS FILE NUMBER : 20

NON-NORMALIZED (PURE) SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
157	2798.0566	0.9252930E+01	209	2396.0046	0.8959961F+01	261	1993.9529	0.1051025E+02			
158	2790.3247	0.9075928E+01	210	2388.2129	0.8886719E+01	262	1986.2212	0.1051025E+02			
159	2782.5930	0.8898926E+01	211	2380.5413	0.8898926E+01	263	1978.4895	0.1033325E+02			
160	2774.8613	0.8862305E+01	212	2372.8006	0.8878271E+01	264	1970.7576	0.1012573E+02			
161	2767.1294	0.8697510E+01	213	2365.0776	0.9020996F+01	265	1963.0259	0.1003418E+02			
162	2759.3977	0.8666992E+01	214	2357.3459	0.8929443E+01	266	1955.2942	0.1007600E+02			
163	2751.6660	0.8630371E+01	215	2349.6143	0.8935547E+01	267	1947.5623	0.1002808E+02			
164	2743.9343	0.8557129E+01	216	2341.8823	0.9020996E+01	268	1939.8306	0.9997559E+01			
165	2736.2024	0.8551025E+01	217	2334.1506	0.9002686E+01	269	1932.0989	0.9985352E+01			
166	2728.4707	0.8563232E+01	218	2326.4189	0.8990479E+01	270	1924.3669	0.9936523E+01			
167	2720.7390	0.9416748E+01	219	2318.6870	0.8947754E+01	271	1916.6353	0.9979248E+01			
168	2713.0071	0.8270264E+01	220	2310.9553	0.8941650F+01	272	1908.9036	0.9997559E+01			
169	2705.2754	0.8325195E+01	221	2303.2236	0.8966064E+01	273	1901.1716	0.1004028E+02			
170	2697.5437	0.8428955E+01	222	2295.4917	0.8984375E+01	274	1893.4369	0.1005249E+02			
171	2689.8119	0.8343506E+01	223	2287.7600	0.8990479E+01	275	1885.7043	0.1001587E+02			
172	2682.0801	0.8264160E+01	224	2280.0283	0.8917236E+01	276	1877.9766	0.1011963E+02			
173	2674.3484	0.8386230E+01	225	2272.2964	0.8984375E+01	277	1870.2446	0.1010132E+02			
174	2666.6165	0.8331299E+01	226	2264.5647	0.9002686E+01	278	1862.5129	0.1014404E+02			
175	2658.8848	0.8264160E+01	227	2256.8330	0.8996592E+01	279	1854.7812	0.1024170E+02			
176	2651.1531	0.8312988E+01	228	2249.1013	0.8966064E+01	280	1847.0493	0.1029053E+02			
177	2643.4211	0.8203125E+01	229	2241.3694	0.9027100F+01	281	1839.3176	0.1050415E+02			
178	2635.6895	0.8239746E+01	230	2233.6377	0.8996592E+01	282	1831.5859	0.1050415E+02			
179	2627.9578	0.8190918E+01	231	2225.9060	0.9002686E+01	283	1823.8540	0.1060181E+02			
180	2620.2261	0.8251953E+01	232	2218.1741	0.9106445E+01	284	1816.1223	0.1063843E+02			
181	2612.4941	0.8306885E+01	233	2210.4424	0.9155273F+01	285	1808.3906	0.1063843E+02			
182	2604.7625	0.8258057E+01	234	2202.7107	0.9173584E+01	286	1800.6587	0.1062622E+02			
183	2597.0308	0.8227539E+01	235	2194.9798	0.9197998F+01	287	1792.9270	0.1070557E+02			
184	2589.2988	0.8203125E+01	236	2187.2471	0.9246826E+01	288	1785.1953	0.1085205E+02			
185	2581.5671	0.8154297E+01	237	2179.5154	0.9277344E+01	289	1777.4634	0.1086426E+02			
186	2573.8354	0.8337402E+01	238	2171.7834	0.9240723F+01	290	1769.7317	0.1119385E+02			
187	2566.1035	0.8325195E+01	239	2164.0518	0.9350586E+01	291	1762.0000	0.1163940E+02			
188	2558.3718	0.8294678E+01	240	2156.3201	0.9411621F+01	292	1754.2683	0.1311177E+02			
189	2550.6401	0.8374023E+01	241	2148.5981	0.9490967F+01	293	1746.5364	0.2147827E+02			
190	2542.9082	0.8374023E+01	242	2140.8564	0.9594727F+01	294	1738.8047	0.3020020F+02			
191	2535.1765	0.8367920E+01	243	2133.1249	0.9600801E+01	295	1731.0730	0.2465210F+02			
192	2527.4448	0.8404541F+01	244	2125.3531	0.9527588F+01	296	1723.3411	0.2027588F+02			
193	2519.7129	0.8508301F+01	245	2117.6611	0.9527588E+01	297	1715.6094	0.1746216E+02			
194	2511.9812	0.8465576E+01	246	2109.9294	0.9643555E+01	298	1707.8777	0.1563110F+02			
195	2504.2495	0.8477830E+01	247	2102.1978	0.9661865F+01	299	1700.1458	0.1473999E+02			
196	2496.5178	0.8471680E+01	248	2094.4659	0.9766143E+01	300	1692.4141	0.1410522E+02			
197	2488.7859	0.8563232E+01	249	2086.7341	0.9924316F+01	301	1684.6824	0.1376343E+02			
198	2481.0542	0.8532715E+01	250	2079.0024	0.9960939F+01	302	1676.9504	0.1346436F+02			
199	2473.3225	0.8593750E+01	251	2071.2705	0.1005859E+02	303	1669.2187	0.1376343E+02			
200	2465.5906	0.8721924E+01	252	2063.5388	0.1015625E+02	304	1661.4871	0.1388505E+02			
201	2457.8589	0.8709717E+01	253	2055.8071	0.1021181F+02	305	1653.7551	0.1413574E+02			
202	2450.1272	0.8752441E+01	254	2048.0752	0.1024780E+02	306	1646.0234	0.1450195E+02			
203	2442.3953	0.8782959E+01	255	2040.3435	0.1024780E+02	307	1638.2917	0.1450195E+02			
204	2434.6636	0.8782959E+01	256	2032.6119	0.1032104F+02	308	1630.5601	0.1376294E+02			
205	2426.9319	0.8868403E+01	257	2024.8799	0.1033936F+02	309	1622.8281	0.1299438F+02			
206	2419.2000	0.8905029E+01	258	2017.1487	0.1041870E+02	310	1615.0964	0.1297231F+02			
207	2411.4683	0.8874512E+01	259	2009.4165	0.1054688E+02	311	1607.3647	0.1279297E+02			
208	2403.7366	0.8850098E+01	260	2001.6848	0.1051025E+02	312	1599.6328	0.1241178F+02			

FTIS FILE NUMBER : 20

NON-NORMALIZED (P.O.E) SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
313	1591.2011	0.1275024E+02	365	1189.8494	0.2153320E+02	417	787.7974	0.3029785E+02
314	1584.1694	0.1282959E+02	366	1182.1174	0.2231201E+02	418	780.0657	0.3051754E+02
315	1576.4375	0.1287842E+02	367	1174.3857	0.2250793E+02	419	772.3340	0.3061511E+02
316	1568.7058	0.1311644E+02	368	1166.6541	0.2278930E+02	420	764.6023	0.2979124E+02
317	1560.9741	0.1316528E+02	369	1158.9221	0.2297058E+02	421	756.8704	0.2932739E+02
318	1553.2422	0.1378734E+02	370	1151.1904	0.2325439E+02	422	749.1387	0.2949829E+02
319	1545.5105	0.1417236E+02	371	1143.4587	0.2380676E+02	423	741.4070	0.3012085E+02
320	1537.7788	0.1473999E+02	372	1135.7271	0.2433960E+02	424	733.6750	0.3139038E+02
321	1530.0469	0.1526489E+02	373	1127.9951	0.2478332E+02	425	725.9434	0.3236694E+02
322	1522.3152	0.1558228E+02	374	1120.2634	0.2520263E+02	426	718.2117	0.3186644E+02
323	1514.5835	0.1603394E+02	375	1112.5317	0.2563517E+02	427	710.4797	0.3145752E+02
324	1506.8518	0.1609342E+02	376	1104.7998	0.2614318E+02	428	702.7480	0.3094482E+02
325	1499.1199	0.1618652E+02	377	1097.0681	0.2665893E+02	429	695.0164	0.3073120E+02
326	1491.3882	0.1655884E+02	378	1089.3364	0.2709594E+02	430	687.2844	0.3013916E+02
327	1483.6565	0.1749268E+02	379	1081.6045	0.2754633E+02	431	679.5527	0.2991333E+02
328	1475.9246	0.1904297E+02	380	1073.8728	0.2815307E+02	432	671.8210	0.2913818E+02
329	1468.1929	0.2283934E+02	381	1066.1411	0.2873718E+02	433	664.0891	0.2799072E+02
330	1460.4612	0.2564087E+02	382	1058.4092	0.2909118E+02	434	656.3574	0.2676392E+02
331	1452.7292	0.2924194E+02	383	1050.6775	0.2930962E+02	435	648.6257	0.2594507E+02
332	1444.9976	0.3341675E+02	384	1042.9458	0.2921936E+02	436	640.8940	0.2517908E+02
333	1437.2659	0.4090576E+02	385	1035.2139	0.2900451E+02	437	633.1621	0.2410596E+02
334	1429.5339	0.4971924E+02	386	1027.4822	0.2877014E+02	438	625.4304	0.2301093E+02
335	1421.8022	0.5692139E+02	387	1019.7505	0.2849914E+02	439	617.6987	0.2182715E+02
336	1414.0706	0.5755155E+02	388	1012.0188	0.2816711E+02	440	609.9668	0.2059049E+02
337	1406.3386	0.5264493E+02	389	1004.2869	0.2781372E+02	441	602.2351	0.2040771E+02
338	1398.6069	0.4588623E+02	390	996.5552	0.2688599E+02	442	594.5034	0.2049561E+02
339	1390.8752	0.4008173E+02	391	988.8235	0.2730894E+02	443	586.7715	0.2055390E+02
340	1383.1436	0.3416275E+02	392	981.0916	0.2692870E+02	444	579.0398	0.2022131E+02
341	1375.4116	0.3010254E+02	393	973.3599	0.2720153E+02	445	571.3081	0.2038815E+02
342	1367.6799	0.2743331E+02	394	965.6282	0.2816101E+02	446	563.5762	0.2096130E+02
343	1359.9482	0.2505687E+02	395	957.8962	0.2741393E+02	447	555.8445	0.2092163E+02
344	1352.2163	0.2400002E+02	396	950.1646	0.2632812E+02	448	548.1128	0.2037500E+02
345	1344.4846	0.2250366E+02	397	942.4329	0.2541503E+02	449	540.3809	0.2059590E+02
346	1336.7529	0.2110566E+02	398	934.7009	0.2508825E+02	450	532.6492	0.2092263E+02
347	1329.0210	0.1996440E+02	399	926.9692	0.2544494E+02	451	524.9175	0.2016217E+02
348	1321.2893	0.1952447E+02	400	919.2375	0.2549072E+02	452	517.1859	0.2053161E+02
349	1313.5576	0.1976527E+02	401	911.5054	0.2606453E+02	453	509.4539	0.2053161E+02
350	1305.8257	0.1959220E+02	402	903.7739	0.2542173E+02	454	501.7222	0.2034234E+02
351	1298.0940	0.1873463E+02	403	896.0422	0.2479515E+02	455	493.9905	0.2059261E+02
352	1290.3624	0.1792067E+02	404	888.3105	0.2441833E+02	456	486.2585	0.2034444E+02
353	1282.6304	0.1682325E+02	405	880.5786	0.2412780E+02	457	478.5263	0.2037104E+02
354	1274.8987	0.1573777E+02	406	872.8459	0.2398437E+02	458	470.7952	0.2059027E+02
355	1267.1670	0.1493594E+02	407	865.1152	0.2381469E+02	459	463.0632	0.2065557E+02
356	1259.4353	0.2022705E+02	408	857.3833	0.2346690E+02	460	455.3315	0.2059027E+02
357	1251.7034	0.2006451E+02	409	849.6516	0.2353908E+02	461	447.6000	0.2037104E+02
358	1243.9717	0.2111711E+02	410	841.9199	0.2342712E+02	462	439.8679	0.2037104E+02
359	1236.2400	0.2133242E+02	411	834.1880	0.2339233E+02	463	432.1362	0.2037353E+02
360	1228.5081	0.2074770E+02	412	826.4563	0.2265381E+02	464	424.4045	0.2037353E+02
361	1220.7764	0.2066940E+02	413	818.7246	0.2176880E+02	465	416.6726	0.2037353E+02
362	1213.0447	0.2033155E+02	414	810.9927	0.2129803E+02	466	408.9409	0.2037353E+02
363	1205.3127	0.1983408E+02	415	803.2610	0.2062134E+02	467	401.2092	0.2037353E+02
364	1197.5811	0.2055078E+02	416	795.5293	0.2040016E+02	468	393.4775	0.2037353E+02

24	468	0	1	4004	13902	393	31295
4004.2119	393.4773	-7.7318					



FTIS FILE NUMBER : 24

NON-NORMALIZED (PURE) SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
1	4004.2114	0.3995361E+02	54	3602.1602	0.4986572E+01	105	3200.1084	0.1311035E+02
2	3996.4800	0.3995361E+02	55	3594.4285	0.4919434E+01	106	3192.3767	0.1226807E+02
3	3989.7483	0.3995361E+02	56	3586.6968	0.4888916E+01	107	3184.6448	0.1194124E+02
4	3981.0166	0.5194092E+01	57	3578.9648	0.4833984E+01	108	3176.9131	0.1173061E+02
5	3973.2847	0.5126953E+01	58	3571.2332	0.5086436E+01	109	3169.1814	0.1086426E+02
6	3965.5530	0.5438232E+01	59	3563.5015	0.5267334E+01	110	3161.4495	0.1055767E+02
7	3957.8213	0.4742432E+01	60	3555.7695	0.4943848E+01	111	3153.7178	0.9973142E+01
8	3950.0894	0.4931641E+01	61	3548.0378	0.4919434E+01	112	3145.9961	0.9716707E+01
9	3942.3577	0.4321289E+01	62	3540.3062	0.5340576E+01	113	3138.2542	0.9295654E+01
10	3934.6260	0.3784180E+01	63	3532.5742	0.5413818E+01	114	3130.5225	0.9405518E+01
11	3926.8940	0.3717041E+01	64	3524.8425	0.4949951E+01	115	3122.7908	0.9143086E+01
12	3919.1624	0.4669218E+01	65	3517.1108	0.5535889E+01	116	3115.0591	0.8929443E+01
13	3911.4307	0.4675293E+01	66	3509.3789	0.6317139E+01	117	3107.3271	0.8495906E+01
14	3903.6987	0.4705911E+01	67	3501.6472	0.6365967E+01	118	3099.5955	0.8092268E+01
15	3895.9670	0.3767285E+01	68	3493.9155	0.5786133E+01	119	3091.8638	0.9112549E+01
16	3888.2354	0.3668213E+01	69	3486.1838	0.5560303E+01	120	3084.1318	0.5964209E+01
17	3880.5037	0.3649902E+01	70	3478.4519	0.6457416E+01	121	3076.4001	0.8338379E+01
18	3872.7717	0.4577637E+01	71	3470.7202	0.7257080E+01	122	3068.6685	0.8764648E+01
19	3865.0400	0.4498291E+01	72	3462.9985	0.76360840E+01	123	3060.9365	0.8551025E+01
20	3857.3083	0.4272461E+01	73	3455.2566	0.7361602E+01	124	3053.2048	0.8526611E+01
21	3849.5764	0.4028320E+01	74	3447.5249	0.7646050E+01	125	3045.4731	0.8996562E+01
22	3841.8447	0.4174905E+01	75	3439.7932	0.7995605E+01	126	3037.7412	0.8994144E+01
23	3834.1130	0.4736329E+01	76	3432.0613	0.8374023E+01	127	3030.0095	0.5649648E+01
24	3826.3811	0.4693609E+01	77	3424.3295	0.8747919E+01	128	3022.2778	0.1030884E+02
25	3818.6494	0.4638672E+01	78	3416.5979	0.8978271E+01	129	3014.5459	0.1130371E+02
26	3810.9177	0.5316162E+01	79	3408.8660	0.9179688E+01	130	3006.8142	0.1132813E+02
27	3803.1858	0.5603027E+01	80	3401.1343	0.9594727E+01	131	2999.0825	0.1083984E+02
28	3795.4541	0.4943848E+01	81	3393.4026	0.1034546E+02	132	2991.3508	0.1165141E+02
29	3787.7224	0.4840088E+01	82	3385.6707	0.1068115E+02	133	2983.6189	0.1274414E+02
30	3779.9905	0.5029297E+01	83	3377.9390	0.1114502E+02	134	2975.8872	0.1474430E+02
31	3772.2588	0.5187989E+01	84	3370.2073	0.1155396E+02	135	2968.1555	0.1764148E+02
32	3764.5271	0.4937764E+01	85	3362.4756	0.1204834E+02	136	2960.4236	0.1948145E+02
33	3756.7954	0.5541992E+01	86	3354.7437	0.1311646E+02	137	2952.6919	0.170011E+02
34	3749.0635	0.5377197E+01	87	3347.0120	0.1372070E+02	138	2944.9602	0.1981505E+02
35	3741.3320	0.5346680E+01	88	3339.2803	0.1446609E+02	139	2937.2283	0.2260723E+02
36	3733.6003	0.5450439E+01	89	3331.5483	0.1608274E+02	140	2929.4966	0.2311225E+02
37	3725.8684	0.5627441E+01	90	3323.8167	0.1751317E+02	141	2921.7649	0.3761045E+02
38	3718.1367	0.5419922E+01	91	3316.0850	0.1877441E+02	142	2914.0330	0.2186248E+02
39	3710.4050	0.5493164E+01	92	3308.3533	0.2014771E+02	143	2906.3013	0.1951904E+02
40	3702.6731	0.5316162E+01	93	3300.6213	0.2187500E+02	144	2898.5696	0.1735830E+02
41	3694.9414	0.5413818E+01	94	3292.8896	0.2223617E+02	145	2890.8376	0.1564344E+02
42	3687.2097	0.5469314E+01	95	3285.1577	0.2273760E+02	146	2883.1060	0.1575745E+02
43	3679.4778	0.5364990E+01	96	3277.4260	0.2216757E+02	147	2875.3743	0.1611264E+02
44	3671.7461	0.5521682E+01	97	3269.6943	0.2144715E+02	148	2867.6426	0.1463767E+02
45	3664.0144	0.5187989E+01	98	3261.9626	0.2033691E+02	149	2859.9106	0.1463767E+02
46	3656.2825	0.5297852E+01	99	3254.2307	0.1898804E+02	150	2852.1790	0.1694711E+02
47	3648.5508	0.5352783E+01	100	3246.4990	0.1774434E+02	151	2844.4473	0.1551514E+02
48	3640.8191	0.5047607E+01	101	3238.7673	0.1715699E+02	152	2836.7153	0.1243286E+02
49	3633.0872	0.5169678E+01	102	3231.0354	0.1611733E+02	153	2828.9836	0.1021729E+02
50	3625.3555	0.5114746E+01	103	3223.3037	0.1538696E+02	154	2821.2520	0.1375000E+02
51	3617.6238	0.5157471E+01	104	3215.5720	0.1420298E+02	155	2813.5200	0.8935567E+01
52	3609.8921	0.5017090E+01	104	3207.8401	0.1360474E+02	156	2805.7883	0.9401270E+01

FTIS FILE NUMBER : 24

NON-NORMALIZED (PURE) SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
151	2798.0566	0.8483887E+01	209	2396.0046	0.8251953E+01	261	1903.9529	0.9729004E+01
152	2790.3247	0.8422952E+01	210	2389.2729	0.8229193E+01	262	1906.2212	0.9625244E+01
153	2782.5930	0.8306885E+01	211	2380.5413	0.8202299E+01	263	1978.4895	0.9527588E+01
154	2774.8613	0.8190918E+01	212	2372.3096	0.8306885E+01	264	1970.7576	0.9417755E+01
155	2767.1294	0.8135986E+01	213	2365.0776	0.8300781E+01	265	1963.0259	0.930.58E+01
156	2759.3977	0.8068348E+01	214	2357.3459	0.8215332E+01	266	1955.7942	0.9295654E+01
157	2751.6660	0.8007813E+01	215	2349.6143	0.8239746E+01	267	1947.5623	0.930178E+01
158	2743.9343	0.7916260E+01	216	2341.8823	0.8264160E+01	268	1939.8306	0.9283447E+01
159	2736.2024	0.7879639E+01	217	2334.1506	0.835195E+01	269	1932.0989	0.9307861E+01
160	2728.4707	0.7899560E+01	218	2326.4189	0.8276367E+01	270	1924.3669	0.934482E+01
161	2720.7390	0.7922363E+01	219	2318.6870	0.8227539E+01	271	1916.6353	0.9393311E+01
162	2713.0071	0.7714844E+01	220	2310.9553	0.8209229E+01	272	1908.9036	0.9490967E+01
163	2705.2754	0.7647705E+01	221	2303.2236	0.8245850E+01	273	1901.1716	0.9478760E+01
164	2697.5437	0.7659912E+01	222	2295.4917	0.8227539E+01	274	1893.4399	0.9509277E+01
165	2689.8118	0.7666016E+01	223	2287.7600	0.8264160E+01	275	1885.7083	0.9598623E+01
166	2682.0801	0.7550049E+01	224	2280.0283	0.8255850E+01	276	1877.9766	0.9637451E+01
167	2674.3484	0.7666016E+01	225	2272.2964	0.8276367E+01	277	1870.2446	0.9649658E+01
168	2666.6165	0.7501221E+01	226	2264.5647	0.8325195E+01	278	1862.5129	0.9735107E+01
169	2658.8848	0.7348633E+01	227	2256.8330	0.8276367E+01	279	1854.7812	0.980246E+01
170	2651.1531	0.7464600E+01	228	2249.1013	0.8239746E+01	280	1847.0493	0.9863281E+01
171	2643.4211	0.7427979E+01	229	2241.3694	0.8367920E+01	281	1839.3176	0.9957041E+01
172	2635.6895	0.7458494E+01	230	2233.6377	0.8312988E+01	282	1831.5859	0.9957041E+01
173	2627.9578	0.7336426E+01	231	2225.9060	0.8398438E+01	283	1823.8540	0.1018677E+02
174	2620.2261	0.7373047E+01	232	2218.1741	0.8453369E+01	284	1816.1223	0.1015015E+02
175	2612.4941	0.7379150E+01	233	2210.4424	0.8477703E+01	285	1808.3906	0.1004639E+02
176	2604.7625	0.7373047E+01	234	2202.7107	0.8453369E+01	286	1800.6587	0.1010742E+02
177	2597.0308	0.7287508E+01	235	2194.9788	0.8447266E+01	287	1792.9270	0.1013184E+02
178	2589.2939	0.7208040E+01	236	2187.2471	0.8587645E+01	288	1785.1953	0.1026001E+02
179	2581.5671	0.7160840E+01	237	2179.5154	0.8526411E+01	289	1777.4634	0.1025391E+02
180	2573.8354	0.7415771E+01	238	2171.7834	0.8599854E+01	290	1769.7317	0.1044312E+02
181	2566.1035	0.7400668E+01	239	2164.0518	0.860889E+01	291	1762.0000	0.1098633E+02
182	2558.3718	0.7427979E+01	240	2156.3201	0.860889E+01	292	1754.2683	0.1279807E+02
183	2550.6401	0.7464600E+01	241	2148.5884	0.8709717E+01	293	1746.5364	0.1215087E+02
184	2542.9082	0.7562255E+01	242	2140.8564	0.8740234E+01	294	1738.8047	0.1304138E+02
185	2535.1765	0.7531738E+01	243	2133.1249	0.8740234E+01	295	1731.0730	0.1304138E+02
186	2527.4448	0.7623291E+01	244	2125.3931	0.8703613E+01	296	1723.3411	0.1304138E+02
187	2519.7129	0.7635498E+01	245	2117.6611	0.8819500E+01	297	1715.6094	0.1304138E+02
188	2511.9812	0.7672705E+01	246	2109.9294	0.8819500E+01	298	1707.8777	0.1304138E+02
189	2504.2495	0.7727051E+01	247	2102.1974	0.8943904E+01	299	1700.1458	0.1304138E+02
190	2496.5178	0.7781987E+01	248	2094.4659	0.8943904E+01	300	1692.4141	0.1304138E+02
191	2488.7859	0.7806336E+01	249	2086.7341	0.9082031E+01	301	1684.6824	0.1304138E+02
192	2481.0542	0.7873535E+01	250	2079.0024	0.9082031E+01	302	1676.9504	0.1304138E+02
193	2473.3225	0.7928467E+01	251	2071.2705	0.9094219E+01	303	1669.2187	0.1304138E+02
194	2465.5906	0.7966777E+01	252	2063.5388	0.9149170E+01	304	1661.4871	0.1304138E+02
195	2457.8589	0.8087154E+01	253	2055.8071	0.9246826E+01	305	1653.7551	0.1304138E+02
196	2450.1272	0.8062764E+01	254	2048.0752	0.9301758E+01	306	1646.0234	0.1304138E+02
197	2442.3953	0.8007813E+01	255	2040.3435	0.9332775E+01	307	1638.2917	0.1304138E+02
198	2434.6636	0.8007813E+01	256	2032.6118	0.9362779E+01	308	1630.5601	0.1304138E+02
199	2426.9319	0.8148193E+01	257	2024.8799	0.9466553E+01	309	1622.8284	0.1304138E+02
200	2419.2000	0.8197021E+01	258	2017.1482	0.9509277E+01	310	1615.0964	0.1304138E+02
201	2411.4683	0.8111572E+01	259	2009.4165	0.9588623E+01	311	1607.3647	0.1304138E+02
202	2403.7366	0.8233643E+01	260	2001.6848	0.9637451E+01	312	1599.6329	0.1304138E+02

FITS FILE NUMBER : 24

NON-NORMALIZED (PURE) SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
313	1591.9011	0.1193237E+02	365	1189.8494	0.2113037E+02	417	797.7974	0.2842407E+02			
314	1594.1694	0.1212769E+02	366	1182.1174	0.2290039E+02	418	780.0657	0.2971704E+02			
315	1576.4375	0.1271093E+02	367	1174.3857	0.2499390E+02	419	772.3340	0.2898560E+02			
316	1569.7058	0.1253052E+02	368	1166.6541	0.2781372E+02	420	764.6023	0.2811597E+02			
317	1560.9741	0.1287842E+02	369	1158.9221	0.2068140E+02	421	756.8704	0.2759789E+02			
318	1553.2422	0.1345825E+02	370	1151.1904	0.3272705E+02	422	749.1387	0.2775269E+02			
319	1545.5105	0.1400757E+02	371	1143.4587	0.3836670E+02	423	741.4070	0.2857173E+02			
320	1537.7788	0.1461192E+02	372	1135.7271	0.4386597E+02	424	733.6750	0.3009644E+02			
321	1530.0469	0.1518555E+02	373	1127.9951	0.4793701E+02	425	725.9434	0.3110352E+02			
322	1522.3157	0.1563721E+02	374	1120.2634	0.5162964E+02	426	718.2117	0.3048706E+02			
323	1514.5835	0.1617432E+02	375	1112.5317	0.5573120E+02	427	710.4797	0.3008423E+02			
324	1506.8518	0.1614380E+02	376	1104.7998	0.6076660E+02	428	702.7480	0.2968140E+02			
325	1499.1199	0.1629639E+02	377	1097.0691	0.6571655E+02	429	695.0164	0.2932129E+02			
326	1491.3882	0.1657104E+02	378	1089.3364	0.6986594E+02	430	687.2844	0.2894897E+02			
327	1483.6565	0.1735340E+02	379	1081.6045	0.7434082E+02	431	679.5527	0.2858042E+02			
328	1475.9246	0.1931152E+02	380	1073.8728	0.8053589E+02	432	671.8210	0.2779231E+02			
329	1468.1929	0.2328491E+02	381	1066.1411	0.8615723E+02	433	664.0891	0.2661133E+02			
330	1460.4612	0.2615356E+02	382	1058.4092	0.8936768E+02	434	656.3574	0.2545156E+02			
331	1452.7292	0.2967529E+02	383	1050.6775	0.9150391E+02	435	648.6257	0.2471875E+02			
332	1444.9975	0.3386230E+02	384	1042.9458	0.9044800E+02	436	640.8940	0.2519551E+02			
333	1437.2659	0.4130859E+02	385	1035.2139	0.8842131E+02	437	633.1621	0.2729614E+02			
334	1429.5339	0.4979248E+02	386	1027.4822	0.8633423E+02	438	625.4304	0.2920532E+02			
335	1421.8022	0.5675659E+02	387	1019.7505	0.8330007E+02	439	617.6987	0.3085593E+02			
336	1414.0706	0.6732677E+02	388	1012.0188	0.8001094E+02	440	609.9668	0.3249839E+02			
337	1406.3386	0.8231934E+02	389	1004.2869	0.7644653E+02	441	602.2351	0.3407170E+02			
338	1398.6069	0.9536133E+02	390	996.5552	0.7481079E+02	442	594.5034	0.3586914E+02			
339	1390.8752	0.9526377E+02	391	988.9235	0.7105713E+02	443	586.7715	0.3786814E+02			
340	1383.1436	0.9469707E+02	392	981.0516	0.6722412E+02	444	579.0398	0.4017061E+02			
341	1375.4116	0.9351818E+02	393	973.3599	0.6988525E+02	445	571.3091	0.4275007E+02			
342	1367.6799	0.9642712E+02	394	965.6282	0.7915649E+02	446	563.5762	0.4531396E+02			
343	1359.9432	0.9700931E+02	395	957.8962	0.7140503E+02	447	555.8445	0.4791003E+02			
344	1352.2163	0.9299805E+02	396	950.1646	0.6033325E+02	448	548.1129	0.5056934E+02			
345	1344.4846	0.9143555E+02	397	942.4329	0.5150757E+02	449	540.3809	0.5307271E+02			
346	1336.7579	0.9040136E+02	398	934.7009	0.5769653E+02	450	532.6492	0.5585027E+02			
347	1329.0210	0.8772599E+02	399	926.9692	0.5177002E+02	451	524.9175	0.5856030E+02			
348	1321.2893	0.841431E+02	400	919.2375	0.5205889E+02	452	517.1858	0.6177710E+02			
349	1313.5575	0.851807E+02	401	911.5056	0.5769653E+02	453	509.4539	0.6492017E+02			
350	1305.8257	0.8181769E+02	402	903.7739	0.5120239E+02	454	501.7222	0.6772339E+02			
351	1298.0940	0.761475E+02	403	896.0422	0.4466472E+02	455	493.9905	0.7055396E+02			
352	1290.3623	0.7175513E+02	404	888.3105	0.4102173E+02	456	486.2586	0.7305396E+02			
353	1282.6304	0.6762085E+02	405	880.5786	0.3870239E+02	457	478.5269	0.7507525E+02			
354	1274.8987	0.641475E+02	406	872.8469	0.3603337E+02	458	470.7952	0.7657373E+02			
355	1267.1670	0.6136106E+02	407	865.1152	0.3561401E+02	459	463.0632	0.7746461E+02			
356	1259.4353	0.5925040E+02	408	857.3833	0.3410685E+02	460	455.3315	0.7777457E+02			
357	1251.7034	0.5750098E+02	409	849.6516	0.3290405E+02	461	447.5999	0.7843617E+02			
358	1243.9717	0.5610161E+02	410	841.9199	0.3209005E+02	462	439.8679	0.7894717E+02			
359	1236.2400	0.5499092E+02	411	834.1880	0.3114674E+02	463	432.1362	0.7935261E+02			
360	1228.5081	0.5422095E+02	412	826.4563	0.3027954E+02	464	424.4045	0.7965351E+02			
361	1220.7764	0.5363784E+02	413	818.7244	0.2977905E+02	465	416.6726	0.7985305E+02			
362	1213.0447	0.5319583E+02	414	810.9927	0.2911377E+02	466	408.9409	0.7995361E+02			
363	1205.3127	0.5284098E+02	415	803.2610	0.2866821E+02	467	401.2092	0.7995361E+02			
364	1197.5811	0.5260073E+02	416	795.5293	0.2839935E+02	468	393.4775	0.7995361E+02			

28	468	0	1	4004	13902	393	31295
4004.2119	393.773	-7.7318					

FTIS FILE NUMBER : 29

NON-NORMALIZED (PURE) SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
1	4004.2119	0.5039063E+02	53	3602.1602	0.4962158E+01	105	3200.1054	0.1465454E+02
2	3996.4800	0.5039063E+02	54	3594.4285	0.4809570E+01	106	3192.3767	0.1389771E+02
3	3988.7483	0.5039063E+02	55	3586.6668	0.5059814E+01	107	3184.6448	0.1307373E+02
4	3981.0166	0.5039063E+02	56	3578.9048	0.5017090E+01	108	3176.9131	0.1249400E+02
5	3973.2847	0.5039063E+02	57	3571.2332	0.5139160E+01	109	3169.1814	0.1190186E+02
6	3965.5530	0.5039063E+02	58	3563.5015	0.5194092E+01	110	3161.4495	0.1135864E+02
7	3957.8213	0.5039063E+02	59	3555.7695	0.5157471E+01	111	3153.7178	0.1082153E+02
8	3950.0894	0.5039063E+02	60	3548.0379	0.5169678E+01	112	3145.9861	0.1024170E+02
9	3942.3577	0.5039063E+02	61	3540.3062	0.5267334E+01	113	3138.2542	0.1004470E+02
10	3934.6260	0.5039063E+02	62	3532.5742	0.5145264E+01	114	3130.5225	0.1029663E+02
11	3926.8940	0.5039063E+02	63	3524.8425	0.5120850E+01	115	3122.7908	0.9906006E+01
12	3919.1624	0.5039063E+02	64	3517.1108	0.5125098E+01	116	3115.0591	0.9698486E+01
13	3911.4307	0.5039063E+02	65	3509.3789	0.6384277E+01	117	3107.3271	0.9564209E+01
14	3903.6987	0.5039063E+02	66	3501.6472	0.6298828E+01	118	3099.5955	0.9729004E+01
15	3895.9670	0.5039063E+01	67	3493.9155	0.5450439E+01	119	3091.8638	0.9475498E+01
16	3888.2354	0.3472900E+01	68	3486.1838	0.5627441E+01	120	3084.1316	0.1021729E+02
17	3880.5037	0.2954102E+01	69	3478.4519	0.6341553E+01	121	3076.4001	0.1003418E+02
18	3872.7717	0.3527832E+01	70	3470.7202	0.7055664E+01	122	3068.6685	0.9399414E+01
19	3865.0400	0.4235840E+01	71	3462.9895	0.7391357E+01	123	3060.9365	0.9252937E+01
20	3857.3083	0.3693730E+01	72	3455.2566	0.7965088E+01	124	3053.2048	0.9283447E+01
21	3849.5764	0.3417969E+01	73	3447.5249	0.8276367E+01	125	3045.4731	0.8527588E+01
22	3841.8447	0.3283691E+01	74	3439.7932	0.8502197E+01	126	3037.7412	0.1019897E+02
23	3834.1130	0.3663.13E+01	75	3432.0613	0.8795166E+01	127	3030.0095	0.1122437E+02
24	3826.3811	0.40344.4E+01	76	3424.3296	0.9155273E+01	128	3022.2778	0.1040649E+02
25	3818.6494	0.4461670E+01	77	3416.5979	0.9771729E+01	129	3014.5459	0.1250610E+02
26	3810.9177	0.4982813E+01	78	3408.8660	0.1015015E+02	130	3006.8142	0.1232910E+02
27	3803.1858	0.4809570E+01	79	3401.1343	0.1046143E+02	131	2999.0825	0.1211548E+02
28	3795.4541	0.4459463E+01	80	3393.4026	0.1117554E+02	132	2991.3508	0.1198219E+02
29	3787.7224	0.4449463E+01	81	3385.6707	0.1175537E+02	133	2983.6185	0.1459351E+02
30	3779.9905	0.4144287E+01	82	3377.9390	0.1236572E+02	134	2975.8872	0.1732798E+02
31	3772.2588	0.4206875E+01	83	3370.2073	0.1280063E+02	135	2968.1555	0.2100375E+02
32	3764.5271	0.4754439E+01	84	3362.4756	0.1373901E+02	136	2960.4238	0.2345059E+02
33	3756.7954	0.4779053E+01	85	3354.7437	0.1461792E+02	137	2952.6919	0.2094727E+02
34	3749.0635	0.4797363E+01	86	3347.0120	0.1583822E+02	138	2944.9602	0.2353992E+02
35	3741.3320	0.4633984E+01	87	3339.2803	0.1694986E+02	139	2937.2283	0.2701416E+02
36	3733.6003	0.5242023E+01	88	3331.5483	0.1860982E+02	140	2929.4964	0.2504983E+02
37	3725.8684	0.5230713E+01	89	3323.8167	0.2047729E+02	141	2921.7649	0.2772877E+02
38	3718.1367	0.5407715E+01	90	3316.0850	0.2222220E+02	142	2914.0340	0.2504983E+02
39	3710.4050	0.5529785E+01	91	3308.3530	0.2406614E+02	143	2906.3013	0.2219238E+02
40	3702.6731	0.5749512E+01	92	3300.6213	0.2554932E+02	144	2898.5696	0.1952515E+02
41	3694.9414	0.5682373E+01	93	3292.8896	0.2655029E+02	145	2890.8376	0.1747437E+02
42	3687.2097	0.5084229E+01	94	3285.1577	0.2684937E+02	146	2883.1060	0.1774733E+02
43	3679.4778	0.5541992E+01	95	3277.4260	0.2599887E+02	147	2875.3743	0.1790161E+02
44	3671.7461	0.5645752E+01	96	3269.6943	0.2523193E+02	148	2867.6426	0.1460352E+02
45	3664.0144	0.5371024E+01	97	3261.9624	0.2366933E+02	149	2859.9106	0.1876831E+02
46	3656.2825	0.5053711E+01	98	3254.2307	0.2247925E+02	150	2852.1790	0.1531479E+02
47	3648.5508	0.4693604E+01	99	3246.4990	0.2105713E+02	151	2844.4473	0.1601895E+02
48	3640.8191	0.5230713E+01	100	3238.7673	0.1966513E+02	152	2836.7153	0.1303204E+02
49	3633.0872	0.5102539E+01	101	3231.0354	0.1864014E+02	153	2828.9836	0.1046143E+02
50	3625.3555	0.5242920E+01	102	3223.3037	0.1746215E+02	154	2821.2520	0.9356489E+01
51	3617.6238	0.5169678E+01	103	3215.5720	0.1600342E+02	155	2813.5200	0.8874512E+01
52	3609.8921	0.5133057E+01	104	3207.8401	0.1553335E+02	156	2805.7893	0.8673096E+01

FTIS FILE NUMBER : 29

NON-NORMALIZED (PURE) SPECTRAL DATA

DATA WORD	WAVE NUMBR	AMPLITUDE	DATA WORD	WAVE NUMBR	AMPLITUDE	DATA WORD	WAVE NUMBR	AMPLITUDE
157	2798.0566	0.8386230E+01	209	2396.0046	0.8221436E+01	261	1993.9529	0.9875488E+01
158	2790.3247	0.8135986E+01	210	2388.2729	0.8154297E+01	262	1986.2212	0.9838867E+01
159	2782.5930	0.8081055E+01	211	2380.5413	0.8050537E+01	263	1978.4895	0.9631348E+01
160	2774.8613	0.8038330E+01	212	2372.8096	0.8123779E+01	264	1970.7576	0.9509277E+01
161	2767.1294	0.7928467E+01	213	2365.0776	0.8190918E+01	265	1963.0259	0.9399414E+01
162	2759.3977	0.7857432E+01	214	2357.3459	0.8105469E+01	266	1955.2942	0.9490967E+01
163	2751.6660	0.7727051E+01	215	2349.6143	0.8074951E+01	267	1947.5623	0.9484863E+01
164	2743.9343	0.7641602E+01	216	2341.8823	0.8129883E+01	268	1939.8306	0.9503174E+01
165	2736.2024	0.7617188E+01	217	2334.1506	0.8233643E+01	269	1932.0989	0.9509277E+01
166	2728.4707	0.7751655E+01	218	2326.4189	0.8099365E+01	270	1924.3669	0.9490967E+01
167	2720.7390	0.7525635E+01	219	2318.6870	0.8044434E+01	271	1916.6353	0.9649659E+01
168	2713.0071	0.7354736E+01	220	2310.9553	0.8050537E+01	272	1908.9036	0.9580176E+01
169	2705.2754	0.7330322E+01	221	2303.2236	0.8062744E+01	273	1901.1716	0.9729004E+01
170	2697.5437	0.7360840E+01	222	2295.4917	0.8142090E+01	274	1893.4399	0.9759521E+01
171	2689.8118	0.7244873E+01	223	2287.7600	0.8166504E+01	275	1885.7083	0.9824660E+01
172	2682.0801	0.7159424E+01	224	2280.0283	0.8087158E+01	276	1877.9766	0.9960938E+01
173	2674.3484	0.7244873E+01	225	2272.2964	0.8099365E+01	277	1870.2446	0.1001587E+02
174	2666.6165	0.7122803E+01	226	2264.5647	0.8166504E+01	278	1862.5129	0.1010742E+02
175	2658.8848	0.7110596E+01	227	2256.8330	0.8117676E+01	279	1854.7812	0.1016946E+02
176	2651.1531	0.7128906E+01	228	2249.1013	0.8074951E+01	280	1847.0493	0.1021118E+02
177	2643.4211	0.7092285E+01	229	2241.3694	0.8178711E+01	281	1839.3176	0.1043091E+02
178	2635.6895	0.7037544E+01	230	2233.6377	0.8221436E+01	282	1831.5859	0.1054577E+02
179	2627.9578	0.6964111E+01	231	2225.9060	0.8251953E+01	283	1823.8540	0.1065674E+02
180	2620.2261	0.7071871E+01	232	2218.1741	0.8337402E+01	284	1816.1223	0.1064453E+02
181	2612.4941	0.7086182E+01	233	2210.4424	0.8377402E+01	285	1808.3906	0.1059570E+02
182	2604.7625	0.7104492E+01	234	2202.7107	0.8312988E+01	286	1800.6587	0.1059570E+02
183	2597.0308	0.6994629E+01	235	2194.9788	0.8398439E+01	287	1792.9270	0.107167E+02
184	2589.2988	0.7031250E+01	236	2187.2471	0.9520508E+01	288	1785.1953	0.1088257E+02
185	2581.5671	0.695008E+01	237	2179.5154	0.8544922E+01	289	1777.4634	0.1094360E+02
186	2573.8354	0.7104492E+01	238	2171.7834	0.8551025E+01	290	1769.7317	0.1130981E+02
187	2566.1035	0.7128906E+01	239	2164.0518	0.8636475E+01	291	1762.0000	0.1201172E+02
188	2558.3719	0.7214355E+01	240	2156.3201	0.855785E+01	292	1754.2683	0.1449585E+02
189	2550.6401	0.7214355E+01	241	2148.5881	0.8691404E+01	293	1746.5364	0.2595825E+02
190	2542.9082	0.7281494E+01	242	2140.8564	0.8770752E+01	294	1738.8047	0.3718977E+02
191	2535.1765	0.7281494E+01	243	2133.1249	0.867306E+01	295	1731.0730	0.2921143E+02
192	2527.4448	0.7409669E+01	244	2125.3931	0.864868E+01	296	1723.3411	0.2297363E+02
193	2519.7129	0.7489014E+01	245	2117.6611	0.8715820E+01	297	1715.6094	0.1916504E+02
194	2511.9812	0.7434082E+01	246	2109.9294	0.8801270E+01	298	1707.8777	0.1665039E+02
195	2504.2495	0.7495117E+01	247	2102.1978	0.8931781E+01	299	1700.1459	0.1536255E+02
196	2496.5178	0.7531738E+01	248	2094.4659	0.8941650E+01	300	1692.4151	0.1447470E+02
197	2488.7859	0.7623291E+01	249	2086.7341	0.9045419E+01	301	1684.6824	0.1361604E+02
198	2481.0542	0.7684326E+01	250	2079.0024	0.9112549E+01	302	1676.9504	0.1361084E+02
199	2473.3225	0.7757568E+01	251	2071.2705	0.9216309E+01	303	1669.2187	0.1400146E+02
200	2465.5906	0.7781982E+01	252	2063.5389	0.9368898E+01	304	1661.4871	0.1404250E+02
201	2457.8589	0.7928467E+01	253	2055.8071	0.9478760E+01	305	1653.7551	0.1449403E+02
202	2450.1272	0.7885742E+01	254	2048.0752	0.9405518E+01	306	1646.0234	0.1497803E+02
203	2442.3953	0.8001709E+01	255	2040.3435	0.9509277E+01	307	1638.2917	0.1459351E+02
204	2434.6636	0.7971191E+01	256	2032.6119	0.9578105E+01	308	1630.5671	0.1334270E+02
205	2426.9319	0.7965088E+01	257	2024.8799	0.9619141E+01	309	1622.8281	0.1372490E+02
206	2419.2000	0.8074951E+01	258	2017.1482	0.9710691E+01	310	1615.0964	0.1301270E+02
207	2411.4683	0.8007813E+01	259	2009.4165	0.9887695E+01	311	1607.3647	0.1299439E+02
208	2403.7366	0.8087158E+01	260	2001.6848	0.9942627E+01	312	1599.6329	0.1330204E+02

FTIS FILE NUMBER : 24

NON-NORMALIZED (PURE) SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
313	1591.9011	0.1311646F+02	365	1189.8494	0.2593388E+02	417	787.7974	0.3139038F+02			
314	1584.1694	0.1336670F+02	366	1182.1174	0.2826538F+02	418	790.0657	0.3188477F+02			
315	1576.4375	0.1346201F+02	367	1174.3857	0.3090879E+02	419	772.3340	0.3223877F+02			
316	1568.7058	0.1397705F+02	368	1166.5541	0.3435669F+02	420	764.6023	0.3139038F+02			
317	1560.9741	0.1431885F+02	369	1158.9221	0.3646240E+02	421	756.8704	0.3066406F+02			
318	1553.2472	0.1512451E+02	370	1151.1904	0.3979492F+02	422	749.1387	0.3104858E+02			
319	1545.5105	0.1586914F+02	371	1143.4587	0.4609985F+02	423	741.4070	0.3214111E+02			
320	1537.7788	0.1673584E+02	372	1135.7271	0.5200195F+02	424	733.6750	0.3393555F+02			
321	1530.0469	0.1743164F+02	373	1127.9951	0.5616455F+02	425	725.9434	0.3521118E+02			
322	1522.3152	0.1804810E+02	374	1120.2634	0.5993047F+02	426	718.2117	0.3447876E+02			
323	1514.5835	0.1866455E+02	375	1112.5317	0.6455078F+02	427	710.4797	0.3411965F+02			
324	1506.8518	0.1868896F+02	376	1104.7998	0.7031860E+02	428	702.7480	0.3353271F+02			
325	1499.1199	0.1898193F+02	377	1097.0681	0.7661133E+02	429	695.0164	0.3323364F+02			
326	1491.3892	0.1943359F+02	378	1089.3364	0.8706055E+02	430	687.2844	0.3273691F+02			
327	1483.6565	0.2074585E+02	379	1081.6045	0.8706055E+02	431	679.5527	0.3291626E+02			
328	1475.9246	0.2301636E+02	380	1073.8728	0.9419556E+02	432	671.8210	0.3178711F+02			
329	1468.1929	0.2823486E+02	381	1066.1411	0.1004578F+03	433	664.0891	0.3051147E+02			
330	1460.4612	0.3175659F+02	382	1058.4092	0.1043823E+03	434	656.3574	0.2948608F+02			
331	1452.7292	0.3595581F+02	383	1050.6775	0.1070068F+03	435	648.6257	0.2853394F+02			
332	1444.9976	0.4088745F+02	384	1042.9458	0.1067078E+03	436	640.8940	0.3001099E+02			
333	1437.2659	0.4932751E+02	385	1035.2139	0.1040466E+03	437	633.1621	0.3444946E+02			
334	1429.5339	0.5913696E+02	386	1027.4822	0.1014832E+03	438	625.4304	0.3002444E+02			
335	1421.8022	0.6744385E+02	387	1019.7505	0.9910889E+02	439	617.6987	0.1712769F+03			
336	1414.0706	0.6826172F+02	388	1012.0188	0.9580688E+02	440	609.9468	0.1152832E+03			
337	1406.3386	0.6273804E+02	389	1004.2869	0.9193726E+02	441	602.2351	0.9993896F+02			
338	1398.6069	0.5458984E+02	390	996.5552	0.9035645E+02	442	594.5034	0.8622437E+02			
339	1390.8752	0.4732666F+02	391	988.8235	0.8608398E+02	443	586.7715	0.7870483F+02			
340	1383.1436	0.4111328F+02	392	981.0916	0.8131109F+02	444	579.0398	0.7435913E+02			
341	1375.4116	0.3436279F+02	393	973.3599	0.8460693F+02	445	571.3081	0.7116089E+02			
342	1367.6799	0.3082886E+02	394	965.6282	0.9476929E+02	446	563.5762	0.7043402F+02			
343	1359.9482	0.2860718E+02	395	957.8962	0.8590088F+02	447	555.8445	0.6442139F+02			
344	1352.2163	0.2651367F+02	396	950.1646	0.7290273F+02	448	548.1128	0.5911621E+02			
345	1344.4846	0.2455444F+02	397	942.4329	0.6131577F+02	449	540.3809	0.6947221F+02			
346	1336.7529	0.2263184E+02	398	934.7009	0.4964111F+02	450	532.6492	0.6878752F+02			
347	1329.0210	0.2127686E+02	399	926.9692	0.6199951F+02	451	524.9175	0.6534424F+02			
348	1321.2893	0.2094116F+02	400	919.2375	0.6254277F+02	452	517.1858	0.6252441F+02			
349	1313.5576	0.2126465F+02	401	911.5056	0.6911011F+02	453	509.4539	0.5992666E+02			
350	1305.8257	0.2101440F+02	402	903.7739	0.6063843F+02	454	501.7222	0.5735674E+02			
351	1298.0940	0.2020264F+02	403	896.0422	0.5233154E+02	455	493.9904	0.5599076F+02			
352	1290.3623	0.2043457F+02	404	888.3105	0.4765015F+02	456	486.2585	0.5745950F+02			
353	1282.6304	0.2036743F+02	405	880.5786	0.4421997F+02	457	478.5269	0.520410E+02			
354	1274.8987	0.2045898F+02	406	872.8469	0.4270957F+02	458	470.7952	0.6061433F+02			
355	1267.1670	0.2145996E+02	407	865.1152	0.4039301E+02	459	463.0632	0.6007080F+02			
356	1259.4353	0.2282102F+02	408	857.3833	0.3849481F+02	460	455.3315	0.5817971F+02			
357	1251.7034	0.2349022F+02	409	849.6516	0.3703003E+02	461	447.5999	0.5333522F+02			
358	1243.9717	0.2470654F+02	410	841.9199	0.3562727F+02	462	439.8679	0.5019362F+02			
359	1236.2400	0.2447651F+02	411	834.1880	0.3468118F+02	463	432.1362	0.5120063F+02			
360	1228.5091	0.2427368F+02	412	826.4563	0.3371592F+02	464	424.4045	0.5079063F+02			
361	1220.7764	0.2401227F+02	413	818.7246	0.3287964F+02	465	416.6726	0.5130063F+02			
362	1213.0447	0.2362061E+02	414	810.9927	0.3220825F+02	466	408.9409	0.5010061F+02			
363	1205.3127	0.2352295F+02	415	803.2610	0.3153697F+02	467	401.2092	0.5039953F+02			
364	1197.5811	0.2441406F+02	416	795.5293	0.3141479F+02	468	393.4775	0.5010063F+02			

SUMMARY OF "AMPLITUDE" NORMALIZATION :

	FILE NO 16	FILE NO 20	FILE NO 24	FILE NO 28
MAXIMUM AMPLITUDE NEAREST 2850.0 W.N. =	0.2190552E+02	0.1731567E+02	0.1694777E+02	0.1491479E+02
WAVE NUMBER AT MAX. AMPLITUDE = (WNMAX) =	2854.6460	2852.1790	2852.1790	2852.1790
BASLINE AMPLITUDE AT (WNMAX) =	0.1203822E+02	0.8775024E+01	0.8116135E+01	0.8191813E+01
PEAK HEIGHT AT (WNMAX) =	0.9867301E+01	0.8540649E+01	0.8851438E+01	0.1072299E+02
NORMALIZING FACTOR FOR THIS FILE =	0.1086719E+01	0.1255123E+01	0.1211435E+01	0.1007000E+01



FTIS FILE NUMBER : 16

NORMALIZED SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
1	4004.2119	0.1052626E+02	53	3603.0212	0.1055942E+02	105	3201.9303	0.1940758E+02			
2	3996.4966	0.1018799E+02	54	3595.3059	0.1039360E+02	106	3194.1150	0.1848561E+02			
3	3988.7815	0.0955838E+01	55	3587.5908	0.1061912E+02	107	3186.3999	0.1799478E+02			
4	3981.0662	0.9962471E+01	56	3579.8755	0.1048646E+02	108	3178.6846	0.1709273E+02			
5	3973.3511	0.1043340E+02	57	3572.1604	0.1073188E+02	109	3170.9695	0.1635649E+02			
6	3965.6357	0.9617565E+01	58	3564.4451	0.1067881E+02	110	3163.2542	0.1608945E+02			
7	3957.9207	0.9949205E+01	59	3556.7300	0.1087780E+02	111	3155.5391	0.1574629E+02			
8	3950.2053	0.9763487E+01	60	3549.0146	0.1067891E+02	112	3147.8237	0.1532178E+02			
9	3942.4902	0.9684909E+01	61	3541.2996	0.1083137E+02	113	3140.1086	0.1487075E+02			
10	3934.7749	0.9697159E+01	62	3533.5842	0.1123557E+02	114	3132.397	0.1477125E+02			
11	3927.0598	0.9358886E+01	63	3525.8691	0.1071198E+02	115	3124.6787	0.1463860E+02			
12	3919.3445	0.9272659E+01	64	3518.1538	0.1103698E+02	116	3116.9629	0.1433349E+02			
13	3911.6294	0.9266026E+01	65	3510.4387	0.1168037E+02	117	3109.2478	0.1423400E+02			
14	3903.9141	0.9989002E+01	66	3502.7234	0.1205180E+02	118	3101.5325	0.1439318E+02			
15	3896.1920	0.91731.7E+01	67	3495.0083	0.1146812E+02	119	3093.8174	0.1462533E+02			
16	3888.4836	0.9305823E+01	68	3487.2930	0.1093749E+02	120	3086.1021	0.1507636E+02			
17	3880.7686	0.9756854E+01	69	3479.5779	0.1211813E+02	121	3078.3870	0.1471154E+02			
18	3873.0532	0.9750221E+01	70	3471.8625	0.1228395E+02	122	3070.6716	0.1437329E+02			
19	3865.3381	0.1013492E+02	71	3464.1475	0.1281458E+02	123	3062.9565	0.1390286E+02			
20	3857.6228	0.9882877E+01	72	3456.4321	0.1338500E+02	124	3055.2412	0.1471154E+02			
21	3849.9077	0.9770120E+01	73	3448.7170	0.1310662E+02	125	3047.5261	0.1433349E+02			
22	3842.1924	0.1008849E+02	74	3441.0017	0.1384929E+02	126	3039.8108	0.1496360E+02			
23	3834.4773	0.1018799E+02	75	3433.2866	0.1345796E+02	127	3032.0957	0.1507636E+02			
24	3826.7620	0.1000890E+02	76	3425.5713	0.1445288E+02	128	3024.3874	0.1609117E+02			
25	3819.0469	0.1002880E+02	77	3417.8562	0.1470493E+02	129	3016.6653	0.173865E+02			
26	3811.3315	0.1047983E+02	78	3410.1409	0.1533504E+02	130	3008.9502	0.1721875E+02			
27	3803.6165	0.1144159E+02	79	3402.4258	0.1536157E+02	131	3001.2349	0.1695345E+02			
28	3795.9011	0.1055942E+02	80	3394.7104	0.1626363E+02	132	2993.5198	0.1770294E+02			
29	3788.1860	0.1042677E+02	81	3386.9954	0.1633659E+02	133	2985.8044	0.1837952E+02			
30	3780.4707	0.1085127E+02	82	3379.2800	0.1748407E+02	134	2978.0894	0.2191935E+02			
31	3772.7556	0.1081147E+02	83	3371.5649	0.1770294E+02	135	2970.3740	0.2444182E+02			
32	3765.0403	0.1012166E+02	84	3363.8496	0.1845245E+02	136	2962.6589	0.2535057E+02			
33	3757.3252	0.1055942E+02	85	3356.1345	0.1939432E+02	137	2954.9436	0.2456126E+02			
34	3749.6099	0.1103035E+02	86	3348.4192	0.2029637E+02	138	2947.2285	0.2694908E+02			
35	3741.8950	0.1043340E+02	87	3340.7041	0.2152344E+02	139	2939.5132	0.3031854E+02			
36	3734.1797	0.1045993E+02	88	3332.9888	0.2294949E+02	140	2931.7981	0.3132679E+02			
37	3726.4646	0.1069209E+02	89	3325.2737	0.2483974E+02	141	2924.0828	0.3145827E+02			
38	3718.7493	0.1060585E+02	90	3317.5583	0.2655110E+02	142	2916.3677	0.2964862E+02			
39	3711.0342	0.1071861E+02	91	3309.8433	0.2789073E+02	143	2908.6523	0.2670652E+02			
40	3703.3188	0.1073189E+02	92	3302.1279	0.2995374E+02	144	2900.9373	0.2418983E+02			
41	3695.6038	0.1085790E+02	93	3294.4128	0.3031854E+02	145	2893.2219	0.2209384E+02			
42	3687.8884	0.1113648E+02	94	3286.6975	0.3103488E+02	146	2885.5068	0.2221324E+02			
43	3680.1733	0.1055942E+02	95	3278.9824	0.3047772E+02	147	2877.7915	0.2203415E+02			
44	3672.4580	0.1074514E+02	96	3271.2671	0.2925729E+02	148	2870.0764	0.2256479E+02			
45	3664.7429	0.1039360E+02	97	3263.5520	0.2818777E+02	149	2862.3611	0.2294286E+02			
46	3657.0276	0.1041350E+02	98	3255.8367	0.2651794E+02	150	2854.6460	0.2380513E+02			
47	3649.3125	0.1040585E+02	99	3248.1216	0.2550754E+02	151	2846.9307	0.2216497E+02			
48	3641.5972	0.1024105E+02	100	3240.4062	0.2399044E+02	152	2839.2156	0.2184759E+02			
49	3633.8821	0.1049310E+02	101	3232.6912	0.2318164E+02	153	2831.5002	0.2156033E+02			
50	3626.1667	0.1045330E+02	102	3224.9759	0.2197140E+02	154	2823.7857	0.2145654E+02			
51	3618.4517	0.1030074E+02	103	3217.2607	0.2114537E+02	155	2816.0698	0.2106918E+02			
52	3610.7363	0.1052626E+02	104	3209.5454	0.1967788E+02	156	2808.3547	0.2137763E+02			

# FITS FILE NUMBER : 16

## NORMALIZED SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
157	2800.6394	0.1348449E+02	209	2799.4495	0.1261559E+02	261	1998.2578	0.13649011E+02
158	2792.9243	0.1329214E+02	210	2391.7334	0.1253600E+02	262	1990.5425	0.1363041E+02
159	2785.2090	0.1317938E+02	211	2384.0181	0.1252937E+02	263	1982.8274	0.1341153E+02
160	2777.4939	0.1299366E+02	212	2376.3030	0.1266065E+02	264	1975.1121	0.1320591E+02
161	2769.7786	0.1295387E+02	213	2368.5876	0.1260896E+02	265	1967.3970	0.1318601E+02
162	2762.0635	0.1292733E+02	214	2360.8726	0.1260896E+02	266	1959.6816	0.1311969E+02
163	2754.3481	0.1278904E+02	215	2353.1572	0.1256253E+02	267	1951.9666	0.1307989E+02
164	2746.6331	0.1270845E+02	216	2345.4421	0.1263549E+02	268	1944.2512	0.1304672E+02
165	2738.9177	0.1271508E+02	217	2337.7268	0.1270845E+02	269	1936.5361	0.1300693E+02
166	2731.2026	0.1272835E+02	218	2330.0117	0.1257580E+02	270	1928.8208	0.1303346E+02
167	2723.4873	0.1267529E+02	219	2322.2964	0.1256253E+02	271	1921.1057	0.1300029E+02
168	2715.7722	0.1239671E+02	220	2314.5813	0.1248294E+02	272	1913.3904	0.1312632E+02
169	2708.0569	0.1252273E+02	221	2306.8660	0.1246967E+02	273	1905.6753	0.1302019E+02
170	2700.3418	0.1244314E+02	222	2299.1509	0.1248957E+02	274	1897.9600	0.1300693E+02
171	2692.6265	0.1238344E+02	223	2291.4355	0.1254263E+02	275	1890.2440	0.1300979E+02
172	2684.9114	0.1231712E+02	224	2283.7205	0.1244977E+02	276	1882.5295	0.1315285E+02
173	2677.1960	0.1244314E+02	225	2276.0051	0.1253600E+02	277	1874.8145	0.1311969E+02
174	2669.4810	0.1229722E+02	226	2268.2900	0.1239008E+02	278	1867.0991	0.1317938E+02
175	2661.7656	0.1216456E+02	227	2260.5747	0.1237018E+02	279	1859.3840	0.1323908E+02
176	2654.0505	0.1219109E+02	228	2252.8596	0.1233701E+02	280	1851.6687	0.1319265E+02
177	2646.3352	0.1215130E+02	229	2245.1443	0.1240334E+02	281	1843.9536	0.1342479E+02
178	2638.6201	0.1211150E+02	230	2237.4292	0.1231712E+02	282	1836.2383	0.1344469E+02
179	2630.9048	0.1200337E+02	231	2229.7139	0.1235691E+02	283	1828.5232	0.1371000E+02
180	2623.1897	0.1203191E+02	232	2222.0088	0.1240998E+02	284	1820.8079	0.1353755E+02
181	2615.4744	0.1203854E+02	233	2214.2934	0.1242987E+02	285	1813.0928	0.1353755E+02
182	2607.7593	0.1202527E+02	234	2206.5684	0.1240998E+02	286	1805.3774	0.1349776E+02
183	2600.0439	0.1193241E+02	235	2198.8530	0.1238344E+02	287	1797.6624	0.1349112E+02
184	2592.3299	0.1197221E+02	236	2191.1379	0.1240998E+02	288	1789.9470	0.1360388E+02
185	2584.6135	0.1186609E+02	237	2183.4226	0.1242987E+02	289	1782.2319	0.1367684E+02
186	2576.8984	0.1186609E+02	238	2175.7075	0.1244314E+02	290	1774.5166	0.1394368E+02
187	2569.1831	0.1188598E+02	239	2167.9922	0.1246304E+02	291	1766.8015	0.1449269E+02
188	2561.4680	0.1187935E+02	240	2160.2771	0.1248957E+02	292	1759.0862	0.1636311E+02
189	2553.7527	0.1198894E+02	241	2152.5618	0.1254263E+02	293	1751.3711	0.1636311E+02
190	2546.0376	0.1193241E+02	242	2144.8467	0.1262222E+02	294	1743.6558	0.1643359E+02
191	2538.3223	0.1183555E+02	243	2137.1313	0.1253600E+02	295	1735.9407	0.1643359E+02
192	2530.6072	0.1194568E+02	244	2129.4163	0.1254926E+02	296	1728.2253	0.1643359E+02
193	2522.8918	0.1197884E+02	245	2121.7009	0.1249620E+02	297	1720.5103	0.1715697E+02
194	2515.1768	0.1197884E+02	246	2113.9854	0.1262886E+02	298	1712.7949	0.1945401E+02
195	2507.4614	0.1204317E+02	247	2106.2705	0.1267529E+02	299	1705.0798	0.1808765E+02
196	2499.7463	0.1207170E+02	248	2098.5554	0.1274162E+02	300	1697.3645	0.1741110E+02
197	2492.0310	0.1193905E+02	249	2090.8401	0.1282121E+02	301	1689.6494	0.1640291E+02
198	2484.3159	0.1207170E+02	250	2083.1250	0.1302019E+02	302	1681.9341	0.1654221E+02
199	2476.6006	0.1218446E+02	251	2075.4097	0.1299366E+02	303	1674.2190	0.1654221E+02
200	2468.8855	0.1219109E+02	252	2067.6946	0.1313295E+02	304	1666.5037	0.1703304E+02
201	2461.1702	0.1226405E+02	253	2059.9792	0.1325234E+02	305	1658.7886	0.1744477E+02
202	2453.4551	0.1223089E+02	254	2052.2642	0.1329214E+02	306	1651.0732	0.1803142E+02
203	2445.7397	0.1232375E+02	255	2044.5488	0.1334452E+02	307	1643.3582	0.1740464E+02
204	2438.0247	0.1231048E+02	256	2036.8337	0.1328551E+02	308	1635.6428	0.1520393E+02
205	2430.3093	0.1240998E+02	257	2029.1184	0.1339500E+02	309	1627.9277	0.1576417E+02
206	2422.5942	0.1241661E+02	258	2021.4033	0.1354418E+02	310	1620.2124	0.1570644E+02
207	2414.8789	0.1241661E+02	259	2013.6882	0.1363701E+02	311	1612.4973	0.1563352E+02
208	2407.1634	0.1251410E+02	260	2005.9729	0.1367684E+02	312	1604.7820	0.1561362E+02

FTJS FILE NUMBER : 16

NORMALIZED SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	AMPLITUDE
313	1597.0669	0.1550749E+02	365	1195.8760	0.2475362E+02	417	794.6851	0.3328340E+02	
314	1599.3516	0.1566005E+02	366	1199.1606	0.2659753E+02	418	786.9700	0.3361504E+02	
315	1591.6365	0.1552739E+02	367	1180.4456	0.2875319E+02	419	779.2546	0.3410587E+02	
316	1573.9211	0.1607177E+02	368	1172.7302	0.3188388E+02	420	771.5106	0.3314412E+02	
317	1566.2061	0.1587893E+02	369	1165.0151	0.3415896E+02	421	763.8242	0.3265991E+02	
318	1558.4907	0.1699324E+02	370	1157.2998	0.3768758E+02	422	756.1091	0.3218593E+02	
319	1550.7756	0.1725191E+02	371	1149.5847	0.4428059E+02	423	748.3938	0.3302729E+02	
320	1543.0603	0.1810091E+02	372	1141.8694	0.5068124E+02	424	740.6787	0.3557172E+02	
321	1535.3452	0.1873766E+02	373	1134.1543	0.5534410E+02	425	732.9634	0.3685185E+02	
322	1527.6299	0.1912900E+02	374	1126.4390	0.5944981E+02	426	725.2483	0.3611562E+02	
323	1519.9148	0.1993820E+02	375	1118.7239	0.6430502E+02	427	717.5330	0.3568444E+02	
324	1512.1995	0.1968616E+02	376	1111.0085	0.7024138E+02	428	709.8179	0.3524670E+02	
325	1504.4844	0.1995810E+02	377	1103.2935	0.7637005E+02	429	702.1025	0.3498853E+02	
326	1496.7690	0.2019638E+02	378	1095.5781	0.8115234E+02	430	694.3875	0.3429158E+02	
327	1489.0540	0.2148334E+02	379	1087.8630	0.8644533E+02	431	686.6721	0.3397984E+02	
328	1481.3386	0.2312195E+02	380	1080.1477	0.9336334E+02	432	678.9570	0.3297166E+02	
329	1473.6235	0.2786440E+02	381	1072.4326	0.937929E+02	433	671.2417	0.3151244E+02	
330	1465.9082	0.3105478E+02	382	1064.7173	0.1024692E+03	434	663.5266	0.3007312E+02	
331	1458.1931	0.3558499E+02	383	1057.0022	0.1041519E+03	435	655.8113	0.2826900E+02	
332	1450.4778	0.2786440E+02	384	1049.2869	0.1031335E+03	436	648.0962	0.2869350E+02	
333	1442.7627	0.4912917E+02	385	1041.5718	0.1010441E+03	437	640.3809	0.5296956E+02	
334	1435.0474	0.5893245E+02	386	1033.8567	0.9947060E+02	438	632.6658	0.4803056E+02	
335	1427.3323	0.6671936E+02	387	1026.1414	0.9499501E+02	439	624.9504	0.1135005E+03	
336	1419.6169	0.6726988E+02	388	1018.4263	0.911977E+02	440	617.2354	0.1036442E+03	
337	1411.9019	0.6157895E+02	389	1010.7109	0.8771219E+02	441	609.5200	0.9066379E+02	
338	1404.1865	0.5383346E+02	390	1002.9958	0.8586163E+02	442	601.8049	0.8004466E+02	
339	1396.4714	0.4721230E+02	391	995.2805	0.8166971E+02	443	594.0896	0.7309349E+02	
340	1388.7561	0.4146165E+02	392	987.5654	0.7743134E+02	444	586.3745	0.7050006E+02	
341	1381.0410	0.3554499E+02	393	979.8501	0.8048906E+02	445	578.6597	0.6777388E+02	
342	1373.3257	0.3231502E+02	394	972.1350	0.9130054E+02	446	570.9441	0.6766887E+02	
343	1365.6106	0.3029201E+02	395	964.4197	0.8284370E+02	447	563.2288	0.6652651E+02	
344	1357.8953	0.2928899E+02	396	956.7046	0.7045336E+02	448	555.5137	0.6557036E+02	
345	1350.1802	0.2648477E+02	397	948.9893	0.6049117E+02	449	547.7983	0.6581731E+02	
346	1342.4648	0.2472046E+02	398	941.2742	0.6696479E+02	450	540.0832	0.6533461E+02	
347	1334.7498	0.2301033E+02	399	933.5589	0.6048453E+02	451	532.3679	0.6262692E+02	
348	1327.0344	0.2268980E+02	400	925.8437	0.6072995E+02	452	524.6529	0.5897224E+02	
349	1319.3193	0.2306226E+02	401	918.1294	0.6731311E+02	453	516.9375	0.5718140E+02	
350	1311.6040	0.2287653E+02	402	910.4134	0.6023912E+02	454	509.2224	0.5385173E+02	
351	1303.8889	0.2186935E+02	403	902.6980	0.5241904E+02	455	501.5071	0.5265772E+02	
352	1296.1736	0.2190814E+02	404	894.9829	0.4820721E+02	456	493.7920	0.5385835E+02	
353	1288.4585	0.2163620E+02	405	887.2676	0.4554082E+02	457	486.0767	0.5463440E+02	
354	1280.7432	0.2162956E+02	406	879.5525	0.4339946E+02	458	478.3614	0.5503923E+02	
355	1273.0281	0.2231274E+02	407	871.8372	0.4140849E+02	459	470.6462	0.5564022E+02	
356	1265.3127	0.2327451E+02	408	864.1221	0.4006213E+02	460	462.9312	0.5371243E+02	
357	1257.5977	0.2385155E+02	409	856.4067	0.3831770E+02	461	455.2159	0.5093393E+02	
358	1249.8823	0.2446440E+02	410	848.6917	0.3737584E+02	462	447.5007	0.4798833E+02	
359	1242.1672	0.2491331E+02	411	840.9763	0.3621510E+02	463	439.7854	0.4537500E+02	
360	1234.4519	0.2418793E+02	412	833.2612	0.3546559E+02	464	432.0703	0.4537500E+02	
361	1226.7368	0.2387403E+02	413	825.5459	0.3466966E+02	465	424.3550	0.4537500E+02	
362	1219.0215	0.2312155E+02	414	817.8308	0.3413240E+02	466	416.6399	0.4537500E+02	
363	1211.3064	0.2292295E+02	415	810.1155	0.3348901E+02	467	408.9246	0.4537500E+02	
364	1203.5911	0.2343576E+02	416	802.4004	0.3214414E+02	468	401.2095	0.4537500E+02	

FTIS FILE NUMBER : 20

NORMALIZED SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
1	4004.2119	0.6746594E+02	53	3402.1602	0.7801038E+01	105	3200.1084	0.1751784E+02			
2	3996.4800	0.6746594E+02	54	3594.4285	0.7977288E+01	106	3152.3767	0.1653697E+02			
3	3988.7483	0.6746594E+02	55	3586.6968	0.8053921E+01	107	3184.6448	0.1610783E+02			
4	3981.0166	0.6746594E+02	56	3578.9648	0.8460065E+01	109	3176.9131	0.1502734E+02			
5	3973.2847	0.6746594E+02	57	3571.2332	0.8260824E+01	109	3169.1814	0.1452158E+02			
6	3965.5530	0.8046257E+01	58	3563.5015	0.8222508E+01	110	3161.4495	0.1394635E+02			
7	3957.8213	0.6429342E+01	59	3555.7695	0.8184193E+01	111	3153.7179	0.1375527E+02			
8	3950.0894	0.6820161E+01	60	3548.0378	0.8199519E+01	112	3145.9861	0.1308092E+02			
9	3942.3577	0.6207112E+01	61	3540.3062	0.8506043E+01	113	3138.2542	0.1267477E+02			
10	3934.6260	0.6283744E+01	62	3532.5742	0.8521370E+01	114	3130.5225	0.1251345E+02			
11	3926.8940	0.6207112E+01	63	3524.8425	0.8007942E+01	115	3122.7908	0.1245254E+02			
12	3919.1624	0.6214776E+01	64	3517.1108	0.8559685E+01	116	3115.0591	0.1226096E+02			
13	3911.4307	0.5747377E+01	65	3509.3789	0.9287680E+01	117	3107.3271	0.1190080E+02			
14	3903.6947	0.6789500E+01	66	3501.6472	0.9640182E+01	118	3099.5955	0.1214502E+02			
15	3895.9670	0.6130482E+01	67	3493.9155	0.8544359E+01	119	3091.8639	0.1220732E+02			
16	3888.2354	0.5333519E+01	68	3486.1838	0.8552022E+01	120	3084.1318	0.1269776E+02			
17	3880.5037	0.5042371E+01	69	3478.4519	0.9563552E+01	121	3076.4071	0.1250618E+02			
18	3872.7717	0.6636247E+01	70	3470.7202	0.9601867E+01	122	3068.6685	0.1196248E+02			
19	3865.0400	0.6659236E+01	71	3462.9885	0.9931380E+01	123	3060.9365	0.1164075E+02			
20	3857.3083	0.6889129E+01	72	3455.2566	0.1016894E+02	124	3053.2048	0.1152531E+02			
21	3849.5764	0.6322060E+01	73	3447.5249	0.1061340E+02	125	3045.4731	0.1194210E+02			
22	3841.8447	0.6459995E+01	74	3439.7932	0.1102720E+02	126	3037.7412	0.1263646E+02			
23	3834.1130	0.5754990E+01	75	3432.0613	0.1157178E+02	127	3030.0095	0.1298896E+02			
24	3826.3811	0.6728204E+01	76	3424.3296	0.1225330E+02	128	3022.2778	0.1364799E+02			
25	3818.6494	0.6812498E+01	77	3416.5979	0.1221498E+02	129	3014.5459	0.1494305E+02			
26	3810.9177	0.7004075E+01	78	3408.8660	0.1225330E+02	130	3006.8142	0.1482810E+02			
27	3803.1859	0.7302937E+01	79	3401.1343	0.1286635E+02	131	2999.0825	0.1459921E+02			
28	3795.4541	0.7149674E+01	80	3393.4026	0.1314498E+02	132	2991.3508	0.1541816E+02			
29	3787.7224	0.7187990E+01	81	3385.6707	0.1359436E+02	133	2983.6189	0.1648947E+02			
30	3779.9905	0.6873802E+01	82	3377.9390	0.1483377E+02	134	2975.8872	0.1941429E+02			
31	3772.2588	0.6758856E+01	83	3370.2073	0.1521126E+02	135	2968.1555	0.2275174E+02			
32	3764.5271	0.7065390E+01	84	3362.4756	0.1617649E+02	136	2960.4236	0.2374794E+02			
33	3756.7954	0.7050054E+01	85	3354.7437	0.1704274E+02	137	2952.6919	0.2282071E+02			
34	3749.0635	0.6781845E+01	86	3347.0120	0.1802736E+02	138	2944.9602	0.230256E+02			
35	3741.3320	0.7287610E+01	87	3339.2803	0.1875160E+02	139	2937.2283	0.2484392E+02			
36	3733.6003	0.7607461E+01	88	3331.5483	0.2079765E+02	140	2929.4966	0.3011598E+02			
37	3725.8684	0.7854680E+01	89	3323.8167	0.2267511E+02	141	2921.7649	0.2486317E+02			
38	3718.1367	0.7686091E+01	90	3316.0850	0.2444530E+02	142	2914.0330	0.2739702E+02			
39	3710.4050	0.8306802E+01	91	3308.3530	0.2643770E+02	143	2906.3013	0.2484144E+02			
40	3702.6731	0.8732140E+01	92	3300.6213	0.2820797E+02	144	2898.5696	0.2220000E+02			
41	3694.9414	0.8237835E+01	93	3292.8896	0.2927707E+02	145	2890.8376	0.2027292E+02			
42	3687.2097	0.7824027E+01	94	3285.1577	0.2910446E+02	146	2883.1049	0.2029189E+02			
43	3679.4778	0.8099899E+01	95	3277.4257	0.2865233E+02	147	2875.3743	0.2030317E+02			
44	3671.7461	0.7992761E+01	96	3269.6938	0.2800079E+02	148	2867.6426	0.2105054E+02			
45	3664.0144	0.7693754E+01	97	3261.9624	0.2626910E+02	149	2859.9106	0.2126640E+02			
46	3656.2825	0.7701418E+01	98	3254.2307	0.2537253E+02	150	2852.1790	0.2174072E+02			
47	3648.5508	0.7977289E+01	99	3246.4990	0.2388588E+02	151	2844.4473	0.1997333E+02			
48	3640.8191	0.8199519E+01	100	3238.7673	0.2292071E+02	152	2836.7153	0.1623911E+02			
49	3633.0872	0.8046257E+01	101	3231.0354	0.2116549E+02	153	2828.9836	0.1455034E+02			
50	3625.3555	0.8245498E+01	102	3223.3037	0.2021526E+02	154	2821.2520	0.1255216E+02			
51	3617.6238	0.8076910E+01	103	3215.5720	0.1905814E+02	155	2813.5200	0.1206172E+02			
52	3609.8921	0.8053921E+01	104	3207.8401	0.1878415E+02	156	2805.7883	0.1180884E+02			

FTIS FILE NUMBER : 20

NORMALIZED SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
157	2799.0566	0.1161726E+02	209	2396.0066	0.1124943E+02	261	1993.9529	0.1210586E+02
158	2790.3247	0.1139503E+02	210	2388.2729	0.1115748E+02	262	1986.2212	0.1319586E+02
159	2782.5930	0.1117280E+02	211	2380.5413	0.1117280E+02	263	1978.4895	0.1297363E+02
160	2774.8613	0.1112682E+02	212	2372.8096	0.1127242E+02	264	1970.7576	0.1271309E+02
161	2767.1294	0.1091992E+02	213	2365.0776	0.1122607E+02	265	1963.0259	0.1259914E+02
162	2759.3977	0.1088161E+02	214	2357.3459	0.1121112E+02	266	1955.2942	0.1265178E+02
163	2751.6660	0.1093563E+02	215	2349.6143	0.1121878E+02	267	1947.5623	0.1259068E+02
164	2743.9343	0.1074367E+02	216	2341.8823	0.1132607E+02	268	1939.8306	0.1255216E+02
165	2736.2024	0.1073601E+02	217	2334.1506	0.1130308E+02	269	1932.0989	0.1253684E+02
166	2728.4707	0.1075133E+02	218	2326.4189	0.1128775E+02	270	1924.3669	0.1247553E+02
167	2720.7390	0.1056742E+02	219	2318.6870	0.1123411E+02	271	1916.6353	0.1252917E+02
168	2713.0071	0.1038350E+02	220	2310.9553	0.1122645E+02	272	1908.9036	0.1255216E+02
169	2705.2754	0.1045247E+02	221	2303.2236	0.1125710E+02	273	1901.1716	0.1260580E+02
170	2697.5437	0.1058274E+02	222	2295.4917	0.1128009E+02	274	1893.4399	0.1262113E+02
171	2689.8113	0.1047546E+02	223	2287.7600	0.1128775E+02	275	1885.7083	0.1257515E+02
172	2682.0801	0.1037586E+02	224	2280.0283	0.1119579E+02	276	1877.9766	0.1270542E+02
173	2674.3484	0.1052919E+02	225	2272.2964	0.1128009E+02	277	1870.2446	0.1268243E+02
174	2666.6165	0.1046013E+02	226	2264.5647	0.1130308E+02	278	1862.5129	0.1273608E+02
175	2658.8848	0.1037584E+02	227	2256.8330	0.1119541E+02	279	1854.7812	0.1285869E+02
176	2651.1531	0.1043715E+02	228	2249.1013	0.1125710E+02	280	1847.0493	0.1291909E+02
177	2643.4211	0.1029921E+02	229	2241.3694	0.1133373E+02	281	1839.3176	0.1310390E+02
178	2635.6895	0.1034519E+02	230	2233.6377	0.1129541E+02	282	1831.5859	0.1318920E+02
179	2627.9578	0.1028389E+02	231	2225.9060	0.1130309E+02	283	1823.8540	0.1331081E+02
180	2620.2261	0.1036051E+02	232	2218.1741	0.1143335E+02	284	1816.1223	0.1335679E+02
181	2612.4941	0.1042948E+02	233	2210.4424	0.1149465E+02	285	1808.3906	0.135679E+02
182	2604.7625	0.1036818E+02	234	2202.7107	0.1151764E+02	286	1800.6587	0.1334146E+02
183	2597.0308	0.1032986E+02	235	2194.9788	0.1149306E+02	287	1792.9270	0.1344108E+02
184	2589.2982	0.1029215E+02	236	2187.2471	0.1160906E+02	288	1785.1953	0.142500E+02
185	2581.5671	0.1023790E+02	237	2179.5154	0.1164792E+02	289	1777.4634	0.1364032E+02
186	2573.8354	0.1046780E+02	238	2171.7834	0.1160194E+02	290	1769.7317	0.1405413E+02
187	2566.1035	0.1045247E+02	239	2164.0518	0.1173987E+02	291	1762.0000	0.1461354E+02
188	2558.3718	0.1041416E+02	240	2156.3201	0.1131650E+02	292	1754.2583	0.1671373E+02
189	2550.6401	0.1051378E+02	241	2148.5881	0.1191617E+02	293	1746.5364	0.2696045E+02
190	2542.9082	0.1046780E+02	242	2140.8564	0.1204640E+02	294	1738.8047	0.3791702E+02
191	2535.1765	0.1050611E+02	243	2133.1248	0.1295406E+02	295	1731.0730	0.3095126E+02
192	2527.4448	0.1055209E+02	244	2125.3931	0.1196210E+02	296	1723.3411	0.254582E+02
193	2519.7129	0.1068236E+02	245	2117.6611	0.1196210E+02	297	1715.6094	0.2192413E+02
194	2511.9812	0.1062872E+02	246	2109.9294	0.1210770E+02	298	1707.8777	0.1962520E+02
195	2504.2495	0.1064405E+02	247	2102.1978	0.1213069E+02	299	1700.1458	0.1850838E+02
196	2496.5178	0.1063639E+02	248	2094.4658	0.1229929E+02	300	1692.4141	0.1770943E+02
197	2488.7859	0.1075133E+02	249	2086.7341	0.1246027E+02	301	1684.6824	0.1695844E+02
198	2481.0542	0.1071302E+02	250	2079.0024	0.1250618E+02	302	1676.9504	0.1690480E+02
199	2473.3225	0.1078965E+02	251	2071.2705	0.1262879E+02	303	1669.2187	0.1728029E+02
200	2465.5906	0.1095057E+02	252	2063.5389	0.1275146E+02	304	1661.4871	0.1743355E+02
201	2457.8589	0.1093525E+02	253	2055.8071	0.1282031E+02	305	1653.7551	0.1774774E+02
202	2450.1272	0.1094849E+02	254	2048.0752	0.1286635E+02	306	1646.0234	0.1820752E+02
203	2442.3953	0.1102729E+02	255	2040.3435	0.1284167E+02	307	1638.2917	0.1793932E+02
204	2434.6636	0.1102720E+02	256	2032.6118	0.1295931E+02	308	1630.5601	0.1665192E+02
205	2426.9319	0.1113449E+02	257	2024.8799	0.1298130E+02	309	1622.8281	0.1621474E+02
206	2419.2000	0.1119047E+02	258	2017.1482	0.1308092E+02	310	1615.0964	0.1616147E+02
207	2411.4683	0.1114215E+02	259	2009.4165	0.1324134E+02	311	1607.3647	0.1606186E+02
208	2403.7366	0.1111150E+02	260	2001.6848	0.1319586E+02	312	1599.6328	0.1608434E+02

FTIS FILE NUMBER : 20

NORMALIZED SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
313	1591.9011	0.1600821E+02	365	1189.8494	0.2703542E+02	417	787.7974	0.3903963E+02
314	1584.1694	0.1610793E+02	366	1192.1174	0.2902783E+02	418	780.0657	0.3831551E+02
315	1576.4375	0.1616914E+02	367	1174.3857	0.3148769E+02	419	772.3340	0.3843811E+02
316	1568.7058	0.1646800E+02	368	1166.6541	0.3502037E+02	420	764.6023	0.3740359E+02
317	1560.9741	0.1652931E+02	369	1158.9221	0.3729631E+02	421	756.8704	0.3682120E+02
318	1553.2422	0.1731094E+02	370	1151.1904	0.4085965E+02	422	749.1387	0.3703577E+02
319	1545.5105	0.1779377E+02	371	1143.4587	0.4779477E+02	423	741.4070	0.3781740E+02
320	1537.7788	0.1850630E+02	372	1135.7271	0.5484846E+02	424	733.6750	0.3941133E+02
321	1530.0469	0.1916541E+02	373	1127.9951	0.6005573E+02	425	725.9434	0.4063742E+02
322	1522.3152	0.1956389E+02	374	1120.2634	0.6532028E+02	426	718.2117	0.4000905E+02
323	1514.5835	0.2013097E+02	375	1112.5317	0.7075342E+02	427	710.4797	0.3949562E+02
324	1506.8518	0.2009265E+02	376	1104.7998	0.7712912E+02	428	702.7480	0.3885193E+02
325	1499.1199	0.2032254E+02	377	1097.0681	0.8360445E+02	429	695.0164	0.3859371E+02
326	1491.3882	0.2074999E+02	378	1089.3364	0.8909122E+02	430	687.2844	0.3784039E+02
327	1483.6565	0.2196245E+02	379	1081.6045	0.9474660E+02	431	679.5527	0.3755685E+02
328	1475.9246	0.2390887E+02	380	1073.8728	0.1023637E+03	432	671.8210	0.3658365E+02
329	1468.1929	0.2867532E+02	381	1066.1411	0.1096973E+03	433	664.0891	0.3514297E+02
330	1460.4612	0.3219269E+02	382	1059.4092	0.1141419E+03	434	656.3574	0.3360269E+02
331	1452.7292	0.3671391E+02	383	1050.6775	0.1169853E+03	435	648.6257	0.3131909E+02
332	1444.9976	0.4195547E+02	384	1042.9458	0.1157512E+03	436	640.8940	0.3160263E+02
333	1437.2659	0.5135811E+02	385	1035.2139	0.1130537E+03	437	633.1621	0.3789707E+02
334	1429.5339	0.6242363E+02	386	1027.4822	0.1101111E+03	438	625.4304	0.3705319E+02
335	1421.8022	0.7146608E+02	387	1019.7505	0.1067087E+03	439	617.6987	0.1296507E+03
336	1414.0706	0.7226305E+02	388	1012.0188	0.1025400E+03	440	609.9668	0.1183336E+03
337	1406.3386	0.6610191E+02	389	1004.2869	0.9810303E+02	441	602.2351	0.1009537E+03
338	1398.6069	0.5761119E+02	390	996.5552	0.9654209E+02	442	594.5034	0.8850883E+02
339	1390.8752	0.5032759E+02	391	988.8235	0.9176564E+02	443	586.7715	0.8165802E+02
340	1383.1426	0.4414713E+02	392	981.0916	0.8699153E+02	444	579.0398	0.7811000E+02
341	1375.4116	0.3779442E+02	393	973.3599	0.8041695E+02	445	571.3081	0.7519269E+02
342	1367.6799	0.3444563E+02	394	965.6282	0.1024633E+03	446	563.5762	0.7484552E+02
343	1359.9482	0.3226165E+02	395	957.8962	0.9304527E+02	447	555.8445	0.7434741E+02
344	1352.2163	0.3020795E+02	396	950.1646	0.7945103E+02	448	548.1128	0.7454666E+02
345	1344.4846	0.2825382E+02	397	942.4329	0.6798705E+02	449	540.3809	0.7503329E+02
346	1336.7529	0.2649901E+02	398	934.7009	0.7643944E+02	450	532.6492	0.7657738E+02
347	1329.0210	0.2506599E+02	399	926.9692	0.6836253E+02	451	524.9175	0.7262712E+02
348	1321.2893	0.2460622E+02	400	919.2375	0.6893726E+02	452	517.1858	0.7125910E+02
349	1313.5576	0.2482079E+02	401	911.5056	0.7614058E+02	453	509.4539	0.6945068E+02
350	1305.8257	0.2459855E+02	402	903.7739	0.6807133E+02	454	501.7222	0.6709811E+02
351	1298.0940	0.2396475E+02	403	896.0422	0.6007872E+02	455	493.9905	0.6519000E+02
352	1290.3623	0.2396475E+02	404	888.3105	0.5547319E+02	456	486.2595	0.6766519E+02
353	1282.6304	0.2343300E+02	405	880.5786	0.5182556E+02	457	478.5269	0.7169598E+02
354	1274.8987	0.2352571E+02	406	872.8469	0.5002472E+02	458	470.7952	0.7503807E+02
355	1267.1670	0.2427670E+02	407	865.1152	0.4789438E+02	459	463.0632	0.7615590E+02
356	1259.4353	0.2539551E+02	408	857.3833	0.4603992E+02	460	455.3315	0.7514438E+02
357	1251.7034	0.2592427E+02	409	849.6516	0.4442300E+02	461	447.5999	0.7419795E+02
358	1243.9717	0.2641395E+02	410	841.9199	0.4302872E+02	462	439.8679	0.7404377E+02
359	1236.2400	0.2685694E+02	411	834.1880	0.4192482E+02	463	432.1352	0.7454594E+02
360	1228.5081	0.2620976E+02	412	826.4563	0.4099750E+02	464	424.4045	0.7465944E+02
361	1220.7764	0.2591959E+02	413	818.7246	0.3988634E+02	465	416.6726	0.7465944E+02
362	1213.0447	0.2540318E+02	414	810.9927	0.3929639E+02	466	408.9409	0.7465944E+02
363	1205.3127	0.2502768E+02	415	803.2610	0.3844579E+02	467	401.2092	0.7465944E+02
364	1197.5811	0.2574402E+02	416	795.5293	0.3816991E+02	468	393.4775	0.7465944E+02

FTIS FILE NUMBER : 24

NORMALIZED SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
1	4004.2119	0.4840137E+02	53	3602.1502	0.6040928E+01	105	3200.1024	0.158R239E+02
2	3996.4800	0.4840137E+02	54	3594.4285	0.5959594E+01	106	3192.3767	0.1486702E+02
3	3988.7483	0.4840137E+02	55	3586.6969	0.5922624E+01	107	3184.6448	0.1451500E+02
4	3981.0166	0.6292325E+01	56	3578.9648	0.5856077E+01	108	3176.9131	0.1360503E+02
5	3973.2847	0.6210991E+01	57	3571.2332	0.6174021E+01	109	3169.1814	0.1316139E+02
6	3965.5530	0.6588087E+01	58	3563.5015	0.6381054E+01	110	3161.4495	0.1254768E+02
7	3957.8213	0.5745167E+01	59	3555.7695	0.5989170E+01	111	3153.7178	0.1208186E+02
8	3950.0894	0.5974382E+01	60	3548.0378	0.5959594E+01	112	3145.9861	0.1177131E+02
9	3942.3577	0.5234978E+01	61	3540.3062	0.5469783E+01	113	3138.2542	0.1126117E+02
10	3934.6260	0.4584303E+01	62	3532.5742	0.6558511E+01	114	3130.5225	0.1139421E+02
11	3926.3940	0.4502969E+01	63	3524.8425	0.5996564E+01	115	3122.7908	0.1107627E+02
12	3919.1624	0.5656438E+01	64	3517.1109	0.6706391E+01	116	3115.0591	0.1081749E+02
13	3911.4307	0.5663333E+01	65	3509.3789	0.7652878E+01	117	3107.3271	0.1085445E+02
14	3903.6987	0.5700803E+01	66	3501.6472	0.7711981E+01	118	3099.5955	0.1090621E+02
15	3895.9670	0.4806124E+01	67	3493.9155	0.7009547E+01	119	3091.8638	0.1103930E+02
16	3888.2354	0.4443816E+01	68	3486.1838	0.6735968E+01	120	3084.1318	0.1158646E+02
17	3880.5037	0.4421634E+01	69	3478.4519	0.7992954E+01	121	3076.4001	0.1131288E+02
18	3872.7717	0.5545527E+01	70	3470.7202	0.7815497E+01	122	3068.6685	0.1061744E+02
19	3865.0400	0.5449406E+01	71	3462.9885	0.8791511E+01	123	3060.9365	0.1035905E+02
20	3857.3083	0.5175826E+01	72	3455.2566	0.8912709E+01	124	3053.2048	0.1032947E+02
21	3849.5764	0.4880064E+01	73	3447.5249	0.9257325E+01	125	3045.4731	0.1049881E+02
22	3841.8447	0.5057521E+01	74	3439.7932	0.9686169E+01	126	3037.7412	0.1138682E+02
23	3834.1130	0.5737773E+01	75	3432.0613	0.1014467E+02	127	3030.0095	0.1168997E+02
24	3826.3811	0.5686024E+01	76	3424.3296	0.1051437E+02	128	3022.2778	0.1248935E+02
25	3818.6494	0.5619468E+01	77	3416.5979	0.1087663E+02	129	3014.5459	0.1369374E+02
26	3810.9177	0.6440207E+01	78	3408.8660	0.1112064E+02	130	3006.8142	0.1372333E+02
27	3803.1858	0.6787726E+01	79	3401.1343	0.1162343E+02	131	2999.0825	0.1313131E+02
28	3795.4541	0.5999170E+01	80	3393.4025	0.1253280E+02	132	2991.3502	0.1411522E+02
29	3787.7224	0.5863471E+01	81	3385.6707	0.1293957E+02	133	2983.6189	0.1543875E+02
30	3779.9905	0.6092697E+01	82	3377.9390	0.1350151E+02	134	2975.8872	0.1787977E+02
31	3772.2588	0.6284731E+01	83	3370.2073	0.1399691E+02	135	2968.1555	0.2147052E+02
32	3764.5271	0.5981776E+01	84	3362.4756	0.1459583E+02	136	2960.4238	0.2733914E+02
33	3756.7954	0.6713786E+01	85	3354.7437	0.1588979E+02	137	2952.6919	0.2169409E+02
34	3749.0635	0.6514147E+01	86	3347.0120	0.1662195E+02	138	2944.9602	0.2407498E+02
35	3741.3320	0.6477177E+01	87	3339.2803	0.1786398E+02	139	2937.2283	0.2762411E+02
36	3733.6003	0.6602875E+01	88	3331.5483	0.1948378E+02	140	2929.4966	0.2873322E+02
37	3725.8684	0.6817303E+01	89	3323.8167	0.2138354E+02	141	2921.7649	0.2970377E+02
38	3718.1367	0.6545905E+01	90	3316.0850	0.2274405E+02	142	2914.0330	0.2648543E+02
39	3710.4050	0.6656634E+01	91	3308.3530	0.2440771E+02	143	2906.3013	0.2364612E+02
40	3702.6731	0.6440207E+01	92	3300.6213	0.2550023E+02	144	2898.5695	0.2103964E+02
41	3694.9414	0.6558511E+01	93	3292.8896	0.2701041E+02	145	2890.8372	0.1901765E+02
42	3687.2097	0.6636434E+01	94	3285.1577	0.2754279E+02	146	2883.1047	0.1901409E+02
43	3679.4779	0.6493359E+01	95	3277.4260	0.2685513E+02	147	2875.3743	0.1905463E+02
44	3671.7461	0.6671604E+01	96	3269.6944	0.2598264E+02	148	2867.6426	0.1961638E+02
45	3664.0144	0.6284931E+01	97	3261.9624	0.2463693E+02	149	2859.9105	0.2011917E+02
46	3656.2825	0.6418024E+01	98	3254.2307	0.2300294E+02	150	2852.1790	0.2005453E+02
47	3648.5508	0.6486571E+01	99	3246.4990	0.2178464E+02	151	2844.4473	0.1879564E+02
48	3640.8191	0.6114850E+01	100	3238.7573	0.2078464E+02	152	2836.7153	0.1534165E+02
49	3633.0872	0.6262750E+01	101	3231.0254	0.2192765E+02	153	2828.9834	0.1277762E+02
50	3625.3555	0.6166203E+01	102	3223.2937	0.1964037E+02	154	2821.2520	0.1135724E+02
51	3617.6238	0.6247961E+01	103	3215.5620	0.1720542E+02	155	2813.5200	0.1087487E+02
52	3609.8921	0.6077899E+01	104	3207.8301	0.16448131E+02	156	2805.7883	0.1067220E+02

F15 FILE NUMBER : 24

NORMALIZED SPECTRAL DATA

DATA WAVE	WAVE NUMBER	AMPLITUDE	DATA WAVE	WAVE NUMBER	AMPLITUDE	DATA WAVE	WAVE NUMBER	AMPLITUDE	DATA WAVE	WAVE NUMBER	AMPLITUDE
157	2798.0566	0.1027771E+02	209	2396.0746	0.9996738E+01	261	1993.9529	0.1178610E+02	311	1607.3647	0.1447013E+02
158	2790.3247	0.1020377E+02	210	2388.2729	0.9996738E+01	262	1986.2212	0.1166040E+02	312	1599.6328	0.1447013E+02
159	2782.5930	0.1006328E+02	211	2380.5413	0.9996738E+01	263	1978.4895	0.1154209E+02			
160	2774.8613	0.9922198E+01	212	2372.8096	0.1006328E+02	264	1970.7576	0.1140900E+02			
161	2767.1294	0.9856252E+01	213	2365.0776	0.1005589E+02	265	1963.0259	0.1126112E+02			
162	2759.3977	0.9774918E+01	214	2357.3459	0.9981951E+01	266	1955.2942	0.1126112E+02			
163	2751.6660	0.9700977E+01	215	2349.6143	0.9981951E+01	267	1947.5623	0.1126112E+02			
164	2743.9343	0.9590066E+01	216	2341.8923	0.1001153E+02	268	1939.8306	0.1126112E+02			
165	2736.2024	0.9545702E+01	217	2334.1506	0.1008547E+02	269	1932.0989	0.1126112E+02			
166	2728.4707	0.9696189E+01	218	2326.4189	0.1002631E+02	270	1924.3669	0.1126112E+02			
167	2720.7390	0.9597461E+01	219	2318.6870	0.9967162E+01	271	1916.6353	0.1137942E+02			
168	2713.0071	0.9346064E+01	220	2310.9553	0.9944981E+01	272	1908.9036	0.1149773E+02			
169	2705.2754	0.9264729E+01	221	2303.2236	0.9989345E+01	273	1901.1716	0.1148294E+02			
170	2697.5437	0.9279517E+01	222	2295.4917	0.9967162E+01	274	1893.4399	0.1151991E+02			
171	2689.8110	0.9246911E+01	223	2287.7600	0.1001153E+02	275	1885.7083	0.1161603E+02			
172	2682.0801	0.9146424E+01	224	2280.0283	0.9989345E+01	276	1877.9756	0.1167519E+02			
173	2674.3484	0.9286911E+01	225	2272.2964	0.1002631E+02	277	1870.2446	0.1168997E+02			
174	2666.6165	0.9082421E+01	226	2264.5647	0.1008547E+02	278	1862.5129	0.1179349E+02			
175	2658.8848	0.8902421E+01	227	2256.8330	0.1002631E+02	279	1854.7812	0.1187492E+02			
176	2651.1531	0.9042908E+01	228	2249.1013	0.9981951E+01	280	1847.0493	0.1194976E+02			
177	2643.4211	0.8998544E+01	229	2241.3694	0.1013723E+02	281	1839.3176	0.1207446E+02			
178	2635.6895	0.9035514E+01	230	2233.6377	0.1007069E+02	282	1831.5859	0.1215031E+02			
179	2627.9578	0.8887632E+01	231	2225.9060	0.1017420E+02	283	1823.8540	0.1229628E+02			
180	2620.2261	0.8931997E+01	232	2218.1741	0.1024074E+02	284	1816.1223	0.1229628E+02			
181	2612.4944	0.8937391E+01	233	2210.4424	0.1027032E+02	285	1808.3906	0.1217959E+02			
182	2604.7625	0.8931997E+01	234	2202.7107	0.1024074E+02	286	1800.6587	0.1224453E+02			
183	2597.0308	0.8828491E+01	235	2194.9788	0.1023335E+02	287	1792.9270	0.1227610E+02			
184	2589.2989	0.8843264E+01	236	2187.2471	0.1043941E+02	288	1785.1953	0.1242938E+02			
185	2581.5671	0.8917209E+01	237	2179.5154	0.1032947E+02	289	1777.4634	0.1242198E+02			
186	2573.8354	0.9983755E+01	238	2171.7834	0.1041820E+02	290	1769.7317	0.1265120E+02			
187	2566.1035	0.9976361E+01	239	2164.0518	0.1049214E+02	291	1762.0000	0.1330927E+02			
188	2558.3718	0.9998544E+01	240	2156.3201	0.1055129E+02	292	1754.2683	0.1559530E+02			
189	2550.6401	0.9019904E+01	241	2148.5881	0.1058826E+02	293	1746.5364	0.2603658E+02			
190	2542.9092	0.9161212E+01	242	2140.8564	0.1058826E+02	294	1738.8047	0.2954657E+02			
191	2535.1765	0.9124242E+01	243	2133.1248	0.1058826E+02	295	1731.0730	0.2954657E+02			
192	2527.4448	0.9235152E+01	244	2125.3931	0.1054390E+02	296	1723.3411	0.2939449E+02			
193	2519.7129	0.9249941E+01	245	2117.6611	0.1051432E+02	297	1715.6094	0.2043712E+02			
194	2511.9812	0.9235152E+01	246	2109.9294	0.1068448E+02	298	1707.8777	0.1817453E+02			
195	2504.2495	0.9294303E+01	247	2102.1978	0.1071397E+02	299	1700.1458	0.1700629E+02			
196	2496.5178	0.9360851E+01	248	2094.4659	0.1083226E+02	300	1692.4141	0.1594894E+02			
197	2488.7859	0.9427398E+01	249	2086.7341	0.1089831E+02	301	1684.6824	0.1517926E+02			
198	2481.0542	0.9456974E+01	250	2079.0024	0.1100233E+02	302	1676.9504	0.1521693E+02			
199	2473.3225	0.9538308E+01	251	2071.2705	0.1101712E+02	303	1669.2187	0.1559402E+02			
200	2465.5905	0.9604855E+01	252	2063.5388	0.1108366E+02	304	1661.4871	0.1557184E+02			
201	2457.8589	0.9627037E+01	253	2055.8071	0.1120197E+02	305	1653.7551	0.1604506E+02			
202	2450.1272	0.9797099E+01	254	2048.0752	0.1126851E+02	306	1646.0234	0.1633907E+02			
203	2442.3953	0.9767524E+01	255	2040.3435	0.1130549E+02	307	1638.2917	0.1626680E+02			
204	2434.6636	0.9709777E+01	256	2032.6119	0.1134245E+02	308	1630.5601	0.1482504E+02			
205	2426.9319	0.9871039E+01	257	2024.8799	0.1146915E+02	309	1622.8281	0.1447013E+02			
206	2419.2000	0.9824019E+01	258	2017.1482	0.1151991E+02	310	1615.0944	0.1439619E+02			
207	2411.4683	0.9826675E+01	259	2009.4165	0.1161603E+02	311	1607.3647	0.1447013E+02			
208	2403.7366	0.9974556E+01	260	2001.6848	0.1167519E+02						



# FITS FILE NUMBER : 24

## NORMALIZED SPECTRAL DATA

DATA WAVE	WAVE NUMBER	AMPLITUDE	DATA WAVE	WAVE NUMBER	AMPLITUDE	DATA WAVE	WAVE NUMBER	AMPLITUDE	DATA WAVE	WAVE NUMBER	AMPLITUDE
313	1591.9011	0.1445534E+02	365	1189.8494	0.2559814E+02	417	787.7974	0.3443402E+02	469	787.7974	0.3443402E+02
314	1584.1694	0.1469195E+02	366	1182.1174	0.2774242E+02	418	780.0657	0.3478894E+02	470	780.0657	0.3478894E+02
315	1576.4375	0.1478068E+02	367	1174.3957	0.3027485E+02	419	772.3340	0.3511427E+02	471	772.3340	0.3511427E+02
316	1568.7058	0.1517996E+02	368	1166.6541	0.3369463E+02	420	764.6023	0.3406433E+02	472	764.6023	0.3406433E+02
317	1560.9741	0.1560142E+02	369	1158.9221	0.3595720E+02	421	756.8704	0.3342104E+02	473	756.8704	0.3342104E+02
318	1553.2422	0.1630385E+02	370	1151.1904	0.3964682E+02	422	749.1387	0.3362064E+02	474	749.1387	0.3362064E+02
319	1545.5105	0.1696931E+02	371	1143.4587	0.4647891E+02	423	741.4070	0.3455234E+02	475	741.4070	0.3455234E+02
320	1537.7788	0.1770132E+02	372	1135.7271	0.5314095E+02	424	733.6750	0.3645099E+02	476	733.6750	0.3645099E+02
321	1530.0469	0.1849636E+02	373	1127.9951	0.5807277E+02	425	725.9434	0.3769001E+02	477	725.9434	0.3769001E+02
322	1522.3152	0.1894351E+02	374	1120.2634	0.6254616E+02	426	718.2117	0.3693321E+02	478	718.2117	0.3693321E+02
323	1514.5835	0.1959419E+02	375	1112.5317	0.6751495E+02	427	710.4797	0.3644571E+02	479	710.4797	0.3644571E+02
324	1506.8518	0.1955722E+02	376	1104.7998	0.7361504E+02	428	702.7430	0.3595720E+02	480	702.7430	0.3595720E+02
325	1499.1199	0.1974207E+02	377	1097.0681	0.7921160E+02	429	695.0164	0.3552095E+02	481	695.0164	0.3552095E+02
326	1491.3982	0.2007460E+02	378	1089.3364	0.8463954E+02	430	687.2844	0.3506992E+02	482	687.2844	0.3506992E+02
327	1483.6565	0.2124308E+02	379	1081.6045	0.9005937E+02	431	679.5527	0.3474457E+02	483	679.5527	0.3474457E+02
328	1475.9246	0.2339473E+02	380	1073.8728	0.9756432E+02	432	671.8210	0.3366504E+02	484	671.8210	0.3366504E+02
329	1468.1929	0.2820825E+02	381	1066.1411	0.1043742E+03	433	664.0891	0.3223799E+02	485	664.0891	0.3223799E+02
330	1460.4612	0.3168344E+02	382	1058.4092	0.1082635E+03	434	656.3574	0.3083313E+02	486	656.3574	0.3083313E+02
331	1452.7292	0.3594980E+02	383	1050.6775	0.1109514E+03	435	648.6257	0.2933054E+02	487	648.6257	0.2933054E+02
332	1444.9976	0.4102211E+02	384	1042.9454	0.1095722E+03	436	640.8940	0.2752588E+02	488	640.8940	0.2752588E+02
333	1437.2659	0.5034285E+02	385	1035.2139	0.1071174E+03	437	633.1621	0.2572963E+02	489	633.1621	0.2572963E+02
334	1429.5339	0.6032056E+02	386	1027.4822	0.1045887E+03	438	625.4304	0.2395247E+02	490	625.4304	0.2395247E+02
335	1421.8022	0.6875714E+02	387	1019.7505	0.1009138E+03	439	617.6987	0.2246673E+03	491	617.6987	0.2246673E+03
336	1414.0706	0.6733383E+02	388	1012.0188	0.9692844E+02	440	609.9668	0.2144002E+03	492	609.9668	0.2144002E+03
337	1406.3386	0.6338148E+02	389	1004.2469	0.9261032E+02	441	602.2351	0.2059358E+02	493	602.2351	0.2059358E+02
338	1398.6069	0.5495244E+02	390	996.5552	0.9063871E+02	442	594.5034	0.1957708E+02	494	594.5034	0.1957708E+02
339	1390.8752	0.4788377E+02	391	988.8235	0.8609138E+02	443	586.7715	0.1859381E+02	495	586.7715	0.1859381E+02
340	1383.1436	0.4179109E+02	392	981.0916	0.8143793E+02	444	579.0398	0.1769016E+02	496	579.0398	0.1769016E+02
341	1375.4116	0.3555792E+02	393	973.3599	0.8466173E+02	445	571.3081	0.1719201E+02	497	571.3081	0.1719201E+02
342	1367.6799	0.3200879E+02	394	965.6282	0.9593266E+02	446	563.5763	0.1719552E+02	498	563.5763	0.1719552E+02
343	1359.9482	0.2992344E+02	395	957.8962	0.9650704E+02	447	555.8445	0.1715447E+02	499	555.8445	0.1715447E+02
344	1352.2163	0.2786073E+02	396	950.1646	0.7309006E+02	448	548.1129	0.1709217E+02	500	548.1129	0.1709217E+02
345	1344.4846	0.2596785E+02	397	942.4329	0.6239827E+02	449	540.3809	0.1714411E+02	501	540.3809	0.1714411E+02
346	1336.7529	0.2409717E+02	398	934.7009	0.6099593E+02	450	532.6492	0.1707718E+02	502	532.6492	0.1707718E+02
347	1329.0210	0.2268491E+02	399	926.9692	0.6271622E+02	451	524.9175	0.1673792E+02	503	524.9175	0.1673792E+02
348	1321.2893	0.2230780E+02	400	919.2375	0.6304374E+02	452	517.1858	0.1635423E+02	504	517.1858	0.1635423E+02
349	1313.5576	0.2243340E+02	401	911.5056	0.5989583E+02	453	509.4530	0.1590211E+02	505	509.4530	0.1590211E+02
350	1305.8257	0.2226344E+02	402	903.7739	0.6207858E+02	454	501.7222	0.1598495E+02	506	501.7222	0.1598495E+02
351	1298.0947	0.2133919E+02	403	896.0422	0.5408737E+02	455	493.9905	0.1571394E+02	507	493.9905	0.1571394E+02
352	1290.3623	0.2150925E+02	404	888.3105	0.4967531E+02	456	486.2585	0.1504253E+02	508	486.2585	0.1504253E+02
353	1282.6304	0.2134657E+02	405	880.5786	0.4688559E+02	457	478.5269	0.1473292E+02	509	478.5269	0.1473292E+02
354	1274.8987	0.2137919E+02	406	872.8469	0.4474132E+02	458	470.7952	0.1426433E+02	510	470.7952	0.1426433E+02
355	1267.1670	0.2221937E+02	407	865.1152	0.4314420E+02	459	463.0632	0.1413187E+02	511	463.0632	0.1413187E+02
356	1259.4353	0.2320799E+02	408	857.3833	0.4131787E+02	460	455.3315	0.1384822E+02	512	455.3315	0.1384822E+02
357	1251.7034	0.2329709E+02	409	849.6516	0.3986125E+02	461	447.5999	0.1330836E+02	513	447.5999	0.1330836E+02
358	1243.9717	0.2472565E+02	410	841.9193	0.3895305E+02	462	439.8679	0.1277124E+02	514	439.8679	0.1277124E+02
359	1236.2400	0.25006579E+02	411	834.1880	0.3773171E+02	463	432.1362	0.1244317E+02	515	432.1362	0.1244317E+02
360	1228.5081	0.2440646E+02	412	826.4553	0.3658181E+02	464	424.4045	0.1244013E+02	516	424.4045	0.1244013E+02
361	1220.7764	0.2427461E+02	413	818.7246	0.3607550E+02	465	416.6726	0.1244413E+02	517	416.6726	0.1244413E+02
362	1213.0447	0.2374225E+02	414	810.9927	0.3526955E+02	466	408.9409	0.1244413E+02	518	408.9409	0.1244413E+02
363	1205.3127	0.2350555E+02	415	803.2610	0.3472974E+02	467	401.2092	0.1244413E+02	519	401.2092	0.1244413E+02
364	1197.5811	0.2442176E+02	416	795.5293	0.3433970E+02	468	393.4775	0.1244413E+02	520	393.4775	0.1244413E+02

## FITS FILE NUMBER : 28

## NORMALIZED SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
1	4004.2119	0.5039063E+02	53	3632.1402	0.4962158E+01	105	3700.1084	0.1465454E+02
2	3996.4800	0.5039063E+02	54	3594.4285	0.4809570E+01	106	3192.3767	0.1389771E+02
3	3988.7483	0.5039063E+02	55	3586.6968	0.5059814E+01	107	3184.6448	0.1307373E+02
4	3981.0166	0.5039063E+02	56	3578.9648	0.5017090E+01	108	3176.9131	0.1249390E+02
5	3973.2847	0.5039063E+02	57	3571.2332	0.5139160E+01	109	3169.1814	0.1190186E+02
6	3965.5530	0.5039063E+02	58	3563.5015	0.5194092E+01	110	3161.4495	0.1135864E+02
7	3957.8213	0.5039063E+02	59	3555.7695	0.5157471E+01	111	3153.7178	0.1082153E+02
8	3950.0894	0.5039063E+02	60	3548.0378	0.5169678E+01	112	3145.9861	0.1024170E+02
9	3942.3577	0.5039063E+02	61	3540.3062	0.5267334E+01	113	3138.2542	0.1029663E+02
10	3934.6260	0.5039063E+02	62	3532.5742	0.5145264E+01	114	3130.5225	0.1026470E+02
11	3926.8940	0.5039063E+02	63	3524.8425	0.5120850E+01	115	3122.7908	0.9906006E+01
12	3919.1624	0.5039063E+02	64	3517.1108	0.5172509E+01	116	3115.0591	0.968486E+01
13	3911.4307	0.5039063E+02	65	3509.3789	0.6384277E+01	117	3107.3271	0.9564209E+01
14	3903.6987	0.5039063E+02	66	3501.6472	0.6399828E+01	118	3099.5955	0.9729004E+01
15	3895.9670	0.5039063E+01	67	3493.9155	0.5650439E+01	119	3091.8638	0.9375488E+01
16	3888.2354	0.5039063E+01	68	3486.1838	0.5627441E+01	120	3084.1318	0.1021729E+02
17	3880.5037	0.5039063E+01	69	3478.4519	0.6341553E+01	121	3076.4001	0.1003418E+02
18	3872.7717	0.5039063E+01	70	3470.7202	0.7055664E+01	122	3068.6685	0.9399414E+01
19	3865.0400	0.5039063E+01	71	3462.9885	0.7391357E+01	123	3060.9365	0.925930E+01
20	3857.3083	0.4698730E+01	72	3455.2566	0.7965088E+01	124	3053.2048	0.9283447E+01
21	3849.5764	0.3417969E+01	73	3447.5249	0.8276367E+01	125	3045.4731	0.9527588E+01
22	3841.8447	0.3284691E+01	74	3439.7932	0.8502107E+01	126	3037.7412	0.1019897E+02
23	3834.1130	0.3668213E+01	75	3432.0613	0.8795166E+01	127	3030.0095	0.1040549E+02
24	3826.3811	0.4034424E+01	76	3424.3296	0.9155273E+01	128	3022.2778	0.1122437E+02
25	3818.6494	0.4461670E+01	77	3416.5979	0.91729E+01	129	3014.5459	0.1232910E+02
26	3810.9177	0.4882813E+01	78	3408.8660	0.1015015E+02	130	3006.8142	0.1232910E+02
27	3803.1859	0.4809570E+01	79	3401.1343	0.1046613E+02	131	2999.0825	0.1211548E+02
28	3795.4541	0.4449463E+01	80	3393.4026	0.1117554E+02	132	2991.3508	0.1299218E+02
29	3787.7224	0.4449463E+01	81	3385.6707	0.1175537E+02	133	2983.6149	0.1659351E+02
30	3779.9905	0.4144247E+01	82	3377.9390	0.1236572E+02	134	2975.8872	0.1732488E+02
31	3772.2584	0.4296875E+01	83	3370.2073	0.1289053E+02	135	2968.1555	0.2109375E+02
32	3764.5271	0.4756630E+01	84	3362.4756	0.1373901E+02	136	2960.4236	0.2185359E+02
33	3756.7954	0.4779053E+01	85	3354.7437	0.1461792E+02	137	2952.6919	0.2094727E+02
34	3749.0635	0.4779053E+01	86	3347.0120	0.1583862E+02	138	2944.9602	0.2363892E+02
35	3741.3320	0.4833984E+01	87	3339.2803	0.1649464E+02	139	2937.2283	0.2701416E+02
36	3733.6003	0.5240023E+01	88	3331.5483	0.1860962E+02	140	2929.4966	0.2916772E+02
37	3725.8684	0.5230713E+01	89	3323.8167	0.2047729E+02	141	2921.7649	0.2772827E+02
38	3718.1367	0.5407715E+01	90	3316.0850	0.2222230E+02	142	2914.0330	0.2504883E+02
39	3710.4050	0.5529785E+01	91	3308.3530	0.2405616E+02	143	2906.3013	0.2210238E+02
40	3702.6731	0.5749512E+01	92	3300.6213	0.2554932E+02	144	2898.5695	0.1757515E+02
41	3694.9414	0.5682373E+01	93	3292.8896	0.2655029E+02	145	2890.8376	0.1747437E+02
42	3687.2097	0.5084224E+01	94	3285.1577	0.2684326E+02	146	2883.1057	0.1776733E+02
43	3679.4778	0.5541992E+01	95	3277.4260	0.2799877E+02	147	2875.3743	0.1793161E+02
44	3671.7461	0.5645752E+01	96	3269.6943	0.2523193E+02	148	2867.6426	0.1860352E+02
45	3664.0144	0.5371094E+01	97	3261.9624	0.2366943E+02	149	2859.9106	0.1876331E+02
46	3656.2825	0.5053711E+01	98	3254.2307	0.2247925E+02	150	2852.1790	0.1891479E+02
47	3648.5508	0.4633606E+01	99	3246.4990	0.2105713E+02	151	2844.4473	0.1691895E+02
48	3640.8191	0.5233713E+01	100	3238.7673	0.1966551E+02	152	2836.7155	0.1309704E+02
49	3633.0872	0.5102539E+01	101	3231.0354	0.1864014E+02	153	2828.9836	0.1046343E+02
50	3625.3555	0.5242920E+01	102	3223.3037	0.1745216E+02	154	2821.2510	0.5356699E+01
51	3617.6239	0.5189678E+01	103	3215.5720	0.1600342E+02	155	2813.5200	0.8874512E+01
52	3609.8921	0.5133057E+01	104	3207.8401	0.1553345E+02	156	2805.7882	0.8673094E+01

FTIS FILE NUMBER : 28

NORMALIZED SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
157	2798.0566	0.8396230F+01	209	2396.0046	0.8221436F+01	261	1993.9529	0.9875499E+01	317	1509.6398	0.1330209E+02
158	2790.3247	0.8135986E+01	210	2388.2729	0.8154297E+01	262	1986.2212	0.9838467E+01	318	1497.4895	0.1304209E+02
159	2782.5930	0.8081055E+01	211	2380.5413	0.8050527E+01	263	1978.4895	0.9631349E+01	319	1485.0253	0.1273012E+02
160	2774.8612	0.8038330E+01	212	2372.8096	0.8123779E+01	264	1970.7576	0.9509277E+01	320	1472.0211	0.1240209E+02
161	2767.1294	0.7924467E+01	213	2365.0776	0.8190918E+01	265	1963.0759	0.9389414E+01	321	1459.5311	0.1204209E+02
162	2759.3977	0.7967432E+01	214	2357.3459	0.8105469E+01	266	1955.2942	0.9490967E+01	322	1446.1450E+02	
163	2751.6660	0.7727051E+01	215	2349.6143	0.8074951E+01	267	1947.5623	0.9444863E+01	323	1433.4279E+02	
164	2743.9343	0.7641602E+01	216	2341.8923	0.8129883E+01	268	1939.8306	0.9509277E+01	324	1420.2000	0.1304209E+02
165	2736.2024	0.7617188E+01	217	2334.1506	0.8233364E+01	269	1932.0989	0.9509277E+01	325	1407.3647	0.1273012E+02
166	2728.4707	0.7751465E+01	218	2326.4189	0.8099365E+01	270	1924.3669	0.9490967E+01	326	1394.2000	0.1240209E+02
167	2720.7390	0.7525635E+01	219	2318.6870	0.8045444E+01	271	1916.6353	0.9631349E+01	327	1381.0253	0.1204209E+02
168	2713.0071	0.7354736E+01	220	2310.9553	0.8050527E+01	272	1908.9036	0.9444863E+01	328	1368.2000	0.1172021E+02
169	2705.2754	0.7330322E+01	221	2303.2236	0.8062744E+01	273	1901.1716	0.9729004E+01	329	1355.0253	0.1140209E+02
170	2697.5437	0.7360340E+01	222	2295.4917	0.8142909E+01	274	1893.4399	0.9759521E+01	330	1342.0211	0.1104209E+02
171	2689.8118	0.7244873E+01	223	2287.7600	0.8166504E+01	275	1885.7083	0.9826660E+01	331	1329.0253	0.1072021E+02
172	2682.0801	0.7159424E+01	224	2280.0283	0.8087158E+01	276	1877.9766	0.9960938E+01	332	1316.0211	0.1040209E+02
173	2674.3484	0.7244873E+01	225	2272.2964	0.8099365E+01	277	1870.2444	0.1001587E+02	333	1303.0253	0.1004209E+02
174	2666.6165	0.7122903E+01	226	2264.5647	0.8166504E+01	278	1862.5129	0.1010742E+02	334	1290.0211	0.0972021E+02
175	2658.8848	0.7110594E+01	227	2256.8330	0.8117676E+01	279	1854.7812	0.1016846E+02	335	1277.0253	0.0940209E+02
176	2651.1531	0.7128906E+01	228	2249.1013	0.8074951E+01	280	1847.0493	0.1021119E+02	336	1264.0211	0.0904209E+02
177	2643.4211	0.7092285E+01	229	2241.3694	0.8178711E+01	281	1839.3176	0.1043091E+02	337	1251.0253	0.0872021E+02
178	2635.6395	0.7037354E+01	230	2233.6377	0.8221436E+01	282	1831.5859	0.1054077E+02	338	1238.0211	0.0840209E+02
179	2627.9578	0.6964111E+01	231	2225.9060	0.8251953E+01	283	1823.8540	0.1065674E+02	339	1225.0253	0.0804209E+02
180	2620.2261	0.7067871E+01	232	2218.1741	0.8337402E+01	284	1816.1223	0.1064453E+02	340	1212.0211	0.0772021E+02
181	2612.4941	0.7086102E+01	233	2210.4424	0.8337402E+01	285	1808.3906	0.1059570E+02	341	1199.0253	0.0740209E+02
182	2604.7625	0.7104492E+01	234	2202.7107	0.8312988E+01	286	1800.6587	0.1059570E+02	342	1186.0211	0.0704209E+02
183	2597.0308	0.6994629E+01	235	2194.9788	0.8398838E+01	287	1792.9270	0.1071167E+02	343	1173.0253	0.0672021E+02
184	2589.2988	0.7031250E+01	236	2187.2471	0.8520508E+01	288	1785.1953	0.1088257E+02	344	1160.0211	0.0640209E+02
185	2581.5671	0.6958008E+01	237	2179.5154	0.8544922E+01	289	1777.4634	0.1094363E+02	345	1147.0253	0.0604209E+02
186	2573.8354	0.7104492E+01	238	2171.7834	0.8551025E+01	290	1769.7317	0.1113958E+02	346	1134.0211	0.0572021E+02
187	2566.1035	0.7128906E+01	239	2164.0518	0.8636475E+01	291	1762.0000	0.1091172E+02	347	1121.0253	0.0540209E+02
188	2558.3718	0.7135010E+01	240	2156.3201	0.8654785E+01	292	1754.2683	0.1094363E+02	348	1108.0211	0.0504209E+02
189	2550.6401	0.7214355E+01	241	2148.5881	0.8691406E+01	293	1746.5364	0.1094363E+02	349	1095.0253	0.0472021E+02
190	2542.9082	0.7291493E+01	242	2140.8564	0.8770752E+01	294	1738.8047	0.1094363E+02	350	1082.0211	0.0440209E+02
191	2535.1765	0.7281494E+01	243	2133.1248	0.8673096E+01	295	1731.0730	0.1094363E+02	351	1069.0253	0.0404209E+02
192	2527.4448	0.7409666E+01	244	2125.3931	0.8648632E+01	296	1723.3411	0.1094363E+02	352	1056.0211	0.0372021E+02
193	2519.7129	0.7487015E+01	245	2117.6511	0.8715820E+01	297	1715.6094	0.1094363E+02	353	1043.0253	0.0340209E+02
194	2511.9812	0.7434082E+01	246	2109.9294	0.8801270E+01	298	1707.8777	0.1094363E+02	354	1030.0211	0.0304209E+02
195	2504.2495	0.7495117E+01	247	2102.1978	0.8831787E+01	299	1700.1458	0.1094363E+02	355	1017.0253	0.0272021E+02
196	2496.5178	0.7531730E+01	248	2094.4658	0.8541650E+01	300	1692.4141	0.1094363E+02	356	1004.0211	0.0240209E+02
197	2488.7859	0.7623291E+01	249	2086.7341	0.9045410E+01	301	1684.6824	0.1094363E+02	357	991.0253	0.0204209E+02
198	2481.0542	0.7694326E+01	250	2079.0024	0.9121254E+01	302	1676.9504	0.1094363E+02	358	978.0211	0.0172021E+02
199	2473.3225	0.7755688E+01	251	2071.2705	0.9215309E+01	303	1669.2187	0.1094363E+02	359	965.0253	0.0140209E+02
200	2465.5906	0.7791982E+01	252	2063.5388	0.9368896E+01	304	1661.4871	0.1094363E+02	360	952.0211	0.0104209E+02
201	2457.8589	0.7928467E+01	253	2055.8071	0.9478760E+01	305	1653.7551	0.1094363E+02	361	939.0253	0.0072021E+02
202	2450.1272	0.7985747E+01	254	2048.0752	0.9405512E+01	306	1646.0234	0.1094363E+02	362	926.0211	0.0040209E+02
203	2442.3953	0.7901704E+01	255	2040.3435	0.9509277E+01	307	1638.2917	0.1094363E+02	363	913.0253	0.0004209E+02
204	2434.6636	0.7971791E+01	256	2032.6118	0.9558105E+01	308	1630.5601	0.1094363E+02	364	900.0211	0.0000000E+02
205	2426.9319	0.7966509E+01	257	2024.8799	0.9619141E+01	309	1622.8284	0.1094363E+02	365	887.0253	0.0000000E+02
206	2419.2000	0.8074951E+01	258	2017.1482	0.9710593E+01	310	1615.0967	0.1094363E+02	366	874.0211	0.0000000E+02
207	2411.4693	0.8007813E+01	259	2009.4165	0.9887595E+01	311	1607.3647	0.1094363E+02	367	861.0253	0.0000000E+02
208	2403.7366	0.8087159E+01	260	2001.6848	0.9946267E+01	312	1600.0000	0.1094363E+02	368	848.0211	0.0000000E+02

# FTIS FILE NUMBER : 28

## NORMALIZED SPECTRAL DATA

DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE	DATA WORD	WAVE NUMBER	AMPLITUDE
313	1591.9011	0.1311546F+02	315	1139.8494	0.2593338E+02	417	797.7974	0.2139038F+02
314	1594.1694	0.1336670F+02	316	1192.1174	0.2826538E+02	419	790.0657	0.3188477F+02
315	1576.4375	0.1356201E+02	317	1174.3957	0.3090820E+02	419	772.3740	0.3223777F+02
316	1568.7058	0.1397705E+02	318	1166.6541	0.3435669E+02	420	764.6023	0.3139038F+02
317	1560.9741	0.1471885E+02	319	1158.9221	0.3646240E+02	421	756.8704	0.3066406F+02
318	1553.2422	0.1512451E+02	320	1151.1904	0.3979492E+02	422	749.1387	0.3104858F+02
319	1545.5105	0.1586914E+02	321	1143.4587	0.4609985E+02	423	741.4070	0.3214111E+02
320	1537.7788	0.1673584E+02	322	1135.7271	0.5200195E+02	424	733.6750	0.3393555F+02
321	1530.0469	0.1743164E+02	323	1127.9951	0.5616455E+02	425	725.9434	0.3521118F+02
322	1522.3152	0.1804810E+02	324	1120.2534	0.5993042E+02	426	718.2117	0.3447876E+02
323	1514.5835	0.1866455E+02	325	1112.5317	0.6455078E+02	427	710.4797	0.3411865E+02
324	1506.8518	0.1868896E+02	326	1104.7998	0.7031860F+02	428	702.7480	0.3353271E+02
325	1499.1199	0.1898193E+02	327	1097.0681	0.7661133E+02	429	695.0164	0.3323364E+02
326	1491.3882	0.1943359E+02	328	1089.3364	0.8208008E+02	430	687.2844	0.3283691E+02
327	1483.6565	0.2074585E+02	329	1081.6045	0.8706053E+02	431	679.5527	0.3291626E+02
328	1475.9244	0.2201636E+02	330	1073.8728	0.9419550E+02	432	671.8210	0.3178711E+02
329	1468.1929	0.2323486E+02	331	1066.1411	0.1004578E+03	433	664.0891	0.3051147F+02
330	1460.4612	0.3175659E+02	332	1058.4032	0.1043823E+03	434	656.3574	0.2948608E+02
331	1452.7292	0.3595581E+02	333	1050.6775	0.1070068E+03	435	648.6257	0.2953394E+02
332	1444.9976	0.4088745E+02	334	1042.9458	0.1067078E+03	436	640.8940	0.2901099E+02
333	1437.2659	0.4922251E+02	335	1035.2139	0.1040466E+03	437	633.1621	0.2844404E+02
334	1429.5339	0.5913696E+02	336	1027.4822	0.1014832E+03	438	625.4304	0.2800244E+02
335	1421.8022	0.6764385E+02	337	1019.7505	0.9910889E+02	439	617.6987	0.276969E+03
336	1414.0706	0.6826172E+02	338	1012.0188	0.9580608E+02	440	609.9668	0.2715283E+03
337	1406.3386	0.6273894E+02	339	1004.2869	0.9193726E+02	441	602.2351	0.2699389E+02
338	1398.6069	0.5458994E+02	340	996.5552	0.9035645E+02	442	594.5034	0.2622343E+02
339	1390.8752	0.4732666E+02	341	988.8235	0.8608398E+02	443	586.7715	0.2797943E+02
340	1383.1436	0.4111328E+02	342	981.0916	0.8131103E+02	444	579.0398	0.2743591E+02
341	1375.4116	0.3436279E+02	343	973.3593	0.8460693E+02	445	571.3081	0.2716089E+02
342	1367.6799	0.3082816E+02	344	965.6282	0.9476929E+02	446	563.5762	0.2704849E+02
343	1359.9432	0.2860715E+02	345	957.8962	0.8590088E+02	447	555.8445	0.2694213E+02
344	1352.2163	0.2651367E+02	346	950.1646	0.7280273E+02	448	548.1128	0.2691162E+02
345	1344.4846	0.2455444E+02	347	942.4329	0.6131592E+02	449	540.3809	0.2694702E+02
346	1336.7529	0.2263194E+02	348	934.7009	0.6964111E+02	450	532.6492	0.2687805E+02
347	1329.0210	0.2127596E+02	349	926.9692	0.6199951E+02	451	524.9175	0.2653424E+02
348	1321.2893	0.2094115E+02	350	919.2375	0.6254273E+02	452	517.1858	0.2625245E+02
349	1313.5576	0.2126465E+02	351	911.5056	0.6411011E+02	453	509.4539	0.2598266E+02
350	1305.8257	0.2101440E+02	352	903.7739	0.6063843E+02	454	501.7222	0.2573547E+02
351	1298.0940	0.2020264E+02	353	896.0422	0.5733154E+02	455	493.9905	0.2559597E+02
352	1290.3623	0.2043577E+02	354	888.3105	0.4765015E+02	456	486.2585	0.2574550E+02
353	1282.6304	0.2075743E+02	355	880.5786	0.4421997E+02	457	478.5269	0.25920410E+02
354	1274.8937	0.2045898E+02	356	872.8469	0.4230957E+02	458	470.7952	0.2604614E+02
355	1267.1670	0.2145996E+02	357	865.1152	0.4039307E+02	459	463.0632	0.2600790E+02
356	1259.4353	0.22222104E+02	358	857.3833	0.3849497E+02	460	455.3315	0.2581787E+02
357	1251.7034	0.2348022E+02	359	849.6516	0.3703703E+02	461	447.5909	0.2533252E+02
358	1243.9717	0.2420654E+02	360	841.9199	0.3562622E+02	462	439.8479	0.2503906E+02
359	1236.2400	0.2476765E+02	361	834.1880	0.3468013E+02	463	432.1267	0.2503906E+02
360	1228.5081	0.2477368E+02	362	826.4563	0.3371597E+02	464	424.4047	0.2503906E+02
361	1220.7764	0.2407227E+02	363	818.7246	0.3287764E+02	465	416.6826	0.2503906E+02
362	1213.0447	0.2362061E+02	364	810.9927	0.3220895E+02	466	408.9609	0.2503906E+02
363	1205.3127	0.2352295E+02	365	803.2613	0.3153697E+02	467	401.2057	0.2503906E+02
364	1197.5811	0.24441406E+02	366	795.5293	0.3141479E+02	468	393.4475	0.2503906E+02
120	0.1507636E+02	0.1423400F+02						
129	0.123865E+02	0.1380286E+02						
136	0.2535057E+02	0.1655345E+02						
141	0.3165837E+02	0.2456126E+02						
144	0.2221326E+02	0.2220396E+02						





11.41	11.29	11.19	11.10	11.02	10.62	10.38	14.90	16.53	16.04
11.74	13.72	21.34	21.85	15.96	66.93	79.39	68.17	52.35	37.43
19.54	22.35	24.31	23.65	24.36	0.0	0.0	0.0	0.0	0.0
35.25	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
124164									
130136141146150									
8.98	8.74	8.55	8.35	8.19					
10.11	9.41	11.25	19.25	19.64	12.95	7.53	11.97	13.28	18.25
20.69	20.29	22.16	62.85	73.31	61.66	56.07	31.03	31.63	0.0
0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
123164									

FOURIER TRANSFORM INFRARED SPECTROSCOPY -- NORMALIZED PEAK HEIGHT INFORMATION : FILE NUMBER 16

\*\*\* - DENOTES A VALID PEAK.  
CB. - DENOTES A PEAK MEASURED FROM A "COMMON BASELINE".

PEAK NUMBER	DATA WORD	WAVE NUMBER	PEAK HEIGHT	BASELINE AMPLITUDE
1	120	3096.1021	1.0579338	14.0184288
*** 2 CB.	130	3008.9502	3.6027365	13.6160135
*** 3 CB.	136	2962.6589	11.8947153	13.4558554
*** 4 CB.	141	2924.0829	18.3359680	13.3223896
*** 5 CB.	146	2885.5068	9.0243320	13.1889248
*** 6 CB.	150	2854.6460	10.7229767	13.0821533
7	166	2731.2026	0.1757608	12.5525799
8	169	2708.0569	0.1459227	12.3768101
*** 9	294	1743.6558	21.6212769	14.8127050
*** 10	306	1651.0732	1.9715034	16.0298157
11	314	1589.3516	0.1426058	15.5174437
12	316	1573.9211	0.3681087	15.7031651
*** 13	323	1510.9148	0.7959442	19.1422577
*** 14	336	1419.6160	45.9810103	21.2878639
15	349	1319.3193	0.5129395	22.5493164
16	352	1296.1736	0.1556339	21.7522583
*** 17	359	1242.1672	2.4651947	22.3481140
*** 18	393	1057.0022	44.3265228	59.8350220
*** 19	394	972.1350	20.6452484	70.6552987
*** 20	398	941.2742	6.4769460	60.4878387
*** 21	401	918.1284	11.9130341	55.5032806
22	419	779.2566	1.2622736	32.8535912
*** 23	425	732.9634	5.4464474	31.4053650



FTIRVIEW TRANSFORM INFRARED SPECTROSCOPY --- NORMALIZED PEAK HEIGHT INFORMATION : FILE NUMBER 20

\*\*\* - DENOTES A VALID PEAK.

CH. - DENOTES A PEAK MEASURED FROM A "COMMON BASELINE".

PEAK NUMBER	DATA WORD	WAVE NUMBER	PEAK HEIGHT	BASFLINE AMPLITUDE
1	120	3084.1318	0.9578886	11.7398720
*** 2 CB.	130	3006.8142	3.4200430	11.4080601
*** 3 CB.	136	2960.4236	12.4571257	11.2908144
*** 4 CB.	141	2921.7649	18.6699829	11.1931105
*** 5 CB.	146	2892.1060	9.1964798	11.0954056
*** 6 CB.	150	2852.1790	10.7229776	11.0172415
7	166	2728.4707	0.2328259	10.6185064
8	170	2697.5437	0.2030725	10.3796721
*** 9	294	1738.8047	23.0211334	14.8958855
10	306	1646.0234	1.6773682	16.5301514
11	312	1590.6328	0.0398047	16.0350342
*** 12	323	1514.5835	0.5473633	19.5836029
*** 13	336	1414.0706	49.9136200	22.3494263
14	349	1311.5576	0.5134430	24.3073425
15	352	1290.3623	0.1609344	23.6483154
*** 16	359	1236.2400	2.5067139	24.3601227
*** 17	343	1050.6775	49.9581604	66.9271393
*** 18	394	965.6287	23.0735474	79.3897858
*** 19	398	934.7009	9.2646484	68.1747894
*** 20	401	911.5054	13.2904053	62.8501740
21	419	772.3340	1.0077057	37.4304047
*** 22	425	725.9434	5.3981179	35.2491309

FT/IR TRANSFORM INFRARED SPECTROSCOPY -- NORMALIZED PEAK HEIGHT INFORMATION : FILE NUMBER 24

\*\*\* - DENOTES A VALID PEAK.

CR. - DENOTES A PEAK MEASURED FROM A "COMMON BASELINE".

PEAK NUMBER	DATA WORD	WAVE NUMBER	PEAK HEIGHT	BASELINE AMPLITUDE
1	114	3130.5225	0.2809744	11.1132374
2	120	3084.1318	1.0309410	10.5555153
*** 3 CR.	130	3006.8142	3.5085621	10.2147713
*** 4 CR.	136	2960.4236	12.2890711	10.1000738
*** 5 CR.	141	2921.7649	18.7878723	10.0044928
*** 6 CR.	146	2883.1060	9.1824494	9.9089108
*** 7 CR.	150	2852.1790	10.7229738	9.8324461
8	166	2723.4707	0.2107296	9.4754591
9	288	1785.1953	0.0813360	12.3480406
*** 10	294	1738.8047	23.2737274	13.5707464
11	303	1669.2187	0.1996393	15.3943853
12	306	1646.0234	1.4019089	15.1311569
13	311	1607.3647	0.0837904	14.4159079
*** 14	323	1514.5835	0.6020966	18.9920959
*** 15	336	1414.0706	48.4013824	20.9324951
16	345	1313.5576	0.4485779	21.9849243
17	352	1290.3623	0.1700592	21.3391876
*** 18	359	1236.2400	2.5230255	22.5427551
*** 19	393	1050.6775	49.0267639	62.6246460
*** 20	394	965.6232	22.0711517	73.8221130
*** 21	398	934.7009	7.3385225	62.5572357
*** 22	401	911.5056	12.4565277	57.4393095
23	419	772.3340	1.3028412	33.9114319
*** 24	425	725.9434	5.4251251	32.2548828

## FOURIER TRANSFORM INFRARED SPECTROSCOPY -- NORMALIZED PEAK HEIGHT INFORMATION : FILE NUMBER 28

\*\*\* - DENOTES A VALID PEAK.  
 CB. - DENOTES A PEAK MEASURED FROM A "COMMON BASELINE".

PEAK NUMBER	DATA WORD	WAVE NUMBER	PEAK HEIGHT	BASELINE AMPLITUDE
1	113	3130.2542	0.1907387	10.1058922
2	120	3084.1318	0.8087158	9.4085693
*** 3 CB.	130	3006.8142	3.3512774	8.9778242
*** 4 CB.	136	2965.4226	13.1085663	8.7420197
*** 5 CB.	141	2921.7549	19.1827545	8.5455170
*** 6 CB.	146	2882.1060	9.4183207	8.3490133
*** 7 CB.	150	2852.1790	10.7229824	8.1918125
8	166	2728.4707	0.2217607	7.5297041
*** 9	294	1738.8047	25.2192535	11.9684595
10	306	1546.0214	1.7920082	13.2760191
*** 11	336	1414.0706	50.0033264	18.2583923
12	349	1313.5170	0.5696716	20.6949768
13	352	1200.3623	0.1495514	20.2850189
*** 14	359	1036.2400	2.5153351	22.1611786
*** 15	383	1050.6775	44.1548920	62.8510440
*** 16	394	965.4282	21.4562531	73.3130341
*** 17	398	924.7002	7.9833994	61.6577149
*** 18	401	911.5056	13.0371399	56.0729675
19	419	772.3340	1.2115479	11.0272217
*** 20	421	755.3424	3.5814200	11.6297607

FILE \*\*\*\*\* PEAK HEIGHT TABULATION FOR ALL THE FILES \*\*\*\*\*

DATA WORD	FILE 16	FILE 20	FILE 24	FILE 28	FILE 0	FILE 0	FILE 0	FILE 0	FILE 0	FILE 0
113	0.0	0.0	0.0	0.191	0.0	0.0	0.0	0.0	0.0	0.0
114	0.0	0.0	0.281	0.0	0.0	0.0	0.0	0.0	0.0	0.0
120	1.058	0.958	1.031	0.909	0.0	0.0	0.0	0.0	0.0	0.0
130	3.603	3.420	3.505	3.351	0.0	0.0	0.0	0.0	0.0	0.0
136	11.895	12.457	12.289	13.109	0.0	0.0	0.0	0.0	0.0	0.0
141	18.336	18.670	18.748	19.183	0.0	0.0	0.0	0.0	0.0	0.0
146	9.024	9.196	9.182	9.418	0.0	0.0	0.0	0.0	0.0	0.0
150	10.723	10.723	10.723	10.723	0.0	0.0	0.0	0.0	0.0	0.0
166	0.176	0.133	0.211	0.222	0.0	0.0	0.0	0.0	0.0	0.0
169	0.146	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
170	0.0	0.203	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
288	0.0	0.0	0.081	0.0	0.0	0.0	0.0	0.0	0.0	0.0
294	21.621	23.021	23.274	25.219	0.0	0.0	0.0	0.0	0.0	0.0
303	0.0	0.0	0.200	0.0	0.0	0.0	0.0	0.0	0.0	0.0
306	1.972	1.677	1.402	1.702	0.0	0.0	0.0	0.0	0.0	0.0
311	0.0	0.0	0.084	0.0	0.0	0.0	0.0	0.0	0.0	0.0
312	0.0	0.050	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
314	0.143	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
316	0.368	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
323	0.796	0.547	0.602	0.0	0.0	0.0	0.0	0.0	0.0	0.0
336	45.982	49.314	48.401	50.003	0.0	0.0	0.0	0.0	0.0	0.0
349	0.513	0.513	0.449	0.570	0.0	0.0	0.0	0.0	0.0	0.0
352	0.154	0.161	0.170	0.150	0.0	0.0	0.0	0.0	0.0	0.0
355	2.465	2.507	2.523	2.515	0.0	0.0	0.0	0.0	0.0	0.0
383	44.327	49.958	48.027	44.155	0.0	0.0	0.0	0.0	0.0	0.0

***** PEAK HEIGHT TABULATION FOR ALL THE FILES *****									
DATA	FILE	FILE	FILE	FILE	FILE	FILE	FILE	FILE	FILE
WORD	16	20	24	28	0	0	0	0	0
394	20.645	23.074	22.071	21.456	0.0	0.0	0.0	0.0	0.0
398	6.477	8.265	7.339	7.983	0.0	0.0	0.0	0.0	0.0
401	11.813	13.290	12.457	13.037	0.0	0.0	0.0	0.0	0.0
419	1.252	1.008	1.303	1.212	0.0	0.0	0.0	0.0	0.0
425	5.446	5.388	5.425	3.581	0.0	0.0	0.0	0.0	0.0
NORMAL END OF PEAK HEIGHT TABLE :									
TOTAL NUMBER OF FILES LISTED = 4.									

# F490 - FTIS REGRESSION ANALYSIS

NO. OF VARIABLES	16	NO. OF DEPENDENT VARIABLES	4
NO. OF OBSERVATIONS	4		
F LEVEL TO ENTER VARIABLE	4.060	F LEVEL TO REMOVE VARIABLE	3.780
TRANSFORMED DATA VALUES.			
OBSERVATION / FTIS FILE NO.			
1/16	3.60274 VAR( 1) 45.98192 VAR( 6) 11.81303 VAR(11) 237.00000 VAR(16)	11.89472 VAR( 2) 2.46519 VAR( 7) 5.44649 VAR(12)	18.33597 VAR( 3) 44.32652 VAR( 8) 1284.00000 VAR(13)
2/20	3.42004 VAR( 1) 42.91362 VAR( 6) 13.79041 VAR(11) 265.00000 VAR(16)	12.45713 VAR( 2) 2.50671 VAR( 7) 5.38926 VAR(12)	18.66998 VAR( 3) 49.95816 VAR( 8) 1467.00000 VAR(13)
3/24	3.50856 VAR( 1) 48.40138 VAR( 6) 12.45653 VAR(11) 272.00000 VAR(16)	12.28907 VAR( 2) 2.52303 VAR( 7) 5.42513 VAR(12)	18.78787 VAR( 3) 48.02676 VAR( 8) 1437.00000 VAR(13)
4/28	3.35128 VAR( 1) 50.00333 VAR( 6) 13.03714 VAR(11) 274.00000 VAR(16)	13.10857 VAR( 2) 2.51534 VAR( 7) 3.59142 VAR(12)	19.18275 VAR( 3) 44.15499 VAR( 8) 1619.00000 VAR(13)
			9.02433 VAR( 4) 20.64525 VAR( 9) 0.21400 VAR(14)
			21.62128 VAR( 5) 6.47694 VAR(10) 0.21100 VAR(15)
			9.19645 VAR( 4) 23.07355 VAR( 9) 0.23100 VAR(14)
			23.02113 VAR( 5) 8.26465 VAR(10) 0.22900 VAR(15)
			9.18249 VAR( 4) 22.07115 VAR( 9) 0.25100 VAR(14)
			23.27373 VAR( 5) 7.33859 VAR(10) 0.25100 VAR(15)
			9.41832 VAR( 4) 21.45625 VAR( 9) 0.21500 VAR(14)
			25.21925 VAR( 5) 7.99340 VAR(10) 0.20900 VAR(15)

MEAN	STD. DEV.	VARIABLE	DATA WORD / PHYS. PROP.
3.4707	0.1090E+00	1	130
12.4374	0.5058E+00	2	136
18.7441	0.3496E+00	3	141
9.2054	0.1624E+00	4	146
23.2838	0.1481E+01	5	294
48.5751	0.1878E+01	6	336
2.5026	0.2581E-01	7	359
46.6166	0.2855E+01	8	383
21.8115	0.1024E+01	9	394
7.5159	0.7931E+00	10	398
12.6493	0.6577E+00	11	401
4.9603	0.9196E+00	12	425
1451.7500	0.1373E+03	13	MODULUS
0.2277	0.1735E-01	14	STRAIN AT BREAK
0.2247	0.1947E-01	15	STRAIN AT MAXIMUM STRESS
242.0000	0.1711E+02	16	MAXIMUM STRESS

SIMPLE CORRELATION COEFFICIENTS. (ROW R' C. J.)

1.0000	R(1,1)	-0.9571	R(1,2)	-0.8816	R(1,3)	-0.9358	R(1,4)	-0.9030	R(1,5)
-0.9563	R(1,6)	-0.6975	R(1,7)	-0.1330	R(1,8)	-0.4932	R(1,9)	-0.9153	R(1,10)
-0.9133	R(1,11)	0.7472	R(1,12)	-0.9676	R(1,13)	0.0744	R(1,14)	0.1405	R(1,15)
-0.9054	R(1,16)								
1.0000	R(2,2)	0.9581	R(2,3)	0.9928	R(2,4)	0.9780	R(2,5)	0.8465	R(2,6)
0.6615	R(2,7)	-0.1275	R(2,8)	0.7453	R(2,9)	0.7627	R(2,10)	0.7568	R(2,11)
-0.8856	R(2,12)	0.9877	R(2,13)	-0.1700	R(2,14)	-0.2348	R(2,15)	0.7779	R(2,16)
1.0000	R(3,3)	0.9819	R(3,4)	0.9959	R(3,5)	0.8014	R(3,6)	0.7924	R(3,7)
-0.1103	R(3,8)	0.2225	R(3,9)	0.6798	R(3,10)	0.6600	R(3,11)	-0.8447	R(3,12)
0.9725	R(3,13)	0.0444	R(3,14)	-0.0187	R(3,15)	0.8738	R(3,16)		
1.0000	R(4,4)	0.9926	R(4,5)	0.8376	R(4,6)	0.7160	R(4,7)	-0.1230	R(4,8)
0.2393	R(4,9)	0.7375	R(4,10)	0.7267	R(4,11)	-0.4841	R(4,12)	0.9897	R(4,13)
-0.9927	R(4,14)	-0.1574	R(4,15)	0.8193	R(4,16)				
1.0000	R(5,5)	0.9062	R(5,6)	0.7459	R(5,7)	-0.1460	R(5,8)	0.2018	R(5,9)
0.6212	R(5,10)	0.6757	R(5,11)	-0.8798	R(5,12)	0.9815	R(5,13)	-0.0274	R(5,14)
-0.1007	R(5,15)	0.9408	R(5,16)						
1.0000	R(6,6)	0.9341	R(6,7)	0.4133	R(6,8)	0.7215	R(6,9)	0.9827	R(6,10)
0.9751	R(6,11)	-0.5295	R(6,12)	0.9007	R(6,13)	0.1816	R(6,14)	0.1181	R(6,15)
0.9740	R(6,16)								
1.0000	R(7,7)	0.4384	R(7,8)	0.6263	R(7,9)	0.7199	R(7,10)	0.6851	R(7,11)
-0.3479	R(7,12)	0.7674	R(7,13)	0.6247	R(7,14)	0.5719	R(7,15)	0.9945	R(7,16)
1.0000	R(8,8)	0.7267	R(8,9)	0.5119	R(8,10)	0.5025	R(8,11)	0.5543	R(8,12)
0.9022	R(8,13)	0.7169	R(8,14)	0.7084	R(8,15)	0.3537	R(8,16)		
1.0000	R(9,9)	0.7972	R(9,10)	0.7999	R(9,11)	0.2066	R(9,12)	0.3595	R(9,13)
0.5762	R(9,14)	0.5443	R(9,15)	0.5350	R(9,16)				
1.0000	R(10,10)	0.9097	R(10,11)	-0.4164	R(10,12)	0.8151	R(10,13)	0.1676	R(10,14)
0.1014	R(10,15)	0.7869	R(10,16)						
1.0000	R(11,11)	-0.0170	R(11,12)	0.8037	R(11,13)	0.1214	R(11,14)	0.0628	R(11,15)
0.7569	R(11,16)								

1.0000	R(12,12)	-0.8250	R(12,13)	0.4826	R(12,14)	0.5330	R(12,15)	-0.4825	R(12,16)
1.0000	R(13,13)	-0.0157	R(13,14)	-0.0830	R(13,15)	0.8655	R(13,16)		
1.0000	R(14,14)	0.9978	R(14,15)	0.4864	R(14,16)				
1.0000	R(15,15)	0.4274	R(15,16)						
1.0000	R(16,16)								



TRIAL NUMBER 1 FOR VARIABLE (13)  
 PURE CONST.  $B(0) = 0.1452E+04$   
 COEFFICIENTS

0.0	B(1)	0.0	B(2)	0.0	B(3)	0.0	B(4)	0.0	B(5)
0.0	B(6)	0.0	B(7)	0.0	B(8)	0.0	B(9)	0.0	B(10)
0.0	B(11)	0.0	B(12)						

STANDARD ERROR OF COEFFICIENTS

0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
0.0	0.0	0.0							

STANDARD ERROR OF ESTIMATE  
 $0.1373E+03$

TRIAL NUMBER 2  
 VARIABLE GOING IN = 4  
 F LEVEL 143.0815  
 PURE CONST.  $B(0) = -0.6752E+04$   
 COEFFICIENTS

0.0	B(1)	0.0	B(2)	0.0	B(3)	0.0	B(4)	0.0	B(5)
0.0	B(6)	0.0	B(7)	0.0	B(8)	0.0	B(9)	0.0	B(10)
0.0	B(11)	0.0	B(12)						

STANDARD ERROR OF COEFFICIENTS

0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
0.0	0.0	0.0							

STANDARD ERROR OF ESTIMATE  
 $0.2410E+02$

TRIAL NUMBER 3  
 VARIABLE GOING IN = 6  
 F LEVEL 10.3992  
 PURE CONST.  $B(0) = -0.5539E+04$   
 COEFFICIENTS

0.0	B(1)	0.0	B(2)	0.0	B(3)	0.0	B(4)	0.0	B(5)
0.1755E+02	B(6)	0.0	B(7)	0.0	B(8)	0.0	B(9)	0.0	B(10)

STANDARD ERROR OF COEFFICIENTS

0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0	0.0
0.0	0.0	0.0							

STANDARD ERROR OF ESTIMATE  
 $0.1755E+02$

70

ACTUAL VS. PREDICTED RESULTS FOR VARIABLE (13)

OBSERVATION	ACTUAL	PREDICTED	DEVIATION
1	1284.0000	1285.4805	-1.4805
2	1467.0000	1469.3125	-2.3125
3	1437.0000	1433.4219	3.5781
4	1619.0000	1618.7930	0.2070

MEAN = 1451.75000      STD. DEV. = 137.311081      COEFFICIENT OF VARIATION = 0.094583

TRIAL NUMBER 1 FOR VARIABLE (14)

PURE CONST. B(0) = 0.2277E+00

COEFFICIENTS

0.0	B( 1)	0.0	B( 2)	0.0	B( 3)	0.0	B( 4)	0.0	B( 5)
0.0	B( 6)	0.0	B( 7)	0.0	B( 8)	0.0	B( 9)	0.0	B(10)
0.0	B(11)	0.0	B(12)						

STANDARD ERROR OF COEFFICIENTS

0.0		0.0		0.0		0.0		0.0	
0.0		0.0		0.0		0.0		0.0	
0.0		0.0							

STANDARD ERROR OF ESTIMATE

2.1735E-01

MULTIPLE CORRELATION COEFFICIENT

0.0

COEFFICIENTS

0.0	B( 1)	0.0	B( 2)	0.0	B( 3)	0.0	B( 4)	0.0	B( 5)
0.0	B( 6)	0.0	B( 7)	0.0	B( 8)	0.0	B( 9)	0.0	B(10)
0.0	B(11)	0.0	B(12)						

ACTUAL VS. PREDICTED RESULTS FOR VARIABLE (14)

OBSERVATION	ACTUAL	PREDICTED	DEVIATION
1	0.2140	0.2277	-0.0137
2	0.2310	0.2277	0.0033
3	0.2510	0.2277	0.0233
4	0.2150	0.2277	-0.0127
MEAN =	0.227750	STD. DEV. = 0.017347	COEFFICIENT OF VARIATION = 0.076165

TRIAL NUMBER 1 FOR VARIABLE (J5)  
 PURE CONST. B(0) = 0.2247E+00  
 COEFFICIENTS

0.0	B( 1)	0.0	B( 2)	0.0	B( 3)	0.0	B( 4)	0.0	B( 5)
0.0	B( 6)	0.0	B( 7)	0.0	B( 8)	0.0	B( 9)	0.0	B(10)
0.0	B(11)	0.0	B(12)						

STANDARD ERROR OF COEFFICIENTS

0.0		0.0	0.0		0.0				
0.0		0.0	0.0		0.0				
0.0		0.0	0.0		0.0				

STANDARD ERROR OF ESTIMATE  
 0.1947E-01

MULTIPLE CORRELATION COEFFICIENT  
 0.0  
 COEFFICIENTS

0.0	B( 1)	0.0	B( 2)	0.0	B( 3)	0.0	B( 4)	0.0	B( 5)
0.0	B( 6)	0.0	B( 7)	0.0	B( 8)	0.0	B( 9)	0.0	B(10)
0.0	B(11)	0.0	B(12)						

# ACTUAL VS. PREDICTED RESULTS FOR VARIABLE (15)

OBSERVATION	ACTUAL	PREDICTED	DEVIATION
1	0.2110	0.2247	-0.0137
2	0.2280	0.2247	0.0033
3	0.2510	0.2247	0.0263
4	0.2090	0.2247	-0.0157
MEAN =	0.224750	STD. DEV. = 0.019466	COEFFICIENT OF VARIATION = 0.086610

TRIAL NUMBER 1 FOR VARIABLE (16)

PURE CONST.  $B(0) = 0.2620F+03$

COEFFICIENTS

0.0	B(1)	0.0	B(2)	0.0	B(3)	0.0	B(4)	0.0	B(5)
0.0	B(6)	0.0	B(7)	0.0	B(8)	0.0	B(9)	0.0	B(10)
0.0	B(11)	0.0	B(12)						

STANDARD ERROR OF COEFFICIENTS

0.0	0.0	0.0	0.0	0.0
0.0	0.0	0.0	0.0	0.0
0.0	0.0			

STANDARD ERROR OF ESTIMATE

0.1711E+02

TRIAL NUMBER 2

VARIABLE GOING IN = 7

F LEVEL

94.3252

PURE CONST.  $B(0) = -0.1371E+04$

COEFFICIENTS

0.0	B(1)	0.0	B(2)	0.0	B(3)	0.0	B(4)	0.0	B(5)
0.0	B(6)	0.4525E+03	B(7)	0.0	B(8)	0.0	B(9)	0.0	B(10)
0.0	B(11)	0.0	B(12)						

STANDARD ERROR OF COEFFICIENTS

0.0	0.0	0.0	0.0	0.0
0.0	0.8228E+02	0.0	0.0	0.0
0.0	0.0			

STANDARD ERROR OF ESTIMATE

0.3679E+01

TRIAL NUMBER 3

VARIABLE GOING IN = 2

F LEVEL

25.2573

PURE CONST.  $B(0) = -0.1219E+04$

COEFFICIENTS

0.0	B(1)	0.7621E+01	B(2)	0.0	B(3)	0.0	B(4)	0.0	B(5)
0.0	B(6)	0.5537E+03	B(7)	0.0	B(8)	0.0	B(9)	0.0	B(10)



0.0 R(11) 0.0 B(12)

STANDARD ERROR OF COEFFICIENTS

0.0 0.2145E+01 0.0 0.0  
0.0 0.4203E+02 0.0 0.0  
0.0 0.0 0.0 0.0

STANDARD ERROR OF ESTIMATE

0.1409E+01

MULTIPLE CORRELATION COEFFICIENT

0.99867

COEFFICIENTS

0.0 B(1) 0.7621E+01 B(2) 0.0 0.0 0.0  
0.0 R(6) 0.553/E+03 R(7) 0.0 0.0 0.0  
0.0 B(11) 0.0 R(12) 0.0 0.0 0.0  
0.0 0.0 0.0 0.0 0.0 0.0

ACTUAL VS. PREDICTED RESULTS FOR VARIABLE (16)

OBSERVATION	ACTUAL	PREDICTED	DEVIATION
1	237.0000	237.1721	-0.1721
2	265.0000	264.4463	0.5537
3	272.0000	272.1968	-0.1968
4	274.0000	274.1846	-0.1846

MEAN = 262.000000 STN. DEV. = 1/.107498 COEFFICIENT OF VARIATION = 0.065296

\*\*\*\*\* FTIS - NONLINEAR REGRESSION MODE FOR \*\*\*\*\*  
 DEPENDENT VARIABLE = VAR(13)

NO. OF VARIABLES	5	NO. OF DEPENDENT VARIABLES	1
NO. OF OBSERVATIONS	4		
F LEVEL TO ENTER VARIABLE	8.530	F LEVEL TO REMOVE VARIABLE	5.540
TRANSFORMED DATA VALUES.			
OBSERVATION / FTIS FILE NO.			
1/16	9.02433 VAR( 1)	81.43857 VAR( 2)	45.98192 VAR( 3)
2/20	9.19648 VAR( 1)	84.57523 VAR( 2)	49.91362 VAR( 3)
3/24	9.18249 VAR( 1)	84.31810 VAR( 2)	49.47138 VAR( 3)
4/28	5.41832 VAR( 1)	88.70476 VAR( 2)	50.00333 VAR( 3)
			2114.33667 VAR( 4)
			2491.36938 VAR( 4)
			2342.69360 VAR( 4)
			2500.33252 VAR( 4)
			1284.00000 VAR( 5)
			1467.00000 VAR( 5)
			1437.00000 VAR( 5)
			1619.00000 VAR( 5)

MEAN	STD.. DEV.	VARIABLE
9.2054	0.1624E+00	1
84.7592	0.2790E+01	2
48.5751	0.1878E+01	3
23.21416	0.1804E+03	4
1451.7500	0.1175E+03	5

SIMPLE CORRELATION COEFFICIENTS. (ROW B, COL.)

1.0000	R( 1, 1)	0.9982	R( 1, 2)	0.8376	R( 1, 3)	0.8379	R( 1, 4)	0.9897	R( 1, 5)
1.0000	R( 2, 2)	0.9352	R( 2, 3)	0.8365	R( 2, 4)	0.9910	R( 2, 5)		
1.0000	R( 3, 3)	1.0003	R( 3, 4)	0.9007	R( 3, 5)				
1.0000	R( 4, 4)	0.9010	R( 4, 5)						
1.0000	R( 5, 5)								

TRIAL NUMBER 1 FOR VARIABLE ( 5 )  
 PURE CONST. B(0) = 0.1452E+04  
 COEFFICIENTS  
 0.0 B( 1) 0.0 B( 2) 0.0 B( 3) 0.0 B( 4)  
 STANDARD ERROR OF COEFFICIENTS  
 0.0 0.0 0.0 0.0  
 STANDARD ERROR OF ESTIMATE  
 0.1373E+03

TRIAL NUMBER 2  
 VARIABLE GOING IN = 2  
 F LEVEL 163.9341  
 PURE CONST. B(0) = -0.2406E+04  
 COEFFICIENTS  
 0.0 B( 1) 0.4551E+02 B( 2) 0.0 B( 3) 0.0 B( 4)  
 STANDARD ERROR OF COEFFICIENTS  
 0.0 0.4353E+01 0.0 0.0  
 STANDARD ERROR OF ESTIMATE  
 0.2254E+02

TRIAL NUMBER 3  
 VARIABLE GOING IN = 3  
 F LEVEL 102.2383  
 PURE CONST. B(0) = -0.2418E+04  
 COEFFICIENTS  
 0.0 B( 1) 0.3625E+02 B( 2) 0.1765E+02 B( 3) 0.0 B( 4)  
 STANDARD ERROR OF COEFFICIENTS  
 0.0 0.1551E+01 0.2469E+01 0.0  
 STANDARD ERROR OF ESTIMATE  
 0.4416E+01

MULTIPLE CORRELATION COEFFICIENT

0.99983

\*\*\*\* CONCLUSION - NON-LINEAR REGRESSION MODEL PROVIDES THE BEST FIT FOR VARIABLE (13) -  
\*\*\*\*\*

# ACTUAL VS. PREDICTED RESULTS FOR VARIABLE ( 5 )

OBSERVATION	ACTUAL	PREDICTED	DEVIATION
1	1284.0000	1295.6121	-1.6121
2	1467.0000	1468.7078	-1.7078
3	1437.0000	1437.6968	4.3032
4	1619.0000	1619.9836	-0.9836
MEAN =	1451.7500	STD. DEV. = 137.311081	COEFFICIENT OF VARIATION = 0.094583

\*\*\*\*\* FTIS - NONLINEAR REGRESSION MODE FOR \*\*\*\*\*  
 DEPENDENT VARIABLE = VAR(16)

NO. OF VARIABLES	5	NO. OF DEPENDENT VARIABLES	1
NO. OF OBSERVATIONS	4		
F LEVEL TO ENTER VARIABLE	9.530	F LEVEL TO REMOVE VARIABLE	5.540
TRANSFORMED DATA VALUES.			
OBSERVATION / FTIS FILE NO.			
1/16	11.87472 VAR( 1)	141.48424 VAR( 2)	2.46519 VAR( 3)
2/20			6.07719 VAR( 4)
3/24	12.45713 VAR( 1)	155.17998 VAR( 2)	2.50671 VAR( 3)
			6.28361 VAR( 4)
4/28	12.28907 VAR( 1)	151.07126 VAR( 2)	2.52303 VAR( 3)
			6.36566 VAR( 4)
	13.10857 VAR( 1)	171.83450 VAR( 2)	2.51534 VAR( 3)
			6.32691 VAR( 4)
			237.00000 VAR( 5)
			265.00000 VAR( 5)
			272.00000 VAR( 5)
			274.00000 VAR( 5)



MEAN	STD. DEV.	VARIABLE	
12.4374	0.5058E+00	1	
154.8799	0.1267E+02	2	
2.5026	0.2581E-01	3	
6.2633	0.1286E+00	4	
262.0000	0.1711E+02	5	
SIMPLE CORRELATION COEFFICIENTS. (ROW BY COL.)			
1.0000 R( 1, 1)	1.0001 R( 1, 2)	0.6615 R( 1, 3)	0.6619 R( 1, 4)
1.0000 R( 2, 2)	0.6515 R( 2, 3)	0.6514 R( 2, 4)	0.7698 R( 2, 5)
1.0000 R( 3, 3)	0.9991 R( 3, 4)	0.9845 R( 3, 5)	
1.0000 R( 4, 4)	0.9849 R( 4, 5)		
1.0000 R( 5, 5)			
			0.7779 R( 1, 5)

TRIAL NUMBER 1 FOR VARIABLE ( 5 )  
 PURE CONST. B(0) = 0.2620E+03  
 COEFFICIENTS  
 0.0 B( 1 ) 0.0 B( 2 ) 0.0 B( 3 ) 0.0 B( 4 )  
 STANDARD ERROR OF COEFFICIENTS  
 0.0 0.0 0.0 0.0  
 STANDARD ERROR OF ESTIMATE  
 0.1711E+02

TOTAL NUMBER 2  
 VARIABLE GOING IN = 4  
 F LEVEL 97.1742  
 PURE CONST. A(0) = -0.5587E+03  
 COEFFICIENTS  
 0.0 A( 1 ) 0.0 B( 2 ) 0.0 B( 3 ) 0.1310E+03 B( 4 )  
 STANDARD ERROR OF COEFFICIENTS  
 0.0 0.0 0.0 0.1628E+02  
 STANDARD ERROR OF ESTIMATE  
 0.3626E+01

TRIAL NUMBER 3  
 VARIABLE GOING IN = 2  
 F LEVEL 41.1636  
 PURE CONST. B(0) = -0.4844E+03  
 COEFFICIENTS  
 0.0 A( 1 ) 0.3006E+00 B( 2 ) 0.0 B( 3 ) 0.1117E+03 B( 4 )  
 STANDARD ERROR OF COEFFICIENTS  
 0.0 0.6627E-01 0.0 0.6531E+01  
 STANDARD ERROR OF ESTIMATE  
 0.1104E+01

MULTIPLE CORRELATION COEFFICIENT

0.99931

\*\*\* CONCLUSION : NON-LINEAR REGRESSION MODEL PROVIDES THE BEST FIT FOR VARIABLE (16) .  
\*\*\*\*\*

ACTUAL VS. PREDICTED RESULTS FOR VARIABLE ( 5 )

OBSERVATION	ACTUAL	PREDICTED	DEVIATION
1	237.0000	237.1731	-0.1731
2	265.0000	264.3555	0.6445
3	272.0000	272.2717	-0.2717
4	274.0000	274.2000	-0.2000
MEAN =	262.000000	STD. DEV. = 17.107498	COEFFICIENT OF VARIATION = 0.065296

## 6.0 E410 - COMPUTER PROGRAM DESCRIPTION

E410 is designed to create and maintain a permanent file for solid propellant mechanical properties. Once established, this file will be utilized for optional input to the E490 computer program.

This FORTRAN code consists of two parts, a main and a subroutine called RENAME. The main program creates the initial file and prints out data contained in the file. RENAME is designed strictly for updating or accessing an existing propellant file.

## 7.0 E410 - INPUT INSTRUCTIONS

E410 is always executed utilizing one of two possible run modes; the CREATE mode for building a new propellant file or the UPDATE mode designed to update or access an existing file.

### 7.0.1 CREATE MODE

The first card in the input deck must have CREATE in the first six columns. This tells the computer to stay in the main program.

The second card is the title of the chart that is printed when the data set is established. The title can be 40 alphanumeric characters long. It is read in on a 10A4 format.

The next sequence of three cards is repeated for each set of data. The first 15 columns on data card type 1 are ignored. Any comments with alphanumeric data can be placed here. The next five columns are allocated to the specimens' age and the next two sets of ten columns are allocated for the aging temperature and test temperature respectively. This card is read with an (A8,7X,I5,2I10) format.

Data cards 2 & 3 are used to input the values for modulus, strain at break, strain at maximum stress, maximum stress, strain energy density and strain endurance respectively. Ten columns are allocated for each of these plus four extra variables that can be allocated later. Five variables are read from each of

these cards. Only the first 50 columns are read on each. All five fields are end to end. The format is (5(F10.3)) for both cards.

#### 7.0.2 UPDATE MODE

In this mode, the first card in the input deck must have UPDATE in the first six columns. Also, the title card is eliminated and the first of the three data cards is changed. Data card 1 in this mode should contain the action to be performed in the first 8 columns, skip the next two columns, and have the line number that the action is to be performed on in the next five columns. These two new inputs replace the columns ignored on the first data card in the CREATE mode. The remainder of this card and the two following data cards are the same as in the CREATE mode.

All actions should be nested together by type, and these groups then input ONLY in the following sequence;

1. All CHANGE cards.
2. All ADD cards.
3. All DELETE cards.

If there are no CHANGE cards in the input, then the proper input sequence should be all the ADD cards first, then any DELETE cards next. This correct sequencing of the input cards is mandatory in order to maintain the integrity of the existing properties file. Care should also be taken to put the proper data cards after the CHANGE and ADD action cards; the other action cards do not require any data other than the action and line number. All strings of actions need to be followed by either ENDLIST or ENDNOLST.

<u>Action</u>	<u>Purpose</u>
CHANGE - -	Changes data on specified line to the data following this action. [Requires data and line number]
ADD - -	Adds a line to data set. It will place the data on the line specified and move everything else down one record number. [Requires data and line number]
DELETE - -	Deletes line from data. [Requires line number only]
ENDLIST - -	Marks the end of all update actions and specifies that a listing of the data is desired. This action can be used as the only member in the string of actions if only a listing is desired. [Requires no other information]
ENDNOLST - -	Marks the end of all update actions and specifies no listing of data desired. [Requires no other information]

TABLE 3: TYPES OF ACTIONS PERFORMED ON PROPELLANT FILE DATA.



## 8.0 E410 - SAMPLE CASE

The card deck shown in Figure 5 is the CDC Job Control Language and input data that was used to create the propellant properties file listed in Table 4.

Card Column 1

E410 SAMPLE CASE

---

```

SEQ,THA01.
THA,T100,P4. 1A,SMITHO,30
ATTACH,LGO,F410GO,IO=SMITHO,WR=1.
MAP,PART.
REQUEST,TAPER8,*PF.
LGO.
CATALOG,TAPER8,TPH8156V613,IO=SMITHO,RP=999.
CREATE
TEST CASE --- SEPT. 27, 1970.
ADD      0      77      77
1119.0    0.0    24.00  165.4    0.0
0.0      0.0    0.0    0.0      0.0
ADD      2      190      77
0843.0    0.0    21.00  126.7    0.0
0.0      0.0    0.0    0.0      0.0
ADD      2      170      77
774.0     0.0    22.6   121.2    0.0
0.0      0.0    0.0    0.0      0.0
ADD      2      150      77
845.0     0.00   21.00   116.1    0.
0.0      0.0    0.0    0.0      0.0
ADD      2      130      77
773.0     0.0    23.3   117.8    0.0
0.0      0.0    0.0    0.0      0.0
ADD      2      110      77
932.0     0.0    24.3   144.9    0.0
0.0      0.0    0.0    0.0      0.0
ADD      2      75      77
1171.0    0.0    21.3   151.8    0.0
0.0      0.0    0.0    0.0      0.0
ADD      4      190      77
1428.     0.0    15.8   155.2    0.0
0.0      0.0    0.0    0.0      0.0
ADD      4      170      77
1041.0    0.0    15.0   134.0    0.0
0.0      0.0    0.0    0.0      0.0

```

Figure 5. E410 - SAMPLE CASE 93

ADD	4	150	77		
767.0	0.0	22.7	120.55	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	4	130	77		
787.0	0.0	23.9	122.3	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	4	110	77		
912.0	0.0	24.9	150.2	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	4	75	77		
1187.0	0.0	24.1	153.4	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	8	190	77		
2939.0	0.0	5.6	214.6	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	8	170	77		
1401.0	0.0	15.1	149.6	0.00	
0	0.0	0.0	0.0	0.0	
ADD	8	150	77		
904.0	0.0	19.0	122.6	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	8	130	77		
764.0	0.0	23.2	121.8	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	8	110	77		
1042.0	0.0	21.7	146.0	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	8	75	77		
985.0	0.0	25.5	152.7	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	13	190	77		
3496.0	0.0	8.3	225.6	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	13	170	77		
1788.0	0.0	13.5	171.0	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	13	150	77		
1017.0	0.0	19.4	132.1	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	13	130	77		
729.0	0.0	25.7	124.4	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	13	110	77		
972.0	0.0	21.8	136.0	0.0	
0.0	0.0	0.0	0.0	0.0	
ADD	13	75	77		
1058.0	0.0	23.6	150.9	0.0	
0.0	0.0	0.0	0.0	0.0	

Figure 5: E410 - SAMPLE CASE (Cont.)

## FOURIER TRANSFORM INFRARED SPECTROSCOPY - E410 PHYSICAL PROPERTIES MASTER TAPE GENERATOR (AFRPL/PCD - THINOKOL/HUNTSVILLE)

INVESTIGATIONS FOR SOLID PROPELLANT

TEST CASE --- 5405. 27, 1519.

NO.	REC AGE		AGE TEST TEMP	MODULUS (PSI)	STRAIN AT BREAK (%)	MAXIMUM STRESS (PSI)	MAX AUM STRESS	STRAIN ENERGY DENSITY (TRN)	STRAIN ENDURANCE (%)
	(WK)	(F)							
1	0	77	77	1119.000	0.0	24.000	165.400	0.0	0.0
2	2	150	77	943.000	0.0	21.000	126.700	0.0	0.0
3	2	170	77	774.000	0.0	22.600	121.200	0.0	0.0
4	2	150	77	945.000	0.0	21.000	116.100	0.0	0.0
5	2	120	77	773.000	0.0	23.300	117.800	0.0	0.0
6	2	110	77	932.000	0.0	24.300	144.900	0.0	0.0
7	2	75	77	1171.000	0.0	21.300	151.800	0.0	0.0
8	4	190	77	1429.000	0.0	15.900	155.200	0.0	0.0
9	4	170	77	1041.000	0.0	15.000	134.000	0.0	0.0
10	4	150	77	767.000	0.0	22.700	129.550	0.0	0.0
11	4	120	77	797.000	0.0	23.800	122.300	0.0	0.0
12	4	110	77	212.000	0.0	24.900	150.200	0.0	0.0
13	4	75	77	1197.000	0.0	24.100	163.400	0.0	0.0
14	8	90	77	2026.000	0.0	9.000	214.600	0.0	0.0
15	8	170	77	1401.000	0.0	15.100	140.600	0.0	0.0
16	8	150	77	904.000	0.0	19.300	123.600	0.0	0.0
17	8	110	77	703.000	0.0	23.200	121.800	0.0	0.0
18	8	110	77	1140.000	0.0	21.700	146.000	0.0	0.0
19	7	77	77	984.000	0.0	25.500	152.700	0.0	0.0
20	13	190	77	2454.000	0.0	9.000	225.900	0.0	0.0
21	13	170	77	1793.000	0.0	13.500	171.000	0.0	0.0
22	13	150	77	1317.000	0.0	16.400	133.100	0.0	0.0
23	13	120	77	120.000	0.0	25.700	124.400	0.0	0.0
24	13	110	77	923.000	0.0	21.800	134.000	0.0	0.0
25	13	75	77	1354.000	0.0	23.600	153.800	0.0	0.0

**TABLE 4: E410 SAMPLE CASE OUTPUT**

APPENDIX A

FORTRAN SOURCE LISTINGS

PROGRAM F490FIN(INPUT,OUTPUT,TAPE1=INPUT,TAPE6=OUTPUT,TAPE1,TAPE2,  
X,TAPE21,TAPE22,TAPE23,TAPE24,TAPE25,TAPE26,TAPE27,TAPE11,TAPE12,TA  
XPE13,TAPE14,TAPE15,TAPE10)

```

C*****
C
C      F490 - FOURIER TRANSFORM INFRARED SPECTROSCOPY PROGRAM.
C
C      THIOKOL CORPORATION / HUNTSVILLE , ALABAMA 35897
C
C      PRINCIPAL INVESTIGATOR    W. W. SCHWARZ
C                                TELEPHONE (205) - 492 - 8388
C
C      SCIENTIST PROGRAMMER    D. C. SMITH
C                                TELEPHONE (205) - 182 - 8215
C
C*****
C
C      READS RPPMS FT/S DATA TAPES.
C
C      SEPTEMBER , 1979.
C
C      FORTRAN IV - R EXTENDED LANGUAGE
C
C      IBM 4600 - (SERPL)
C
C*****
C
C      COMMON /COPPEL/ IPEFF(40),MPFILE(30),
C                      NUNSYS(10),NUMBER,NUMORD,NUMORDS,IREG
C      COMMON /SPELCP/ IFILE(30),NPT29(30)
C      COMMON /POINTS/ NWORDDS(30)
C      COMMON /PTINEF/ IMAX,EMIN,NORM
C      COMMON /TWINPK/ RLIMIT(6),RLIMIT(6),RFILE(6),NUMORD,NUMORDS
C
C*****
C      DIMENSION AMPL(468),WAVNOR(468),PRNAME(60),FAC(130)
C      DIMENSION PMORD(30),ATL(30),BNA(30),HW(30),DIFF(30)
C      DIMENSION RLIMIT(6),INFILE(30),MPRECS(30),NPHOTO(6),NPHYSP(10)
C      DIMENSION WSA(132),IHEAD(101),IDATA(471)
C      REAL RLIMIT(6)
C      AMPL(1) /ZFILE/ IPEFF,IOTAL,INFILE,MULTPK,NPHOTO,RLIMIT,RLIMIT,
C      VALID,WAVNOR,MPDISK,MPRECS,NPHYSP,IPLDT,NMPLDT
C      DATA INFILE /30*0/,
C      NPHOTO /6*0/,
C      IPLDT /0/,
C      RLIMIT /6*0.0/,
C      MPDISK /0/,
C      MPRECS /30*0/,

```



	GO TO 25	0000091
20	WRITE(N2,725) MULTPK	0000092
	DO 22 KK=1,MULTPK	0000093
	WRITE (N2,730) KK,NPHOTO(KK),LLIMIT(KK),RLIMIT(KK)	0000100
22	CONTINUE	0000101
25	WRITE(N2,981) WAVNOR	0000102
	V = VALID * 100.	0000021
	WRITE (N2,995) V	0000024
	IF (MPDISK.NE.0) GO TO 26	0000103
	WRITE (N2,735)	0000104
	GO TO 27	0000105
26	MPDISK = MPDISK * 10	0000106
	WRITE (N2,740) MPDISK	0000107
	WRITE (N2,745) MPREFS	0000108
	WRITE (N2,750)	0000109
		0000110
	ISP = 1	0000111
	DO 28 IS1=1,10	0000112
	IF (NPHYSP(IS1).EQ.0) GO TO 28	0000113
	IS2 = NPHYSP(IS1)	0000114
	IS3 = IS2 * 6	0000115
	IS4 = IS3 - 5	0000116
	IS7 = IS8 + 5	0000117
	ID = 0	0000118
	DO 29 IS6 = IS8,IS7	0000119
	PRREF(IS6) = PRNAME(IS4+ID)	0000120
	ID = ID + 1	0000121
28	CONTINUE	0000122
	ISP = ISP + 5	0000123
	WRITE (N2,760) IS* (PRNAME(IS5),IS5=IS4,IS3)	0000124
29	CONTINUE	0000125
		0000126
		0000127
27	IRIG = MPDISK	0000128
	IF (IRIG.EQ.1) WRITE(N2,982)	0000129
	IF (NMPLT.EQ.1) WRITE(N2,983)	0000130
	IF (IRIG.EQ.0.AND.NMPLT.EQ.0) WRITE(N2,984)	0000131
	WRITE(N2,490)	
	FORMAT(1H1)	
		0000141
	READ HEADER RECORD FROM FTIS FILE.	0000132
		0000133
	CALL XREAD(WSA,IHEAD,"(10I16I)",STATUS)	AUG. 79
		0000135
	TEST FOR MATCHING FILE NUMBER.	0000136
		0000137
	IF (IHEAD(1).EQ.INFILE(LLLOC+1)) GO TO 9	
	CALL XREAD(WSA,TDATA,"(47I16I)",STATUS)	0000459
	CALL XREAD(WSA,IHEAD,"(10I16I)",STATUS)	AUG. 79
	IF (STATUS.NE.-2) STOP "EOF MARKER NOT FOUND."	
	WRITE(6,985)	

```

995 FORMAT(7,2X,"EOF ENCOUNTERED",/)
GO TO 5
9 LLOC = LLOC + 1
WRITE(6,1000) IHEAD(1),STATUS,IHEAD
1000 FORMAT(1H1,3X,"FILE ",13,27X,"STATUS = ",F4.1,/,26(1X,4(120,5X),/)
X)
C
C      CALCULATION OF CONSTANTS.
C
IFILE(LLOC) = IHEAD(1)
NWORDS(LLOC) = IHEAD(5)
NSPR = IHEAD(6)
NEXP = IHEAD(7)
NW = NWORDS(LLOC)
C
TEP1 = IHEAD(10)
TEP2 = IHEAD(11)
TEP3 = IHEAD(12)
TEP4 = IHEAD(13)
F = (2.*NEXP)/324.6250
X1 = TEP1 + (TEP2/PKLOC)
X2 = TEP3 + (TEP4/PKLOC)
DELTA = (X2 - X1)/(NW - 1)
WRITE(N2,900) X1,X2,DELTA
C
CALL XHEAD(56,10DATA,"(4718161)",STATUS)
C
C      CALCULATION OF WAVE NUMBER AND AMPLITUDE.
C      (X) (Y)
C
DO 30 M=1,10
WAVENO(M) = X1 + (M-1)*DELTA
AMPLTD(M) = F * 10DATA(M+2)
40 CONTINUE
DO 10 K=1,940 IF1(LLOC)
K=0
50 K=K+1
IF K=52
X=X1+52
Y=Y1+52
CALL I(X2,950, K,WAVENO(K),AMPLTD(K),K1,WAVENO(K1),AMPLTD(K1),K2,WAVENO(K2),AMPLTD(K2))
IF (K.EQ.52) GO TO 55
IF (K.EQ.208) GO TO 51
IF (K.EQ.364) GO TO 70
GO TO 50
65 WRITE(12,940) IFILE(LLOC)
K=K+104
GO TO 50
70 IF (AMPLTD(K.EQ.0)) GO TO 65

```



DO 68 M=1,6	00001910
ALIMIT(M6) = LLIMIT(M6)	00001920
ZLIMIT(M6) = RLIMIT(M6)	00001930
MEILF(M6) = NPHOTO(M6)	00001940
68 CONTINUE	00001950
C	00001960
DO 69 K=1,MNTPK	00001970
IF (ALIMIT(K).GT.WAVNOR.AND.ZLIMIT(K).LT.WAVNOR) GO TO 69	00001980
GO TO 60	00001990
69 JAX = K	00002000
NORMC0 = 1	00002010
60 CONTINUE	00002020
C	00002030
C.....WRITE INFORMATION FOR THIS FILE (IFILE(I)) TO PRIMARY FILE 20.	00002040
C	00002050
C.....WRITE INFORMATION FOR THIS FILE (IFILE(I)) TO PRIMARY FILE 20.	00002060
C	00002070
C.....WRITE INFORMATION FOR THIS FILE (IFILE(I)) TO PRIMARY FILE 20.	00002080
C	00002090
67 WRITE(20) LLOC,AMPLTD,WAVEND	00002100
IF (IPLN.NE.1) GO TO 400	00002110
C PLWAVE SUBROUTINE RESIDES HERE...	
CALL PLWAVE(TTAP,IFILE(LLOC),NW,WAVEND,AMPLTD)	
C	00002140
400 CALL APALC(AMPLTD,WAVEND,WAVNDR,LLOC,JAX,AN,W,BA,AP22)	00002150
C	00002160
IF (LLOC.EQ.0) = AK	00002170
WRITE(20) = W	00002180
RNA(LLOC) = BA	00002190
APNORM(LLOC) = AP	00002200
IF (LLOC.EQ.ITOTAL) GO TO 420	00002210
CALL XREAD(MSC,THEAD,"(1D18I6I)",STATUS)	AUG. 79
IF (STATUS.NE.-2) STOP "EOF MARKER NOT FOUND."	
GO TO 5	
420 APX = APNORM(1)	00002230
DO 430 I=2,ITOTAL	00002240
IF (APNORM(I).LT.APX) GO TO 430	00002250
APX = APNORM(I)	00002260
430 CONTINUE	00002270
C	00002280
C CALCULATE THE NORMALIZING FACTOR FOR EVERY SPECTRUM.	00002290
C	00002300
DO 440 MEM=1,ITOTAL	00002310
FACT(MEM) = APX / APNORM(MEM)	00002320
440 CONTINUE	00002330
C	00002340
REWIND 20	00002350
C	00002360
C NORMALIZE ALL SPECTRAL FILES FROM EYES TAP.	00002370
C	00002380

445 READ(20) L1,AMPLTD,NAVENO	0000239
XMIN = AMPLTD(1) * FACT(L1)	0000240
XMAX = XMIN	0000241
NWWW = NWORDS(L1)	0000242
DO 460 L2=1,NWWW	0000243
AMPLTD(L2) = AMPLTD(L2) * FACT(L1)	0000244
IF (AMPLTD(L2).GE.XMAX) XMAX = AMPLTD(L2)	0000245
IF (AMPLTD(L2).LE.XMIN) XMIN = AMPLTD(L2)	0000246
460 CONTINUE	0000247
DIFF(L1) = XMAX - XMIN	0000248
IF (L1.EQ.1) FMIN = XMIN	0000249
IF (L1.EQ.1) FMAX = XMAX	0000250
IF (XMAX.GE.FMAX) FMAX = XMAX	0000251
IF (XMIN.LE.FMIN) FMIN = XMIN	0000252
C	0000253
C.....WRITE NEW DATA TO SECONDARY DISK FILE 21.	0000254
C.....(NORMALIZED AMPLITUDE = AMPLTD).	0000255
C	0000256
WRITE(21) L1,AMPLTD,NAVENO	0000257
C	0000258
IF (L1.GE.ITOTAL) GO TO 450	0000259
GO TO 445	0000260
C	0000261
C    PRINT OUT NORMALIZATION INFORMATION.	0000262
C	0000263
450 LH = 1	0000264
L1=4	0000265
REWIND 20	0000266
REWIND 21	0000267
C	0000268
460 WRITE(N2,860)	0000269
WRITE(N2,880) (FILE(LW),LW=LH,L1)	0000270
WRITE(N2,882) (VNOR.(ATL(LW),LW=LH,L1)	0000271
WRITE(N2,883) (WW(LW),LW=LH,L1)	0000272
WRITE(N2,884) (BNA(LW),LW=LH,L1)	0000273
WRITE(N2,885) (APNORM(LW),LW=LH,L1)	0000274
WRITE(N2,886) (FACT(LW),LW=LH,L1)	0000275
IF (L1.EQ.4) GO TO 470	0000276
L1 = L1 + 1	0000277
LH = LH + 1	0000278
GO TO 460	0000279
470 NORH = 1	0000280
C	0000281
C    TESTING OF NORMALIZED DATA.	0000282
C	0000283
475 READ(21) M,AMPLTD,NAVENO	0000284
WRITE(N2,941) IFLE(M)	0000285
K=0	0000286
480 K=K+1	0000287
	0000288



C	NUMDEP = 0	0000349
	NUMIND = KOUNT	0000350
	DO 625 I3=1,10	0000351
	IF (NPHYS(I3).NE.0) NUMDEP = NUMDEP + 1	0000352
	NPHYS(I3) = NPHYS(I3)	0000353
625	CONTINUE	0000354
	NUMORS = ITOAL	0000355
	DO 630 LL=1,NUMORS	0000356
	MPEILE(LL) = MPRFCS(LL)	0000357
630	CONTINUE	0000358
		0000359
		0000360
C	.....CALL THE MULTIPLE LINEAR REGRESSION ROUTINE;	0000361
C	CALL REGRES	0000362
C		0000363
C	*****	0000364
C	*****	0000365
C	650 IF (IPLOT.EQ.1.OR.NMPLOT.EQ.1) CALL PLOT(0.,0.,999)	0000366
C		0000369
C	.....FORMAT STATEMENTS.	0000370
C		0000371
	715 FORMAT(A2,6X,3A2,4X,4A2)	0000372
	720 FORMAT(/5X,"NUMBER OF COMMON BASELINE AREAS SPECIFIED 0.")	0000373
	725 FORMAT(/5X,"NUMBER OF COMMON BASELINE AREAS SPECIFIED ",I1,"")	0000374
	730 FORMAT(10X,"AREA ",I1," IS BEST PICTURED IN FILE ",I2," :",/10X,"	0000375
	LEFT LIMIT ESTIMATE = ",F7.2,/10X,"RIGHT LIMIT ESTIMATE = ",F7.2)	
	735 FORMAT(/5X,"NO STATISTICAL CORRELATION WITH PHYSICAL PROPERTIES WILL	0000378
	BE PERFORMED. (NODISK=0).")	0000379
	740 FORMAT(/5X,"STATISTICAL CORRELATION WITH PHYSICAL PROPERTIES HAS	0000380
	BEEN REQUESTED //5X,"PHYSICAL PROPERTY INPUT DISK FT",I2,"FOOL."	0000381
	750 FORMAT(/5X,"PHYSICAL PROP. RECORD NUMBERS ",10(I3,5X)/36X,10(I3,5000382	
	760 FORMAT(/5X,"THE FOLLOWING PHYSICAL PROPERTIES WILL BE USED"/5X,"AS	0000383
	DEPENDENT VARIABLES ")	0000384
	765 FORMAT(10X,"PROP. NO. ",I1," - ",6A4)	0000385
	770 FORMAT(10X,"SUMMARY OF AMPLITUDE NORMALIZATION "//	0000386
	15X,"FILE",15X),"FILE"/55X,3("NO",18X),"NO")	0000387
	775 FORMAT(54X,3(13,17X),13/)	0000388
	780 FORMAT(/7X,"MAXIMUM AMPLITUDE NEAREST ",F6.1," WN. = ",4(F15.7,5000389	
	1X)/)	0000390
	785 FORMAT(7X,"AVE NUMBER AT MAX. AMPLITUDE = (WNMAX) = ",F10.4,3(10000391	
	1X,F10.4)/)	0000392
	790 FORMAT(7X,"BASELINE AMPLITUDE AT (WNMAX) = ",8X,4(E15.7,5X)/)	0000393
	795 FORMAT(7X,"PEAK HEIGHT AT (WNMAX) = ",15X,4(E15.7,5X)/)	0000394
	800 FORMAT(1H1,/(1X8I10))	0000395
	900 FORMAT(1X,3(5X,F10.4))	0000396
	930 FORMAT(6X,13,4X,F10.4,2X,E15.7,2(7X,13,4X,F10.4,2X,E15.7))	0000397
		0000398
		0000399

```

940 FORMAT(1H1,4X,"FTIS FILE NUMBER  ",I3,20X,"NON-NORMALIZED (PURE) S00004
SPECTRAL DATA"/5X,3("DATA",8X,"WAVE",8X,"AMPLITUDE",8X)/ 00004
25X,3("WORD",7X,"NUMBER",24X)/) 00004
941 FORMAT(1H1,4X,"FTIS FILE NUMBER  ",I3,20X,"NORMALIZED SPECTRAL DAT00004
1A"/5X,3("DATA",8X,"WAVE",8X,"AMPLITUDE",8X)/ 00004
25X,3("WORD",7X,"NUMBER",24X)/) 00004
950 FORMAT(7X,"NORMALIZING FACTOR FOR THIS FILE = ",5X,4(E15.7,5X)) 00004
970 FORMAT(1H1,4X,"FOURIER TRANSFORM INFRARED SPECTROSCOPY COMPUTER PRO0704
LOGRAM"/5X,"E490 - THICKOL CORP. / HUNTSVILLE , ALABAMA 35807"/5X00004
2,"INPUT DESCRIPTION "/) 00004
972 FORMAT(/5X,"NAME OF FTIS INPUT TAPE  ",A4) 00004
975 FORMAT(/5X,"TOTAL NUMBER OF FILES (THIS RUN)  ",I3,"") 00004
980 FORMAT(/5X,"LIST OF FILES TO BE ANALYZED  ",10(I3,5X)/36X,10(I3,500004
1X)/36X,10(I3,5X)) 00004
981 FORMAT(/5X,"NORMALIZING WAVE NUMBER  ",F6.1,"") 00004
982 FORMAT(/5X,"NON-NORMALIZED SPECTRAL PLOT HAS BEEN REQUESTED.") 00004
983 FORMAT(/5X,"NORMALIZED SPECTRAL PLOT HAS BEEN REQUESTED.") 00004
984 FORMAT(/5X,"NO SPECTRAL PLOTS HAVE BEEN REQUESTED.") 00004
995 FORMAT(/5X,"PEAK HEIGHT VALIDITY =",F4.1," PER CENT.") 00000
STOP 0 00004
END 00004

```

```

SUBROUTINE APCAL(AMPLTD, WAVEND, WAVNOR, N, JAX, APMAX, WE, DZ, APZ)
C
C-----SEPTEMBER 4, 1972 - CDC REVISION.
C
C-----LOGIC WILL CALCULATE THE NORMALIZATION AMPLITUDE VALUE (APNORM)
C-----FOR THE PEAK NEAREST WAVE NUMBER (WAVNOR) CM-1.
C
COMMON /SPECTR/ IFILE(30), NPT29(30)
COMMON /TWINPK/ ALIMIT(6), ZLIMIT(6), MFILE(6), NORMCO, NUMBER
DIMENSION AMPLTD(468), WAVEND(468), A(25)
I=0
J=0
DO 1 J=1,25
  A(J) = 0.0E0
C
C-----SEARCH FOR MAXIMUM AMPLITUDE NEAREST TO WAVNOR.
C
5 I=1
  IF (WAVEND(I).GT.(WAVNOR+20.0)) GO TO 5
  IF (WAVEND(I).LT.(WAVNOR-20.0)) GO TO 10
  IF (J.EQ.0) II=I-1
  J=J+1
  A(J) = AMPLTD(I)
  GO TO 5
C
10 APMAX = A(I)
  DO 15 JJ=2,J
    IF (A(JJ).LT.APMAX) GO TO 15
    APMAX = A(JJ)
  IREF = J
C
15 CONTINUE
  MK = II+IREF
  NPT29(IREF) = MK
C
C-----CHECK TO SEE IF A "COMMON BASELINE" IS TO BE USED
C-----FOR THE NORMALIZATION AMPLITUDE.
C
  IF (NORMCO.EQ.1) GO TO 45
C
  NOW LOCATE "MINIMUM" POINTS ON BOTH SIDES OF APMAX PEAK.
C
  K=MK-1
  20 IF (AMPLTD(K).GT.AMPLTD(K+1)) GO TO 30
  K=K-1
  GO TO 20
  30 KMIN = K+1
  KMIN = K+1
C
  K = MK
  40 IF (AMPLTD(K+1).GT.AMPLTD(K)) GO TO 40

```

```

      K = K+1
      GO TO 35
40  RMAX = AMPLTD(K)
      KMAX = K
      GO TO 100
C
C.....I = THE I-TH DATA WORD IN THE SPECTRUM FILE.
C
      45  I=0
      200  I=I+1
           IF (WAVENO(I).GT.ALIMIT(JAX))GO TO 200
           IS = I
      210  IF (AMPLTD(I-1).GT.AMPLTD(I)) GO TO 220
           BMIN = AMPLTD(I-1)
           KMIN = I-1
           I=I-1
           GO TO 210
      220  I=IS
      230  IF (WAVENO(I+1).LE.ZLIMIT(JAX))GO TO 240
           I=I+1
           GO TO 230
      240  I=I+1
      245  IF (AMPLTD(I+1).GT.AMPLTD(I)) GO TO 250
           I=I+1
           GO TO 245
      250  RMAX = AMPLTD(I)
           KMAX = I
C
      100  CHGX = WAVENO(KMAX) - WAVENO(KMIN)
           CHGY = BMAX - BMIN
C
C.....CHECK FOR CHGY = 0.0.
C
           IF (CHGY.EQ.0.0E0) GO TO 110
C
C.....'MULTIPLE' TANGENTIAL CHECK.
C
           IF (BMIN.GT.BMAX) GO TO 80
C
C.....RIGHT = SHORT SIDE OF THE PEAK(S).
C
      75  SLOPE1 = (CHGY/CHGX)
           CY = RMAX - AMPLTD(KMIN+1)
           CX = WAVENO(KMAX) - WAVENO(KMIN+1)
           SLOPE2 = (CY/CX)
           SIGN = SLOPE1 * SLOPE2
           ZZ1 = ABS(SLOPE1)
           ZZ2 = ABS(SLOPE2)
           IF (SIGN.LT.0.0E0.OR.ZZ2.LT.ZZ1) GO TO 110
           BMIN = AMPLTD(KMIN+1)

```

```

      CHGX = CX
      CHGY = CY
      KMIN = KMIN + 1
      GO TO 75
C.....LEFT = SHORT SIDE OF THE PEAK(S).
C
      80 SLOPE1 = (CHGY/CHGX)
      CY = AMPLTD(KMAX-1) - BMIN
      CX = WAVENO(KMAX-1) - WAVENO(KMIN)
      SLOPE2 = (CY/CX)
      SIGN = SLOPE1 * SLOPE2
      ZZA = ABS(SLOPE1)
      ZZB = ABS(SLOPE2)
      IF (SIGN.LT.0.OEQ.OR.ZZB.LT.ZZA) GO TO 110
      BMAX = AMPLTD(KMAX-1)
      CHGX = CX
      CHGY = CY
      KMAX = KMAX - 1
      GO TO 80
C
      110 D1 = (CHGY*(WAVENO(MK)-WAVENO(KMIN)))/CHGX
      D2 = D1 + BMIN
      AP2 = APMAX - D2
      ME = WAVENO(MK)
C.....SEND THE VALUES FOR (APMAX,ME,D2,AP2) BACK TO MAIN.
C
      RETURN
      END

```



```

SUBROUTINE BASELN(DIFF,VALUE)                                00000006
C   ...SEPTEMBER 13,1979 - REVISION FOR CDC 6600.          00000007
C   ...                                     00000008
C   ...                                     00000009
COMMON /POINTS/ NWORDS(30)                                00000010
DIMENSION AMPLTD(468),WAVENO(468),IDWORD(468)              00000011
DIMENSION CMIN1(50),CMIN2(50),PEAK(50),AB(50),HT(50),DIFF(30) 00000012
C   *****                                00000013
C   ...THIS SUBROUTINE SEARCHES FOR AND RECORDS VALID PEAKS 00000014
C   ...WITHIN EACH INDIVIDUAL SPECTRUM.                    00000015
C   ...                                     00000016
C   ...                                     00000017
C   ... *****                                00000018
C   ... ALLOW FOR MAX. 50 PEAKS PER INDIVIDUAL SPECTRUM.   00000019
C   ...                                     00000020
C   ... I = THE I-TH DATA WORD IN THE SPECTRUM FILE.     00000021
C   ...                                     00000022
C   ... IP = THE IP-TH VALID PEAK IN THE SPECTRUM (LEFT-RIGHT). 00000023
C   ...                                     00000024
C   ... IDWORD(IPK) = 0 NO PEAK AT DATA WORD "IPK".        00000025
C   ...                                     00000026
C   ...                                     00000027
C   ...                                     00000028
C   ...                                     00000029
C   ...                                     00000030
C   ...                                     00000031
C   ...                                     00000032
C   ...                                     00000033
C   ... *****                                00000034
C   ... INITIALIZE SOME VARIABLES.                          00000035
C   ...                                     00000036
C   ...                                     00000037
C   ...                                     00000038
C   ...                                     00000039
C   ...                                     00000040
C   ... IDWORD(1) = 0                                       00000041
C   ... CONTINUE                                           00000042
C   ... IDWORD(1) = 0                                       00000043
C   ... AB(1) = 0.0                                         00000044
C   ... HT(1) = 0.0                                         00000045
C   ... CONTINUE                                           00000046
C   ... I = 0                                               00000047
C   ...                                     00000048
C   ...                                     00000049
C   ... *****                                00000050
C   ...                                     00000051
C   ... READ(21) IF,AMPLTD,WAVENO                          00000052
C   ... NW = NWORDS(IF)                                     00000053
C   ...                                     00000054
C   ... I = I+1                                             00000055

```

IF (WAVENO(I).GT.3200.E0) GO TO 10	00000560
15 IF (WAVENO(I).LT.2700.E0) GO TO 90	00000570
IF (J.NE.0) GO TO 20	00000580
CMIN1(IP) = AMPLTD(I)	00000590
N1=I	00000600
20 IF (WAVENO(I+1).LT.700.E0.OR.(I+1).GE.NW) GO TO 130	00000610
IF (AMPLTD(I+1).GT.AMPLTD(I)) GO TO 30	00000620
IF (J.NE.0) GO TO 40	00000630
CMIN1(IP) = AMPLTD(I+1)	00000640
N1=I+1	00000650
I=I+1	00000660
GO TO 20	00000670
30 J=1	00000680
I=I+1	00000690
GO TO 20	00000700
40 IOWEN=I	00000710
WRITE(6,400) IOWEN,AMPLTD(I),CMIN1(IP)	00000720
400 FORMAT(3X,I5,3X,F15.7,3X,E15.7)	00000730
A = AMPLTD(I) - CMIN1(IP)	00000740
PEAK(IP) = AMPLTD(I)	00000750
IPK = I	00000760
50 IF (AMPLTD(I+1).GT.AMPLTD(I)) GO TO 60	00000770
I=I+1	00000780
GO TO 50	00000790
60 CMIN2(IP)= AMPLTD(I)	00000800
N2=I	00000810
IMIN2 = N2	00000820
AMP2 = CMIN2(IP)	00000830
B = PEAK(IP) - AMP2	00000840
C = AMAX1(A,B)	00000850
DIFF = DIFF(I-1) * VALUE	00000860
IF (C.GT.022) IWORD(IPK) = 2	00000870
***VALID PEAK BELOW ***	00000880
CALCULATE EQUATION FOR LINE CONNECTING TWO MINIMUM POINTS.	00000890
ADJUST BASELINE FOR THE TANGENTIAL CONDITION;	00000900
THUS BASELINE MINIMUM POINTS WILL LIE ON THE TANGENT.	00000910
CHGX = WAVENO(N2) - WAVENO(N1)	00000920
CHGY = CMIN2(IP) - CMIN1(IP)	00000930
CHECK FOR A = B = C.	00000940
IF (CHGY.EQ.0.0E0) GO TO 110	00000950
IF (C.EQ.B) GO TO 80	00000960
B = SHORT SIDE OF THE IP-TH PEAK.	00000970
75 SLOPE1 = (CHGY/CHGX)	00000980
	00000990
	00001000
	00001010
	00001020
	00001030
	00001040
	00001050

```

CY = CMIN2(IP) - AMPLTD(N1+1)
CX = WAVENO(N2) - WAVENO(N1+1)
SLOPE2 = (CY/CX)
SIGN = SLOPE1 * SLOPE2
ZZ1 = ABS(SLOPE1)
ZZ2 = ABS(SLOPE2)
IF (SIGN.LT.0.0EQ.OR.ZZ2.LT.ZZ1) GO TO 110
CMIN1(IP) = AMPLTD(N1+1)
CHGX = CX
CHGY = CY
N1 = N1 + 1
GO TO 75

```

C  
C  
C

A = SHORT SIDE OF THE IP-TH PEAK.

```

80 SLOPE1 = (CHGY/CHGX)
CY = AMPLTD(N2-1) - CMIN1(IP)
CX = WAVENO(N2-1) - WAVENO(N1)
SLOPE2 = (CY/CX)
SIGN = SLOPE1 * SLOPE2
ZZ1 = ABS(SLOPE1)
ZZ2 = ABS(SLOPE2)
IF (SIGN.LT.0.0EQ.OR.ZZ2.LT.ZZ1) GO TO 110
CMIN2(IP) = AMPLTD(N2-1)
CHGX = CX
CHGY = CY
N1 = N2 - 1
GO TO 80

```

C

90 (CHGY\*(WAVENO(IPK1)-WAVENO(N1)))/CHGX

C

AB(IP) = AMPLITUDE OF THE BASELINE FOR THE IP TH PEAK.

HT(IP) = HEIGHT OF THE IP-TH PEAK.

C

```

C(IP) = PD+CMIN1(IP)
H(IP) = PEAK(IP) - AB(IP)
IF (IDWORD(IPK).EQ.2) GOTO 120
IF WORD(IPK) = 1
IP = IP+1

```

C

(AMP2,IMIN2) THE TRUE MINIMUM POINT ON THE RIGHT SIDE.

C

```

CMIN1(IP) = AMP2
N1=IMIN2
I=I+1
GO TO 15

```

C

```

90 IF (WAVENO(I).LT.1800.E0) GO TO 100
I=I+1
GO TO 90

```

00001060  
00001070  
00001080  
00001090  
00001100  
00001110  
00001120  
00001130  
00001140  
00001150  
00001160  
00001170  
00001180  
00001190  
00001200  
00001210  
00001220  
00001230  
00001240  
00001250  
00001260  
00001270  
00001280  
00001290  
00001300  
00001310  
00001320  
00001330  
00001340  
00001350  
00001360  
00001370  
00001380  
00001390  
00001400  
00001410  
00001420  
00001430  
00001440  
00001450  
00001460  
00001470  
00001480  
00001490  
00001500  
00001510  
00001520  
00001530  
00001540  
00001550

100 IF (WAVENO(I).LT.700.EQ.OR.I.GE.NW) GO TO 130	00001560
IF (K.NE.0) GO TO 20	00001570
J=0	00001580
CMIN1(IP) = AMPLD(I)	00001590
N1=I	00001600
K=1	00001610
GO TO 20	00001620
C	00001630
C*****	00001640
C	00001650
C.....WRITE IDWORD INFORMATION TO THE DISK FILE 22.	00001660
C	00001670
C*****	00001680
C	00001690
130 WRITE(6,5000) IF	00001700
WRITE(6,5001) IDWORD	00001710
5000 FORMAT (5X,I5)	00001720
5001 FORMAT (5(75X,100I1))	00001730
WRITE (22) IF,IDWORD	00001740
WRITE (24) IF,AB,H,	00001750
RETURN	00001760
END	00001770

SUBROUTINE COBASE(DIFF,VALUE,JX)	000001
C	000001
C.....SEPTEMBER 28, 1979. (CHANGED AT AF/RPL : OCT. 24, 1979).	
C.....	000001
C.....CDC REVISION.	000001
C	000001
COMMON /SPECTR/ IFILE(30),NPT29(30)	000001
COMMON /TWINPK/ ALIMIT(6),ZLIMIT(6),MFILE(6),NORMCD,NUMBER	000001
DIMENSION AMPLTD(468),HAVENO(468),IDWORD(468),AB(50),HT(50)	000001
DIMENSION EDGE1(25),EDGE2(25),ISPEC(29),JWORD(25),MASTER(25)	000001
DIMENSION A(25),Y(50),Z(50),DIFF(30)	000001
-----	000001
C	000001
C.....COMMON BASELINE ROUTINE / MARCH 19, 1979. (OCS)	000001
C	000002
C.....THIS SUBROUTINE IS A USER REQUESTED OPTION.	000002
C.....THIS LOGIC IS UTILIZED TO ESTABLISH A COMMON BASELINE	000002
C.....FOR ADJACENT SPECTRAL PEAKS.	000002
C	000002
-----	000002
C	000002
C.....I - THE I-TH DATA WORD IN THE SPECTRUM FILE.	000002
C	000002
C.....IP = THE IP-TH PEAK IN THE COMMON BASELINE (LEFT-TO-RIGHT).	000002
C	000003
C.....MAXIMUM NUMBER OF PEAKS PER BASELINE = 25.	000003
C	000003
C.....MAXIMUM NUMBER OF COMMON BASELINES PER COMPUTER RUN = 6.	000003
C	000003
-----	000003
C	000003
DO 2 J2=1,29	000003
ISPEC(J2) = 0	000003
2 CONTINUE	000003
C	000004
C.....BEGIN BY SEARCHING FOR THE FILE (MFILE(JX)) WITH THE OPTIMUM	000004
C.....PICTURE OF THE MULTIPLE PEAKS.	000004
C	000004
JB = 1	000004
DO 5 I = 1,NUMBER	000004
IF (IFILE(I).EQ.MFILE(JX)) GO TO 4	000004
ISPEC(IB) = I	000004
IB = IB + 1	000004
GO TO 5	000004
4 IF = I	000005
5 CONTINUE	000005
MAXX=NUMBER - 1	000005
LOOP = 1	000005
C	000005
C READ FROM THE NORMALIZED FILE 21.	000005

C		00000
	6 READ (21) IX,AMPLTD,WAVENO	00000
	READ (22) IX,IDWORD	00000
	READ (24) IX,AB,HT	00000
	IF (IX.EQ.IF) GOTO 7	00000
C		00000
C	WRITE TO UPDATED IDWORD FILE 23.	00000
C		00000
	WRITE (23) IX,IDWORD	00000
	WRITE (25) IX,AB,HT	00000
	GO TO 6	00000
C		00000
	7 I=0	00000
	IP=1	00000
	IVALID = 0	00000
	10 I=I+1	00000
	IF (WAVENO(I).GT.ALIMIT(JX)) GO TO 10	00000
	IS = I	00000
	20 IF (AMPLTD(I-1).GT.AMPLTD(I)) GO TO 30	00000
	EDGE1(IP) = AMPLTD(I-1)	00000
	NI=I-1	00000
	I=I-1	00000
	GO TO 20	00000
	30 I=IS	00000
	45 IF (AMPLTD(I+1).LT.AMPLTD(I)) GO TO 40	00000
	I=I+1	00000
	GO TO 35	00000
	40 IDWORD(IP) = I	00000
	ALFFT = AMPLTD(I) - EDGE1(IP)	00000
	55 IF (AMPLTD(I+1).GT.AMPLTD(I)) GO TO 60	00000
	I=I+1	00000
	GO TO 55	00000
C		00000
C	TAKE INTO ACCOUNT THE POSSIBILITY OF A DISAPPEARING 2-ND PEAK.	00000
C	JULY , 1979.	00000
C		00000
	60 IF (WAVENO(I).LE.ZLIMIT(JX)) GO TO 600	00000
	EDGE2(IP)= AMPLTD(I)	00000
C		00000
C	FINISHED DEFINING THE FIRST PEAK.	00000
C		00000
	62 IP=IP+1	00000
	EDGE1(IP)=EDGE2(IP-1)	00000
C		00000
	65 IF (AMPLTD(I+1).LT.AMPLTD(I)) GO TO 70	00001
	I=I+1	00001
	GO TO 65	00001
	70 IDWORD(IP) = I	00001
	71 IF (AMPLTD(I+1).GT.AMPLTD(I)) GO TO 75	00001
	I=I+1	00001



CHGY = CY	000019
N1 = N1 + 1	000019
GO TO 110	000019
C	000019
C.....ALEFT = SHORT SIDE OF THE PEAK(S).	000019
C	000019
170 SLOPE1 = (CHGY/CHGX)	000019
CY = AMPLD(N2-1) - EDGE1(1)	000019
CX = WAVEND(N2-1) - WAVEND(N1)	000019
SLOPE2 = (CY/CX)	000019
SIGN = SLOPE1 * SLOPE2	000019
ZZ1 = ABS(SLOPE1)	000019
ZZ2 = ABS(SLOPE2)	000019
IF (SIGN.LT.0.OR.0EQ.0R.ZZ2.LT.ZZ1) GO TO 110	000019
EDGE2(IP) = AMPLD(N2-1)	000019
CHGX = CX	000019
CHGY = CY	000019
N2 = N2 - 1	000019
GO TO 170	000019
C	000019
C.....ZERO-OUT PREVIOUS PEAK HEIGHT INFORMATION THAT HAS	000019
C.....CALCULATED IN THE "BASELN" SUBROUTINE.	000019
C	000019
C.....ISUM = NO. OF PEAKS FOUND IN BASELINE INTERVAL BETWEEN	000019
C.....N1 + N2 BY BASELN SUBROUTINE.	000019
C	000019
110 ISUM = 0	000019
IF (LOOP.NE.1) GO TO 111	000019
MASN1=M1	000019
MASN2=N2	000019
111 DO 105 LAX=MASN1,MASN2	000019
IF (IDWORD(LAX).NE.0) ISUM = ISUM + 1	000019
IDWORD(LAX) = 0	000019
105 CONTINUE	000019
IF (LOOP.NE.1) GO TO 115	000019
DO 116 N2=1,IP	000019
MASTER(N2) = JWORD(N2)	000019
CONTINUE	000019
NPEAK = IP	000019
116 ISHIFT=NPEAK-ISUM	000019
C	000019
DO 150 J=1,NPEAK	000019
IAJ = MASTER(J)	000019
DEN = (CHGY*(WAVEND(IAJ)-WAVEND(N1)))/CHGX	000019
C	000020
C.....AB(IP) = AMPLITUDE OF THE BASELINE FOR THE IP-TH PEAK.	000020
C	000020
ALJ = DEN + EDGE1(1)	000020



IDWORD(IAD) = 3	000020
IF (IVALID.EQ.1) IDWORD(IAD) = 4	000020
150 CONTINUE	000020
C	000020
C.....MASTER(X) IS THE ARRAY FOR THE REFERENCE DATA WORDS.	000020
C	000020
C.....CHANGE THE FILE NUMBER(S); THEN RETURN TO CALCULATE THE	000021
C....."COMMON BASELINE" PEAK HEIGHTS FOR THE OTHER FILE(S).	000021
C	000021
200 WRITE(6,666) (MASTER(L),L=1,NPEAK)	000021
WRITE(6,700) (A(JJ),JJ=1,NPEAK)	000021
700 FORMAT (10(5X,F6.2))	000021
666 FORMAT (1H0,5X,5I3)	000021
C	000021
K=0	000021
K1 = 0	000022
IC = 0	000022
DO 220 IR = 1,468	000022
IF (IDWORD(IR).EQ.0) GOTO 220	000022
K = K + 1	000022
I = K	000022
IF (IR.GE.MASN1.AND.IR.LE.MASN2) GO TO 225	
IF (IC.EQ.1) L = K-ISHIFT	
Y(K) = AB(L)	000022
Z(K) = HI(L)	000022
GO TO 220	000023
225 K1 = K1 + 1	000023
Y(K) = A(K1)	000023
Z(K)=AMPLTD(IR) - Y(K)	000023
IC = 1	000023
220 CONTINUE	000023
C*****	000023
WRITE (25) IF,Y,Z	000023
WRITE (23) IF,IDWORD	000023
IF (IF.GE.NUMBER) GOTO 230	000023
IH = IF + 1	000024
DO 300 L1 = IH,NUMBER	000024
READ (22) IY,IDWORD	000024
WRITE (23) IY,IDWORD	000024
READ (24) IZ,AB,HT	000024
WRITE (25) IZ,AB,HT	000024
300 CONTINUE	000024
230 CONTINUE	000024
REWIND 21	000024
REWIND 22	000024
REWIND 23	000025
REWIND 24	000025
REWIND 25	000025
C	000025
DO 310 I1=1,NUMBER	000025

HEAD (23) IY, IDWORD	00002
WRITE (22) IY, IDWORD	00002
READ (25) III, AB, HT	00002
WRITE (24) III, AB, HT	00002
310 CONTINUE	00002
C	00002
REWIND 22	00002
REWIND 23	00002
REWIND 24	00002
REWIND 25	00002
IF (LOOP.GT.MAXX) GO TO 400	00002
IF = ISPEC (LOOP)	00002
LOOP = LOOP + 1	00002
GO TO 6	00002
C.....END OF THE CYCLE.	00002
400 CONTINUE	00002
RETURN	00002
END	00002



```

REWIND 21
REWIND 22
REWIND 24

C
C
C***** PRINT-OUT PEAK HEIGHT TABLE.
C
C
5  ILOUN. = 0
   ILCOUNT = 0
   DO 500 I8=1,30
   IPRINT(I8) = 0
500 CONTINUE

   DO 10 I1=1,MXDNH
   DO 35 M=1,30
   LEA(M) = 0
   T(M) = 0.0
35 CONTINUE
   N = 0
   ISTOP = 0
45 READ (22) JJ, IDWORD
   READ (24) JJ, AR, HT
   IF (IDWORD(I1).NE.0) GO TO 50
   IF (JJ.GE.ITOTAL) GO TO 60
   GO TO 45
50 N = 1

   IPRINT(JJ) = IPRINT(JJ) + 1
   JPL = IPRINT(JJ)
   T(JJ) = HT(JPL)
   IF (IDWORD(I1).EQ.2.OR.IDWORD(I1).EQ.4) LEA(JJ) = 1

   IF (JJ.GE.ITOTAL) GO TO 60
   GO TO 45
60 IF (N.EQ.0) GO TO 65

C.....WRITE PEAK HEIGHTS TO FILE 26 , FOR LATER PRINTING.

   WRITE(26) I1,T
   ILCOUNT = ILCOUNT + 1

C.....CHECK IF REGRES WILL BE CALLED.

   IF (IREG.EQ.0) GOTO 65
C.....TO PREPARE FOR CORRELATION OF PEAK HEIGHTS AND MECHANICAL
C.....PROPERTIES ; SEARCH FOR ONLY THOSE PEAKS THAT ARE VALID IN
C.....ALL THE FILES.

```

```

C.....DO NOT INCLUDE THE WAVNOR PEAK IN THE REGRESSION ANALYSIS. 000036
C.....(THE REASON FOR THE NPT29 CHECK BELOW SD(II)-SD(KK) = 0.0). 000034
C
      IF (II.EQ.NPT29(1)) GO TO 65
      DO 210 KAY = 1,ITOTAL
      IF (LFA(KAY).EQ.1) GOTO 210
      ISTOP = 1
210 CONTINUE
      IF (ISTOP.EQ.1) GOTO 65
C      FILE 27 ++++++ VALID PEAKS IN ALL THE FILES.
      WRITE (27) II,T
      IKOUNT = IKOUNT + 1
      65 CONTINUE
      REWIND 22
      REWIND 24
      10 CONTINUE
C****
      REWIND 26
      IF (IREG.EQ.10) REWIND 27
C
      N1 = 1
      N2 = 10
C
      250 WRITE(6,410) (IFILE(I),I=N1,N2)
      DO 240 N3 = 1,ICOUNT
      READ (26) N4,I
      WRITE (6,400) N4,(T(K),K=N1,N2)
      IF (N3.EQ.25.OR.N3.EQ.50) WRITE(6,410) (IFILE(I),I=N1,N2)
      240 CONTINUE
C
      REWIND 26
      IF (N2.GE.ITOTAL) GO TO 70
      IF (N2.EQ.20) GO TO 20
      N1 = 11
      N2 = 20
      GO TO 250
      20 N1 = 21
      N2 = 30
      GO TO 250
C
C*****
C
      70 WRITE(6,420) ITOTAL
C
C*****
C
C***** FORMATS.
C
C*****

```

```

000014
200 FORMAT(1H1,4X,"FOURIER TRANSFORM INFRARED SPECTROSCOPY -- NORMAL 000014
11ZED PEAK HEIGHT INFORMATION FILE NUMBER",I3,7/5X 000014
2,"*** - DENOTES A VALID PEAK."/5X,"CB. - DENOTES A PEAK MEASURED FROM A 000014
COMMON BASELINE ."/5X 000014
4,"PEAK NUMBER DATA WORD WAVE NUMBER PEAK HEIGHT BA000014
55FLINE AMPLITUDE"/5X,"-----" 000014
6 -----") 000015
300 FORMAT(1H0,8X,I2,I3X,I3,8X,F10.4,5X,F12.7,5X,F12.7) 000015
310 FORMAT(1H0,4X,"*** ",I2,I3X,I3,8X,F10.4,5X,F12.7,5X,F12.7) 000015
320 FORMAT(1H0,8X,I2," CR.",9X,I3,8X,F10.4,5X,F12.7,5X,F12.7) 000015
330 FORMAT(1H0,4X,"*** ",I2," CB.",9X,I3,8X,F10.4,5X,F12.7,5X,F12.7) 000015
400 FORMAT(1/8X,I3,7X,10(F7.3,4X)) 000015
410 FORMAT(1H1,7X,"FTIS ***** PEAK HEIGHT TABULATION FOR ALL THE FILES *****" 000015
1ATION FOR ALL THE FILES *****" 000015
27/8X,"DATA",8X, 000015
310("FILE",7X),/8X,"WORD",9X,10(I2,9X)/8X,"-----" 000015
4-----" 000015
5-----"/) 000015
420 FORMAT(1H0,4X,"NORMAL END OF PEAK HEIGHT TABLE "/5X,"TOTAL NUMBER 000015
1 OF FILES LISTED = ",I3,".") 000015
RETURN 000015
END 000015

```



WN(NPTS+1) = 4200.0F0	00006490
WN(NPTS+2) = -200.0F0	00006500
IF (NORM.F0.0) GO TO 40	00006510
CALL AXIS (0.,0., "NORMALIZED AMPLITUDE", 20, 10.0, 90., AP(NPTS+1), AP(NPTS+2))	00006520
GO TO 50	00006530
40 CALL AXIS (0.,0., "AMPLITUDE", 9, 10.0, 90., AP(NPTS+1), AP(NPTS+2))	00006540
50 CALL AXIS (0.,0., "WAVE NUMBER", -11, 20.0, 0., WN(NPTS+1), WN(NPTS+2))	00006550
CALL LINE (WN(1), AP(1), NPTS, 1, 0, 2)	00006560
CALL PLOT (23.0, 0.0, -3)	00006570
RETURN	00006580
END	00006590



```

SUBROUTINE REGRES                                00000060
C                                                    00000070
C          MULTIPLE LINEAR REGRESSION              00000080
C          BY THE STEPWISE METHOD                   00000090
C          DEVELOPED FROM THIOKOL PROGRAM NUMBER - E023. 00000100
C                                                    00000110
C*****00000120
C                                                    00000130
C                                                    00000140
C          N    - TOTAL NUMBER OF VARIABLES ON INPUT (MAX. = 60). 00000150
C          I    - TOTAL NUMBER OF DEPENDENT VARIABLES 00000160
C          M    - TOTAL NUMBER OF OBSERVATIONS        00000170
C          F1   - THE F SIGNIFICANCE LEVEL TO ENTER A VARIABLE INTO THE 00000180
C                  REGRESSION                                00000190
C          F2   - THE F LEVEL TO REMOVE A VARIABLE FROM THE REGRESSION 00000200
C          ITP  - NUMBER OF TRANSFORMED VARIABLES      00000210
C          J    - VARIABLE NUMBER BEING TRANSFORMED   00000220
C          NTRAN - TRANSFORMATION TYPE CODE            00000230
C                  = 1, LOG OF X - LOG10(X(I))        00000240
C                  = 2, POWER OF X - (X(I)+A)**P      00000250
C                  = 3, SQUARE ROOT OF X - Sqrt(X(I)) 00000260
C                  = 4, NATURAL LOG OF X - LOG(X(I))  00000270
C                  = 5, X(I-1) * (X(J)**POWER(I))     00000280
C          CONS - TRANSFORMATION CONSTANT IF ANY       00000290
C          POWER - POWER TO WHICH TRANSFORMED VARIABLE MAY BE TAKEN, 00000300
C                  1.0 WHEN NOT INPUT                 00000310
C          WT   - WEIGHT APPLIED TO THE SPECIFIC OBSERVATION, 00000320
C          X     - DATA FOR VARIABLES 1 UP THRU N.    00000330
C                                                    00000340
C          INCORPORATED AT THIOKOL MAY 26, 1978.      00000350
C                                                    00000360
C          SUBROUTINE WRITTEN ON MAY 22, 1979.        00000370
C                                                    00000380
C          BY OWEN C. SMITH.                          00000390
C                                                    00000400
C*****00000410
C          COMMON /CORREL/ PRREF(60),MPFILE(30),      00000420
C                  NPHYSI(10),NUMDEP,NUMIND,NUMOBS,IRES 00000430
C          COMMON /SELECTR/ IFILE(30),NPT29(30)       00000440
C          DIMENSION PHYS(100,10),XHT(30),IVAR(50)    00000450
C          DIMENSION FDIST(30),SUMX(60),SUMXX(60,60),X(60),B(60),SB(60), 00000460
C          CORNS(60),SD(60),R(60,60),XM(60),POWER(60),B1(60),YY(60),CORRCF(10) 00000470
C          NTRAN(60),IFLAG(60),ISUB(60)              00000480
C          C.....F - DISTRIBUTION PERCENTAGE POINTS (ALPHA = .10). 00000490
C          C.....FROM THE TEXT , PAGE 482.          00000500
C          DATA FDIST / 39.86,8.53,5.54,4.54,4.06,3.78,3.59,3.46,3.36,3.28, 00000510
C          A          3.23,3.13,3.14,3.10,3.07,3.05,3.03,3.01,2.99,2.97, 00000520
C          00000530
C          00000540
C          00000550

```

```

      2.96,2.95,2.94,2.93,2.92,2.91,2.90,2.89,2.89,2.88 / 00000560
C
C *****
C
C      READ MECHANICAL PROPERTIES TAPE ; GENERATED BY E410 PROGRAM.
C
C *****
C
C
C
C      READ(IREG) TITLE
C      4 READ (IREG) N1,N2,N3,N4,(PHYS(N1,N5),N5=1,10)
C      IF (FOF(IREG)) 5,4
C      5 N = NUMIND + NUMDEP
C      NN = N
C      L = NUMDEP
C      M = NUMOBS
C      ICURVE = LINEAR OR NON-LINEAR REGRESSION MODE.
C      ICURVE = 0
C      JCURVE = 0
C      REWIND IREG
C      10 ITR = 0
C      LL = L + 1
C      IF (LL.GT.29) STOP
C      F1 = FDISP(LL)
C      F2 = FDISP(LL+1)
C      IF (ICURVE.EQ.0) GO TO 12
C      NV = NUM1 + JCURVE
C      WRITE (6,705) NV
C      GO TO 11
C      12 WRITE (6,710)
C      11 WRITE (6,720) N,L,M,F1,F2
C
C      INITIALIZE SUMS.
C
C      N = 1 + 1
C      N2=N
C      IOBS=0
C      DO 20 I=1,N
C      SUMX(I)=0.0
C      STRAN(I)=0
C      FLAG(I)=0
C      ISUB(I)=0
C      DO 20 J=1,N
C      20 SUMXX(I,J)=0.0
C
C      IOBS=M
C      30 IF (ITR) 50,50,40
C
C      CALCULATE SUMS, SUMS OF SQUARES AND CROSS PRODUCTS.
C      READ OBSERVATIONS FOR EACH VARIABLE X(I).

```

C	WRITE (6,750)	00000980
C		00000990
C	70 N=N2	00001000
C		00001010
C	C.....SETUP INDEPENDENT VARIABLES.	00001020
C		00001030
C	IF (ICURVE.EQ.1) GO TO 76	00001040
C		00001050
C	DO 3000 MCN = 1,NUMIND	00001060
	READ (27) IAX,XHT	00001070
	X(MCN) = XHT(I0BS+1)	00001080
	IVAR(MCN) = IAX	00001090
3000	CONTINUE	00001100
	REWIND 27	00001110
C		00001120
C	FIND WHICH RECORD IS NEEDED FOR PHYSICAL PROPERTY INPUT.	00001130
C		00001140
C	FIND WHICH RECORD IS NEEDED FOR PHYSICAL PROPERTY INPUT.	00001150
C		00001160
C	LIB = MPFILE(I0BS+1)	00001170
C		00001180
C	C.....SETUP DEPENDENT VARIABLES.	00001190
C		00001200
	DO 75 MGM=1,NUMDEP	00001210
	KIM = NPHYS1(MGM)	00001220
	X(NUMIND+MGM) = PHYS(LIB,KIM)	00001230
75	CONTINUE	00001240
	GO TO 78	00001250
75	READ (14) (X(JB),JB=1,N)	00001260
78	WT = 0.0	00001270
C		00001280
C	C.....TRANSFORMATION POSSIBILITIES.	00001290
C		00001300
	DO 130 J=1,N	00001310
	IF (IFIA(I).NE.0) GO TO 178	00001320
	IF (NTRAN(I)) 180,180,120	00001330
120	NNTR=NTRAN(I)	00001340
	GO TO (130,140,150,160,170,175),NNTR	00001350
	WRITE (6,780)	
	STOP	
C		00001370
C	FOR OF X	00001380
C		00001390
130	X(I)=ALOG10(X(I))	00001400
	GO TO 180	00001410
C		00001420
C	POWER OF X	00001430
C		00001440
140	X(I)=(X(I)+CONS(I))**POWER(I)	00001450

GO TO 180	00001460
C	00001470
C SQUARE ROOT OF X	00001480
C	00001490
150 X(I)=(X(I))**.5	00001500
GO TO 180	00001510
C	00001520
C NATURAL LOG OF X	00001530
C	00001540
160 X(I)=ALOG(X(I))	00001550
GO TO 180	00001560
	00001570
C X(I-1) MULTIPLIED BY (X(J)**POWER(I)). TRANSFORMATION CARD	00001580
C CONTAINS (A) RESULTANT VARIABLE NUMBER, (B) CODE 5, (C) NUMBER	00001590
C OF VARIABLE MULTIPLIER (X(J)), AND (D) POWER OF X(J). (J NOT = I)	00001600
C	00001610
170 JJ=CONS(I)	00001620
X(I)=X(I-1)*(X(JJ)**POWER(I))	00001630
GO TO 180	00001640
C	00001650
C FACTOR TIMES X	00001660
C	00001670
175 X(I)=X(I)*CONS(I)	00001680
GO TO 180	00001690
C	00001700
C VARIABLE NOT CONSIDERED IN THIS PROBLEM.	00001710
C	00001720
178 X(I)=0.950	00001730
180 CONTINUE	00001740
DO 179 I=1,N	00001750
179 ISUB(I)=I	00001760
	00001770
C CHECK TO SEE IF ALL VARIABLES ARE USED IN THIS PROBLEM. (IFLAG)	00001780
C	00001790
ICOUNT=0	00001800
DO 181 I=1,N	00001810
IF (IFLAG(I).EQ.0) GO TO 181	00001820
ICOUNT=ICOUNT+1	00001830
GO 183 K2=1,K3	00001840
K5=I+K2	00001850
X(K5-ICOUNT) = X(K5)	00001860
ISUB(K5-ICOUNT)=K5	00001870
183 CONTINUE	00001880
181 CONTINUE	00001890
C	00001900
N=N-ICOUNT	00001910
C	00001920
IF (WT) 200,190,200	00001930
190 WT=1.0	00001940
	00001950

C		00001960
C	SUM OF OBSERVATIONS	00001970
C		00001980
	200 DO 210 I=1,N	00001990
	SUMX(I)=SUMX(I)+(WT*X(I))	00002000
C		00002010
C	SUMS OF SQUARES AND CROSS PRODUCTS STORED BY COLUMNS.	00002020
C		00002030
	DO 210 J=1,N	00002040
	210 SUMXX(I,J)=SUMXX(I,J)+(WT*X(I)*X(J))	00002050
	IF (WT-1.0) 220,230,220	00002060
	220 OBS=OBS+WT-1.0	00002070
	230 IOBS=IOBS+1	00002080
C		00002090
C	WRITE OBSERVATIONS ON UNIT 11 FOR USE WITH PREDICTIONS.	00002100
C		00002110
	IF (ICURVE.EQ.1) GO TO 235	00002120
	WRITE (11) (X(K),K=1,N)	00002130
	235 WRITE (6,790) IOBS,IFILE(IOBS)	00002140
	WRITE (6,800) (X(I),ISUB(I),I=1,N)	00002150
C		00002160
	IF (M-IOBS) 260,260,70	00002170
	260 IF (ICURVE.EQ.1) GO TO 262	00002180
	WRITE (6,810)	00002190
	GO TO 261	00002200
	262 WRITE (6,815)	00002210
	261 CONTINUE	00002220
	REWIND 11	00002230
	ICT = 1	00002240
	DO 270 I=1,N	00002250
C		00002260
C	MEANS.	00002270
C		00002280
	X(I)=SUMX(I)/OBS	00002290
C		00002300
C	STANDARD DEVIATIONS.	00002310
C		00002320
	SD(I)=(SUMXX(I,I)-(OBS*X(I)*X(I)))/(OBS-1.0)**0.5	00002330
	IF (ICURVE.EQ.1) GO TO 263	00002340
	IF (I.GT.NUMIND) GO TO 265	00002350
		00002360
	WRITE (6,820) XM(I),SD(I),ISUB(I),IVAR(I)	00002370
	GO TO 270	00002380
	263 WRITE (6,819) XM(I),SD(I),ISUB(I)	00002390
	GO TO 270	
	264 ICTT= ICT + 5	00002410
	WRITE (6,821) XM(I),SD(I),ISUB(I),(PRREF(LA),LA=ICT,ICTT)	00002420
	ICT = ICT + 6	00002430
	270 CONTINUE	00002440
C		00002450

DO 300 I=1,N	00002460
R(I,I)=1.0	00002470
K=I+1	00002480
IF (N-K) 300,280,280	00002490
280 DO 290 J=K,N	00002500
R(I,J)=((SUMXX(I,J)-(OBS*XM(I)*XM(J)))/(OBS-1.0))/(SD(I)*SD(J))	00002510
290 CONTINUE	00002520
300 CONTINUE	00002530
WRITE (6,830)	00002540
DO 310 I=1,N	00002550
WRITE (6,840)	00002560
410 WRITE (6,850) (R(I,J),ISUB(I),ISUB(J),J=I,N)	00002570
WRITE (6,860)	00002580
C	00002590
CALCULATE FOR NEW DEPENDENT VARIABLE.	00002600
C	00002610
NEND=0	00002620
320 IF (NEND) 330,350,330	00002630
330 K=N-L+1	00002640
IF (ICURVE.EQ.1) GO TO 335	00002650
READ (12) ((R(I,J),I=1,N),J=1,N)	00002660
REWIND 12	00002670
GO TO 337	00002680
335 READ (15) ((R(I,J),J=1,N),J=1,N)	00002690
REWIND 15	00002700
337 KNEND=K+NEND	00002710
DO 340 I=1,N	00002720
R(K,I)=R(KNEND,I)	00002730
340 R(I,K)=R(I,KNEND)	00002740
R(K,K)=1.0	00002750
SD(K)=SD(KNEND)	00002760
XM(K)=XM(KNEND)	00002770
GO TO 370	00002780
C	00002790
LOWER HALF OF COEFFICIENT MATRIX.	00002800
C	00002810
350 DO 360 I=1,N	00002820
J=I+1	00002830
DO 360 J=J,N	00002840
R(I,J)=R(I,J)	00002850
360 CONTINUE	00002860
C	00002870
WRITE ORIGINAL CORRELATION MATRIX ON UNIT 12 .	00002880
TO BE USED WITH SUCCEEDING DEPENDENT VARIABLES IF ANY.	00002890
C	00002900
TO BE USED WITH SUCCEEDING DEPENDENT VARIABLES IF ANY.	00002890
C	00002900
C	00002910
IF (ICURVE.NE.1) GO TO 364	00002920
WRITE (15) ((R(I,J),I=1,N),J=1,N)	00002930

REWIND 15	00002840
GO TO 359	00002850
364 WRITE (12) ((R(I,J),I=1,N),J=1,N)	00002860
REWIND 12	00002870
C	00002880
C REGRESSION ANALYSIS. CALCULATES BEST REGRESSION EQUATION WITH	00002890
C STANDARD ERROR OF EACH COEFFICIENT IN THE EQUATION.	00003000
C	00003010
C DEGREES OF FREEDOM	00003020
C	00003030
369 CONTINUE	00003040
370 DF=(OBS-1.0)	00003050
K=N-L+1	00003060
L	00003070
C STANDARD ERROR OF DEP. VAR.	00003080
C	00003090
NVIN=0	00003100
NVOUT=0	00003110
NTRY=0	00003120
380 SE=SD(K)*((OBS-1.0)*R(K,K)/DF)**.5	00003130
NTRY=NTRY+1	00003140
K=K-1	00003150
VMIN=1.0E35	00003160
VMAX=0.0	00003170
NMIN=0	00003180
NMAX=0	00003190
DO 470 I=1,K	00003200
IF (R(I,I)) 390,470,390	00003210
390 IF (R(I,I)-0.001) 470,470,400	00003220
400 VRNC=R(I,K+1)*R(K+1,I)/R(I,I)	00003230
IF (VRNC) 440,470,410	00003240
410 IF (VRNC-VMAX) 430,430,420	00003250
420 VMAX=VRNC	00003260
NMAX=I	00003270
430 SE(I)=0.0	00003280
SR(I)=0.0	00003290
GO TO 470	00003300
C	00003310
C CALCULATE REGRESSION COEFFICIENTS AND STANDARD ERRORS.	00003320
C	00003330
440 R(I)=R(I,K+1)*(SD(K+1)/SD(I))	00003340
SR(I)=(SE/SD(I))*(R(I,I)/(OBS-1.0))**.5	00003350
IF (VMIN) 460,450,450	00003360
450 VMIN=VRNC	00003370
NMIN=I	00003380
GO TO 470	00003390
460 IF (VRNC-VMIN) 470,470,450	00003400
470 CONTINUE	00003410
SUMBT=0.0	00003420
DO 480 I=1,K	00003430

480	SUMBI=SUMBI+B(I)*XM(I)	00003440
C		00003450
C	REGRESSION (PURE) CONSTANT	00003460
C		00003470
	RO=XM(K+1)-SUMBI	00003480
C		00003490
C	OUTPUT REGRESSION COEFFICIENTS	00003500
C	STD. ERROR OF COEFFICIENTS	00003510
C	STD. ERROR OF DEPENDENT VAR.	00003520
C		00003530
	IF (NTRY.NE.1) GO TO 481	00003540
	WRITE (6,871) NTRY,KSC	00003550
	GO TO 482	00003560
481	WRITE (6,870) NTRY	00003570
482	IF (NVOUT) 490,500,490	00003580
490	WRITE (6,880) ISUB(NVOUT)	00003590
	WRITE (6,890) FOUT	00003600
	GO TO 520	00003610
500	IF (NVIN) 510,520,510	00003620
510	WRITE (6,900) ISUB(NVIN)	00003630
	WRITE (6,890) FIN	00003640
520	WRITE (6,910) RO	00003650
	WRITE (6,920) (B(I),ISUB(I),I=1,K)	00003660
	WRITE (6,930) (SB(I),I=1,K)	00003670
	WRITE (6,940) SE	00003680
	IF (VMIN) 530,550,550	00003690
530	IF (VMIN*DF/R(K+1,K+1)+F2) 550,550,540	00003700
C		00003710
C	SELECT PIVOT ELEMENT FOR MATRIX INVERSION. (REMOVE VARIABLE)	00003720
C		00003730
540	KP=NMIN	00003740
	NVOUT=NMIN	00003750
	NVIN=0	00003760
	FOUT=(VMIN*DF/R(K+1,K+1))*(-1.0)	00003770
	DF=DF+1.0	00003780
	GO TO 570	00003790
550	IF (VMAX*DF/(R(K+1,K+1)-VMAX)-F1) 660,660,560	00003800
C		00003810
C	(ADD VARIABLE)	00003820
C		00003830
560	KP=NMAX	00003840
	NVIN=NMAX	00003850
	NVOUT=0	00003860
	FIN=VMAX*DF/(R(K+1,K+1)-VMAX)	00003870
	DF=DF-1.0	00003880
570	K=K+1	00003890
C		00003900
C	CALCULATE NEXT MATRIX. INVERT ONE VECTOR, STORE ON ORIGINAL MAT-	00003910
C	RIX. MUST READ BACK IN ORIGINAL MATRIX WHEN SOLVING NEW REGRESS-	00003920
C	ION FOR NEW DEPENDENT VARIABLE.	00003930



	DO 610 I=1,K	00003740
	IF (I-KP) 580,610,580	00003950
580	DO 600 J=1,K	00003960
	IF (J-KP) 590,600,590	00003970
C		00003980
C	I NOT =KP, J NOT =KP	00003990
C		00004000
	590 R(I,J)=(R(KP,KP)*R(I,J)-R(I,KP)*R(KP,J))/R(KP,KP)	00004010
	600 CONTINUE	00004020
	610 CONTINUE	00004030
	DO 620 I=1,K	00004040
	IF (I-KP) 620,630,620	00004050
C		00004060
C	I NOT =KP, J =KP	00004070
C		00004080
	620 R(I,KP)=(-R(I,KP))/R(KP,KP)	00004090
	630 CONTINUE	00004100
	DO 650 J=1,K	00004110
	IF (J-KP) 640,650,640	00004120
C		00004130
C	I =KP, J NOT =KP	00004140
C		00004150
	640 R(KP,J)=R(KP,J)/R(KP,KP)	00004160
	650 CONTINUE	00004170
C		00004180
C	I =KP, J =KP	00004190
C		00004200
	R(KP,KP)=1.0/R(KP,KP)	00004210
	GO TO 380	00004220
		00004230
	RESOLVE FOR NEW DEPENDENT VAR.	00004240
C		00004250
	660 NEND=NEND+1	00004260
C		00004270
C	MULTIPLE CORRELATION COEFFICIENT	00004280
C		00004290
	K=N-1+1	00004300
	RMULT=(1.0*(DE*SF**2)/((OBS-1.0)*SD(K)**2))**.5	00004310
	WRITE (6,950) RMULT	00004320
		00004330
	IF (ICURVE.EQ.1) GO TO 661	00004340
C		00004350
C	PRINT LAST SET OF COEFFICIENTS FOR THIS RMULT.	00004360
C		00004370
	WRITE (13) (B(IAD),IAD=1,NUMIND)	00004380
	WRITE(6,920) (B(IAD),ISUB(IAD),IAD=1,NUMIND)	00004390
	CORRCF(NEND) = RMULT	00004400
	GO TO 666	00004410
	661 IF (RMULT.LT.CORRCF(JCURVE)) GO TO 665	00004420
		00004430

WRITE(6,1000) NV	00004440
GO TO 666	00004450
665 WRITE(6,1010)	00004460
C	00004470
C PREDICTION OF DEPENDENT VARIABLE.	00004480
C	00004490
666 WRITE(6,960) KSC	00004500
YPRED=0.0	00004510
SUMD = 0.	00004520
SUMD2 = 0.	00004530
KPRED=N-L+NEND	00004540
KYY=N-L	00004550
C	00004560
REWIND 14	00004570
DO 680 I=1,M	00004580
IF (ICURVE.EQ.1) GO TO 685	00004590
READ (14) (X(K),K=1,N)	00004600
GO TO 686	00004610
685 READ (11) (X(K),K=1,N)	00004620
686 CONTINUE	00004630
DO 670 J=1,KYY	00004640
YPRED=YPRED+(R(J)*X(J))	00004650
670 CONTINUE	00004660
YPRED=YPRED+BO	00004670
AXX=X(KPRED)	00004680
DELTA=AXX-YPRED	00004690
SUMD = SUMD + AXX	00004700
SUMD2 = SUMD2 +(AXX)**2	00004710
WRITE (6,970) I,X(KPRED),YPRED,DELTA	00004720
680 YPRED=0.0	00004730
REWIND 11	00004740
REWIND 14	00004750
C	00004760
C DEVV = VARIANCE (FROM PAGE 210 IN TEXT).	00004770
C	00004780
SNUMR = SUMD2 - ((SUMD)**2/M)	00004790
DEVV = SNUMR / (M-1)	00004800
C SIGMA = STANDARD DEVIATION.	00004810
SIGMA = SQRT(DEVV)	00004820
XBARR = (SUMD/M)	00004830
CVV = (SIGMA/XBARR)	00004840
WRITE (6,980) XBARR,SIGMA,CVV	00004850
KSC = KSC + 1	00004860
IF (I-NEND) 1200,1200,1100	00004870
1100 WRITE (6,860)	00004880
GO TO 320	00004890
C	00004900
C START CHECK FOR NON-LINEAR REGRESSION ATTEMPTS.	00004910
C	00004920
1200 IF (ICURVE.EQ.1) GO TO 1210	00004930

REWIND 13	00004940
C ENTERING NONLINEAR REGRESSION MODE. (ICURVE=1).	00004950
ICURVE = 1	00004960
MIAMI = L	00004970
NUM1 = NUMIND	00004980
1210 JCURVE = JCURVE + 1	00004990
IF (JCURVE.GT.MIAMI) GO TO 1500	00005000
IF (CORRCF(JCURVE).LT..750E0) GO TO 1350	00005010
READ (13) (B1(J3),J3=1,NUM1)	00005020
1250 READ (11) (YY(II),II=1,NN)	
IF (EOF(11)) 1400,9999	
9999 IA = 1	
DO 1310 J4 = 1, NUM1	00005050
IF (B1(J4).EQ.0.0E0) GO TO 1310	00005060
X(IA) = YY(J4)	00005070
X(IA+1) = X(IA)**2	00005080
IA = IA + 2	00005090
1310 CONTINUE	00005100
X(IA) = YY(NUM1-JCURVE)	00005110
WRITE (14) (X(KA),KA=1,IA)	00005120
GO TO 1250	00005130
1350 READ (13)	00005140
GO TO 1210	00005150
C	00005160
1400 CONTINUE	00005170
REWIND 11	00005180
REWIND 14	00005190
N = IA	00005200
L = 1	00005210
GO TO 10	00005220
C	00005230
C FORMATS.	00005240
C	00005250
705 FORMAT (1H1,44X, ***** FTIS - NONLINEAR REGRESSION MODE FOR *****	00005260
1/58X,"DEPENDENT VARIABLE = VAR(",I2,"")//1H )	00005270
710 FORMAT (1H1,44X,"E490 - FTIS REGRESSION ANALYSIS"//1H )	00005280
720 FORMAT (19H NO. OF VARIABLES ,I3,34H NO. OF DEPENDENT VARIABLE	00005290
IFES ,I3//22H NO. OF OBSERVATIONS ,I5//28H F LEVEL TO ENTER VARIABLE	00005300
BLE ,F10.3,34H F LEVEL TO REMOVE VARIABLE ,F10.3/1H )	00005310
730 FORMAT (12,I2,F10.4,F10.4,5X,I1)	00005320
740 FORMAT (9H VARIABLE,I4,20H TRANSFORMED. TYPE,I4/1H )	00005330
745 FORMAT (9H VARIABLE,I4,20H TRANSFORMED. TYPE,I4,3X," * IFLAG=100005340	
X; THIS VARIABLE IS NOT CONSIDERED PART OF PROBLEM ",I9/1H )	00005350
750 FORMAT (1H ,24HTRANSFORMED DATA VALUES.//1X,"OBSERVATION / FTIS F100005360	
IF NO."/)	00005370
760 FORMAT (44HCOMPUTED GO TO INDEX OUTSIDE ALLOWABLE RANGE)	00005380
790 FORMAT (1H ,I3,"/",I2)	00005390
800 FORMAT (7H ,F10.5,5H VAR(,I2,3H) ,1H ,F10.5,5H VAR(,I2,3H) 00005400	
1 ,1H ,F10.5,5H VAR(,I2,3H) ,1H ,F10.5,5H VAR(,I2,3H) ,1H ,F10.5,5H VAR(,I2,3H) 00005410	
25H VAR(,I2,3H) )	00005420

```

810 FORMAT (1H1/1X,"          MEAN          STD. DEV.          VARIABLE 000054
1      DATA WORD / PHYS. PROP."/) 000054
815 FORMAT (1H1/1X,"          MEAN          STD. DEV.          VARIABLE 000054
1"/) 000054
819 FORMAT(9XF10.4,4XE12.4,7XI6) 000054
820 FORMAT(9XF10.4,4XE12.4,7XI6,16X,I3) 000054
821 FORMAT(9XF10.4,4XE12.4,7XI6,16X,6A4) 000054
830 FORMAT (1H /47H SIMPLE CORRELATION COEFFICIENTS. (ROW BY COL.)/1H 000055
1) 000055
840 FORMAT (1H ) 000055
850 FORMAT (1H F7.4,4H R(,I2,1H,,I2,6H) ,F7.4,4H R(,I2,1H,,I2,6H) 000055
1) ,F7.4,4H R(,I2,1H,,I2,6H) ,F7.4,4H R(,I2,1H,,I2,6H) 000055
2 ,F7.4,4H R(,I2,1H,,I2,6H) ) 000055
860 FORMAT (1H1) 000055
870 FORMAT (1H /14H TRIAL NUMBER ,I5/1H ) 000055
871 FORMAT (1H /14H TRIAL NUMBER ,I5,3X,"FOR VARIABLE (",I2,")"/1H ) 000055
880 FORMAT (23H VARIABLE GOING OUT = ,I4/1H ) 000055
890 FORMAT (20H F LEVEL ,F12.4/1H ) 000056
900 FORMAT (23H VARIABLE GOING IN = ,I4/1H ) 000056
910 FORMAT (20H PURE CONST. B(0) = ,E12.4/1H ) 000056
920 FORMAT (1X"COEFFICIENTS"/(1XE12.4," B(",I2,") ",E12.4," B(",I2, 000056
1") ",F12.4," B(",I2,") ",E12.4," B(",I2,") ",E12.4," B(",I2, 000056
2")")) 000056
930 FORMAT (1H /31H STANDARD ERROR OF COEFFICIENTS// 000056
1(1XE12.4,7XE12.4,7XE12.4,7XE12.4,7XE12.4)) 000056
940 FORMAT (1H /28H STANDARD ERROR OF ESTIMATE //7H ,E12.4//1H ) 000056
950 FORMAT (1H //33H MULTIPLE CORRELATION COEFFICIENT//7H ,F10.500056
1/) 000057
960 FORMAT(1H1,19X,"ACTUAL VS. PREDICTED RESULTS FOR VARIABLE (",I2,") 000057
1"/1X,"OBSERVATION",8X, 000057
2"ACTUAL PREDICTED DEVIATION "/) 000057
970 FORMAT (17,F20.4,5F16.4/1H ) 000057
980 FORMAT(/1X,"MEAN = ",F12.6,4X,"STD. DEV. = ",F12.6,4X,"COEFFICIENT 000057
1 OF VARIATION = ",F12.6) 000057
1000 FORMAT(1X,"**** CONCLUSION NON-LINEAR REGRESSION MODEL PROVIDES 000057
1THE BEST FIT FOR VARIABLE (",I2,") ."/1X,"***** 000057
2***** 000057
3") 000058
1010 FORMAT(1X,"**** CONCLUSION LINEAR REGRESSION MODEL PROVIDES THE 000058
1T FIT FOR THIS VARIABLE. ****"/1X,"***** 000058
2***** 000058
3") 000058
1500 CONTINUE 000058
RETURN 000058
END 000058

```

```

PROGRAM F0410F101INPUT,OUTPUT,TAPF5=INPUT,TAPF6=OUTPUT,TAPF8,TAPF9)
C*****
C
C      E410 / FOURIER TRANSFORM INFRARED SPECTROSCOPY.
C
C      PHYSICAL PROPERTIES MASTER TAPE GENERATOR.
C*****
C
C      THICKOL CORPORATION / HUNTSVILLE , ALABAMA 35807
C
C      PRINCIPAL INVESTIGATOR   W. W. SCHWARZ
C                                TELEPHONE (205) - 882 - 2388
C
C      SCIENTIFIC PROGRAMMER   O. C. SMITH
C                                TELEPHONE (205) - 882 - 8215
C*****
C
C      SEPTEMBER , 1970.
C
C      FORTRAN IV - A LEVEL 21 LANGUAGE
C
C      CDC 6600 -- (AFRPI)
C
C      PREPARED IN FULFILLMENT OF CONTRACT F04611 - 78 - C - 0027
C*****
C
C      DIMENSION PH(10,10),TITLE(10),IAGEWK(100),ITEMPA(100),ITESTT(100)
C      DIMENSION IREC(100)
C      COMMON N1,N2
C      DATA CREATE/'CREATE' ,%/UPDATE/'UPDATE' %/
C
C
C      MASTER = 3
C      N1 = 5
C      N2 = 6
C      NREC = 0
C
C      READ (N1,400) RUNTYP
C
C      IF (RUNTYP.EQ.UPDATE) GO TO 650
C      IF (RUNTYP.NE.CREATE) GO TO 20
C      GO TO 30
C20  CONTINUE
C      WRITE(N2,470) RUNTYP
C      STOP 11
C30  CONTINUE
C

```

```

      READ (N1,480) TITLE
      WRITE (N2,520) TITLE
      WRITE(N2,560)
      WRITE (N2,530)
      WRITE (N2,535)
      WRITE (N2,540)
      WRITE (N2,550)
      WRITE (N2,545)
      WRITE (N2,550)
C
      ICOUNT=11
10  READ (N1,500) FLAG,K1,K2,K3
      IF (FLAG(N1)) 200,999
999 ICOUNT=ICOUNT+2
      NREC=NREC+1
      IAGEWK(NREC) = K1
      ITEMPA(NREC) = K2
      ITESTT(NREC) = K3
      READ (N1,510) (PH(NREC,M),M=1,10)
      WRITE (N2,600) NREC,IAGEWK(NREC),ITEMPA(NREC),ITESTT(NREC),(PH(NREC,
      -C,JJ),JJ=1,10)
      IREC(NREC)=NREC
      IF (ICOUNT.GE.60) GO TO 40
      GO TO 10
40  CONTINUE
      WRITE(N2,570)
      WRITE(N2,540)
      WRITE (N2,530)
      WRITE (N2,535)
      WRITE (N2,540)
      WRITE (N2,530)
      WRITE (N2,545)
      WRITE (N2,550)
      ICOUNT=7
      GO TO 10
C
650 CONTINUE
      CALL PENAME
C
      GO TO 250
C
200 CONTINUE
      WRITE(MASTER) TITLE
      DO 240 N=1,NREC
      WRITE(MASTER) IREC(N),IAGEWK(N),ITEMPA(N),ITESTT(N),(PH(N,JJ),
      1JJ=1,10)
240 CONTINUE
C
C *****
C

```

```

470  FORMAT(10I,1X,"-(1 ",AR," IS AN INVALID MODE. NO CHANGES APPLIED T
      1) DATA-SET. EXECUTION TERMINATED N G . . .",/,10I)
480  FORMAT (10A6)
490  FORMAT (AR)
500  FORMAT(AR,7X,15,2I10)
510  FORMAT (5(F10.2))
520  FORMAT ((1I/1X,"FOURIER TRANSFORM INFRARED SPECTROSCOPY - F410 ON
      -YSICAL PROPERTIES MASTER TAPE GENERATOR"
      7/1X,"HIAEPL/PCO - INOXOL/INTSVILLE)"7/1X,"PROPERTIES FOR SOLID "
      ARPELLANT",13X,10A6)
530  FORMAT (1X,"REF AGE AGE TEST MODULUS STRAIN AT STRAIN AT
      -MAXIMUM STRAIN STRAIN")
535  FORMAT (1X,"HC. TEMP TEMP BREAK MAXIMUM
      - STRESS ENERGY ENDURANCE")
540  FORMAT (45X,"STRESS",15X,"DENSITY")
545  FORMAT (5X,"(WK) (F) (F) (PSI) (I) (I) (PS
      -I) (TPO) (I) ")
550  FORMAT (1X,"-----
      -----")
560  FORMAT (/1X,14,1X,14,1X,14,7X,15,1X,10(F9.3,2X))
560  FORMAT (1X,"-----
      -----
      -----")
570  FORMAT(10I)
C
C*****
C
250  CONTINUE
      REMIND MASTER
      STOP
      END

```

# SUBROUTINE RENAME

```

C*****
C
C.....SEPTEMBER , 1979.
C
C
C      CHANGES THE DATA SET ESTABLISHED IN C410 MAIN.
C
C*****
C      COMMON N1,N2
C      DIMENSION JAGENK(100),JTEMPA(100),JTESTT(100),NREC(100),
C      IPHI(100,10),SAVGAT(100,10),IAGENK(100),ITEMPA(100),ITESTT(100)
C      DIMENSION COUNT(30),V(10),TITLE(10)
C      DATA ENDLIST/ENDLIST*/,CHANGE/"CHANGE  "/*DELETE/"DELETE  "/*
C      IADD/"ADD      "/*ENDNIL/"ENDNILST*/
C      MASTER=8
C
C      WRITE(N2,10)
C      LINTOT=0
C      I=1
C      READ (MASTER) TITLE
C12  READ(MASTER) NREC(I),JAGENK(I),JTEMPA(I),JTESTT(I),(PHI(I,JJ),JJ=1
C      - ,10)
C      IF (I.EQ.(MASTER)) 14,999
C888 IF(NREC(I).GT.LINTOT) LINTOT=NREC(I)
C      I=I+1
C      GO TO 12
C14  CONTINUE
C      REMIND MASTER
C*****
C
C      CHECK THE INPUT SEQUENCE FOR ALL THE ACTION CARDS.
C
C*****
C      IFIRST = 0
C      ISHIFT = 0
C      LAST = 0
C      LOOP = 0
C      I1 = 0
C      MIDDLE = 0
C5  READ(N1,20)ACTION,NP,K1,K2,K3
C      IF (EOF(N1)) 30,277
C77  LOOP = LOOP + 1
C      IF (ACTION.NE.CHANGE.AND.(LOOP.EQ.1)) IFIRST = 1
C      IF (ACTION.EQ.CHANGE) GO TO 70
C      IF (ACTION.EQ.AND) GO TO 60
C      IF (ACTION.EQ.DELETE) LAST = 1
C      GO TO 25
C70  IF(IFIRST.EQ.1.AND.LOOP.GT.1) GO TO 55

```



```

      IF (MIDDLE.NF.1) GO TO 25
      IF (LAST.NF.1) GO TO 25
      WRITE (6,710) LOOP,ACTION
      STOP 60
55  WRITE (6,710) LOOP,ACTION
      STOP 62
60  MIDDLE = 1
      LI=LI+1
      KOUNT(LI) = NR
      IF (LAST.EQ.0) GO TO 25
      WRITE (6,710) LOOP,ACTION
      STOP 64
C
C      WRITE TO TEMP. DISK FILE 9.
C
25  WRITE (9) ACTION,NR,K1,K2,K3
      IF (ACTION.NF.AND.AND.ACTION.NF.CHANGE) GO TO 5
      READ(NI,125) (V(K),K=1,10)
      WRITE (9) (V(K),K=1,10)
      GO TO 5
C
30  CONTINUE
      RETURN
C
15  READ(9) ACTION,NR,K1,K2,K3
      IF (END(9)) 650,666
666  IF(ACTION.EQ.ENDNOI) GO TO 600
      IF(ACTION.EQ.ENDIST) GO TO 500
      IF(NR.LE.0) GO TO 450
      IF(NP.GT.100) GO TO 450
      IF(ACTION.EQ.CHANGE) GO TO 100
      IF(ACTION.EQ.DELETE) GO TO 200
      IF(ACTION.EQ.ADD) GO TO 300
      GO TO 400
C *****
C ACTION = CHANGE ROUTINE
C *****
100  CONTINUE
      JAGFWK(NP)=K1
      JTEMPA(NR)=K2
      JTESTT(NR)=K3
C
      READ(9) (PHI(NP,I),I=1,10)
C
      WRITE(N2,125) NR
      WRITE(N2,170) NR, JAGFWK(NR), JTEMPA(NR), JTESTT(NR), (PHI(NR, JJ)
      -,JJ=1,10)
      IF(NR.GT.LINTOT) LINTOT=NR
      GO TO 15
C *****

```

C ACTION = DELETE ROUTINE

\*\*\*\*\*

```
200 CONTINUE
    DO 205 LX=1,L1
    IF (MP.GT.KOUNT(LX)) ISHIFT=ISHIFT+1
205 CONTINUE
    MI = NR+ISHIFT
    DO 230 J=MI,LINTOT
    M=J+1
    IF(M.GE.100) STOP 72
    NREC(J)=NREC(M)
    JAGFWK(J)=JAGFWK(M)
    JTEMPA(J)=JTEMPA(M)
    JTESTT(J)=JTESTT(M)
```

C

```
    DO 220 I=1,10
    PH1(J,I)=PH1(M,I)
220 CONTINUE
230 CONTINUE
    LINTOT=LINTOT-1
    WRITE(N2,240) MP
    GO TO 15
```

\*\*\*\*\*

C ACTION = ADD ROUTINES

\*\*\*\*\*

```
300 CONTINUE
```

C

```
    IF(MP.GT.LINTOT) GO TO 340
    DO 330 J=1,LINTOT
    IF(MP.EQ.NREC(J)) GO TO 310
```

C

```
    GO TO 230
```

C

```
310 CONTINUE
    I=LINTOT+1
    DO 320 K=J,L
    DO 314 I=1,10
    SAVDAT(K,I)=PH1(K,I)
    CONTINUE
    JAGFWK(K)=JAGFWK(K)
    JTEMPA(K)=JTEMPA(K)
    JTESTT(K)=JTESTT(K)
```

```
320 CONTINUE
```

C

```
    DO 315 K=J,L
    M=K+1
    DO 317 I=1,10
    PH1(M,I)=SAVDAT(K,I)
317 CONTINUE
    JAGFWK(M)=JAGFWK(K)
```

```

      JTEMPA(M)=JTEMPA(K)
      JTESTT(M)=JTESTT(K)
315  CONTINUE
330  CONTINUE
C
      LINTOT=L
      GO TO 350
C
340  CONTINUE
C
      LINTOT=LINTOT+1
350  CONTINUE
C
      JAGEWK(NR)=K1
      JTEMPA(NR)=K2
      JTESTT(NR)=K3
C
      READ(7) (PHI(NR,I),I=1,10)
      IF (NR.GT.LINTOT) GO TO 380
      WRITE(N2,360) NR
      GO TO 15
380  WRITE(N2,365) NR
      GO TO 15
C
C*****
C ERROR MESSAGES FOR INVALID UPDATE ACTION AND BAD LINE NUMBER
C*****
400  CONTINUE
      WRITE(N2,410) ACTION
      WRITE(N2,420)
      GO TO 500
430  CONTINUE
      WRITE(N2,460) NR
      STOP 70
C****
C ACTION=ENDLIST (LINES RENUMBERED, DATA STORED AND PRINTED)
C*****
500  WRITE(MASTER) TITLE
C
      WRITE(N2,590) TITLE
      WRITE(N2,570)
      WRITE(N2,530)
      WRITE(N2,540)
      WRITE(N2,550)
      WRITE(N2,560)
      WRITE(N2,565)
      WRITE(N2,560)
C
      ICOUNT=11
      DO 580 J=1,LINTOT

```

```

        ICOUNT=ICOUNT+2
        NPEC(J)=J
        WRITE(N2,510) NPEC(J),JAGEWK(J),JTEMPA(J),JTESTT(J),(PH1(J,I),I=
11,10)
        WRITE(MASTER) NPEC(J),JAGEWK(J),JTEMPA(J),JTESTT(J),(PH1(J,I),I=
11,10)
C
        IF(ICOUNT.GE.61) GO TO 575
        GO TO 590
575  CONTINUE
        WRITE(N2,501)
        WRITE(N2,570)
        WRITE(N2,520)
        WRITE(N2,540)
        WRITE(N2,550)
        WRITE(N2,560)
        WRITE(N2,565)
        WRITE(N2,560)
        ICOUNT=7
590  CONTINUE
        GO TO 650
C*****
C ACTION = FNUMPL (LISTING AND RENUMBERING SUPPRESSED)
C*****
600  WRITE(MASTER) TITLE
        GO 640 J=1,LINTOT
        WRITE(MASTER) NPEC(J),JAGEWK(J),JTEMPA(J),JTESTT(J),(PH1(J,I),I=
11,10)
640  CONTINUE
C
C*****  FORMAT.
C
10  FORMAT(1H1," **UPDATE MODE ASSUMED.....NEW DATA FOLLOWS**",//)
20  FORMAT(A8,2X,215,2(10)
125  FORMAT(5F10.3)
125  FORMAT(//," (( CORRECTED DATA FROM LINE ",I3," ))")
170  FORMAT(//,I3,1X,I4,1X,I4,2X,I4,1X,10(F9.3,2X))
240  FORMAT(//," -(( LINE ",I3," DELETED FROM DATA SET")
        FORMAT(//," (( A NEW LINE WILL BE ADDED PRECEDING LINE ",I3)
145  FORMAT(//," (( NEW LINE ",I3," ADDED TO DATA SET ")
145  FORMAT(1H1," -(( ",A9," )) IS AN ILLEGAL UPDATE FUNCTION",/,
1" *-----**")
420  FORMAT(//," ALL CHANGES PRECEDING ABOVE FUNCTION WERE APPLIED",
1/,1H1)
440  FORMAT(1H1," -(( ",I3," )) IS AN INVALID LINE NUMBER. NO CHANGES
APPLIED TO DATA. EXECUTION TERMINAT I N  S . . . ",/,1H1)
501  FORMAT(1H1)
510  FORMAT(//,1X,I3,1X,I4,1X,I4,2X,I4,1X,10(F9.3,2X))
530  FORMAT(1X,"NPEC AGE TEST MODULUS STRAIN AT STRAIN AT "
1AXIMUM STRAIN STRAIN")

```

```

540  FORMAT(1X,"NO.      TEMP  TEMP      BREAK      MAXIMUM
      -STRESS  ENERGY  ENDURANCE")
550  FORMAT(45X,"STRESS",15X,"DENSITY")
560  FORMAT(1X,"-----
      -----")
565  FORMAT (5X,"(4K)  (F)  (F)  (PSI)      (I)      (I)      (PS
      -I)      (TRD)      (I)")
570  FORMAT(" =====
      1=====
      2=====")
570  FORMAT (11H1/1X,"SECURITE TRANSFORM INFRARED SPECTROSCOPY - E410 PH
      -YSICAL PROPERTIES MASTER TAPE GENERATOR"
      2/1X,"(APPL/PCO - THICKOL/HUNTSVILLE)"//1X,"PROPERTIES FOR SOLID P
      3ROPELLANT ",9X,10A4)
710  FORMAT (1H0," ((( ACTION CARD NUMBER -",13," IS AN OUT OF SEQUENC
      IE ",A6," COMMAND. JOB IS TERMINATED....")

```

C

```

650  CONTINUE
      REWIND MASTER
      DEWIND 2
      RETURN
      END

```

**APPENDIX B**

**PROGRAM FLOWCHARTS**

FORTRAN MODULE E490 - MAIN PROGRAM

CHART TITLE - INTRODUCTORY COMMENTS

CHART TITLE - PROCEDURES

(000001)	2.06	(000002)	2.07		
(000002)	2.07 3				
(000004)	2.09	(000090)	2.11		
(000090)	2.11 8				
(000093)	2.14	(000091)	2.12		
(000102)	2.23 20	(000099)	2.21		
(000104)	2.26	(000105)	2.20		
(000105)	2.28 22				
(000106)	2.29 23	(000101)	2.22		
(000112)	2.36 26	(000109)	2.34		
(000119)	2.44	(000131)	3.02		
(000124)	2.47	(000120)	2.48		
(000128)	2.48 29				
(000131)	3.02 28	(000119)	2.44		
(000133)	3.03 27	(000111)	2.35		
(000135)	3.06	(000134)	3.04		
(000136)	3.08	(000135)	3.06		
(000137)	3.10	(000136)	3.08		
(000142)	3.11 5	(000152)	3.18	(000234)	5.14
(000150)	3.18	(000149)	3.16		
(000153)	4.01 9	(000146)	3.13		
(000182)	4.12	(000184)	4.13		
(000184)	4.13 40				
(000187)	4.17 50	(000195)	4.22	(000198)	4.25
(000196)	4.23 55	(000192)	4.20	(000193)	4.21
(000200)	4.26 70	(000194)	4.22		
(000202)	4.28	(000205)	4.29		
(000205)	4.29 68				
(000208)	4.31	(000212)	5.02		
(000210)	5.01 69	(000208)	4.32		
(000212)	5.02 60	(000209)	4.32		
(000220)	5.03 65	(000200)	4.24		
(000225)	5.07 400	(000211)	5.05		
(000234)	5.14	(000233)	5.12		
(000235)	5.15 470	(000231)	5.09		
(000237)	5.17	(000239)	5.19		
(000239)	5.19 430	(000237)	5.17		
(000244)	5.21	(000245)	5.22		
(000245)	5.22 440				
(000251)	5.24 445	(000272)	6.19		
(000254)	6.02	(000250)	6.07		
(000258)	6.05	(000257)	6.03		
(000259)	6.07 460				
(000259)	6.07	(000258)	6.05		
(000262)	6.11	(000251)	6.09		
(000263)	6.13	(000262)	6.11		
(000264)	6.15	(000263)	6.13		
(000266)	6.17	(000264)	6.15		
(000276)	7.01 450	(000271)	6.16		
(000282)	7.04 490	(000282)	7.19		
(000293)	7.19 470	(000283)	7.17		
(000297)	7.20 475	(000314)	7.34		
(000300)	7.25 52	(000300)	7.30	(000311)	7.33
(000309)	7.31 551	(000305)	7.28	(000316)	7.29
(000313)	7.34 800	(000307)	7.30		
(000315)	8.01 300	(000313)	7.34		
(000318)	8.04 510	(000324)	8.08		
(000326)	8.09 500	(000317)	8.03	(000323)	8.09
(000329)	8.12	(000330)	8.14		
(000330)	8.13 600				
(000340)	8.19	(000341)	7.20		
(000341)	8.20 710				
(000343)	8.21 700	(000337)	8.17		
(000344)	8.25	(000342)	8.28		
(000346)	8.27	(000345)	8.25		
(000347)	8.28 625				
(000348)	8.31	(000341)	7.32		
(000361)	8.32 650				
(000369)	8.34 650	(000347)	8.22		
(000371)	8.36	(000360)	8.34		





12/11/79

AUTOFLOW CHART SET - E490 - MAIN PROGRAM

PAGE 01

CHART TITLE - INTRODUCTORY COMMENTS

\*\*\*\*\*  
\*\*\*\*\*  
\*\*\*\*\*

E490 - FOURIER TRANSFORM INFRARED SPECTROSCOPY PROGRAM.

THIOL CORPORATION / HUNTSVILLE, ALABAMA 35807

PRINCIPAL INVESTIGATOR W. W. SCHWARZ  
TELEPHONE (205) - 892 - 8388

SCIENTIFIC PROGRAMMER C. C. SMITH  
TELEPHONE (205) - 982 - 8215

\*\*\*\*\*  
\*\*\*\*\*  
\*\*\*\*\*

READS POLYS F115 DATA TAPES.

SEPTEMBER, 1970.


FORTRAN IV -- H EXTENDED LANGUAGE

CDC 6600 -- (INFRPL)


\*\*\*\*\*  
\*\*\*\*\*  
\*\*\*\*\*

\*\*\*\*\*  
\*\*\*\*\*  
\*\*\*\*\*

CALL XFILE (PPPL PROGRAM).

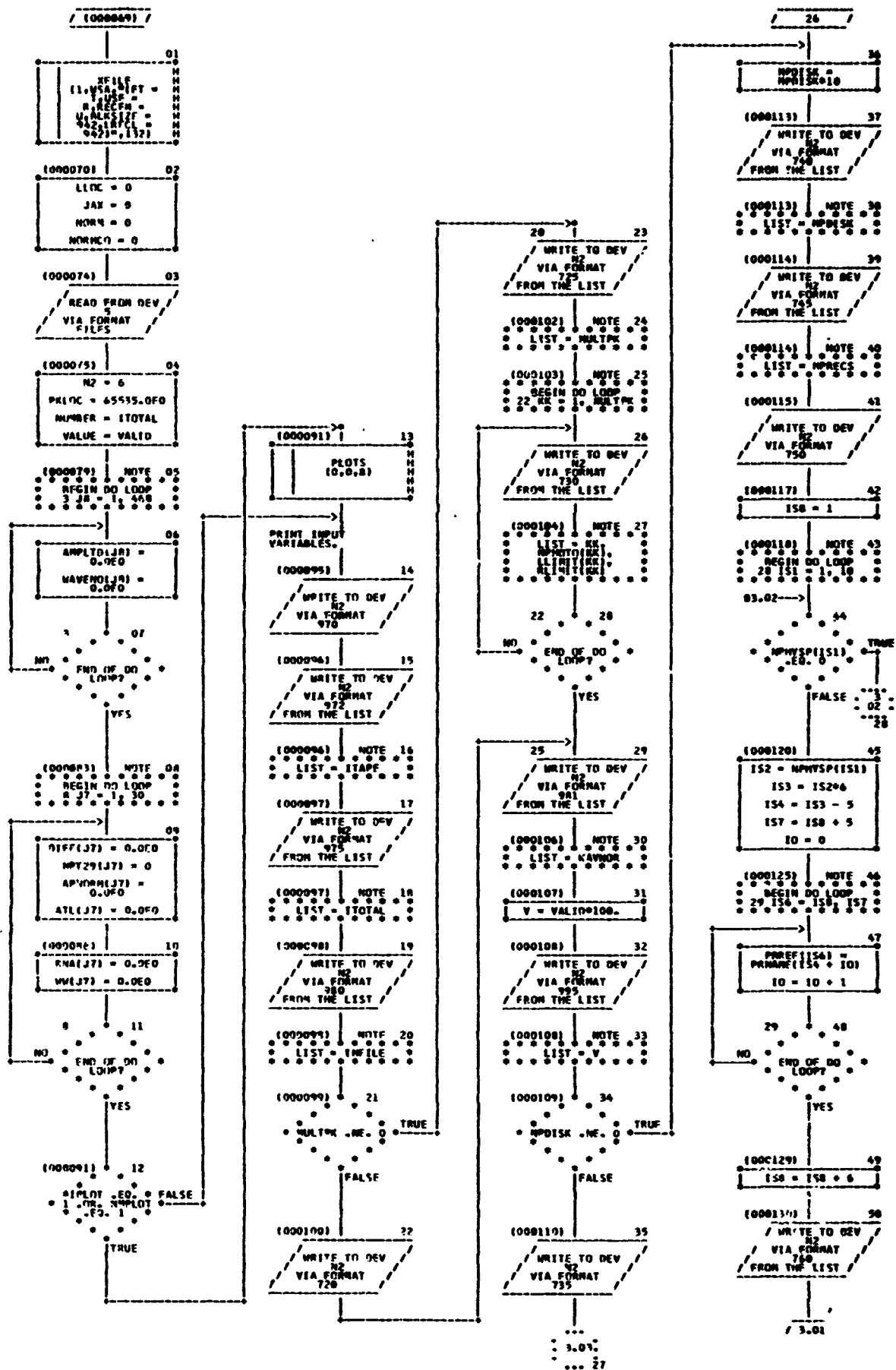**CHAPT TITLE - PROCEDURES**

**CHART TITLE - PROCEDURES**

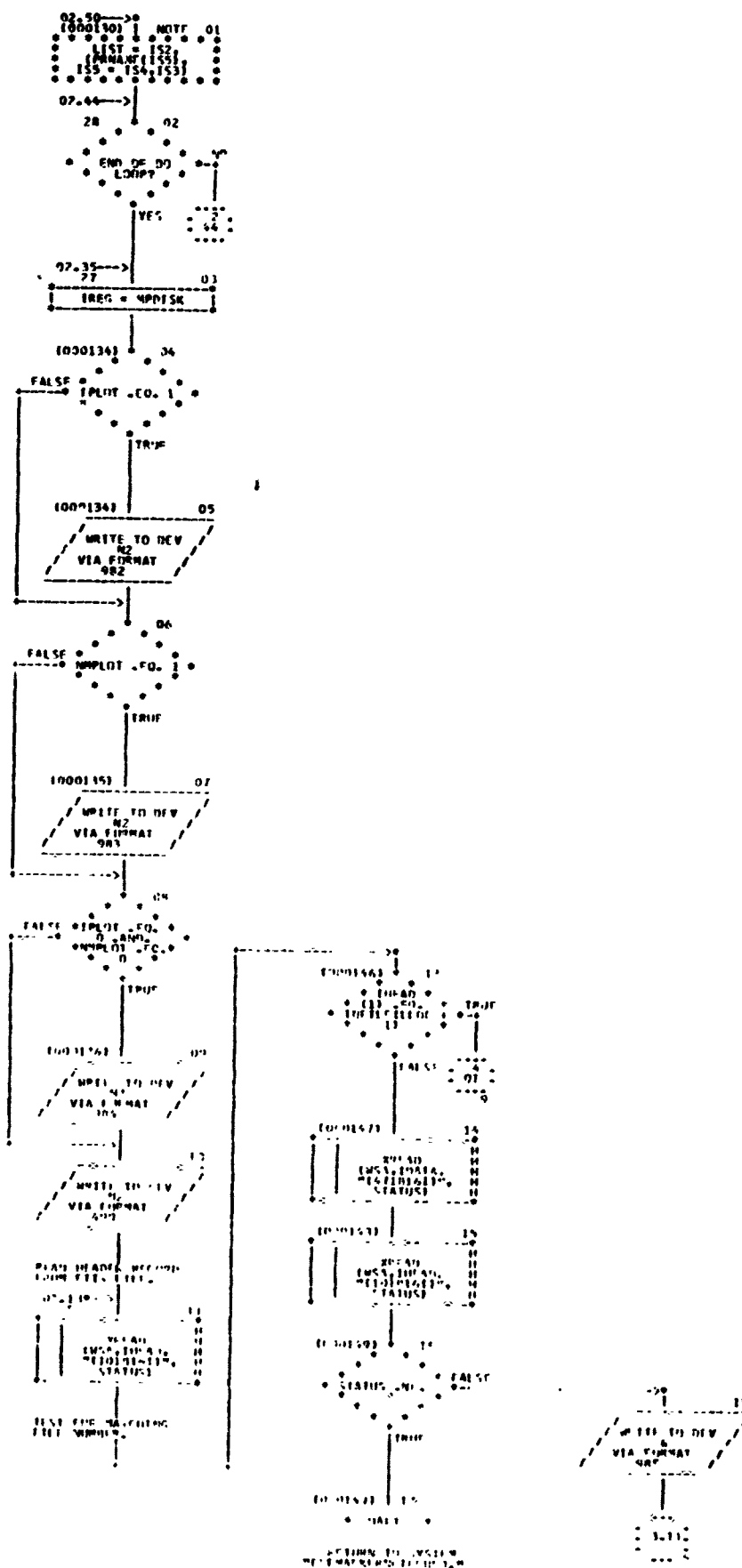

```

      /-----/
      01.17-->|
      -----|----- 01
      LLF = LLF + 1
      -----|-----
      (000154)
      /-----/
      WRITE TO DFV
      VIA FORMAT
      1000
      FROM THE LIST
      -----|-----
      (000154)  NOTE 03
      LIST 04
      STATUS 05
      06
  
```

### CALCULATION OF CONSTANTS

```

(000160) 24
FILE(117) =
  INFAD(1)
NUMRDS(1100) =
  INFAD(5)
VSP4 = INFAD(6)
VSP7 = INFAD(7)

```

(J00164) | 05  
| TM = 460705(1100) |

TEP1	=	INHEAD(10)
TEP2	=	INHEAD(11)
TEP3	=	INHEAD(12)
TEP4	=	INHEAD(13)

(170170) | 3

r w  
170047X2)  
/27.68C)

X1 = TFP1 +  
(T2-27.68) \* r

X2 = TFP2 +  
(T2-27.68) \* w

DELTA = X2 -  
X1 / (W1 - 1)

(003174) 30  
/ WRITE TO DFV /  
/ VIA FLIPAT /  
/ ECU /  
/ END THE LIST /

(000174) 107E 70  
LIST = 11, 17,  
7FLVA

(C09176) 1

RECAP  
(MSG) 11/11/54  
-1471/11/11/54  
STATIST

CALCULATION OF ABOVE  
NUMBER AND ABOVE TIME  
(X) (Y)

10071811 10071811  
• • • • •  
• 40-1 OF 100  
• 40-1 - 1. MM  
• • • • •

WAVENET(41) = K1 +  
(41 - 1909-214  
AMPLTD(41) =  
F0128TA(01) = 73

(J00145) 14  
WRITE TO DEV  
N2  
VIA FORMAT  
340  
FROM THE LIST

```

(000145)  NOTE 14
* * * * *
* LIST =
* IFILE(ILLOC)
* * * * *

```

(0001'6)	16
$\kappa = 0$	

50	1	17
$K = K + 1$ $V1 = K + 57$ $V7 = K1 + 52$		

1P  
WRITE TO NEW  
VIA FORMAT  
FROM THE LIST

```

( ) 0215 )      WCTF 10
* * * * *
*   LIST = K,
*   MAXV=10(K),
*   AMP1=2(K), K1,
*   MAXW=11,
*   AMP1=7(K1), K2,
*   MAXV=10(K2),
*   AMP1=2(K2)
* * * * *

```

20

FBI SF

1991 71

15 JULY 2004

• K.F.O. 364 •

FALSE

4.17.30

```

graph TD
    Start[04-27-29] --> Write[/WRITE TO DEV  
VIA FORMAT  
FROM THE LIST/]
    Write --> WriteList["(000196) WRITE  
LIST  
IF L < LLOC"]
    WriteList --> IncrementK["(000177) K = K + 104"]
    IncrementK --> End[25]
  
```

4.17

\_\_\_\_\_

70 76

MULTIPLY.EQ. 0  
FALSE

10032011 0076 27  
REGTN ON LGCP  
69 47 = 1.6

```
ALIMIT(46) =  
LLIMIT(46)  
  
ZLIMIT(46) =  
PLIMIT(46)  
  
%ELF(46) =  
VOMT(46)
```

49 22

10

11111111

AD R 1, W 12K

\* \* \* \* \*  
 \* ZLIMITAKI \* ST. \*  
 \* VAVIR \*  
 \* ZLIMITAKI \* LT. \*


10017041 32

SECRET


FALSE

1  
...  
: 9.27:  
... 99

## CHART TITLE - PROCEDURES



















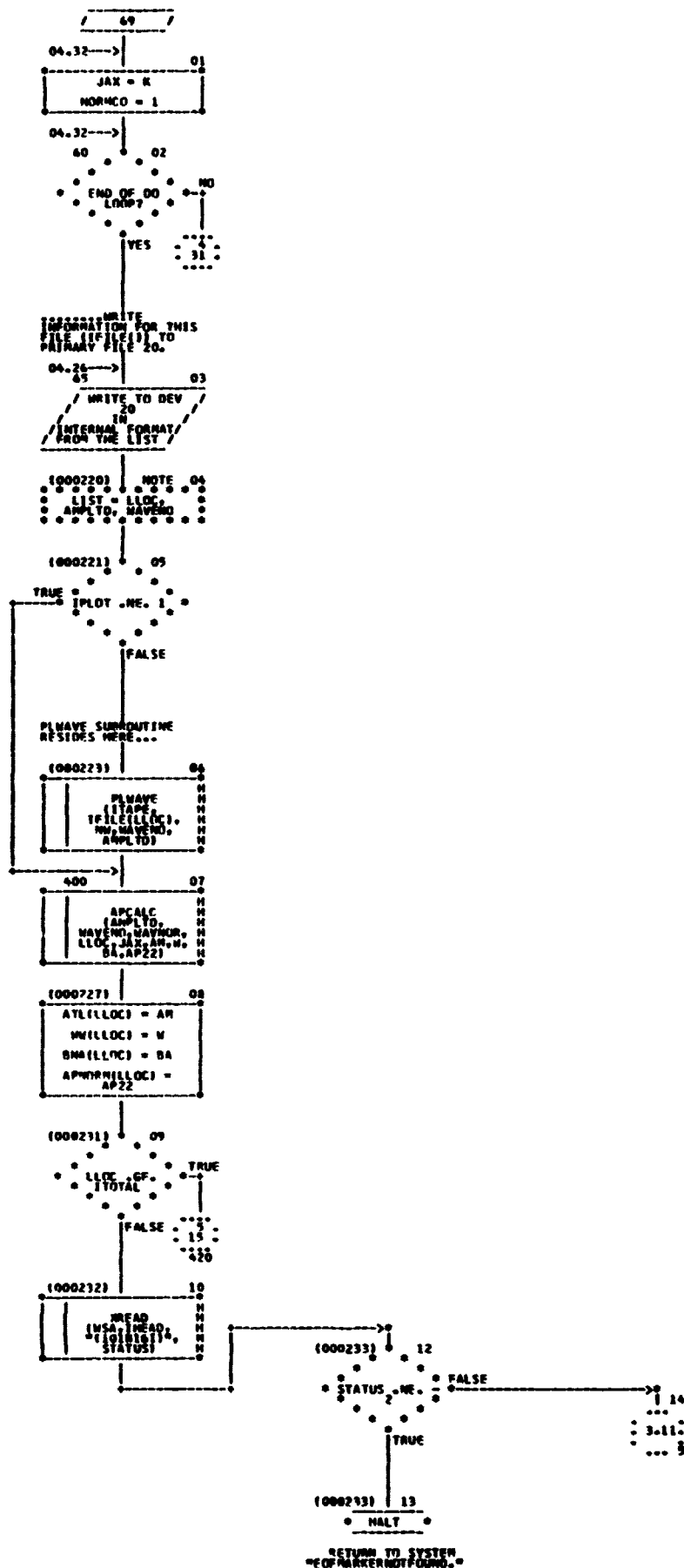


## CHART TITLE - PROCEDURES















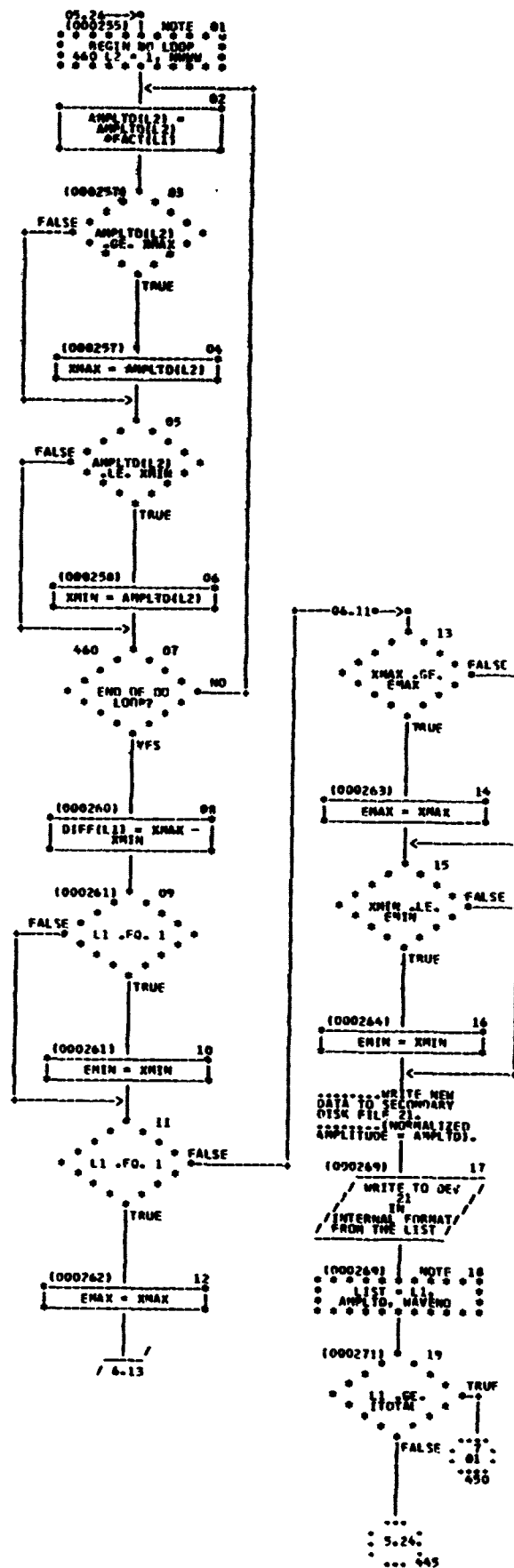


## CHART TITLE - PROCEDURES

PRINT OUT  
NORMALIZATION  
INFORMATION.

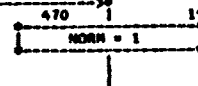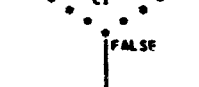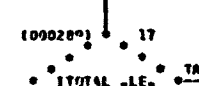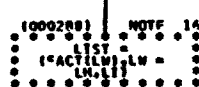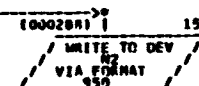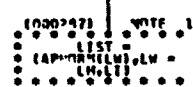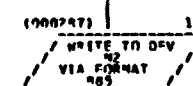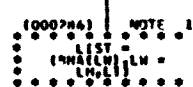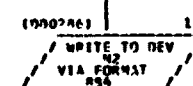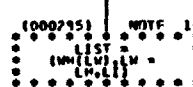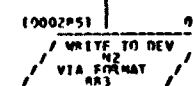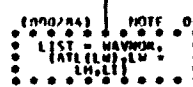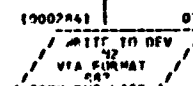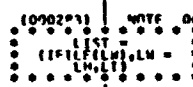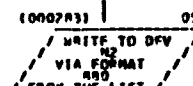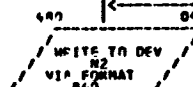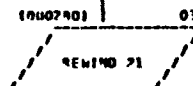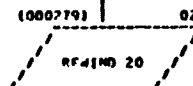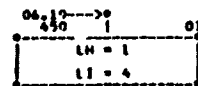

LISTING OF NORMALIZED  
DATA.

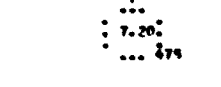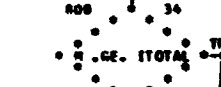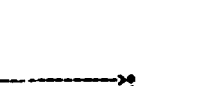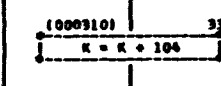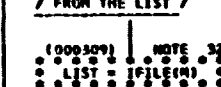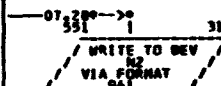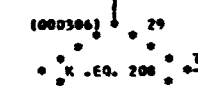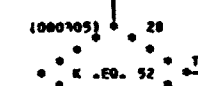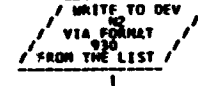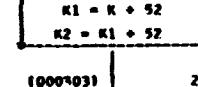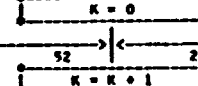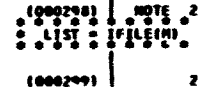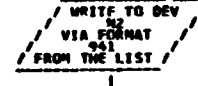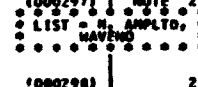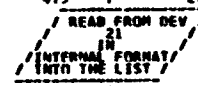

## CHART TITLE - PROCEDURES

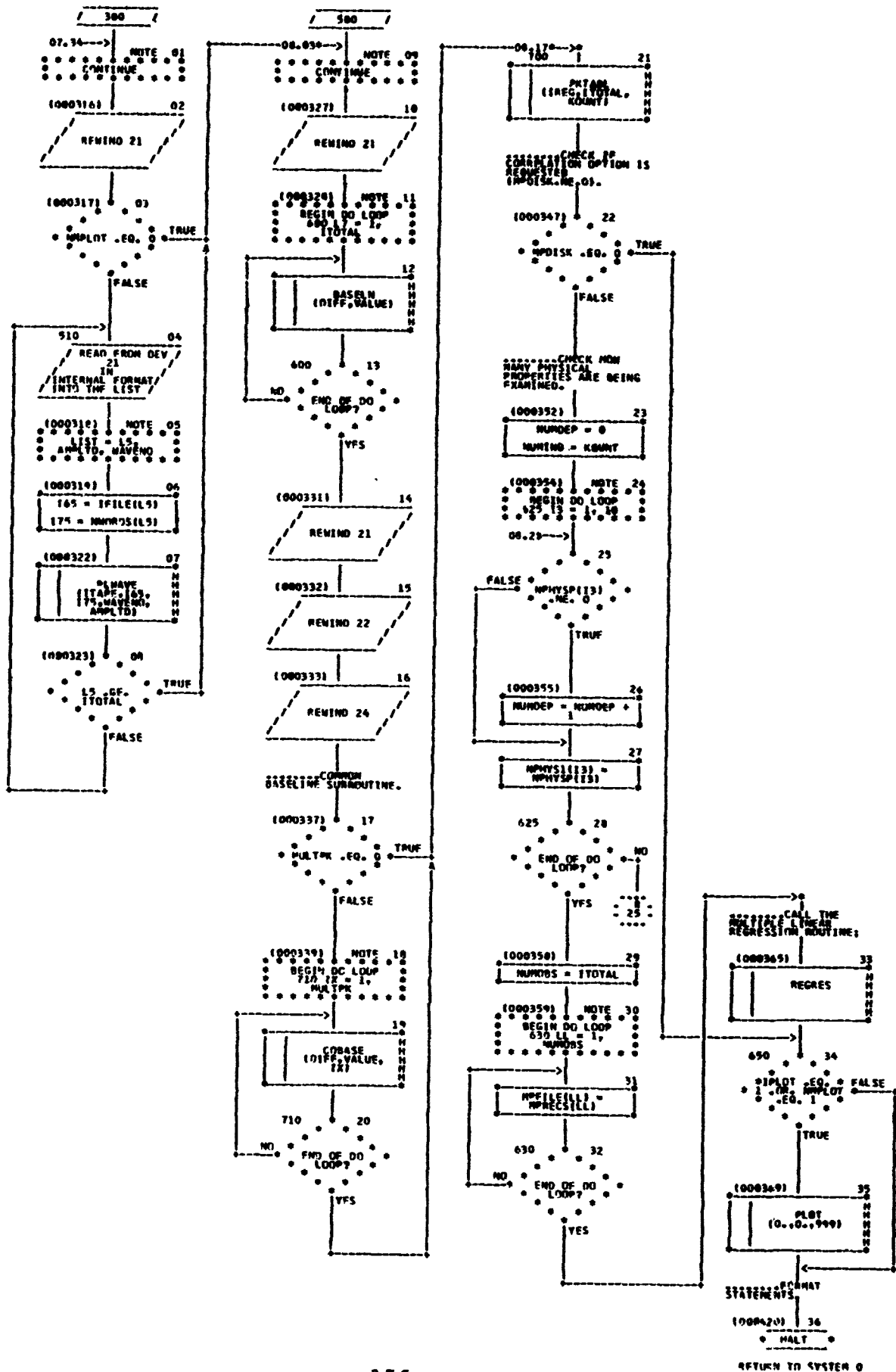




CHART TITLE - NON-PROCEDURAL STATEMENTS

```

PROGRAM F490FTN(INPUT,OUTPUT,TAPE5=INPUT,TAPE6=OUTPUT,TAPE1,TAPE20,
,TAPE21,TAPE22,TAPE23,TAPE24,TAPE25,TAPE26,TAPE27,TAPE11,TAPE12,TA
PE13,TAPE14,TAPE15,TAPE10)
COMMON /CORREL/ PRPF(40),MPFILE(30),
      NPHYS(10),NUMDP,NUMIND,NUMMS,IPEG
COMMON /SPECTR/ IFILE(30),NPT29(30)
COMMON /POINTS/ NMORS(30)
COMMON /PTINF/ ENAX,EMIN,NORM
COMMON /TWINPK/ ALIMIT(4),ZLIMIT(4),MFILE(6),MORCO,NUMBER
      DIMENSION AMPLTD(464),WAVENM(468),PRNAME(60),FACT(30)
      DIMENSION APNORM(30),ATL(30),NMA(30),NM(30),DIFF(30)
      DIMENSION RLIMIT(6),INFILE(30),NPRECS(30),NPHOTO(6),NPHYS(10)
      DIMENSION WSA(17),INFAD(10),IDATA(471)
      REAL LLIMIT(4)
      NAMELIST /FILES/ ITAPE,ITOTAL,INFILE,MULTPK,NPHOTO,LLIMIT,RLIMIT,
      VALID,NAVNOR,NPOISK,NPRECS,NPHYS,IPLOT,NMPLT
      DATA INFILE /30*0/,
      NPHOTO /6*0/,
      IPLOT /0/,
      LLIMIT /6*0.0/,
      NPOISK /0/,
      NPRECS /30*0/,
      MULTPK /0/,
      NMPLT /0/,
      NPHYS /10*0/,
      RLIMIT /6*0.0/,
      VALID /-07/
      DATA PRNAME /"MODU","LUS "," "," "," "," ",
      "STRA","IN A","T DR","FAK "," "," ",
      "STRA","IN A","Y RA","XIMU","N ST","RESS",
      "MAXI","NUM ","STAE","SS "," "," ",
      "STRA","IN E","NERG","Y DE","NSIT","Y ",
      "STRA","IN C","NUMM","ANCE"," "," "," ",
      "PHYS","ICAL"," PROP","PERT","Y NO",". 7 ",
      "PHYS","ICAL"," PROP","PERT","Y NO",". 4 ",
      "PHYS","ICAL"," PROP","PERT","Y NO",". 9 ",
      "PHYS","ICAL"," PROP","PERT","Y NO",". 10"/
400  FORMAT(1H1)
905  FORMAT(/,2X,"END ENCOUNTERED",/)
1000 FORMAT(1H1,7X,"FILE ",I3,2X,"STATUS = ",F4.1,/,26(1X,4(120,5X),/))
      )
715  FORMAT(A2,1X,1A2,4X,4A2)
720  FORMAT(/5X,"NUMBER OF COMMON BASELINE AREAS SPECIFIED 0.")
725  FORMAT(/5X,"NUMBER OF COMMON BASELINE AREAS SPECIFIED ",I1,".")
      )
730  FORMAT(10X,"AREA ",I1," IS BEST PICTURED IN FILE ",I2," : ",F10.4,10X,"
LEFT LIMIT ESTIMATE = ",F7.2,/,10X,"RIGHT LIMIT ESTIMATE = ",F7.2)
735  FORMAT(/5X,"NO STATISTICAL CORRELATION WITH PHYSICAL PROPERTIES W
ILL BE PERFORMED. -- (NPOISK=0).")
740  FORMAT(/5X,"STATISTICAL CORRELATION WITH PHYSICAL PROPERTIES HAS B
EEN REQUESTED //5X,"PHYSICAL PROPERTY INPUT DISK FT=",I2,"F001."
      )
745  FORMAT(/5X,"PHYSICAL PROP. RECORD NUMBERS ",10(13,5X)/36X,10(13,5
X)/36X,10(13,5X))
740  FORMAT(/5X,"THE FOLLOWING PHYSICAL PROPERTIES WILL BE USED//5X,"AS
DEPENDENT VARIABLES ")
760  FORMAT(10X,"PROP. NO. ",I1," - ",I4A4)
840  FORMAT(1H1,4X,"SUMMARY OF AMPLITUDE NORMALIZATION //"
54X,3("FILE",16X),"FILE"/5X,3("NO",18X),"NO")
880  FORMAT(54X,1(13,17X),17/)
882  FORMAT(/7X,"MAXIMUM AMPLITUDE NEAREST ",F4.1," UN. = ",4(E15.7,5
X)/)
883  FORMAT(7X,"WAVE NUMBER AT MAX. AMPLITUDE = (NMMAX) = ",F10.4,3(10
X,F10.4)/)
884  FORMAT(7X,"BASELINE AMPLITUDE AT (NMMAX) = ",8X,4(E15.7,5X)/)
885  FORMAT(7X,"PEAK HEIGHT AT (NMMAX) = ",19X,4(E15.7,5X)/)
890  FORMAT(1H1,/(120(10))
900  FORMAT(1X,3(5X,F10.4))
930  FORMAT(6X,13,4X,F10.4,2X,E15.7,2(7X,13,4X,F10.4,2X,E15.7))
940  FORMAT(1H1,4X,"THIS FILE NUMBER ",I3,20X,"NON-NORMALIZED (PURE) S

```

12/05/79  
CARD ID PAGE/BOX NAME  
TABLE OF CONTENTS AND REFERENCES  
AUTOFLG= CHART SET -  
REFERENCES (SOURCE SEQUENCE NO. AND PAGE/BOX)

FORTRAN MODULE E450 - SUBROUTINE APCALC

CHART TITLE - SUBROUTINE APCALC (AMPLTD, WAVEC, WAVNCR, IFN, JAX, APMAX, ME, O2, AP2)

(000003)	1-01	APCALC		
(000014)	1-03	1		
(000014)	1-03	(000014)	1-04	
(000018)	1-05	5	(000019)	1-06
(000022)	1-10		(000021)	1-08
(000026)	1-11	10	(000020)	1-07
(000028)	1-13		(000031)	1-15
(000031)	1-15	15	(000028)	1-13
(000043)	1-19	20	(000045)	1-20
(000059)	1-21	45	(000038)	1-17
(000060)	1-22	200	(000061)	1-23
(000063)	1-25	210	(000067)	1-26
(000066)	1-27	30	(000043)	1-19
(000050)	1-29	35	(000052)	1-30
(000068)	1-31	220	(000063)	1-25
(000069)	1-32	230	(000071)	1-33
(000053)	1-34	40	(000050)	1-29
(000072)	2-01	240	(000069)	1-32
(000073)	2-02	245	(000075)	2-03
(000076)	2-04	250	(000073)	2-02
(000079)	2-05	100	(000055)	1-34
(000052)	2-08	75	(000104)	2-11
(000122)	2-12	110	(000084)	2-06
(000108)	2-15	80	(000088)	2-07
			(000099)	2-10
			(000120)	2-16
			(000115)	2-17

CHART TITLE - NON-PROCEDURAL STATEMENTS

12/05/79	PROCEDURAL STATEMENT LABEL INDEX		AUTOFLOW CHART SET -		E490 - SUBROUTINE APCALC		PAGE	1
	PG.BX	NAME	PG.BX	NAME	PG.BX	NAME		
	1.01	PCALC						
	1.03	1	1.21	45	2.12	110	1.32	230
	1.05	5	2.08	75	1.22	200	2.01	240
	1.11	10	2.15	80	1.25	210	2.02	245
	1.15	15	2.05	100	1.31	220	2.04	250

CHART TITLE - SUBROUTINE APCALC(AMPLTD,NAVENO,NAVNO,IPN,JAX,APHAX,WE,D2,AP2)

/ APCALC /

1999\*\*\*SEPTEMBER 4:  
1999\*\*\*CDC REVISION:

\*\*\*\*\* LOGIC WILL  
NORMALIZATION  
APPLTD VALUE  
\*\*\*\*\* FOR THE PEAK  
NAVENO HAVE NUMBER  
(NAVENO) CH-1.


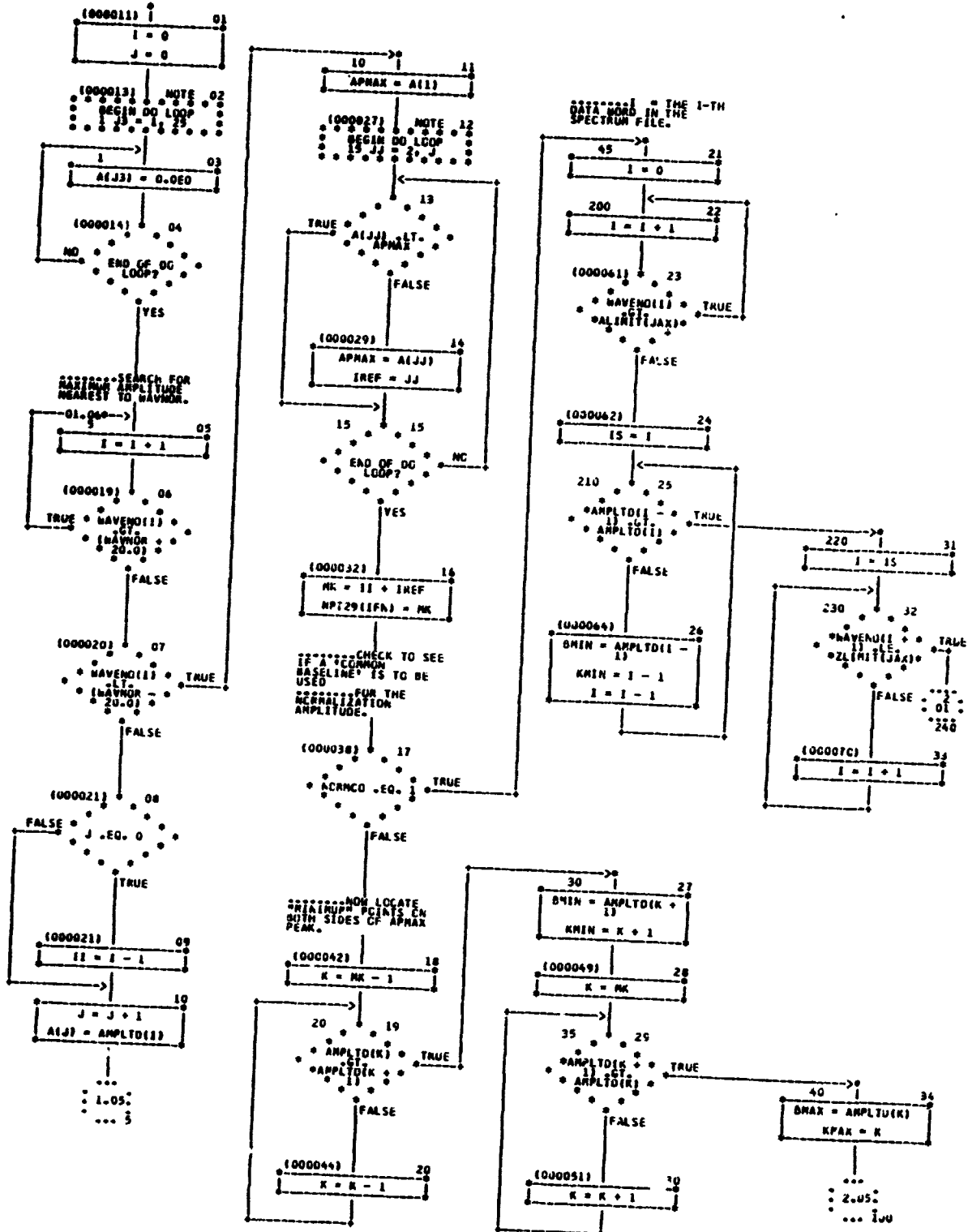


CHART TITLE - SUBROUTINE APCALC(AMPLTD, HAVEND, HAVNDR, IFN, JAX, APHAX, ME, D2, AP2)

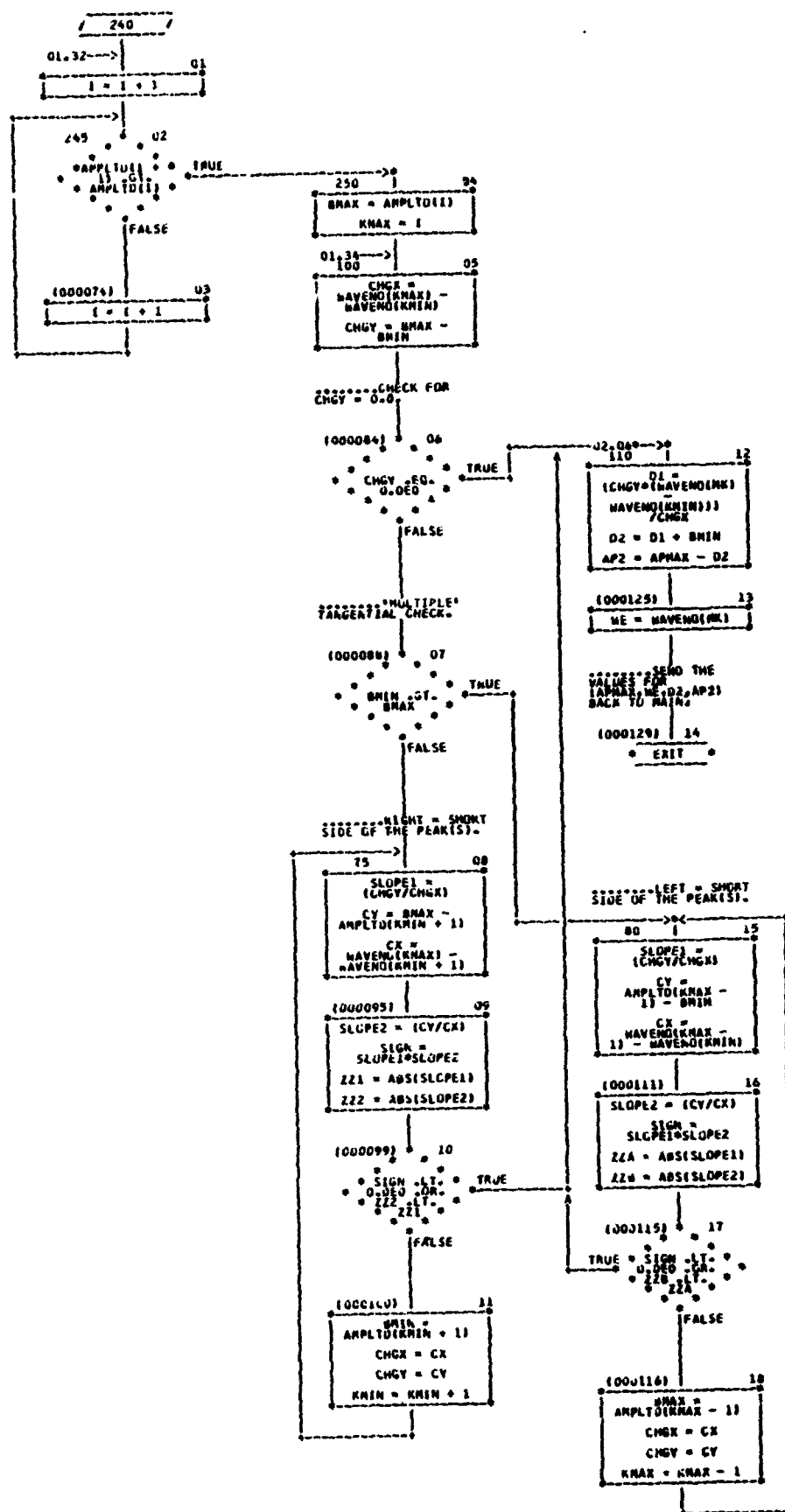

12/05/79

CHART TITLE - MCN-PROCEDURAL STATEMENTS

AUTOFLOW CHART SET -

F490 - SUBROUTINE APCALC

PAGE 03

```
COMMON /SPECTR/ IFILE(30),MPT25(30)
COMMON /TJNPK/ ALIMIT(6),ZLIMIT(6),MFILE(6),NOMMCO,NUMBER
DIMENSION AMPLTD(468),MAVENO(468),A(25)
```

12/20/79	TABLE OF CONTENTS AND REFERENCES	
CARD ID	PAGE/BGA	NAME
FORTRAN MODULE		
BASEN - SUBROUTINE BASELN		
CHART TITLE - SUBROUTINE BASELN(DIFF,VALUE)		
(000003)	1-01	BASELN
(000037)	1-03	(000030) 1-04
(000038)	1-04	1000
(000040)	1-06	
(000042)	1-07	1100
(000050)	2-06	10
(000052)	2-08	15
(000056)	2-11	20
(000148)	2-14	50
(001065)	2-16	130
(000063)	2-25	30
(000151)	2-26	100
(000066)	3-01	40
(000072)	3-05	50
(000116)	3-07	10
(000075)	3-11	60
(000084)	3-15	
(000100)	3-18	75
(000130)	3-19	110
(000155)	3-23	120

CHART TITLE - NON-PROCEDURAL STATEMENTS

12/05/79

PG-BX	NAME	PG-BX	NAME	PG-BX	NAME	PG-BX	NAME	PG-BX	NAME
1.01	15E1A	2.25	3C	3.18	75	2.26	106	4.16	13C
2.06	1U	3.01	4C	3.07	80	3.19	113	1.04	1000
2.08	15	3.05	50	2.14	50	3.23	120	1.07	1100
2.11	2C	3.11	60						



12/05/74

AUTOFLOW CHART SET -

E490 - SUMMUTINE BASELN

PAGE 01

CHART TITLE - SUMMUTINE BASELN(DIFF,VALUE)

/ BASELN /

\*\*\*\*\*SEPTEMBER  
12, 1974 - REVISION  
FOR CDC 660C.

\*\*\*\*\*THIS  
SUBROUTINE SEARCHES  
FOR AND RECORDS VALID  
PEAKS  
\*\*\*\*\*WITHIN EACH  
INDIVIDUAL SPECTRUM.


\*\*\*\*\*ALLOP FOR  
PAR. 50 PEAKS PER  
INDIVIDUAL SPECTRUM.

\*\*\*\*\*I = THE I-TH  
DATA WORD IN THE  
SPECTRUM FILE.

\*\*\*\*\*IP = THE  
I-TH VALID PEAK IN  
THE SPECTRUM  
(LEFT-RIGHT).

\*\*\*\*\*10=CH01IP1 =  
C: 1 NO PEAK AT DATA  
WORD 'IP1'.  
1: 1 INVALID PEAK AT  
DATA WORD 'IP1'.  
2: 1 VALID PEAK AT  
DATA WORD 'IP1'.  
3: 1 INVALID PEAK AT  
DATA WORD 'IP1'.  
4: 1 VALID PEAK AT  
DATA WORD 'IP1'.  
MEASURED  
BY THE "COMASE"  
SUBROUTINE.  
5: 1 VALID PEAK AT  
DATA WORD 'IP1'.  
MEASURED  
BY THE "COMASE"  
SUBROUTINE.

\*\*\*\*\*INITIALIZE  
SCOPE VARIABLES.

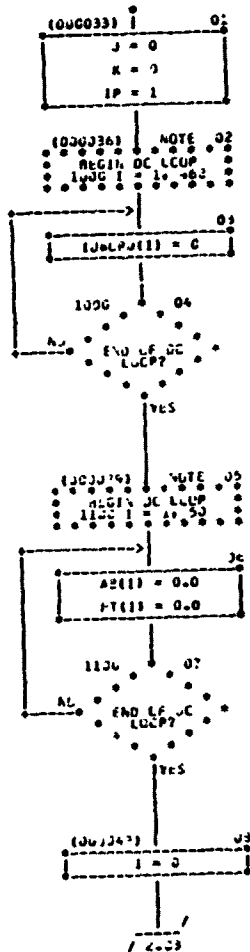




12/05/76

AUTOFLOW CHART SGT -

7490 - SUBROUTINE BASELINE

PAGE

CHART TITLE = SUBROUTINE BASELINE DIFF. VALUE

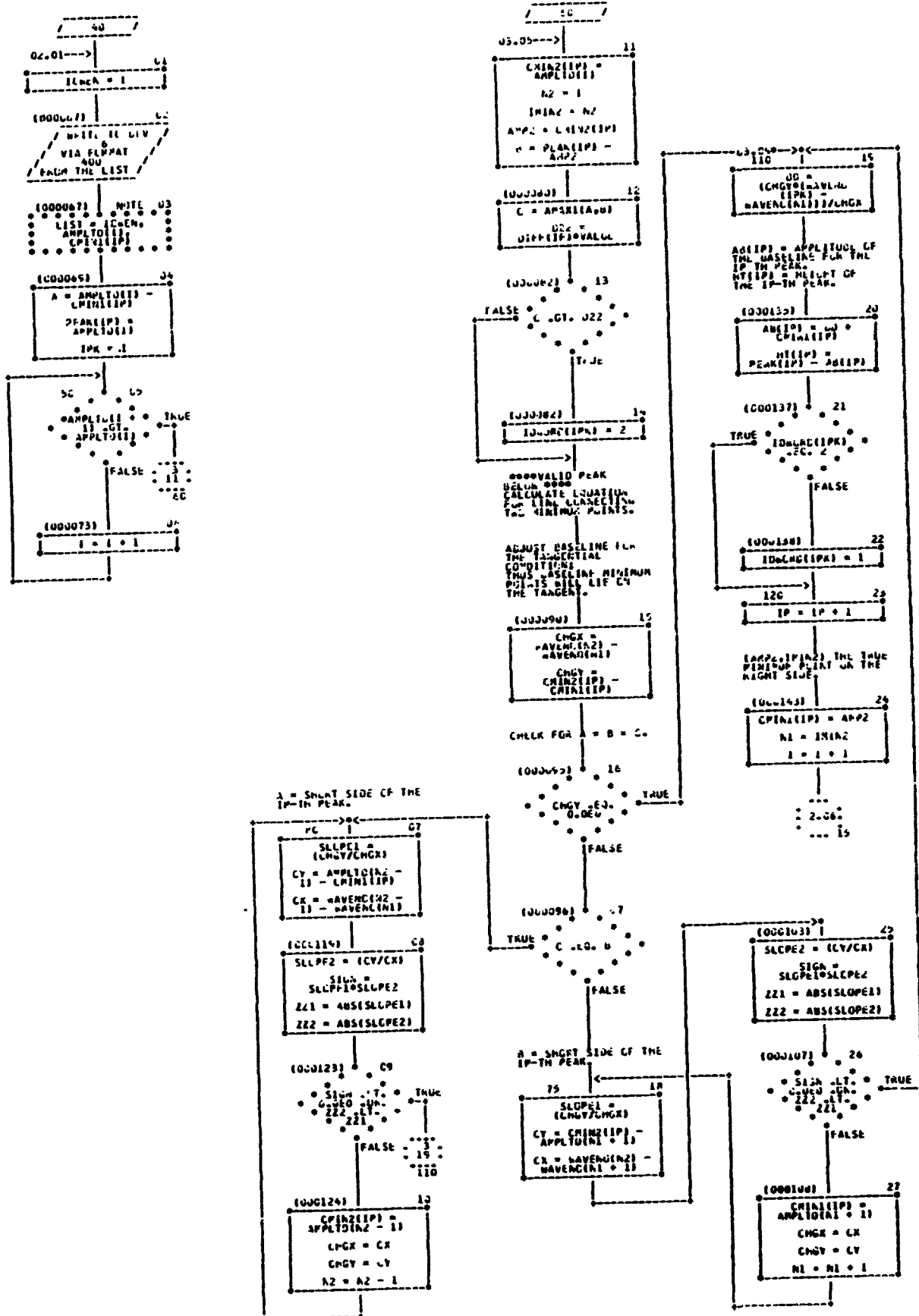

12/05/79

CHANT TITLE - NON-PROCEDURAL STATEMENTS

```
COMMON /POINTS/ NMCPTS(30)
DIMENSION AMPLTD(468),#AVENC(468),IDCJNC(468)
DIMENSION CHIN1(50),CHIN2(50),PEAK(50),AB(50),MT(50),DIFF(30)
FORMAT(3X,15.3X,E15.7,3X,E15.7)
400  FORMAT (3X,15)
5000  FORMAT (5(15X,10011))
5001
```

FORTHMAN "LEONE" 1490 - SUBROUTINE CHARTS

CHART TITLE - SUBROUTINE (CLASSIFIED, VALUE, JAX)

(000003)	1.01	0.45			
(000012)	1.02		(000014)	1.03	
(000034)	1.03	2			
(000041)	1.04		(000046)	1.05	
(000045)	1.05	4	(000048)	1.06	
(000047)	1.06	7	(000049)	1.07	
(000052)	1.11	6	(000051)	1.21	(000051) 1.25
(000063)	1.22	7	(000055)	1.17	
(000067)	1.23	16	(000067)	1.24	
(000068)	1.24	20	(000073)	1.27	
(000074)	1.26	30	(000085)	1.28	
(000075)	1.29	35	(000077)	1.30	
(000074)	2.01	40	(000075)	1.29	
(000081)	2.02	55	(000082)	2.05	
(000087)	2.04	65	(000083)	2.02	
(000092)	2.06	77	(000105)	2.14	
(000095)	2.07	85	(000097)	2.08	
(000107)	2.09	900	(000097)	2.04	
(000098)	2.10	70	(000095)	2.07	
(000094)	2.11	71	(000101)	2.12	
(000102)	2.13	75	(000095)	2.11	
(000115)	2.21	87	(000104)	2.14	
(000116)	2.22	45	(000111)	2.05	
(000128)	2.05	100	(000117)	2.03	
(000141)	2.08	100	(000115)	2.15	
(000157)	2.09	104	(000117)	2.07	(000105) 2.12
(000177)	2.15	110	(000113)	2.06	(000104) 2.11 (000103) 2.14
(000181)	2.16	111	(000170)	2.17	
(000192)	2.20		(000164)	2.73	
(000183)	2.22		(000152)	2.20	
(000184)	2.23	105			
(000186)	2.26		(000135)	4.01	
(000185)	4.01	116			
(000192)	4.03	117	(000175)	4.24	
(000195)	4.05		(000203)	4.05	
(000203)	4.04	1.0			
(000203)	4.05		(000203)	4.07	
(000210)	4.10	200			
(000219)	4.16		(000211)	4.23	
(000224)	4.21		(000227)	4.19	
(000227)	4.22	225	(000222)	4.18	
(000231)	4.23	220	(000219)	4.16	(000224) 4.21
(000236)	4.24		(000242)	4.23	
(000242)	4.25	220			
(000243)	4.26	2.3	(000235)	4.24	
(000251)	4.27		(000255)	4.15	
(000255)	4.28	210			
(000262)	4.29	205	(000261)	4.24	

CHART TITLE - RUN-PRODUCTION STATEMENTS

PROCEDURAL STATEMENT LABEL INDEX			AUTOFLOW CHART SET -			E490 - SUBROUTINE COBASE		
PG-BX	NAME	PG-BX NAME	PG-BX NAME	PG-BX NAME	PG-BX NAME	PG-BX NAME	PG-BX NAME	PG-BX NAME
1-01	CBASE	1-28 3C	2-11 71	3-19 111	4-23 220			
1-03	2	1-29 35	2-13 75	4-03 115	4-22 225			
1-08	4	2-01 40	3-01 80	4-01 116	5-04 230			
1-09	5	2-02 55	3-02 85	3-09 120	5-03 300			
1-11	6	2-04 60	3-05 100	3-08 130	5-19 310			
1-22	7	2-06 62	3-23 105	4-09 150	5-26 400			
1-23	10	2-07 65	3-16 110	4-10 200	2-04 600			
1-26	20	2-10 70						

CHART TITLE - SUBROUTINE CUPASE(UPP,VALUE,JA)

/ CUPASE /

1976... SUPERMAN JR.  
 1976... IL DANGED AT  
 AF/PL 1 OCT. 24,  
 1976.

.....LOC REVISION.

.....COMMON  
 BASELINE MONITOR /  
 MANH 19, 1976.  
 (CCS)

.....THIS  
 SUBROUTINE IS A USLM  
 REQUESTED LPTLM.  
 .....THIS LPTLM IS  
 USED TO ESTABLISH  
 A COMMON BASELINE  
 .....FOR ADJACENT  
 SPECTRAL PEAKS.

.....I = THE 1-TH  
 DATA BLK IN THE  
 SPECTRUM FILE.

.....IP = THE  
 18-TH PEAK IN THE  
 COMMON BASELINE  
 (LEFT-TO-RIGHT).

.....MAXIMUM  
 NUMBER OF PEAKS PER  
 BASELINE = 25.

.....MAXIMUM  
 NUMBER OF COMMON  
 BASELINES PER  
 COMPUTER RUN = 6.

(000032) NOTE 01  
 BEGIN OF LOOP  
 2 2 = 1 25

(000032) ISPELJ21 = U

2 03  
 END OF LOOP?

YES

SEARCHING FOR THE  
 FILE (INFILLJX) WITH  
 THE OPTIMUM

.....PICTURE OF  
 THE MULTIPLE PEAKS.

(000039) 04  
 18 = 1

(000040) NOTE 05  
 BEGIN OF LOOP  
 5 1 = 1 100000

01.05--> 06

(000041) IFILLJX = EQ  
 IFILLJX

FALSE

(000042) 07  
 ISPELJ21 = 1

18 = 18 + 1

08  
 17 = 1

09  
 END OF LOOP?

YES

(000047) 10  
 MAX = NUMBER - 1

LCCP = 1

READ FROM THE  
 NORMALIZED FILE 11

CL-210--> 11

READ FROM DEV  
 INTERNAL FORMAT  
 INTO THE LIST

(000052) NOTE 12  
 LIST = 1X  
 AMPLTD, AVERAGE

(000052) 13  
 READ FROM DEV

INTERNAL FORMAT  
 INTO THE LIST

(000053) NOTE 14  
 LIST = 1X, 15-18

(000054) 15  
 READ FROM DEV

INTERNAL FORMAT  
 INTO THE LIST

(000054) NOTE 16  
 LIST = 1X, 15-18

(000055) 17  
 1X = EQ. IF TRUE

FALSE

WRITE TO UPDATED  
 TONE-0 FILE 23

(000056) 18  
 WRITE TO DEV

INTERNAL FORMAT  
 INTO THE LIST

(000059) NOTE 19  
 LIST = 1X, 15-18

(000060) 20  
 WRITE TO DEV

INTERNAL FORMAT  
 INTO THE LIST

(000060) NOTE 21  
 LIST = 1X, 15-18

1-11

1-11

1-11

1-11

1-11

1-11

1-11

1-11

1-11

1-11

[illegible]



CHART TITLE - SUBROUTINE CBASE(DIFF,VALUE,JX)

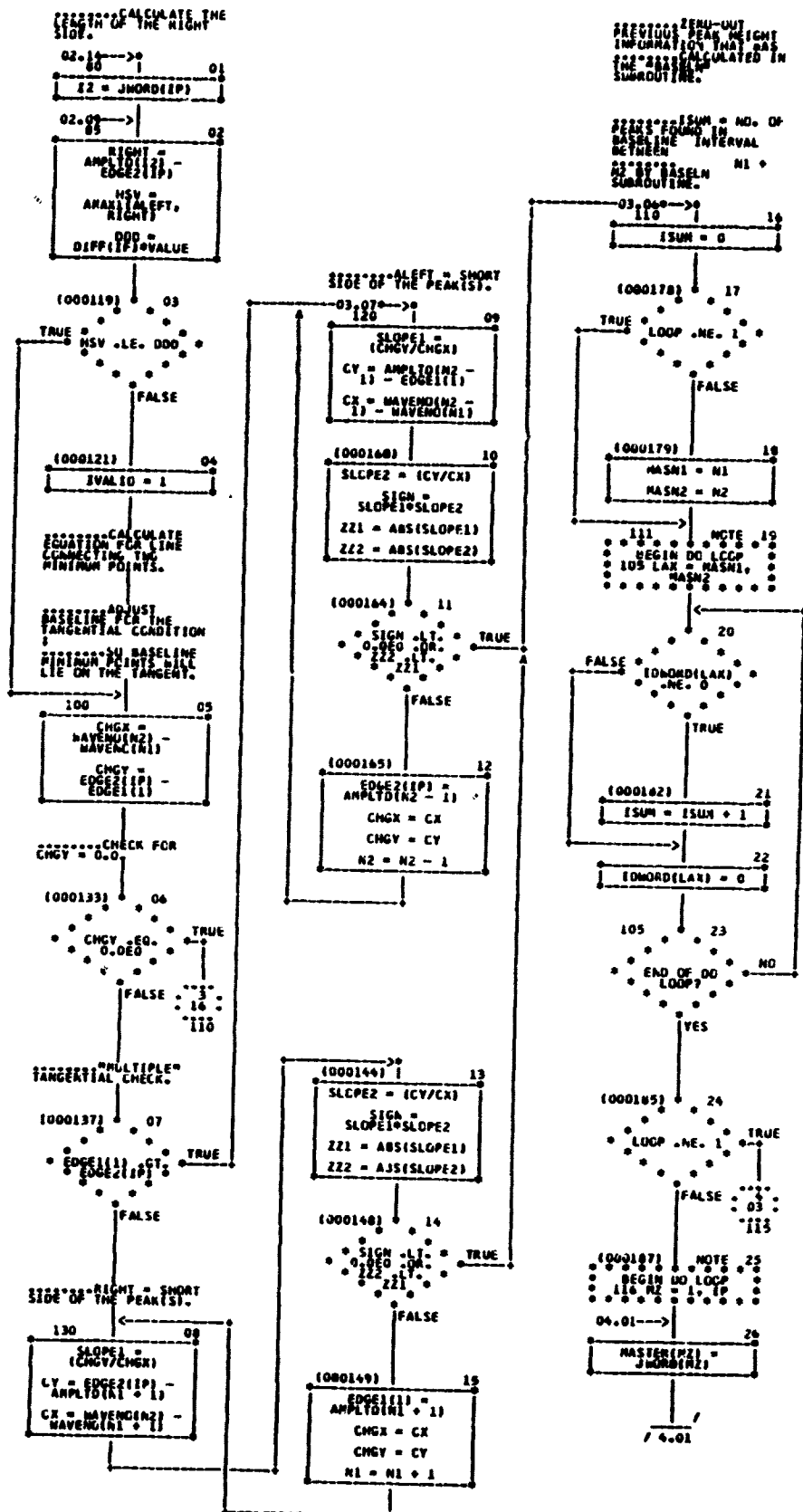


CHART TITLE - SUBROUTINE C0BASE(DIFF,VALUE,JN)


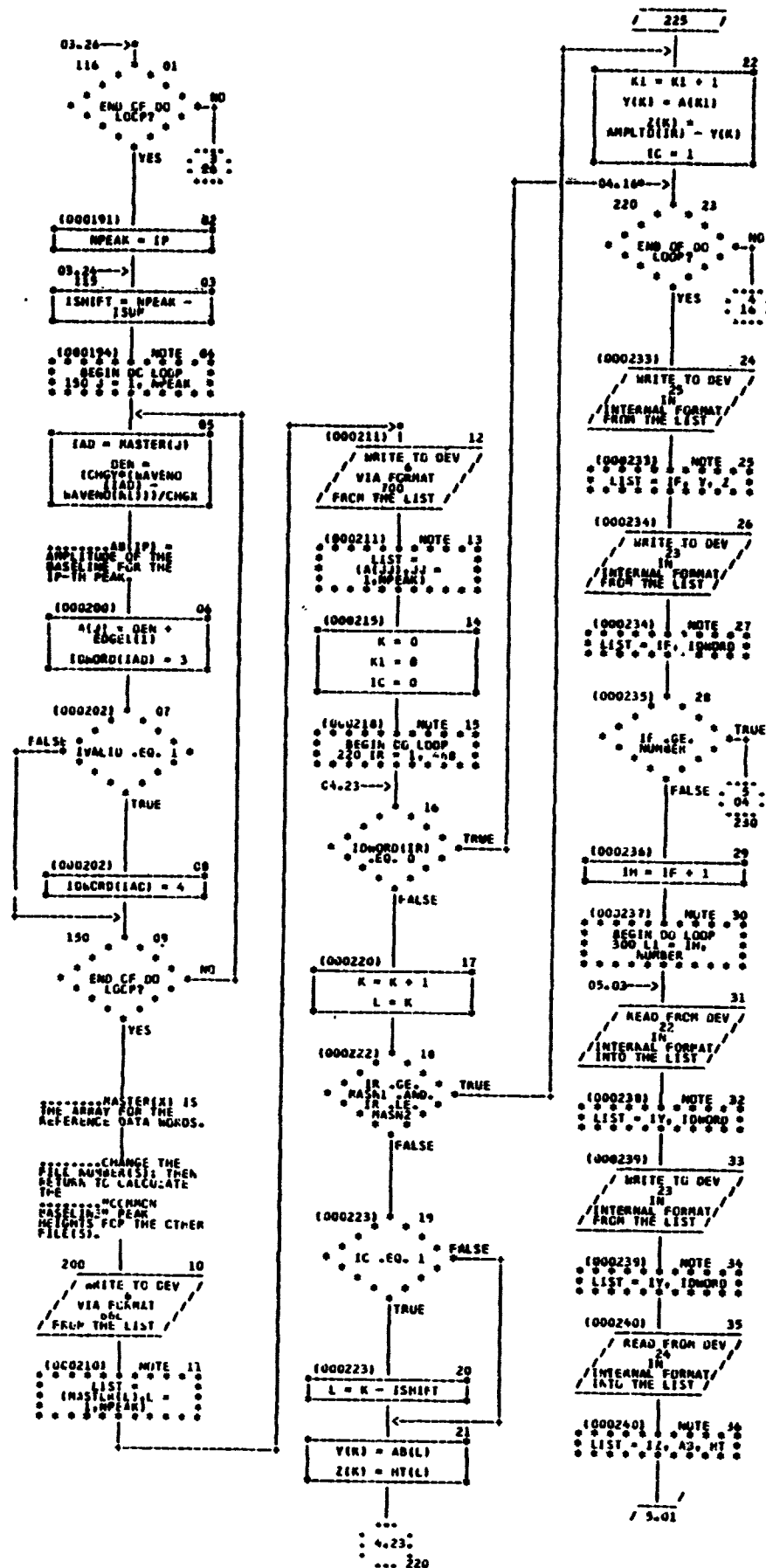


CHART TITLE - SUBROUTINE CDBASE(DIFF,VALUE,JX)

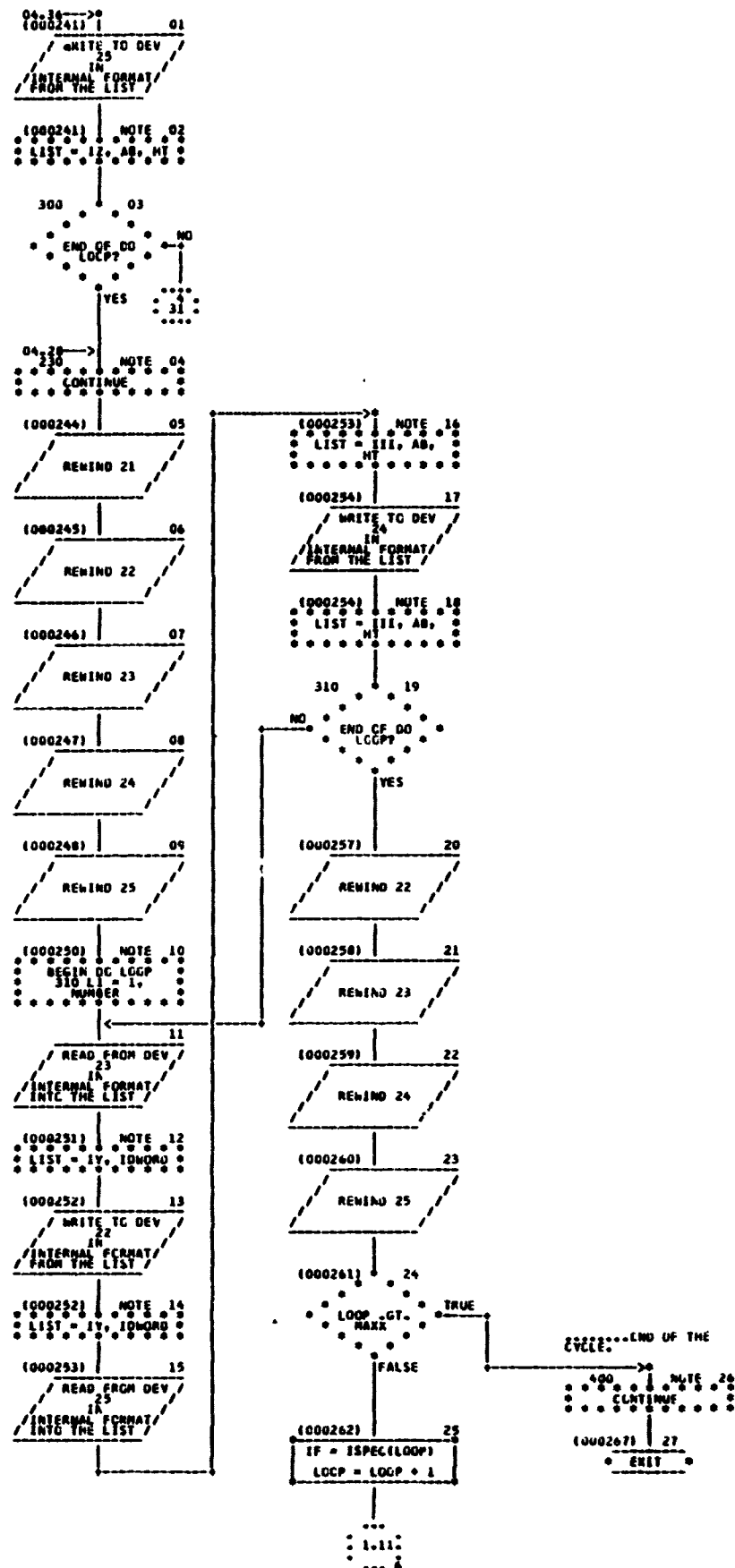

12/05/79

AUTOFLOW CHART SET - E490 - SUBROUTINE COBASE PAGE 06

CHART TITLE - NON-PROCEDURAL STATEMENTS

```
COMMON /SPECTR/ IFILE(30),NPT25(30)
COMMON /TWINPK/ ALIMIT(6),ZLIMIT(6),MFILE(6),NORM(G,NUMBER
DIMENSION AMPLTD(468),WAVENO(468),IDWORD(468),AB(50),MT(50,
DIMENSION EDGE1(25),EDGE2(25),ISPEC(29),JMKRD(25),MASTER(25)
DIMENSION A(25),Y(50),Z(50),DIFF(30)
FORMAT (10(5X,F6.2))
FCRMAT(1140,5X,513)
```

700  
666

## CHART TITLE - SUBROUTINE PKTABL(IRFG,ITOTAL,IKOUNT)

177

PG.BX	NAME	PG.BX	NAME	PG.BX	NAME	PG.BX	NAME		
1.01	KTABL	2.09	5	2.16	35	3.06	65	2.05	230
1.01	1	3.09	10	2.18	45	3.29	70	3.24	240
1.20	2	1.21	15	2.24	50	1.25	110	3.14	250
2.01	3	3.32	20	2.29	60	3.01	210	2.12	500
1.27	4	1.23	25						

CHART TITLE - SUBROUTINE PKTABL (IREG, ITOTAL, INOUNT)

PKTABL /

1979-----SEPTEMBER 24,  
1979. (UPDATED AT  
AF/RPL 1 OCT. 24,  
1979).

THIS SUBROUTINE  
CROSS-REFERENCES ALL  
THE PEAK HEIGHTS  
WITH THEIR RESPECTIVE  
DATA WORD LOCATIONS.



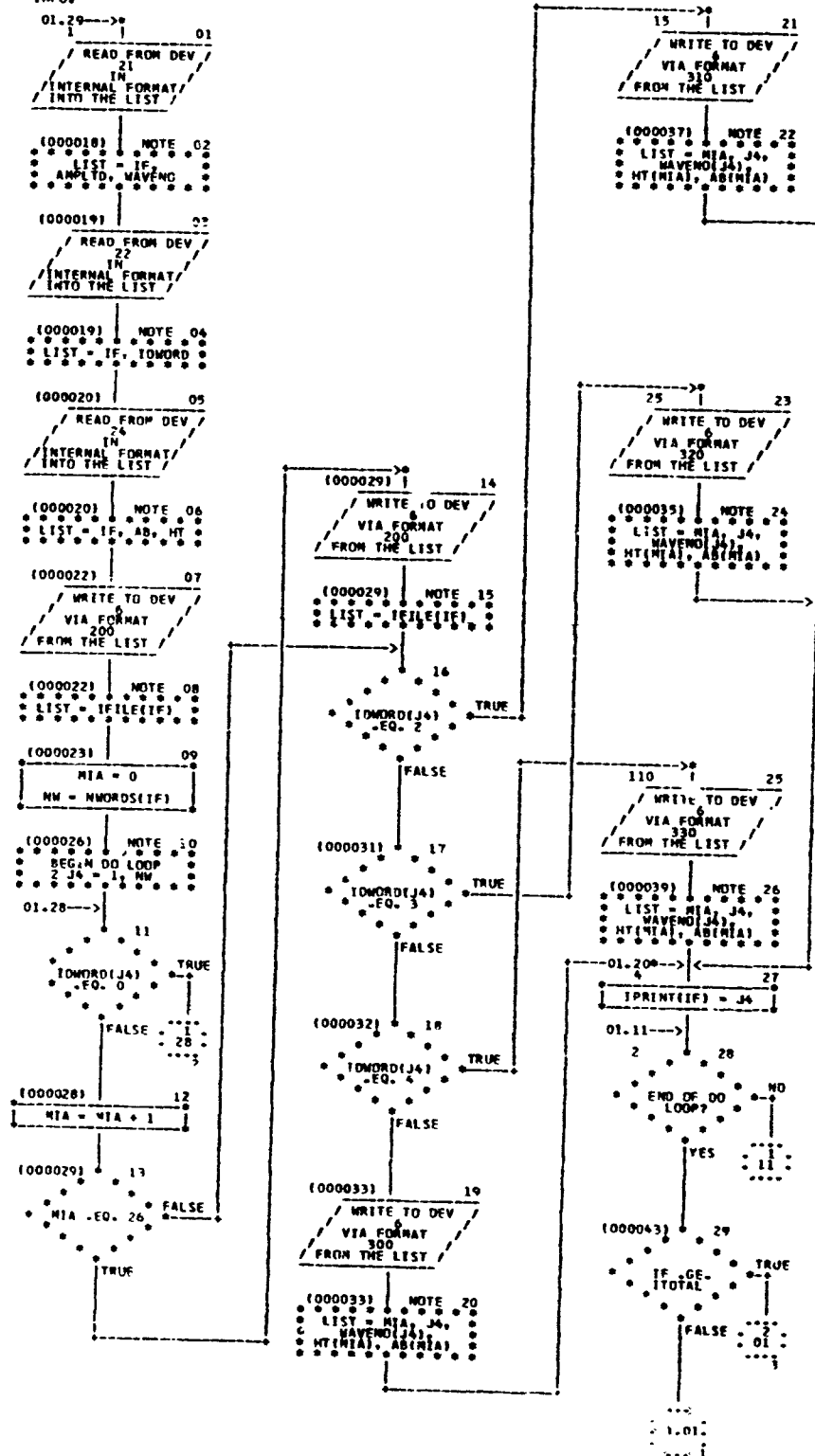
WRITE PEAK HEIGHT  
INFO.

CHART TITLE - SUBROUTINE PKTABL(IREG,ITOTAL,IKOUNT)



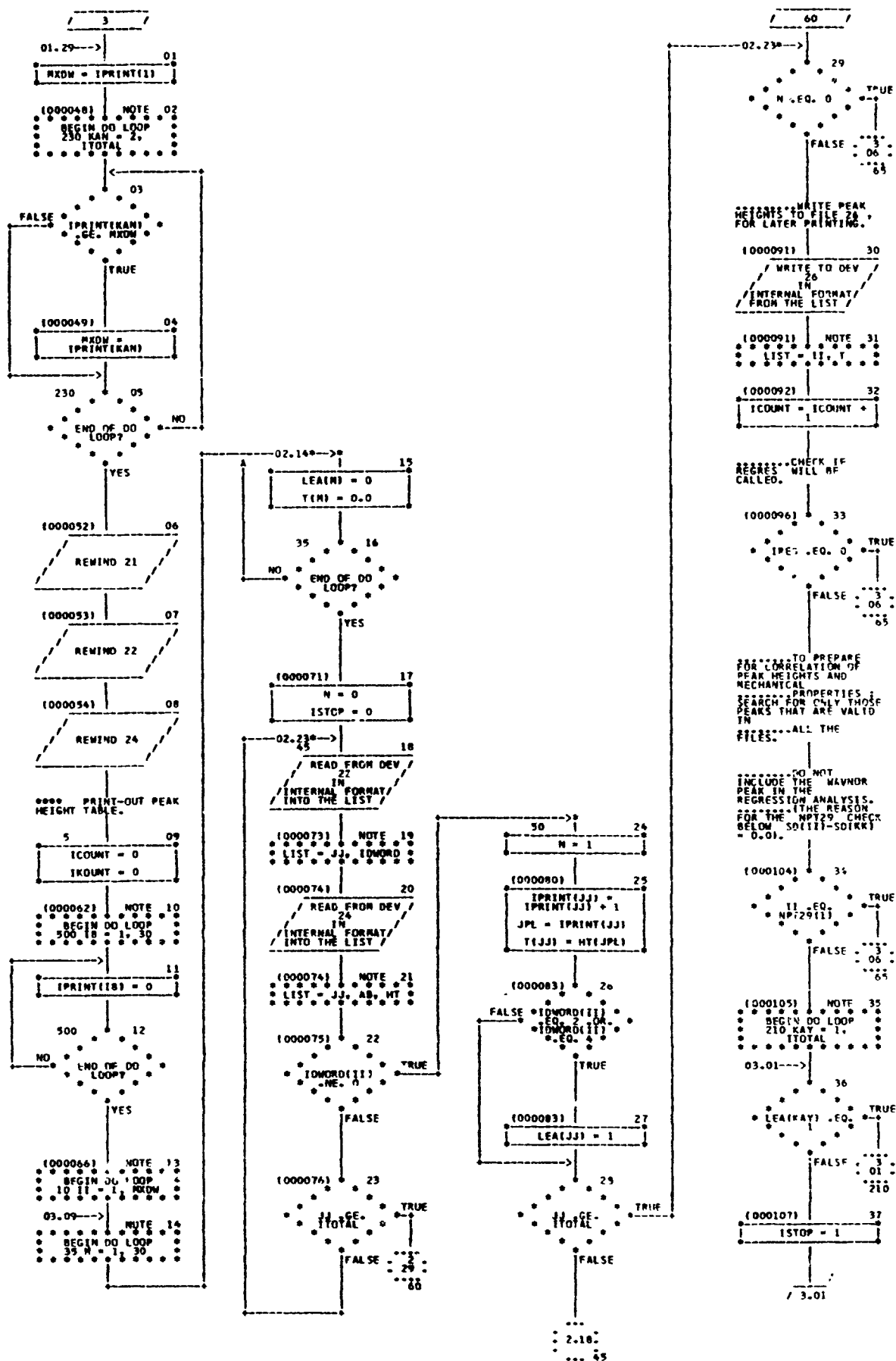

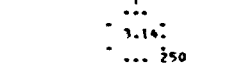

## CHART TITLE - NON-PROCEDURAL STATEMENTS

```

COMMON /SPECTR/ IFILE(30),MPT2(30)
COMMON /POINTS/ NWORDS(30)
DIMENSION AMPLD(468),WAVEMO(468),IDWGRD(468),AB(50),HT(50)
DIMENSION T(30),IPRINT(30),LFA(30)
FORMAT(1H1,4X,"FOURIER TRANSFORM INFRARED SPECTROSCOPY --- NORMAL
IZED PEAK HEIGHT INFORMATION FILE NUMBER",I3,/,5X
,*** - DENOTES A VALID PEAK."/5X,"CB. - DENOTES A PEAK MEASURED F
ROM A COMMON BASELINE."/5X
,PEAK NUMBER DATA WORD WAVE NUMBER PEAK HEIGHT BA
SELINE AMPLITUDE//5X,"-----")
-----
300 FORMAT(1H0,8X,I2,I3X,I3,8X,F10.4,5X,F12.7,5X,F12.7)
310 FORMAT(1H0,4X,*** "I2,I3X,I3,8X,F10.4,5X,F12.7,5X,F12.7)
320 FORMAT(1H0,8X,I2," CB,"9X,I3,8X,F10.4,5X,F12.7,5X,F12.7)
330 FORMAT(1H0,4X,*** "I2," CB,"9X,I3,8X,F10.4,5X,F12.7,5X,F12.7)
400 FORMAT(1H1,7X,10(F7.3,4X))
410 FORMAT(1H1,7X,"FTTS ***** PEAK HEIGHT TABL
ATION FOR ALL THE FILES *****
//8X,"DATA",8X,
10("FILE",7X),/8X,"WORD",9X,10(I2,9X)/8X,"-----")
-----
420 FORMAT(1H0,4X,"NORMAL END OF PEAK HEIGHT TABLE "/5X,"TOTAL NUMBER
OF FILES LISTED = ",I3,")

```

12/11/79      TABLE OF CONTENTS AND REFERENCES      AUTOFLOW CHART SET -  
 CARD ID    PAGE/BOX    NAME      REFERENCES (SOURCE SEQUENCE NO. AND PAGE/BOX)

FORTRAN MODULE    F400 - SUBROUTINE PLWAVE

CHART TITLE - SUBROUTINE PLWAVE (ITAPC, IF, NPTS, WN, ...)

(000004)	1.01	PLWAVE	
(000032)	1.15	(000030)	1.13
(000039)	1.17	20	(000032) 1.15
(000042)	1.19		(000044) 1.23
(000043)	1.21		(000042) 1.19
(000044)	1.23	10	
(000044)	1.23		(000043) 1.21
(000045)	1.24	30	(000035) 1.16
(000057)	1.29	40	(000053) 1.26
(000058)	1.30	50	(000056) 1.28

CHART TITLE - NON-PROCEDURAL STATEMENTS

12/11/79	PROCEDURAL STATEMENT LAPEL INDEX	AUTOFLOW CHART SET -	5450 - SUPERCHARGE PLWAVE	PAGE 1
PG.8X NAME	PG.3X NAME	PG.8X NAME	PG.8X NAME	
1.01 LWAVE	1.17 20	1.24 10	1.23 40	
1.23 10			1.30 50	



CHART TITLE - SUBROUTINE PLWAVE(ITAPE,IF,NPTS,WN,AP)

/ PLWAVE /

.....PLOTS AN  
INDIVIDUAL FILE  
SPECTRUM FROM THE  
FTIS TAPE.

.....THICKOL /  
DUPN C. SMITH  
AUGUST 4, 1978.

.....UPDATED  
AUGUST 21, 1979.

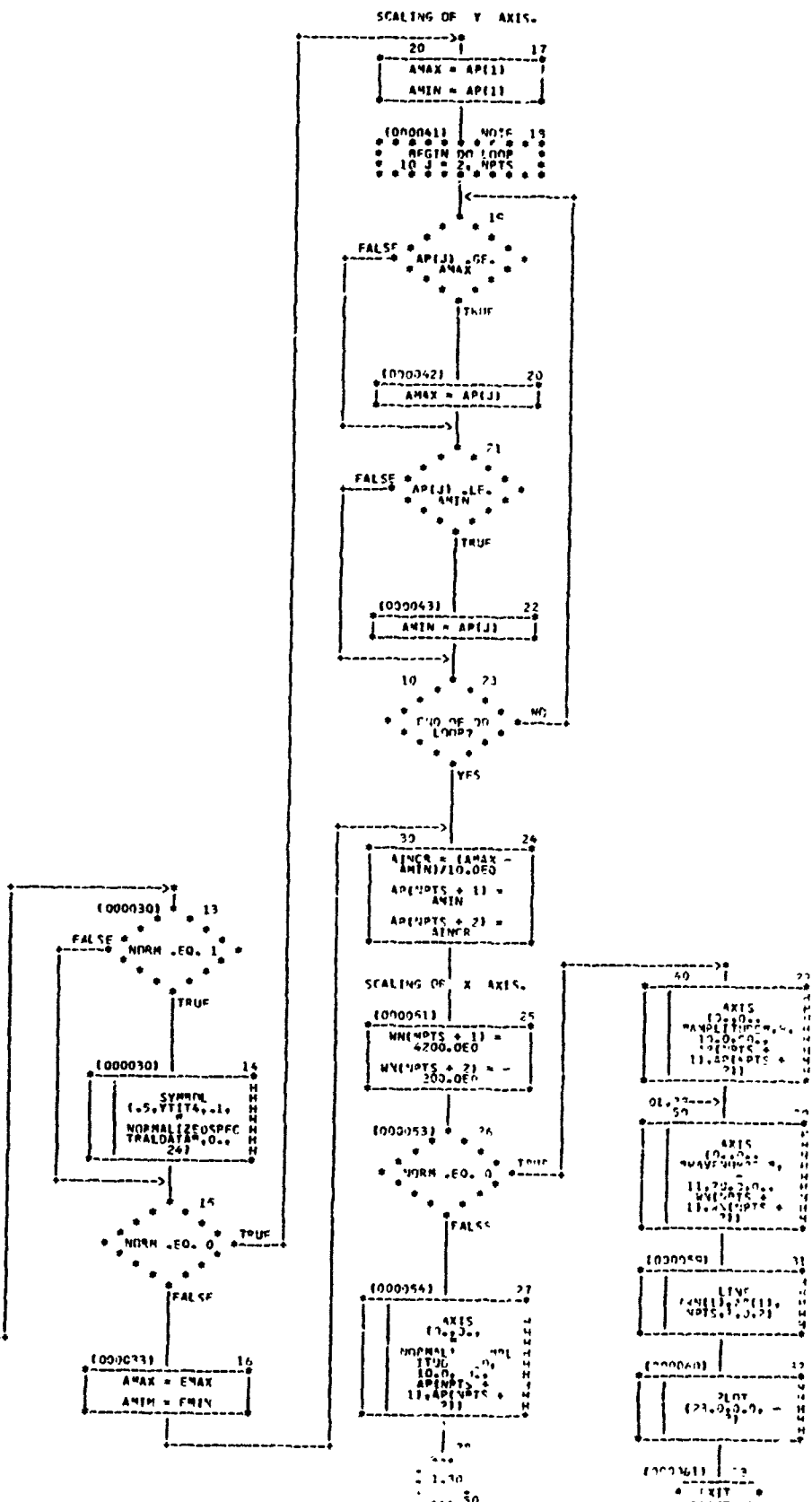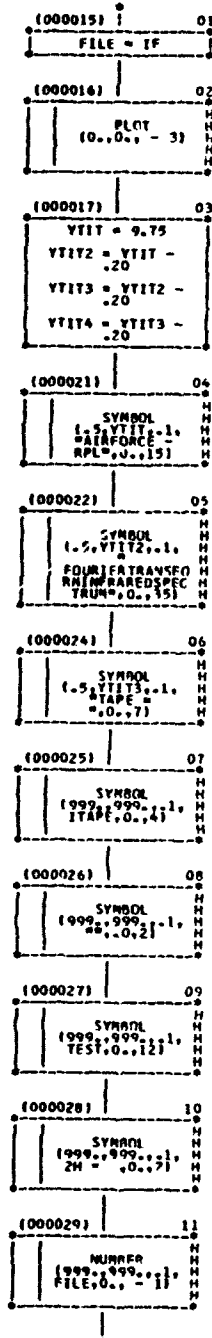

PAGE 02

5490 - SURROUTINE PLWAVE

AUTOFLOW CHART SET -

12/11/79

CHART TITLE - NON-PROCEDURAL STATEMENTS

```
COMMON /OTINF / MAX.CHIN,NDPM  
DIMENSION AP(470),MN(470),TEST(.)  
DATA TEST/'FILE'," NUM",BER " /
```

[illegible]

12/11/77  
TABLE OF CONTENTS AND REFERENCES

AUTOFLOW CHART SET -  
REFERENCES (SOURCE SEQUENCE NO. AND PAGE/NOX)

PAGE 2

CARD ID	PAGE/NOX	NAME	REFERENCES (SOURCE SEQUENCE NO. AND PAGE/NOX)
(000331)	6.15	420	
(000332)	6.16	430	(000332) 6.14
(000341)	6.17	440	(000331) 6.13
(000344)	6.18	450	(000347) 6.20
(000347)	6.20	460	(000341) 6.18
(000348)	6.21	470	(000328) 6.10
(000351)	6.24	490	(000329) 6.11 (000331) 6.13 (000337) 6.16 (000346) 6.19
(000351)	6.24		(000351) 6.25
(000364)	7.01	481	(000361) 6.27
(000365)	7.03	412	(000363) 6.29
(000366)	7.04	490	
(000369)	8.01	511	(000365) 7.03
(000370)	8.02	510	
(000372)	8.06	520	(000369) 7.07 (000369) 8.01
(000377)	8.15	550	(000376) 8.14 (000377) 8.16
(000377)	8.16	570	
(000384)	8.17	642	(000387) 8.15
(000400)	8.27	641	(000442) 8.21
(000381)	9.01	540	(000377) 8.16
(000391)	9.03	560	(000397) 8.15
(000396)	9.05	570	(000396) 9.02
(000403)	9.07		(000411) 9.12
(000404)	9.08	587	
(000405)	9.09		(000410) 9.11
(000409)	9.10	590	(000405) 9.09
(000410)	9.11	600	(000405) 9.09
(000411)	9.12	610	(000403) 9.07
(000413)	9.14		(000410) 9.16
(000417)	9.15	620	(000413) 9.14
(000419)	9.16	630	(000413) 9.14
(000420)	9.18		(000425) 9.20
(000424)	9.19	640	(000420) 9.18
(000425)	9.20	650	(000427) 9.18
(000451)	10.01	665	(000450) 8.27
(000457)	10.02	688	(000449) 8.26 (000452) 8.29
(000466)	10.07		(000487) 10.21
(000469)	10.10	645	(000466) 10.07
(000470)	10.12	686	(000468) 10.09
(000472)	10.14		(000473) 10.15
(000473)	10.15	670	
(000480)	10.20	680	
(000495)	10.31	1100	
(000512)	11.01	4999	
(000514)	11.03		(000518) 11.05
(000518)	11.05	1310	(000514) 11.03
(000519)	11.06	1230	(000494) 10.30
(000506)	11.12	1310	(000500) 11.09 (000523) 11.22
(000510)	11.17	1250	(000521) 11.09
(000592)	11.20	1500	(000507) 11.13
(000522)	11.22	1350	(000508) 11.14
(000525)	11.23	1400	(000511) 11.19

CHART TITLE - HON-PROCEDURAL STATEMENTS



PG.8X	NAME	PG.8X	NAME	PG.8X	NAME	PG.8X	NAME
1.01	EGRES	3.03	180	5.22	330	6.24	480
1.03	4	4.01	181	5.37	335	7.01	481
1.06	5	3.13	182	5.40	337	7.03	482
1.09	10	4.04	190	5.43	340	7.04	490
2.02	11	4.05	200	5.27	350	8.01	500
2.01	12	4.08	210	5.31	360	8.02	510
2.00	20	4.12	220	5.01	364	8.06	520
2.14	70	4.13	230	6.04	369	8.16	530
2.24	75	4.17	235	6.05	370	9.01	540
2.25	76	4.22	260	6.07	380	8.15	550
2.27	78	4.25	261	6.11	390	9.03	560
2.31	120	4.24	252	6.12	400	9.05	570
2.35	130	4.35	263	6.14	410	9.08	580
2.36	140	5.01	265	6.15	420	9.10	590
2.37	150	5.05	270	6.16	430	9.11	600
2.38	160	5.09	280	6.17	440	9.12	610
2.39	170	5.11	290	6.19	450	9.15	620
3.01	175	5.17	300	6.20	460	9.16	630
3.02	178	5.16	310	6.21	470	9.19	640
3.05	179	5.21	320				

**/ REGRES /**

190






[illegible]


CHART TITLE - SIMROUTINE PROGRAMS



[illegible]



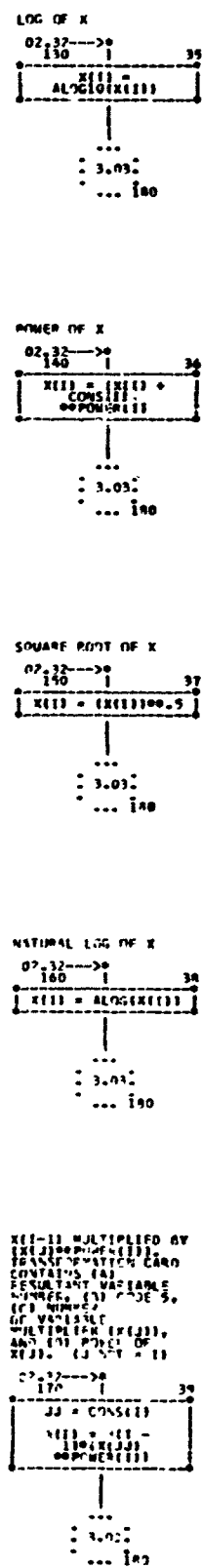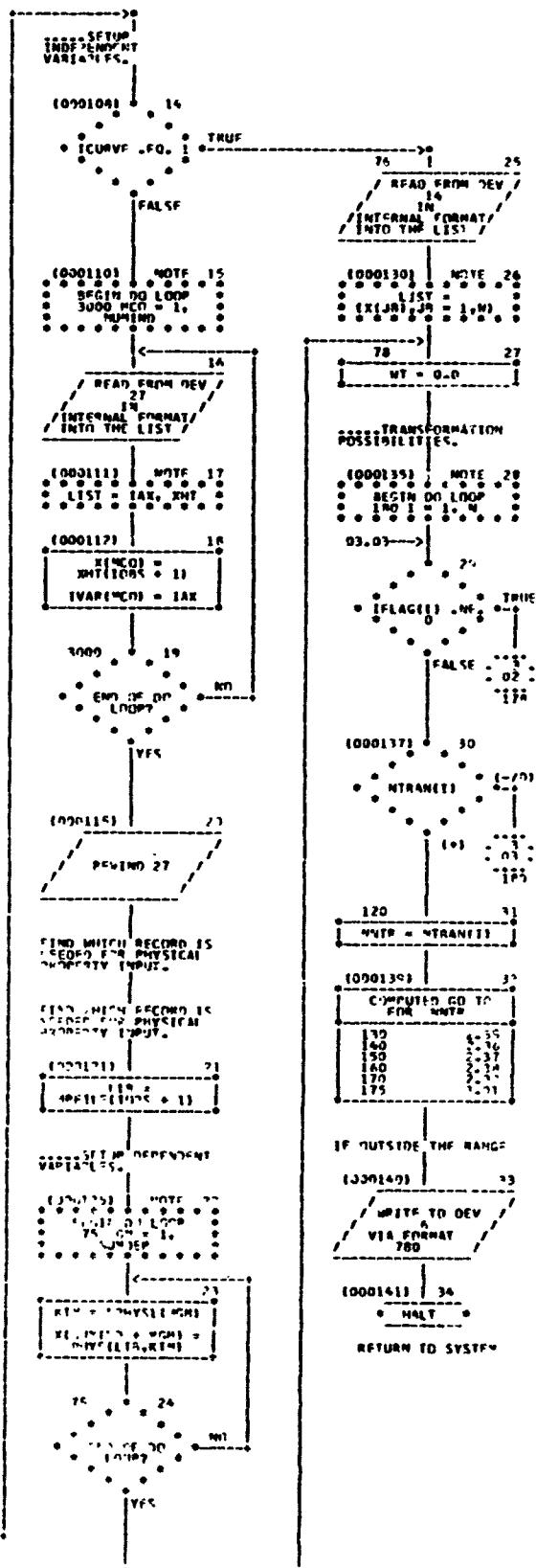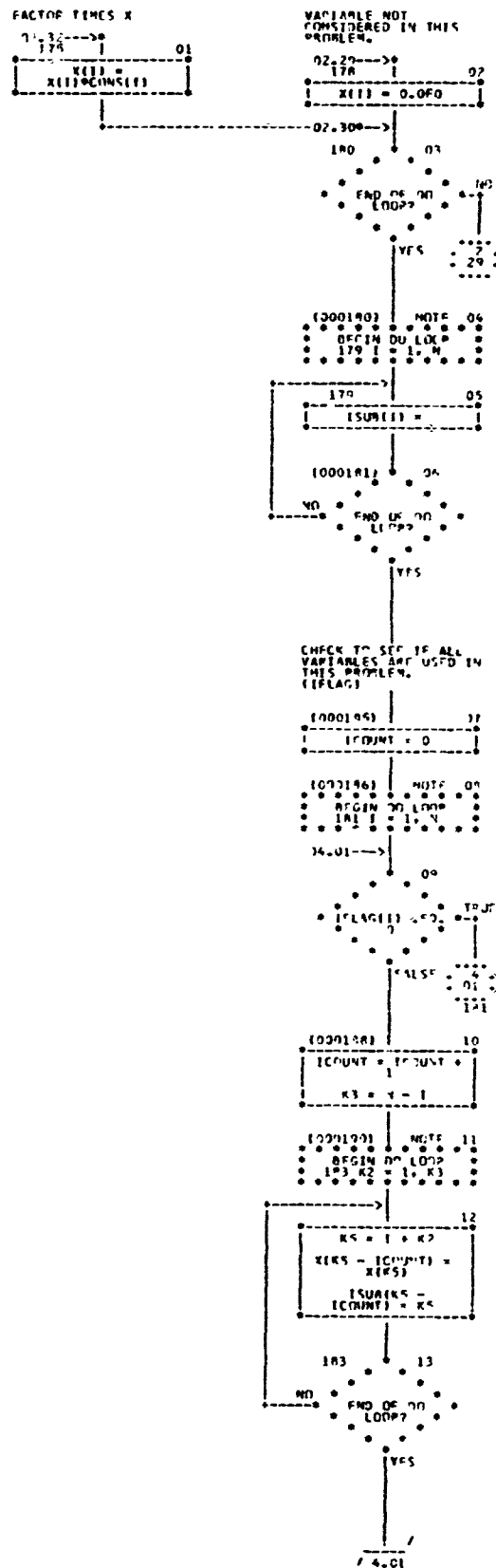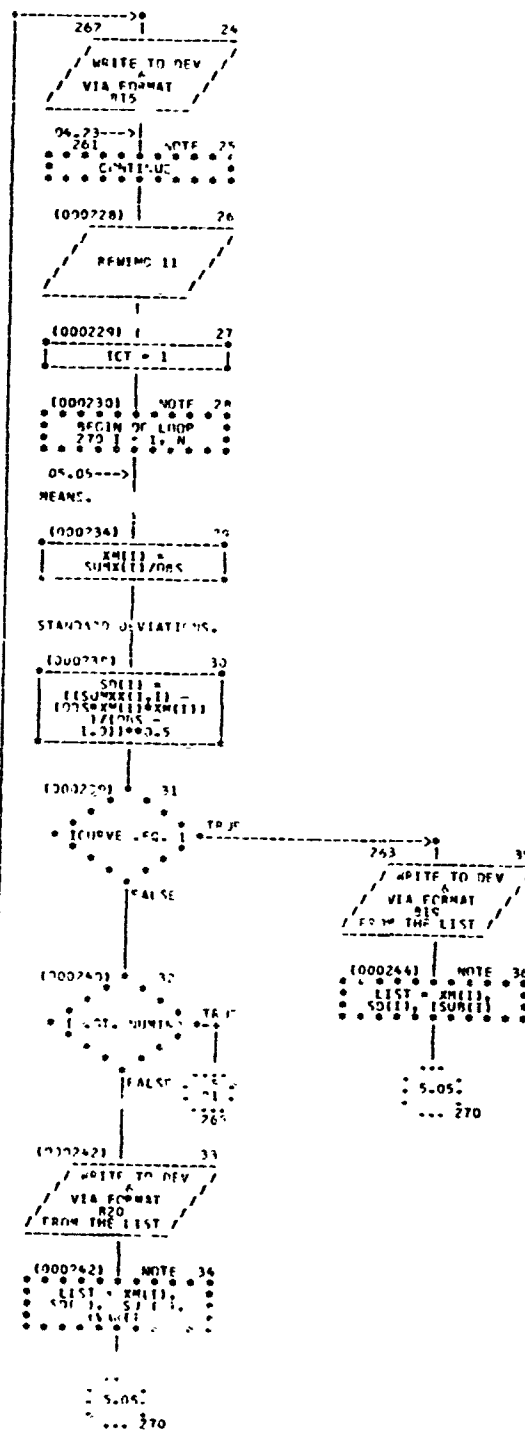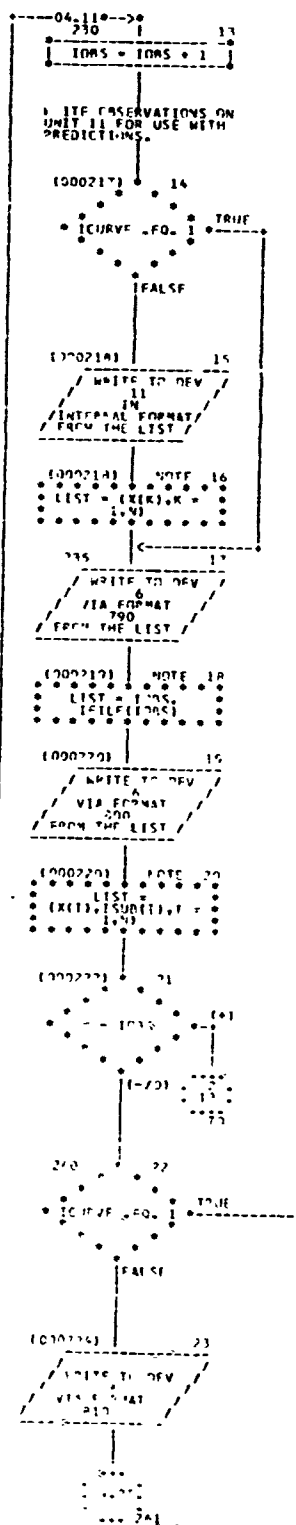

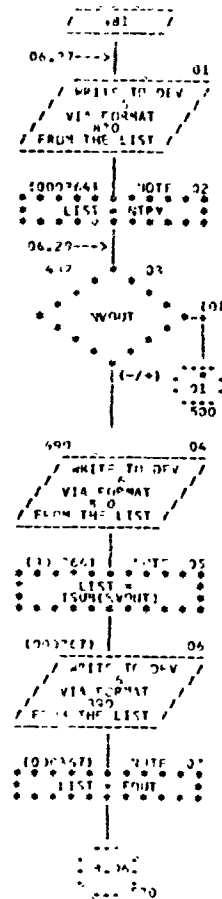




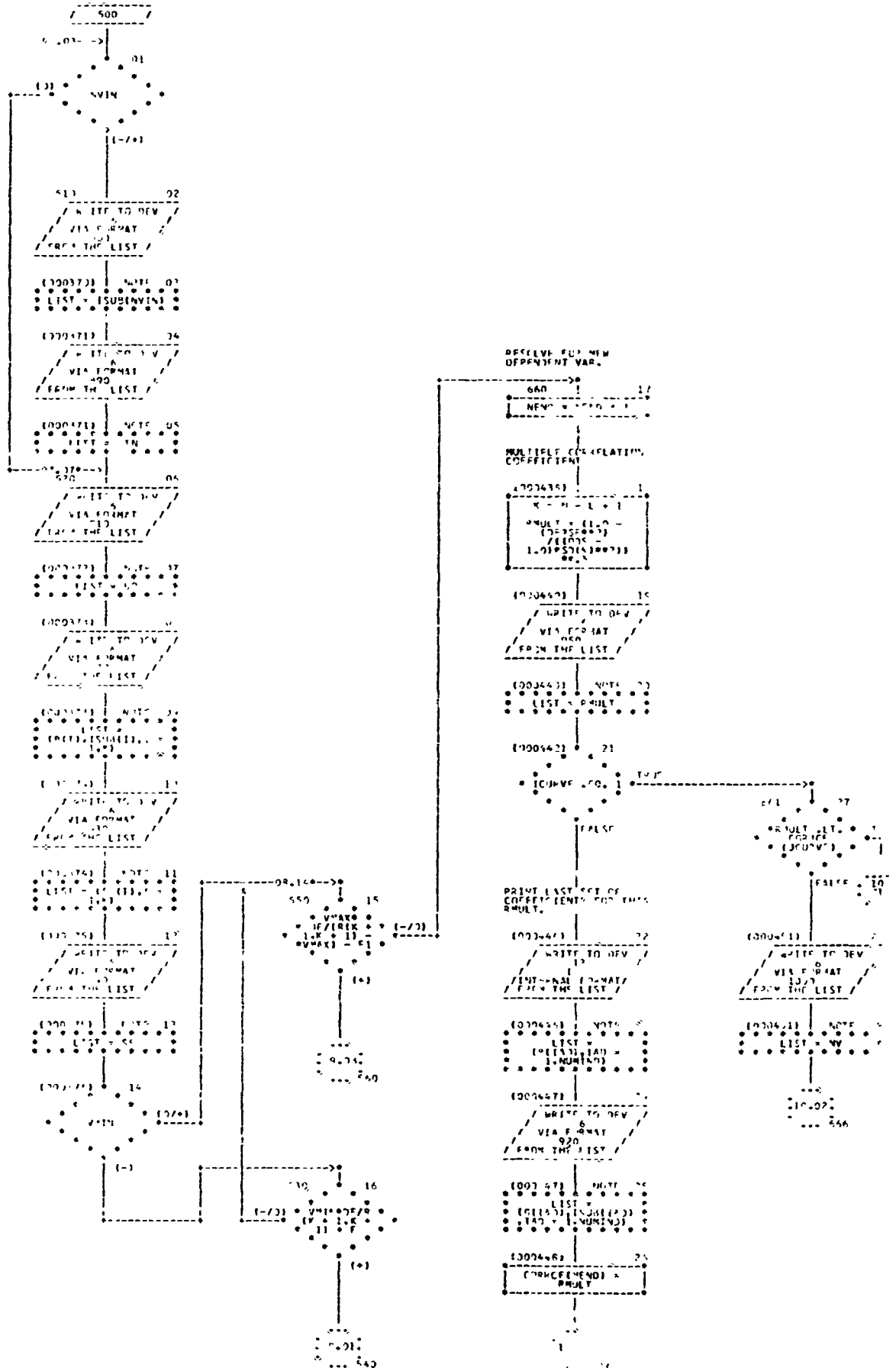


CHART TITLE - SUMMUTINE REGRES

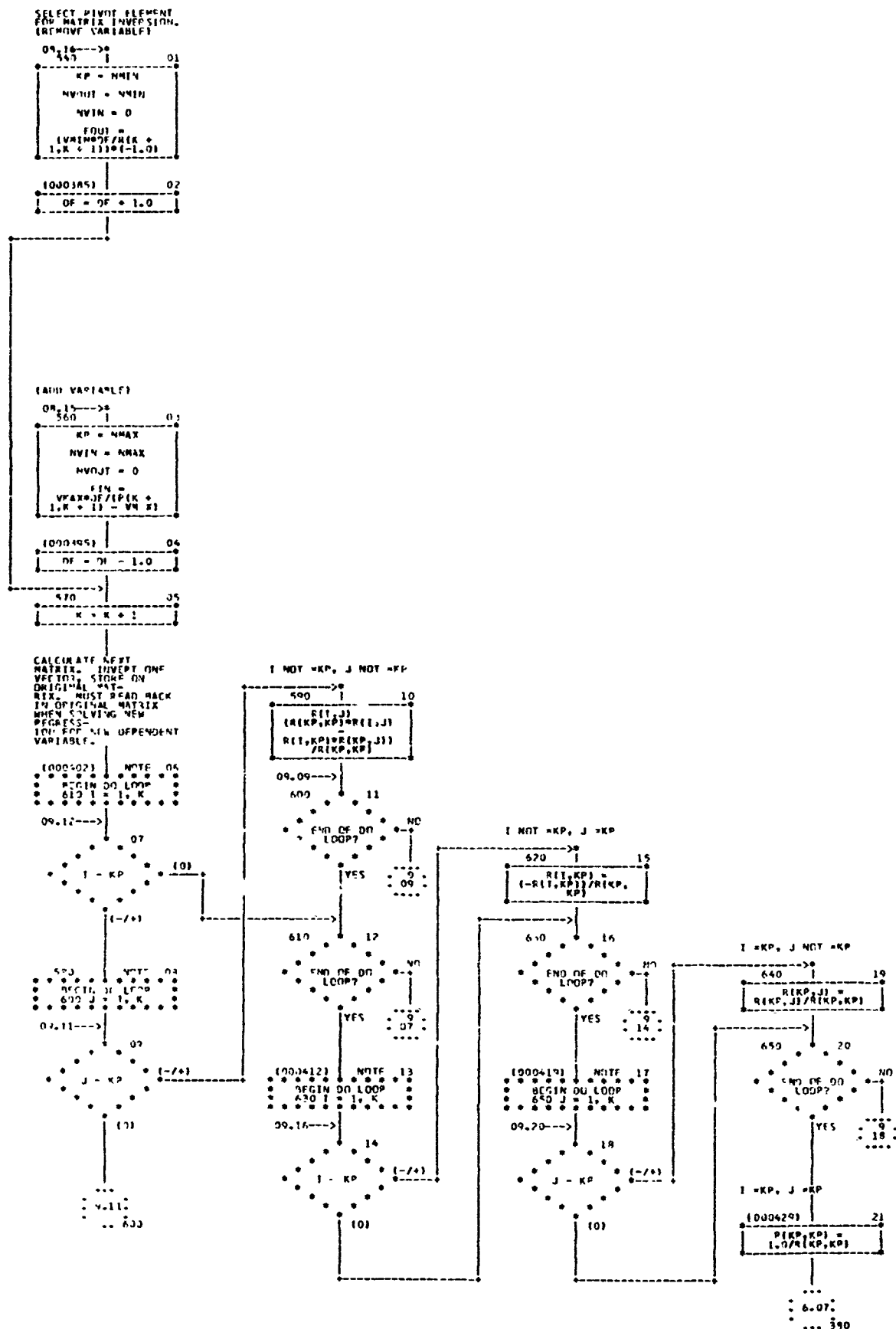

## CHART TITLE - SUMMATION REGRES


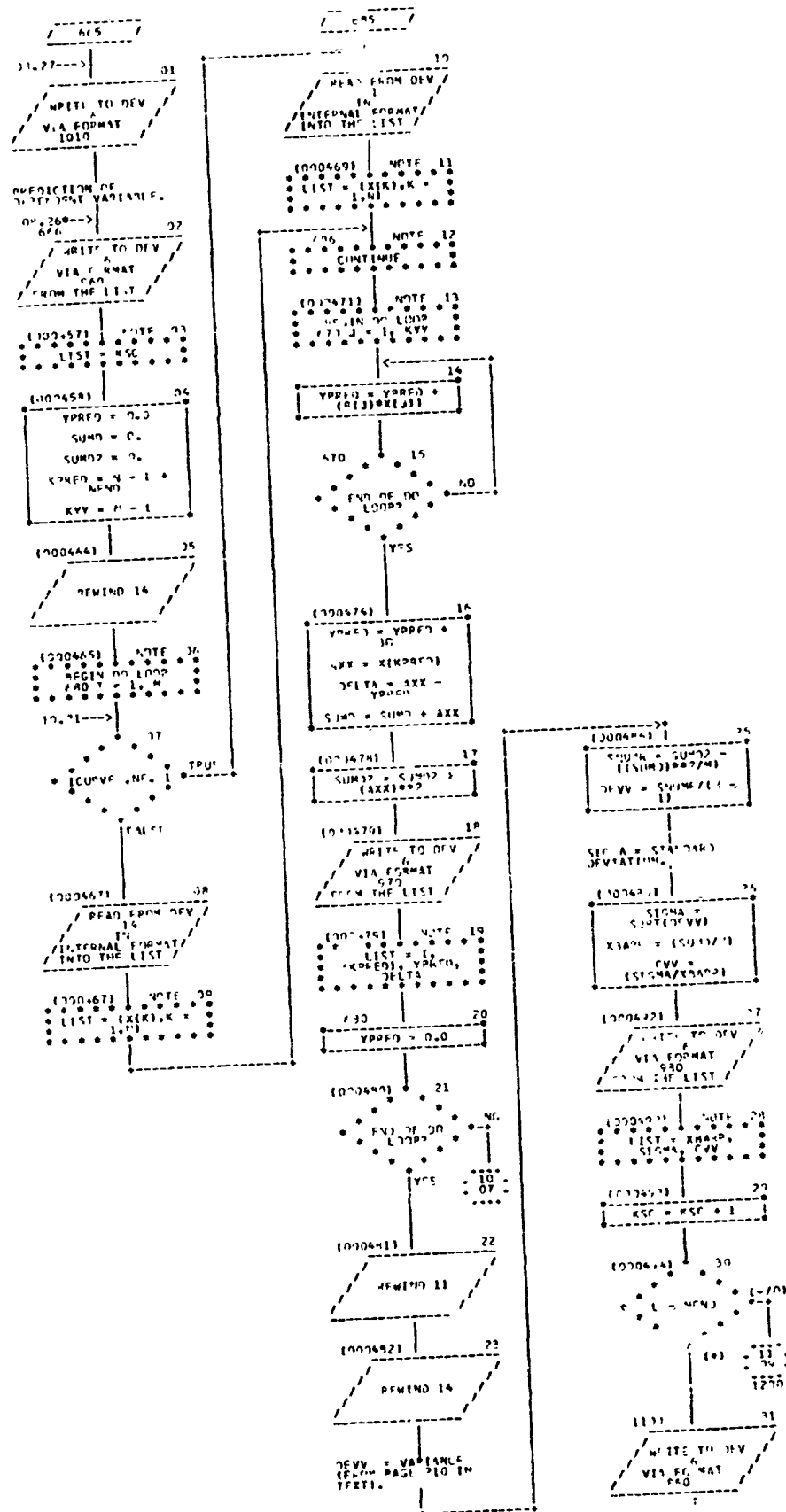


CHART TITLE - SUBROUTINE REGRES



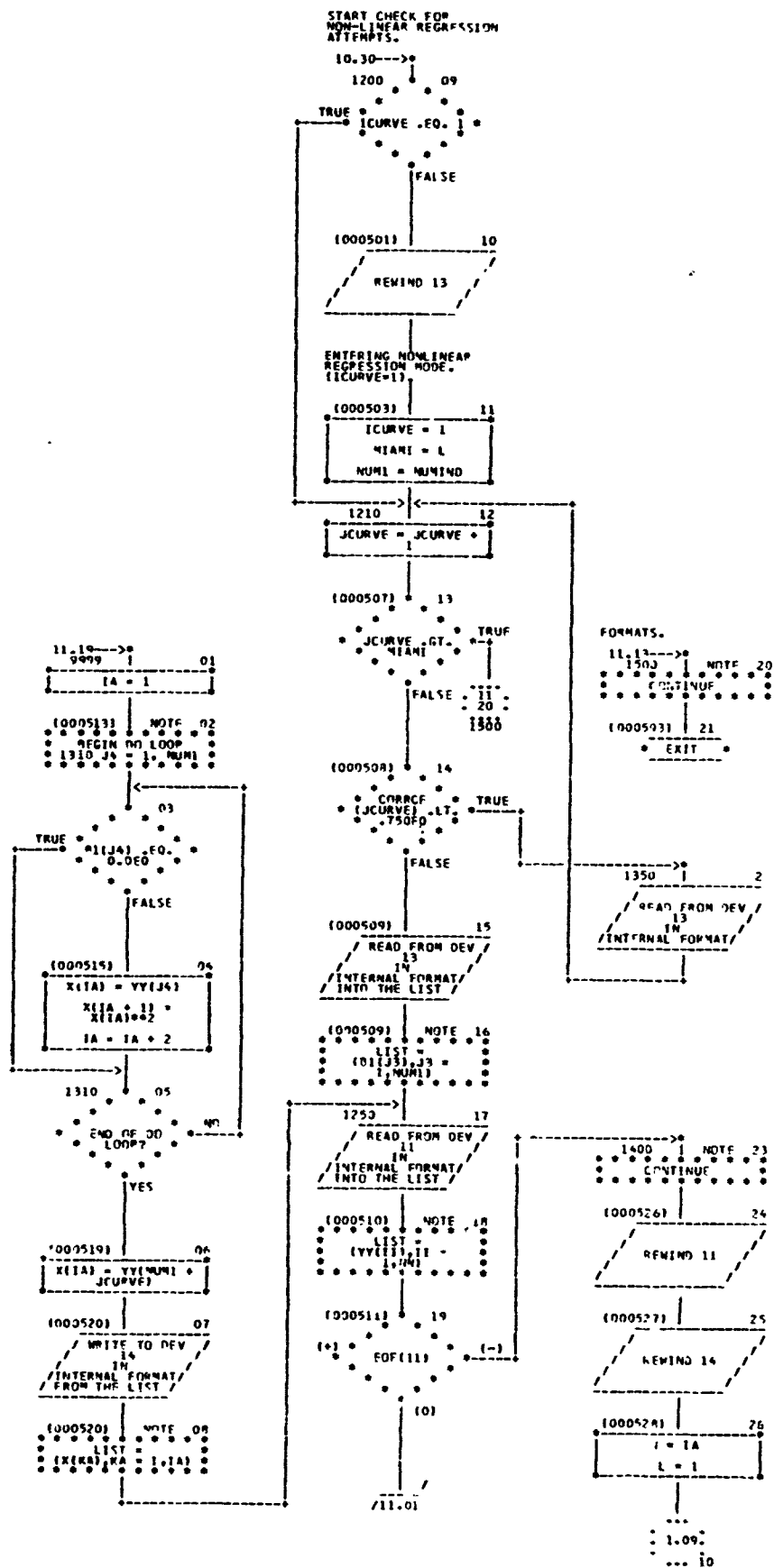

12/12/79 TABLE OF CONTENTS AND REFERENCES AUTOFLOW CHART SET -  
 CARD 10 PAGE/PCF NAME REFERENCES (SOURCE SEQUENCE NO. AND PAGE/BOX)

FORTRAN MODULE E41C - MAIN PROGRAM

CHART TITLE - INSTRUCTIONS

CHART TITLE - PROCEDURES

10000871	2.06	650	(000043)	2.04
(000046)	2.09	70	(000044)	2.05
(000049)	3.01	30	(000045)	2.05
(000062)	3.14	10	(000074)	3.24
(000064)	3.17	999	(000085)	3.44
(000092)	3.25	200	(000063)	3.16
(000095)	3.29		(000097)	3.31
(000097)	3.31	240		
(000129)	3.32	250	(000090)	2.09
(000075)	3.35	40	(000073)	3.24

CHART TITLE - NON-PROCEDURAL STATEMENTS

12/13/79

REFLECTAL DIFFERENCE LEVEL INDEX

TOFLOW CHART SET -

E410 - MAIN PROGRAM

PAGE 1

PG.BX	NAME	PG.BX	NAME	PG.BX	NAME	PG.BX	NAME
3.11	10	3.21	30	3.25	200	3.32	250
2.09	20	3.35	40	3.31	240	2.06	650

12/13/79

CHART TITLE - INTRODUCTORY COMMENTS

AUTOFLOW CHART SET -

E410 - MAIN PROGRAM

PAGE 01

\*\*\*\*\*

E410 / FOURIER TRANSFORM IMPAIRED SPECTROSCOPY.

PHYSICAL PROPERTIES MASTER TAPE GENERATOR.

\*\*\*\*\*

THICKOL CORPORATION / HUNTSVILLE , ALABAMA 35807

PRINCIPAL INVESTIGATOR W. W. SCHWARZ

TELEPHONE (205) - 992 - 8300

SCIENTIFIC PROGRAMMER D. C. SMITH

TELEPHONE (205) - 992 - 8215

\*\*\*\*\*

SEPTEMBER , 1979.

FORMAT IV - C LEVEL 21 LANGUAGE


CDC 4600 -- (AEDOL)

ORDERED IN FULFILLMENT OF CONTRACT F04411 - 79 - 1 - 0007

\*\*\*\*\*



## CHART TITLE - PROCEDURES

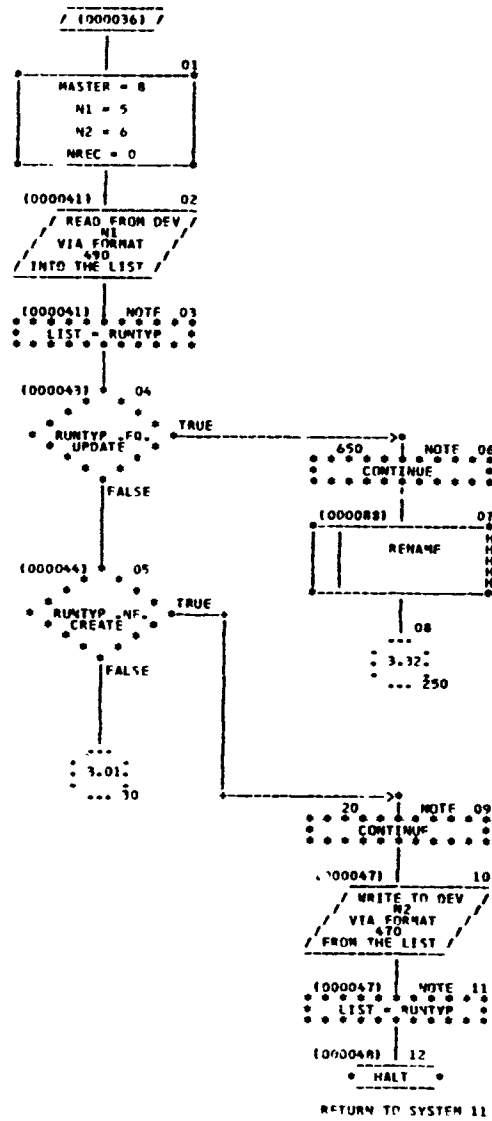





FORTRAN MODULE E410 - SUBROUTINE RENAME

CHART TITLE - SUBROUTINE RENAME

(000004)	1.01	RNAME								
(000023)	1.06	12	(000028)	1.11						
(000026)	1.09	888								
(000027)	1.11		(000026)	1.09						
(000049)	2.03		(000048)	2.01						
(000029)	2.04	14	(000075)	1.08						
(000042)	2.07	5	(000047)	2.31	(000070)	2.35				
(000044)	2.10	777								
(000046)	2.13		(000045)	2.11						
(000050)	2.15	70	(000046)	2.13						
(000057)	2.21	60	(000047)	2.14						
(000055)	2.24	55	(000050)	2.15						
(000046)	2.29	25	(000049)	2.03	(000051)	2.16	(000052)	2.17	(000060)	2.22
(000088)	3.01	100	(000051)	3.21						
(000099)	3.11		(000098)	3.09						
(000072)	3.12	30	(000043)	2.09						
(000075)	3.14	15	(000099)	3.11	(000172)	4.74	(000171)	5.26	(000173)	5.28
(000077)	3.17	666								
(000225)	3.22	600	(000077)	3.17						
(000227)	3.25		(000229)	3.27						
(000229)	3.27	640								
(000269)	3.28	650	(000076)	3.16	(000221)	6.34				
(000182)	3.32	450	(000079)	3.19	(000080)	3.20				
(000126)	4.01	300	(000083)	4.06						
(000130)	4.04		(000154)	5.16						
(000103)	4.07	200	(000082)	4.05						
(000105)	4.09		(000106)	4.11						
(000106)	4.11	205								
(000106)	4.11		(000105)	4.09						
(000109)	4.14		(000119)	4.21						
(000111)	4.17		(000110)	4.15						
(000117)	4.19		(000114)	4.20						
(000118)	4.20	220								
(000119)	4.21	230								
(000134)	5.01	310	(000130)	4.04						
(000137)	5.04		(000143)	5.08						
(000138)	5.05		(000139)	5.06						
(000139)	5.06	314								
(000143)	5.08	320								
(000146)	5.10		(000153)	5.15						
(000148)	5.12		(000149)	5.13						
(000149)	5.13	317								
(000153)	5.15	315								
(000154)	5.16	330	(000132)	4.04						
(000159)	5.18	340	(000178)	4.02						
(000162)	5.20	350	(000157)	5.17						
(000172)	5.27	380	(000169)	5.24						
(000178)	6.01	490	(000084)	4.06						
(000198)	6.05	500	(000078)	3.18	(000181)	6.04				
(000201)	6.18		(000220)	6.34						
(000210)	6.24	575	(000208)	6.23						
(000220)	6.34	580	(000209)	6.23						

CHART TITLE - NON-PROCEDURAL STATEMENTS

PAGE 1

**AUTOLOG CHART -**

ENGL - SURROUTINE REM' ME


PAGE 1

PG. RX	NAME	PG. RX	NAME	PG. RX	NAME	PG. RX	NAME
1.01	NAME	2.21	60	4.01	300	5.18	340
2.07	5	3.15	70	5.01	310	5.20	350
1.06	12	2.01	100	5.06	314	5.27	380
2.04	14	3.07	200	5.15	315	6.01	400
3.14	15	4.11	205	5.13	317	3.32	450
2.09	25	4.20	220	5.09	320	4.05	500
3.12	30	4.21	230	5.14	330	6.24	575
2.24	55						





CHART TITLE - SUBROUTINE REFERENCE

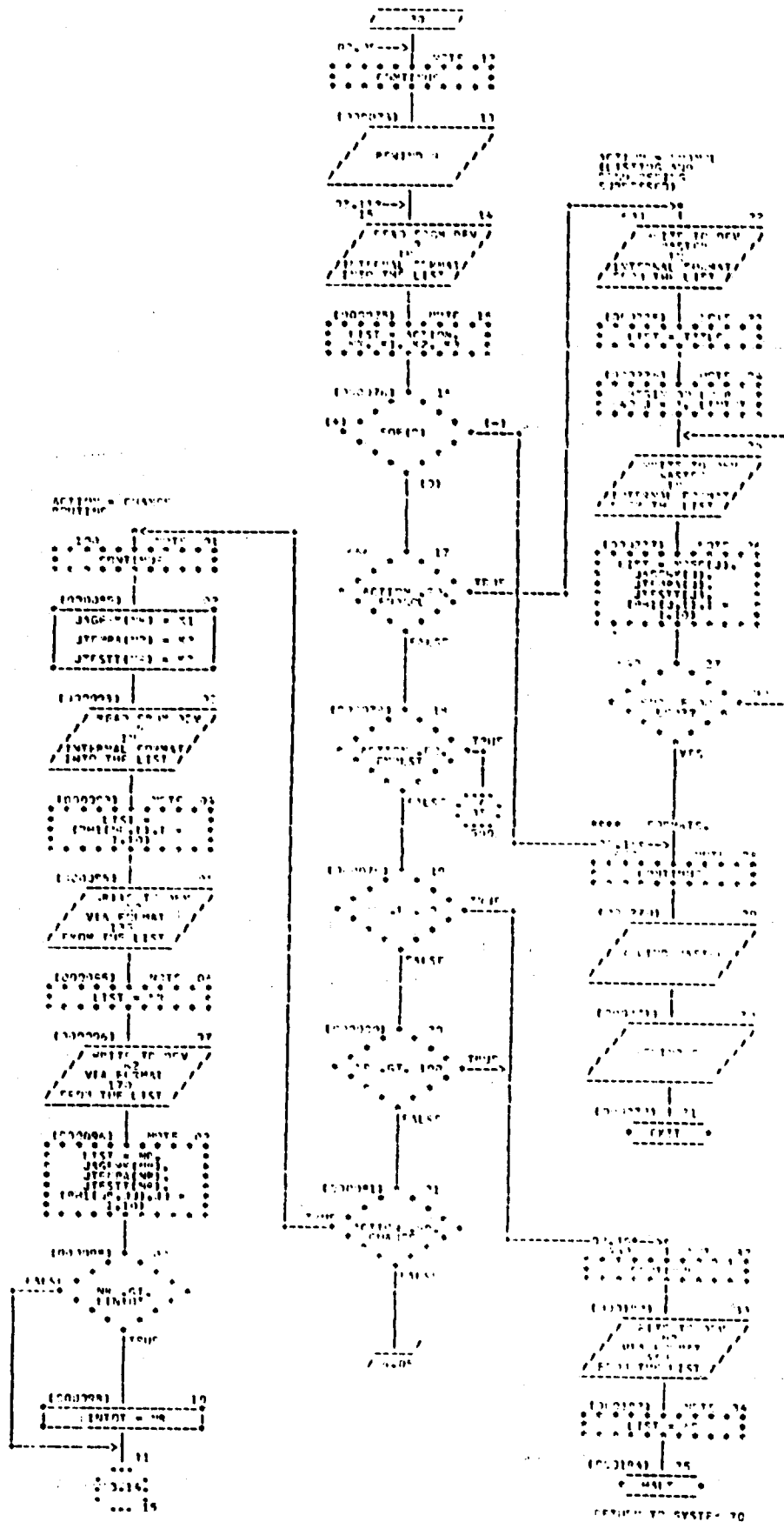



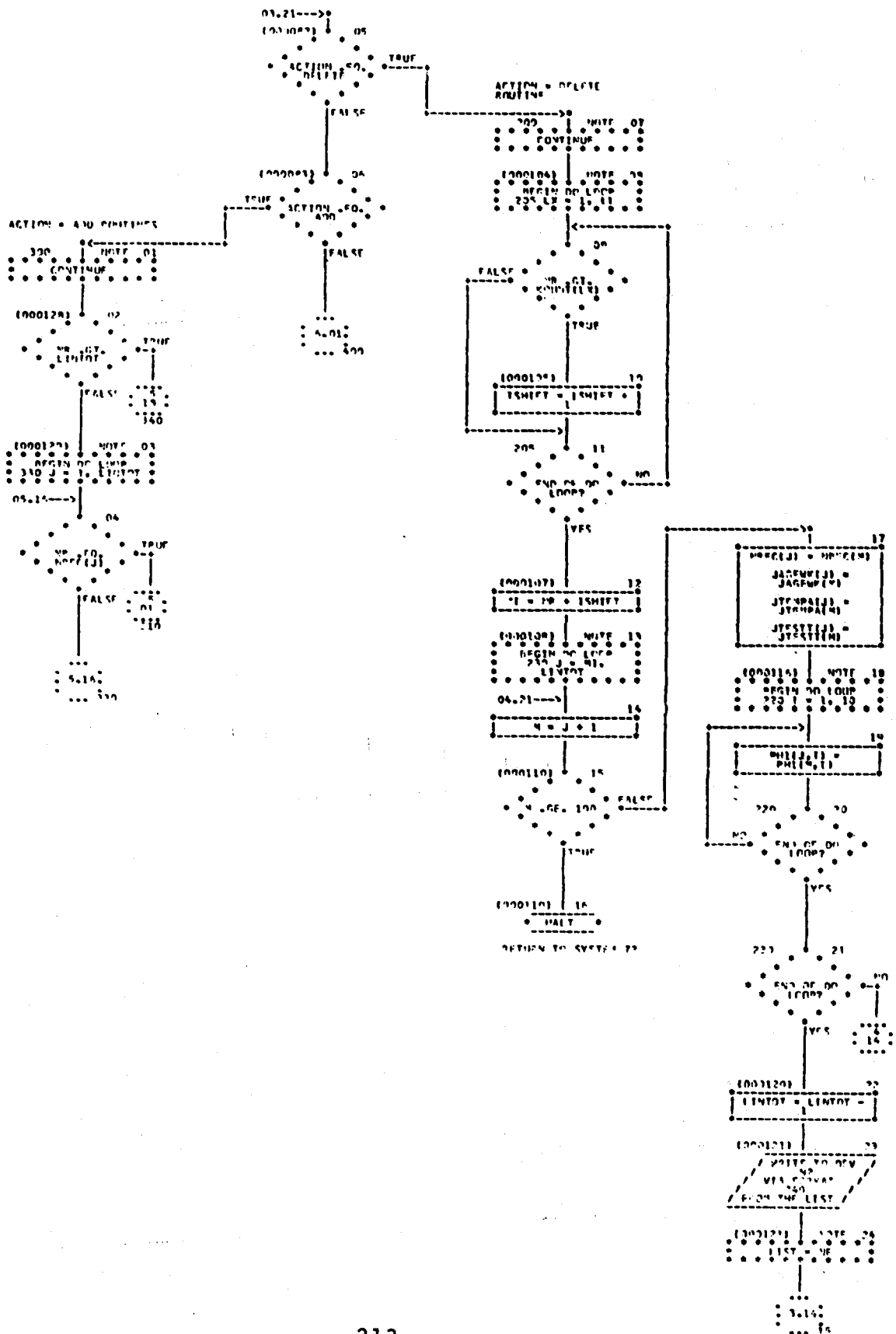


CHART TITLE - SUBROUTINE RENAME

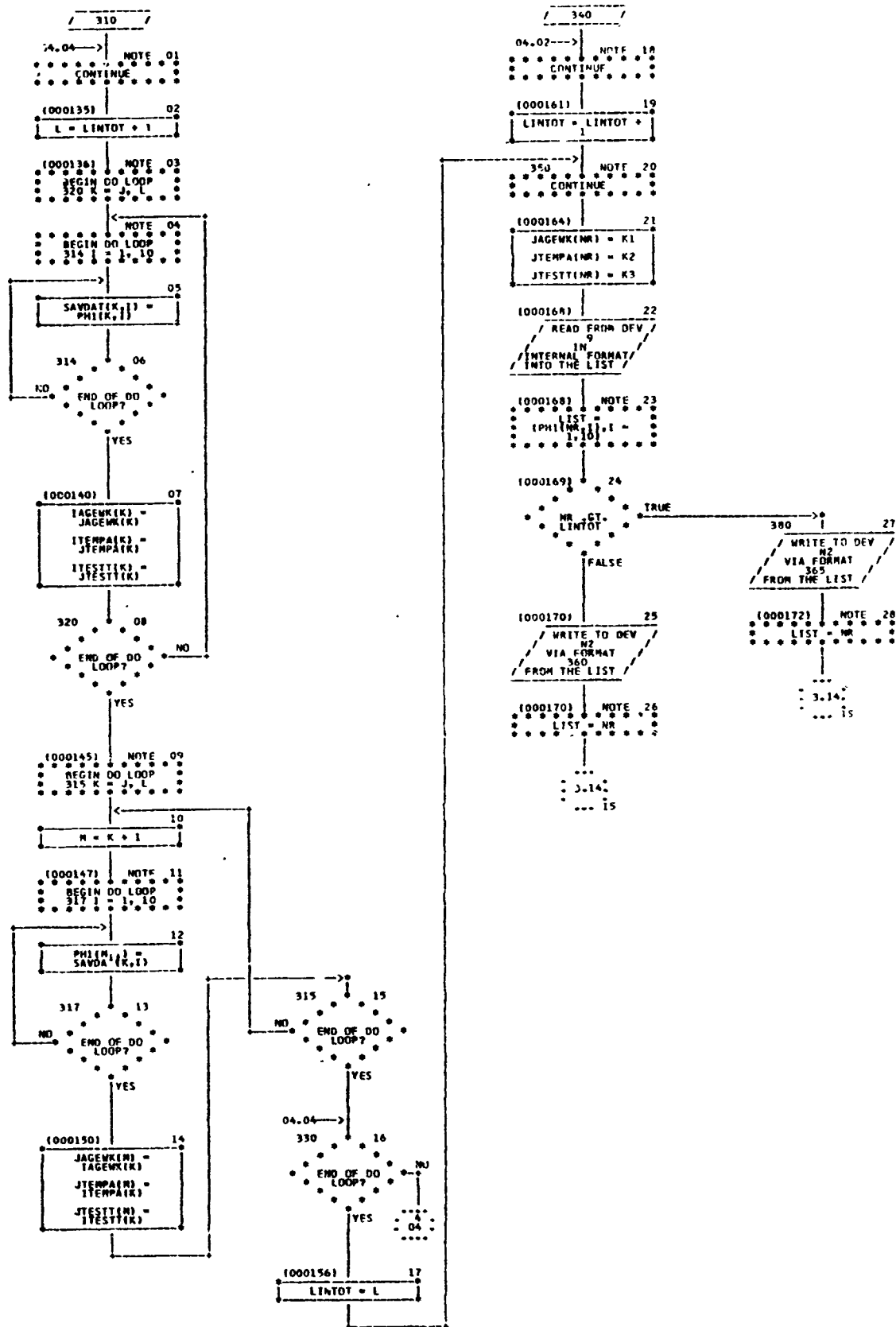


CHART TITLE - SUBROUTINE RENAME

ERROR MESSAGES FOR  
IMBAL TO UPDATE ACTION  
AND BAD LINE NUMBER

04.01--> NOTE 01  
CONTINUE

(000179) 02  
WRITE TO DEV  
VIA FORMAT  
FROM THE LIST

(000179) 03  
NOTE  
LIST - ACTION

(000180) 04  
WRITE TO DEV  
VIA FORMAT  
420

ACTION=ENDLIST (LINES  
RENUMBERED, DATA  
STORED AND PRINTED)

03.18--> 05  
WRITE TO DEV  
VIA FORMAT  
INTERNAL FORMAT  
FROM THE LIST

(000188) 06  
NOTE  
LIST - TITLE

(000190) 07  
WRITE TO DEV  
VIA FORMAT  
FROM THE LIST

(000190) 08  
NOTE  
LIST - TITLE

(000191) 09  
WRITE TO DEV  
VIA FORMAT  
570

(000192) 10  
WRITE TO DEV  
VIA FORMAT  
530

(000193) 11  
WRITE TO DEV  
VIA FORMAT  
540

(000194) 12  
WRITE TO DEV  
VIA FORMAT  
550

(000195) 13  
WRITE TO DEV  
VIA FORMAT  
560

(000196) 14  
WRITE TO DEV  
VIA FORMAT  
565

(000197) 15  
WRITE TO DEV  
VIA FORMAT  
560

(000199) 16  
ICOUNT = 11

(000200) 17  
NOTE  
IFC.N ON LOOP  
580 J = 1 LINTOT

06.34--> 18  
ICOUNT = ICOUNT +  
NREC(J) = J

(000203) 19  
WRITE TO DEV  
VIA FORMAT  
FROM THE LIST

(000203) 20  
NOTE  
LIST = NREC(J),  
JAGEWK(J),  
JTERPA(J),  
JTESTY(J),  
(PHI(J),1) =  
1,10

(000205) 21  
WRITE TO DEV  
VIA FORMAT  
INTERNAL FORMAT  
FROM THE LIST

(000205) 22  
NOTE  
LIST = NREC(J),  
JAGEWK(J),  
JTERPA(J),  
JTESTY(J),  
(PHI(J),1) =  
1,10

(000208) 23  
ICOUNT = 61

TRUE  
FALSE

575--> NOTE 24  
CONTINUE

(000211) 25  
WRITE TO DEV  
VIA FORMAT  
501

(000212) 26  
WRITE TO DEV  
VIA FORMAT  
570

(000213) 27  
WRITE TO DEV  
VIA FORMAT  
520

(000214) 28  
WRITE TO DEV  
VIA FORMAT  
550

(000215) 29  
WRITE TO DEV  
VIA FORMAT  
560

(000216) 30  
WRITE TO DEV  
VIA FORMAT  
565

(000217) 31  
WRITE TO DEV  
VIA FORMAT  
560

(000218) 32  
WRITE TO DEV  
VIA FORMAT  
560

(000219) 33  
ICOUNT = 7

END OF DO LOOP

YES

NO

19

20

295 18

[illegible]